

BeneView T5/ BeneView T8

Patient Monitor

Operator's Manual



© Copyright 2006-2012 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved.
For this Operator's Manual, the issue date is 2012-08.

 **WARNING**


- **Federal Law (USA) restricts this device to sale by or on the order of a physician.**
-

Intellectual Property Statement

SHENZHEN MINDRAY BIO-MEDICAL ELECTRONICS CO., LTD. (hereinafter called Mindray) owns the intellectual property rights to this Mindray product and this manual. This manual may refer to information protected by copyrights or patents and does not convey any license under the patent rights of Mindray, nor the rights of others.

Mindray intends to maintain the contents of this manual as confidential information. Disclosure of the information in this manual in any manner whatsoever without the written permission of Mindray is strictly forbidden.

Release, amendment, reproduction, distribution, rental, adaption and translation of this manual in any manner whatsoever without the written permission of Mindray is strictly forbidden.

mindray ,  **MINDRAY** and **BeneView** are the registered trademarks or trademarks owned by Mindray in China and other countries. All other trademarks that appear in this manual are used only for editorial purposes without the intention of improperly using them. They are the property of their respective owners.

Responsibility on the Manufacturer Party

Contents of this manual are subject to changes without prior notice.

All information contained in this manual is believed to be correct. Mindray shall not be liable for errors contained herein nor for incidental or consequential damages in connection with the furnishing, performance, or use of this manual.

Mindray is responsible for the effects on safety, reliability and performance of this product, only if:

- all installation operations, expansions, changes, modifications and repairs of this product are conducted by Mindray authorized personnel;
- the electrical installation of the relevant room complies with the applicable national and local requirements;
- the product is used in accordance with the instructions for use.



WARNING

- **This equipment must be operated by skilled/trained clinical professionals.**
 - **It is important for the hospital or organization that employs this equipment to carry out a reasonable service/maintenance plan. Neglect of this may result in machine breakdown or personal injury.**
-

Warranty

THIS WARRANTY IS EXCLUSIVE AND IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE.

Exemptions

Mindray's obligation or liability under this warranty does not include any transportation or other charges or liability for direct, indirect or consequential damages or delay resulting from the improper use or application of the product or the use of parts or accessories not approved by Mindray or repairs by people other than Mindray authorized personnel.

This warranty shall not extend to

- Malfunction or damage caused by improper use or man-made failure.
- Malfunction or damage caused by unstable or out-of-range power input.
- Malfunction or damage caused by force majeure such as fire and earthquake.
- Malfunction or damage caused by improper operation or repair by unqualified or unauthorized service people.
- Malfunction of the instrument or part whose serial number is not legible enough.
- Others not caused by instrument or part itself.

Company Contact

Manufacturer:	Shenzhen Mindray Bio-Medical Electronics Co., Ltd.
Address:	Mindray Building, Keji 12th Road South, High-tech industrial park, Nanshan, Shenzhen 518057, P.R. China
Website:	www.mindray.com
E-mail Address:	service@mindray.com
Tel:	+86 755 81888998
Fax:	+86 755 26582680
EC-Representative:	Shanghai International Holding Corp. GmbH(Europe)
Address:	Eiffestraße 80, 20537 Hamburg, Germany
Tel:	0049-40-2513175
Fax:	0049-40-255726

Preface

Manual Purpose

This manual contains the instructions necessary to operate the product safely and in accordance with its function and intended use. Observance of this manual is a prerequisite for proper product performance and correct operation and ensures patient and operator safety.

This manual is based on the maximum configuration and therefore some contents may not apply to your product. If you have any question, please contact us.

This manual is an integral part of the product. It should always be kept close to the equipment so that it can be obtained conveniently when needed.

Intended Audience

This manual is geared for clinical professionals who are expected to have a working knowledge of medical procedures, practices and terminology as required for monitoring of critically ill patients.

Illustrations

All illustrations in this manual serve as examples only. They may not necessarily reflect the setup or data displayed on your patient monitor.

Conventions

- ***Italic text*** is used in this manual to quote the referenced chapters or sections.
- [] is used to enclose screen texts.
- → is used to indicate operational procedures.

Contents

1 Safety	1-1
1.1 Safety Information.....	1-1
1.1.1 Dangers.....	1-1
1.1.2 Warnings.....	1-2
1.1.3 Cautions.....	1-2
1.1.4 Notes.....	1-3
1.2 Equipment Symbols.....	1-3
2 The Basics.....	2-1
2.1 Monitor Description.....	2-1
2.1.1 Intended Use.....	2-1
2.2 Main Unit.....	2-1
2.2.1 Front View	2-1
2.2.2 Side View	2-3
2.2.3 Rear View	2-4
2.2.4 Bottom View (BeneView T8)	2-6
2.3 Satellite Module Rack.....	2-6
2.4 Modules.....	2-7
2.4.1 Multi-Parameter Module	2-8
2.4.2 BeneView T1	2-9
2.5 Display Screen.....	2-10
2.6 QuickKeys	2-12
3 Basic Operations	3-1
3.1 Installation.....	3-1
3.1.1 Unpacking and Checking.....	3-1
3.1.2 Environmental Requirements	3-2
3.2 Getting Started	3-2
3.2.1 Turning Power On.....	3-2
3.2.2 Starting Monitoring.....	3-2
3.3 Disconnecting from Power.....	3-3
3.4 Using a Mouse.....	3-3
3.5 Using Keys	3-4
3.6 Using Keyboards	3-4
3.7 Using the Touchscreen.....	3-5
3.8 Setting the Screen	3-5
3.9 Using the Main Menu	3-6
3.10 Setting Parameters.....	3-6
3.10.1 Switching the Parameters On/Off.....	3-6
3.10.2 Accessing the Parameters Menu	3-7
3.10.3 Removing a Module Conflict	3-7
3.10.4 Removing a Label Conflict.....	3-8

3.11 Using a CF Storage Card.....	3-8
3.12 Changing General Settings.....	3-9
3.12.1 Setting up a Monitor.....	3-9
3.12.2 Changing Language.....	3-9
3.12.3 Adjusting the Screen Brightness.....	3-9
3.12.4 Showing/Hiding the Help.....	3-9
3.12.5 Setting the Date and Time.....	3-10
3.12.6 Adjusting Volume.....	3-10
4 Managing Patients.....	4-1
4.1 Admitting a Patient.....	4-1
4.2 Quick Admitting a Patient.....	4-1
4.3 Querying and Obtaining Patient Information.....	4-2
4.4 Associating Patient Information.....	4-2
4.5 Editing Patient Information.....	4-2
4.6 Discharging a Patient.....	4-3
4.7 Transferring a Patient.....	4-3
4.7.1 Transferring Patient Data via MPM/T1.....	4-4
4.7.2 Transferring Patient Data via Storage Medium.....	4-5
4.8 Connecting to a Central Monitoring System.....	4-6
5 Managing Configurations.....	5-1
5.1 Introduction.....	5-1
5.2 Entering the [Manage Configuration] Menu.....	5-2
5.3 Changing Department.....	5-2
5.4 Setting Default Configuration.....	5-3
5.5 Saving Current Settings.....	5-3
5.6 Editing Configuration.....	5-4
5.7 Deleting a Configuration.....	5-5
5.8 Transferring a Configuration.....	5-5
5.9 Loading a Configuration.....	5-5
5.10 Restoring the Latest Configuration Automatically.....	5-6
5.11 Modifying Password.....	5-6
6 User Screens.....	6-1
6.1 Tailoring Your Screens.....	6-1
6.1.1 Setting the Waveform Sweep Mode.....	6-1
6.1.2 Changing the Wave Line Size.....	6-1
6.1.3 Changing Measurement Colors.....	6-1
6.1.4 Changing Screen Layout.....	6-1
6.2 Viewing Minitrends.....	6-2
6.2.1 Having a Split-Screen View of Minitrends.....	6-2
6.2.2 Setting Minitrends.....	6-2
6.3 Viewing OxyCRG.....	6-3
6.4 Viewing Other Patients.....	6-4
6.4.1 Care Group.....	6-4

6.4.2 Viewing the Care Group Overview Bar	6-4
6.4.3 Understanding the View Other Patient Window	6-5
6.5 Understanding the Big Numerics Screen	6-6
7 Alarms	7-1
7.1 Alarm Categories.....	7-1
7.2 Alarm Levels.....	7-2
7.3 Alarm Indicators	7-2
7.3.1 Alarm Lamp	7-2
7.3.2 Alarm Message.....	7-3
7.3.3 Flashing Numeric.....	7-3
7.3.4 Audible Alarm Tones	7-3
7.3.5 Alarm Status Symbols.....	7-4
7.4 Alarm Tone Configuration	7-4
7.4.1 Setting the Minimum Alarm Volume	7-4
7.4.2 Changing the Alarm Volume.....	7-4
7.4.3 Setting the Interval between Alarm Sounds.....	7-5
7.4.4 Changing the Alarm Tone Pattern.....	7-5
7.4.5 Setting the Reminder Tones	7-6
7.5 Understanding the Alarm Setup Menu	7-6
7.5.1 Setting Alarm Properties for All Parameters.....	7-7
7.5.2 Adjusting Alarm Limits Automatically	7-7
7.5.3 Setting Alarm Delay Time.....	7-11
7.5.4 Setting Recording Length	7-11
7.5.5 Entering CPB Mode.....	7-11
7.6 Pausing Alarms	7-12
7.7 Switching Off All Alarms	7-13
7.8 Silencing the Alarm Sound.....	7-13
7.9 Latching Alarms.....	7-14
7.10 Silencing Technical Alarms.....	7-14
7.11 Testing Alarms	7-15
7.12 When an Alarm Occurs	7-15
7.13 Using Care Group Alarms	7-15
7.13.1 Care Group Auto Alarms	7-15
7.13.2 Setting Care Group Alert Tone	7-16
7.13.3 Silencing Care Group Alarms.....	7-16
8 Monitoring ECG	8-1
8.1 Introduction.....	8-1
8.2 Safety.....	8-1
8.3 Preparing to Monitor ECG.....	8-2
8.3.1 Preparing the Patient and Placing the Electrodes.....	8-2
8.3.2 Choosing AHA or IEC Lead Placement.....	8-2
8.3.3 ECG Lead Placements	8-2
8.3.4 Checking Paced Status	8-4
8.4 Understanding the ECG Display	8-5

8.5 Changing ECG Settings	8-6
8.5.1 Accessing ECG Menus	8-6
8.5.2 Setting Pacemaker Rate (For Mortara only)	8-6
8.5.3 Choosing the Alarm Source	8-6
8.5.4 Setting the ECG Lead Set	8-6
8.5.5 Choosing an ECG Display Screen	8-6
8.5.6 Changing the ECG Filter Settings	8-7
8.5.7 Setting the Notch Filter	8-7
8.5.8 Changing the Pacer Reject Settings	8-7
8.5.9 About the Defibrillator Synchronization	8-8
8.5.10 Changing ECG Wave Settings	8-8
8.5.11 Enabling Smart Lead Off	8-8
8.5.12 Setting the Alarm Level for ECG Lead Off Alarms	8-8
8.5.13 Adjusting QRS Volume	8-9
8.6 About ST Monitoring	8-9
8.6.1 Switching ST On and Off	8-9
8.6.2 Changing ST Filter Settings	8-9
8.6.3 Understanding the ST Display	8-10
8.6.4 Saving the Current ST Segment as Reference	8-10
8.6.5 Changing the Reference Segment	8-11
8.6.6 Deleting a Reference Segment	8-11
8.6.7 Recording the ST Segment	8-11
8.6.8 Changing the ST Alarm Limits	8-11
8.6.9 Setting the ST Alarm Delay Time	8-11
8.6.10 Adjusting ST Measurement Points	8-11
8.7 About Arrhythmia Monitoring	8-12
8.7.1 Understanding the Arrhythmia Events	8-12
8.7.2 Changing Arrhythmia Alarm Settings	8-14
8.7.3 Changing Arrhythmia Threshold Settings	8-14
8.7.4 Setting the Extended Arrh. (For Mindray Algorithm Only)	8-15
8.7.5 Reviewing Arrhythmia Events	8-15
8.8 ECG Relearning	8-16
8.8.1 Initiating an ECG Relearning Manually	8-16
8.8.2 Automatic ECG Relearning	8-16
8.9 12-Lead ECG Monitoring	8-17
8.9.1 Entering the 12-lead ECG Monitoring Screen	8-17
8.9.2 Interpretation of resting 12-lead ECG	8-18
8.9.3 Reviewing Interpretation of resting 12-lead ECG Results	8-19
9 Monitoring Respiration (Resp)	9-1
9.1 Introduction	9-1
9.2 Safety Information	9-1
9.3 Understanding the Resp Display	9-1
9.4 Placing Resp Electrodes	9-2
9.4.1 Optimizing Lead Placement for Resp	9-2
9.4.2 Cardiac Overlay	9-2

9.4.3 Abdominal Breathing	9-3
9.4.4 Lateral Chest Expansion	9-3
9.5 Choosing the Respiration Lead	9-3
9.6 Changing the Apnea Alarm Delay	9-3
9.7 Changing Resp Detection Mode	9-3
9.8 Changing Resp Wave Settings	9-4
9.9 Setting RR Source	9-4
9.10 Setting alarm properties	9-5
10 Monitoring PR	10-1
10.1 Introduction	10-1
10.2 Setting the PR Source	10-1
10.3 Selecting the Active Alarm Source	10-2
10.4 QRS Tone	10-2
11 Monitoring SpO₂	11-1
11.1 Introduction	11-1
11.2 Safety	11-2
11.3 Identifying SpO ₂ Connectors	11-2
11.4 Applying the Sensor	11-2
11.5 Changing SpO ₂ Settings	11-3
11.5.1 Accessing SpO ₂ Menus	11-3
11.5.2 Adjusting the Desat Alarm	11-3
11.5.3 Setting SpO ₂ Sensitivity	11-3
11.5.4 Changing Averaging Time	11-3
11.5.5 Monitoring SpO ₂ and NIBP Simultaneously	11-3
11.5.6 Sat-Seconds Alarm Management	11-4
11.5.7 Changing the Speed of the Pleth/Plethb Wave	11-5
11.5.8 Setting the Alarm Level for SpO ₂ Sensor Off Alarm	11-5
11.5.9 Setting the SpO ₂ Tone Mode	11-5
11.6 Measurement Limitations	11-5
11.7 Masimo Information	11-6
11.8 Nellcor Information	11-6
12 Monitoring NIBP	12-1
12.1 Introduction	12-1
12.2 Safety	12-2
12.3 Measurement Limitations	12-2
12.4 Measurement Methods	12-2
12.5 Setting Up the NIBP Measurement	12-3
12.5.1 Preparing to Measure NIBP	12-3
12.5.2 Starting and Stopping Measurements	12-3
12.5.3 Correcting the Measurement if Limb is not at Heart Level	12-3
12.5.4 Enabling NIBP Auto Cycling and Setting the Interval	12-3
12.5.5 Starting a STAT Measurement	12-4
12.6 Understanding the NIBP Numerics	12-4

12.7 Changing NIBP Settings	12-4
12.7.1 Setting the Initial Cuff Inflation Pressure.....	12-4
12.7.2 Setting NIBP Alarm Properties	12-5
12.7.3 Switching On NIBP End Tone	12-5
12.7.4 Displaying NIBP List.....	12-5
12.7.5 Setting the Pressure Unit	12-5
12.8 Assisting Venous Puncture	12-5
13 Monitoring Temp	13-1
13.1 Introduction.....	13-1
13.2 Safety	13-1
13.3 Making a Temp Measurement	13-1
13.4 Understanding the Temp Display	13-1
13.5 Setting the Temperature Unit	13-2
14 Monitoring IBP.....	14-1
14.1 Introduction.....	14-1
14.2 Safety	14-1
14.3 Setting Up the Pressure Measurement	14-2
14.4 Understanding the IBP Display.....	14-3
14.5 Changing IBP Settings	14-3
14.5.1 Changing a Pressure for Monitoring.....	14-3
14.5.2 Setting Alarm Properties	14-4
14.5.3 Changing Averaging Time	14-4
14.5.4 Setting the Pressure Unit	14-4
14.5.5 Setting Up the IBP Wave.....	14-4
14.6 Measuring PAWP	14-5
14.6.1 Preparing to Measure PAWP.....	14-5
14.6.2 Setting Up the PAWP Measurement.....	14-6
14.6.3 Understanding the PAWP Setup Menu	14-6
14.7 Zeroing the Transducer	14-7
15 Monitoring Cardiac Output	15-1
15.1 Introduction.....	15-1
15.2 Understanding the C.O. Display	15-2
15.3 Influencing Factors.....	15-2
15.4 Setting Up the C.O. Measurement.....	15-3
15.5 Measuring the Blood Temperature	15-6
15.6 Changing C.O. Settings.....	15-6
15.6.1 Setting the Temperature Unit.....	15-6
15.6.2 Setting Alarm Properties.....	15-6
16 Monitoring CCO/SvO₂.....	16-1
16.1 Introduction.....	16-1
16.2 Safety	16-1
16.3 Automatic Communication Detection	16-2

16.4 Connecting the Device.....	16-2
16.4.1 Connecting the Vigilance II® Monitor	16-2
16.4.2 Connecting the Vigileo™ Monitor.....	16-3
16.5 Understanding CCO Parameters.....	16-5
16.6 Understanding the CCO Display	16-6
16.7 Changing CCO Settings.....	16-6
16.7.1 Selecting Vascular Resistance Unit	16-6
16.7.2 Selecting the Displayed Parameters	16-6
16.7.3 Checking the C.O. Measurements.....	16-6
16.7.4 Setting Signal Output.....	16-7
16.7.5 Selecting Alarm Properties.....	16-7
16.8 Understanding SvO ₂ Parameters	16-7
16.9 Understanding the SvO ₂ Display	16-8
16.10 Changing SvO ₂ Settings	16-8
16.10.1 Setting Signal Output	16-8
16.10.2 Selecting Alarm Properties	16-8
17 Monitoring PiCCO	17-1
17.1 Introduction	17-1
17.2 Safety Information	17-1
17.3 Zeroing the Transducer	17-2
17.4 Setting up the PiCCO Measurements.....	17-3
17.5 Preparation for PiCCO Measurements	17-3
17.6 Performing PiCCO Measurements and CCO Calibration	17-6
17.7 Understanding the Displayed PiCCO Parameters.....	17-8
17.7.1 Understanding the CCO Display.....	17-8
17.7.2 Understanding the pArt Display	17-8
17.7.3 Understanding the pCVP Display	17-9
17.8 Understanding PiCCO Parameters	17-9
17.8.1 Spider Vision	17-9
17.8.2 Hemodynamic Parameters.....	17-11
17.8.3 Normal Range Setup	17-12
17.9 Changing PiCCO Settings.....	17-12
17.9.1 Selecting the Displayed Parameters	17-12
17.9.2 Selecting Alarm Properties.....	17-12
18 Monitoring ScvO₂.....	18-1
18.1 Introduction	18-1
18.2 Safety Information	18-1
18.3 Performing ScvO ₂ Measurements.....	18-2
18.4 ScvO ₂ Calibration.....	18-3
18.5 Understanding the ScvO ₂ Display	18-4
18.6 Understanding ScvO ₂ Parameters.....	18-4
18.7 Changing ScvO ₂ Settings.....	18-4
18.7.1 Selecting Hb/Hct.....	18-4
18.7.2 Selecting Alarm Properties.....	18-4

19 Monitoring Carbon Dioxide	19-1
19.1 Introduction	19-1
19.2 Identifying CO ₂ Modules	19-2
19.3 Preparing to Measure CO ₂	19-2
19.3.1 Using a Sidestream CO ₂ Module	19-2
19.3.2 Using a Microstream CO ₂ Module	19-3
19.3.3 Using a Mainstream CO ₂ Module	19-4
19.4 Changing CO ₂ Settings	19-4
19.4.1 Accessing CO ₂ Menus	19-4
19.4.2 Entering the Standby Mode	19-4
19.4.3 Setting the CO ₂ Unit	19-5
19.4.4 Setting up Gas Compensations	19-5
19.4.5 Setting up Humidity Compensation	19-6
19.4.6 Setting the Apnea Alarm Delay	19-6
19.4.7 Choosing a Time Interval for Peak-Picking	19-6
19.4.8 Setting the Flow Rate	19-6
19.4.9 Setting up the CO ₂ Wave	19-7
19.5 Setting RR Source	19-7
19.6 Setting Barometric Pressure Compensation	19-7
19.7 Measurement Limitations	19-8
19.8 Leakage test	19-8
19.9 Troubleshooting the Sidestream CO ₂ Sampling System	19-8
19.10 Removing Exhaust Gases from the System	19-8
19.11 Zeroing the Sensor	19-9
19.11.1 For Sidestream and Microstream CO ₂ Modules	19-9
19.11.2 For Mainstream CO ₂ Modules	19-9
19.12 Calibrating the Sensor	19-9
19.13 Oridion Information	19-10
20 Monitoring tcGas	20-1
20.1 Introduction	20-1
20.2 Safety	20-1
20.3 Connecting a TCM monitor	20-2
20.4 tcGas Parameters	20-2
20.5 Displaying tcGas Parameters	20-3
20.6 Enter the tcGas Setup menu	20-3
20.7 Setting tcpCO ₂ /tcpO ₂ Unit	20-4
20.8 tcGas Display	20-4
21 Monitoring AG	21-1
21.1 Introduction	21-1
21.2 Identifying AG Modules	21-1
21.3 Understanding the AG Display	21-2
21.4 MAC Values	21-2
21.5 Preparing to Measure AG	21-3
21.6 Changing AG Settings	21-4

21.6.1 Setting Gas Unit.....	21-4
21.6.2 Setting the Apnea Alarm Delay	21-4
21.6.3 Changing the Sample Flow Rate.....	21-4
21.6.4 Setting up the O ₂ Compensation.....	21-5
21.6.5 Entering the Standby Mode.....	21-5
21.6.6 Setting up the AG Wave.....	21-5
21.6.7 Setting RR Source	21-5
21.7 Changing the Anesthetic Agent	21-6
21.8 Measurement Limitations	21-6
21.9 Troubleshooting.....	21-6
21.9.1 When the Gas Inlet is Blocked	21-6
21.9.2 When an Internal Occlusion Occurs.....	21-6
21.10 Removing Exhaust Gases from the System	21-7
22 Monitoring ICG.....	22-1
22.1 Introduction	22-1
22.2 Safety	22-1
22.3 Understanding ICG Parameters.....	22-2
22.3.1 Measured Parameters.....	22-2
22.3.2 Calculated Parameters	22-2
22.4 Understanding the ICG Display.....	22-3
22.5 ICG Limitations	22-3
22.6 Preparing to Monitor ICG.....	22-3
22.6.1 Preparing the Patient.....	22-4
22.6.2 Placing ICG Sensors.....	22-4
22.6.3 Setting up the Patient Information	22-4
22.7 Changing ICG Settings.....	22-5
22.7.1 ICG Averaging	22-5
22.7.2 Selecting Secondary Parameters	22-5
22.7.3 Checking Sensors.....	22-5
22.7.4 Changing the ICG Wave Speed	22-5
23 Monitoring BIS.....	23-1
23.1 Introduction	23-1
23.2 Safety Information	23-2
23.3 Understanding the BIS Display.....	23-2
23.3.1 BIS Parameter Area.....	23-2
23.3.2 BIS Waveform Area	23-5
23.3.3 BIS Expand View	23-5
23.4 Setting up the BIS Measurement.....	23-7
23.5 Auto Impedance Check.....	23-7
23.6 Sensor Check.....	23-8
23.7 BIS Sensor Check Window	23-8
23.8 Choosing the BIS Smoothing Rate	23-9
23.9 Changing the Secondary Parameters	23-10
23.10 Changing the EEG Wave Size	23-10

23.11 Changing the Speed of the EEG Wave.....	23-10
23.12 Setting the Trend Length	23-10
23.13 Switching the Filter On or Off.....	23-10
24 Monitoring NMT	24-1
24.1 Introduction.....	24-1
24.2 Safty.....	24-1
24.3 Connecting a TOF-Watch® SX monitor.....	24-1
24.4 NMT Parameters.....	24-2
24.5 Accessing the NMT Setup menu	24-2
24.6 NMT Display.....	24-3
25 Monitoring RM.....	25-1
25.1 Introduction.....	25-1
25.2 Safety Information.....	25-2
25.3 Preparing to Monitor RM	25-3
25.4 Understanding the RM Display.....	25-4
25.5 Changing RM Settings	25-4
25.5.1 Accessing RM Menus	25-4
25.5.2 Setting the Apnea Alarm Delay	25-4
25.5.3 Selecting TV or MV for Display	25-5
25.5.4 Selecting Flow or Vol Waveform for Display.....	25-5
25.5.5 Changing the Wave Sweep Speed	25-5
25.5.6 Changing the Wave Scale.....	25-5
25.5.7 Setting RR Source.....	25-5
25.6 Understanding the Respiratory Loops.....	25-6
25.7 Zeroing the RM Module	25-7
25.8 Calibrating the Flow Sensor.....	25-7
26 Device Integration.....	26-1
26.1 Introduction.....	26-1
26.2 Safety Information.....	26-1
26.3 Supported Devices.....	26-2
26.4 Differences in Displayed Values.....	26-2
26.5 Connecting an External Device	26-3
26.6 Devices Integrated Window.....	26-5
26.7 System Functions of Patient Monitor	26-6
26.7.1 Alarms	26-6
26.7.2 Data Storage	26-6
26.7.3 Recording and Printing.....	26-6
26.8 Integrating the Anesthesia Machine	26-7
26.8.1 Wato 20/30/55/65	26-7
26.8.2 Maquet Flow-i	26-14
26.8.3 Draeger Fabius GS/Fabius Trio/Fabius Plus.....	26-17
26.8.4 Draeger Primus	26-21
26.8.5 GE Datex-Ohmeda Aestiva 7900/Aestiva 7100	26-27

26.8.6 GE Datex-Ohmeda Avance/Aisys	26-30
26.9 Integrating Ventilator.....	26-37
26.9.1 Newport E360	26-37
26.9.2 Puritan Bennett 840	26-41
26.9.3 Maquet Servo-i/Servo-s.....	26-44
26.9.4 Draeger Evita 2.....	26-47
26.9.5 Draeger Evita 4/ Evita2 dura /Evita XL.....	26-50
26.9.6 Hamilton G5.....	26-54
26.9.7 Hamilton C2 /Galileo	26-60
27 Freezing Waveforms	27-1
27.1 Freezing Waveforms	27-1
27.2 Viewing Frozen Waveforms.....	27-1
27.3 Unfreezing Waveforms	27-2
27.4 Recording Frozen Waveforms	27-2
28 Review.....	28-1
28.1 Accessing Respective Review Windows	28-1
28.2 Reviewing Graphic Trends.....	28-1
28.3 Reviewing Tabular Trends.....	28-2
28.4 Reviewing Events	28-3
28.5 Reviewing Waveforms	28-5
29 Calculations	29-1
29.1 Introduction	29-1
29.2 Dose Calculations.....	29-2
29.2.1 Performing Calculations.....	29-2
29.2.2 Selecting the Proper Drug Unit	29-2
29.2.3 Titration Table.....	29-3
29.3 Oxygenation Calculations	29-3
29.3.1 Performing Calculations.....	29-3
29.3.2 Entered Parameters.....	29-4
29.3.3 Calculated Parameters	29-4
29.4 Ventilation Calculations	29-5
29.4.1 Performing Calculations.....	29-5
29.4.2 Entered Parameters.....	29-5
29.4.3 Calculated Parameters	29-5
29.5 Hemodynamic Calculations.....	29-6
29.5.1 Performing Calculations.....	29-6
29.5.2 Entered Parameters.....	29-6
29.5.3 Calculated Parameters	29-7
29.6 Renal Calculations	29-7
29.6.1 Performing Calculations.....	29-7
29.6.2 Entered Parameters.....	29-8
29.6.3 Calculated Parameters	29-8
29.7 Understanding the Review Window.....	29-8

30 Recording	30-1
30.1 Using a Recorder	30-1
30.2 Overview of Recording Types	30-1
30.3 Starting and Stopping Recordings.....	30-2
30.4 Setting up the Recorder	30-3
30.4.1 Accessing the Record Setup Menu	30-3
30.4.2 Selecting Waveforms for Recording	30-3
30.4.3 Setting the Realtime Recording Length	30-3
30.4.4 Setting the Interval between Timed Recordings.....	30-3
30.4.5 Changing the Recording Speed	30-3
30.4.6 Clearing Recording Tasks	30-3
30.5 Loading Paper.....	30-4
30.6 Removing Paper Jam.....	30-4
30.7 Cleaning the Recorder Printhead	30-5
31 Printing	31-1
31.1 Printer	31-1
31.2 Connecting a printer	31-1
31.3 Setting Up the Printer	31-2
31.4 Starting Report Printouts.....	31-2
31.5 Stopping Reports Printouts	31-2
31.6 Setting Up Reports.....	31-3
31.6.1 Setting Up ECG Reports.....	31-3
31.6.2 Setting Up Tabular Trends Reports	31-3
31.6.3 Setting Up Graphic Trends Reports	31-3
31.6.4 Setting Up Realtime Reports	31-4
31.7 End Case Reports.....	31-4
31.8 Printer Statuses.....	31-4
31.8.1 Printer Out of Paper.....	31-4
31.8.2 Printer Status Messages.....	31-4
32 Other Functions	32-1
32.1 Marking Events	32-1
32.2 Privacy Mode.....	32-1
32.3 Night Mode.....	32-2
32.4 Analog Output.....	32-2
32.5 Transferring Data	32-3
32.5.1 Data Export System.....	32-3
32.5.2 Transferring Data by Different Means.....	32-3
32.6 Nurse Call.....	32-4
32.7 Remote Display	32-5
32.8 Built-in PC (only applicable to BeneView T8 patient monitor).....	32-5
32.8.1 Configuring Application Program ShortCuts	32-5
32.8.2 Using PC Software	32-7
32.8.3 Switching over Users	32-7
32.8.4 Remote Login	32-7

32.8.5 Using McAfee Solidifier.....	32-8
32.9 Wireless Network.....	32-8
32.10 Using DVI-VGA Adapter Box.....	32-9
33 Batteries	33-1
33.1 Overview	33-1
33.2 Installing or Replacing a Battery.....	33-2
33.3 Conditioning a Battery	33-3
33.4 Checking a Battery.....	33-4
33.5 Recycling a Battery	33-4
34 Care and Cleaning.....	34-1
34.1 General Points.....	34-1
34.2 Cleaning.....	34-2
34.3 Disinfecting	34-2
35 Maintenance.....	35-1
35.1 Safety Checks.....	35-1
35.2 Maintenance and Testing Schedule.....	35-2
35.3 Checking Monitor and Module Information	35-4
35.4 Calibrating ECG	35-4
35.5 NIBP Leakage Test.....	35-4
35.6 NIBP Accuracy Test.....	35-5
35.7 Calibrating NIBP	35-6
35.8 Calibrating CO ₂	35-6
35.9 Calibrating AG.....	35-8
35.10 Calibrating the Touchscreen	35-9
35.11 Electrical Safty Tests	35-9
35.12 Setting up IP Address	35-9
35.13 Entering/Exiting Demo Mode	35-10
36 Accessories	36-1
36.1 ECG Accessories	36-1
36.2 SpO ₂ Accessories.....	36-3
36.3 NIBP Accessories	36-4
36.4 Temp Accessories	36-5
36.5 IBP/ICP Accessories.....	36-6
36.6 C.O. Accessories.....	36-7
36.7 CCO/SvO ₂ Accessories.....	36-7
36.8 CO ₂ Accessories.....	36-7
36.9 AG Accessories.....	36-8
36.10 ICG Accessories	36-9
36.11 BIS Accessories.....	36-9
36.12 RM Accessories.....	36-9
36.13 PiCCO Accessories.....	36-9
36.14 ScvO ₂ Accessories.....	36-9

36.15 BeneLink Accessories.....	36-10
36.16 Others.....	36-10
A Product Specifications	A-1
A.1 Monitor Safety Specifications	A-1
A.2 Physical Specifications.....	A-3
A.3 Hardware Specifications	A-4
A.4 Data Storage.....	A-7
A.5 Wireless Network.....	A-7
A.6 Measurement Specifications.....	A-8
B EMC.....	B-1
C Default Configurations	C-1
C.1 Parameters Configuration	C-1
C.2 Routine Configuration	C-16
C.3 User Maintenance Items	C-20
D Alarm Messages	D-1
D.1 Physiological Alarm Messages.....	D-1
D.2 Technical Alarm Messages.....	D-2
E Electrical Safety Inspection.....	E-1
E.1 Power Cord Plug	E-1
E.2 Device Enclosure and Accessories	E-2
E.3 Device Labelling	E-2
E.4 Protective Earth Resistance	E-2
E.5 Earth Leakage Test	E-3
E.6 Patient Leakage Current	E-3
E.7 Mains on Applied Part Leakage	E-4
E.8 Patient Auxiliary Current	E-4
F Symbols and Abbreviations.....	F-1
F.1 Symbols.....	F-1
F.2 Abbreviations.....	F-2

1 Safety

1.1 Safety Information

DANGER

- Indicates an imminent hazard that, if not avoided, will result in death or serious injury.
-
-

WARNING

- Indicates a potential hazard or unsafe practice that, if not avoided, could result in death or serious injury.
-
-

CAUTION

- Indicates a potential hazard or unsafe practice that, if not avoided, could result in minor personal injury or product/property damage.
-
-

NOTE

- Provides application tips or other useful information to ensure that you get the most from your product.
-
-

1.1.1 Dangers

There are no dangers that refer to the product in general. Specific “Danger” statements may be given in the respective sections of this manual.

1.1.2 Warnings



WARNINGS

- Before putting the system into operation, the operator must verify that the equipment, connecting cables and accessories are in correct working order and operating condition.
 - The equipment must be connected to a properly installed power outlet with protective earth contacts only. If the installation does not provide for a protective earth conductor, disconnect it from the power line and operate it on battery power, if possible.
 - To avoid explosion hazard, do not use the equipment in the presence of oxygen-rich atmospheres, flammable anesthetics, or other flammable agents (such as gasoline).
 - Do not open the equipment housings. All servicing and future upgrades must be carried out by the personnel trained and authorized by our company only.
 - Do not come into contact with patients during defibrillation. Otherwise serious injury or death could result.
 - Do not rely exclusively on the audible alarm system for patient monitoring. Adjustment of alarm volume to a low level or off may result in a hazard to the patient. Remember that alarm settings should be customized according to different patient situations and always keeping the patient under close surveillance is the most reliable way for safe patient monitoring.
 - The physiological data and alarm messages displayed on the equipment are for reference only and cannot be directly used for diagnostic interpretation.
 - To avoid inadvertent disconnection, route all cables in a way to prevent a stumbling hazard. Wrap and secure excess cabling to reduce risk of entanglement or strangulation by patients or personnel.
 - Dispose of the package material, observing the applicable waste control regulations and keeping it out of children's reach.
-

1.1.3 Cautions



CAUTIONS

- To ensure patient safety, use only parts and accessories specified in this manual.
 - At the end of its service life, the equipment, as well as its accessories, must be disposed of in compliance with the guidelines regulating the disposal of such products. If you have any questions concerning disposal of the equipment, please contact us.
 - Magnetic and electrical fields are capable of interfering with the proper performance of the equipment. For this reason make sure that all external devices operated in the vicinity of the equipment comply with the relevant EMC requirements. Mobile phone, X-ray equipment or MRI devices are a possible source of interference as they may emit higher levels of electromagnetic radiation.
 - Before connecting the equipment to the power line, check that the voltage and frequency ratings of the power line are the same as those indicated on the equipment's label or in this manual.
 - Always install or carry the equipment properly to avoid damage caused by drop, impact, strong vibration or other mechanical force.
-

1.1.4 Notes



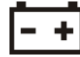







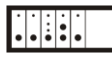














NOTES







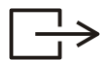

- Put the equipment in a location where you can easily see the screen and access the operating controls.
- Keep this manual in the vicinity of the equipment so that it can be obtained conveniently when needed.
- The software was developed in compliance with IEC60601-1-4. The possibility of hazards arising from software errors is minimized.
- This manual describes all features and options. Your equipment may not have all of them.

1.2 Equipment Symbols

NOTE

- Some symbols may not appear on your equipment.

	Attention: Consult accompanying documents (this manual).		
	Power ON/OFF (for a part of the equipment)		Battery indicator
	Alternating current (AC)		Alarms paused
	Alarm silenced		Record
	Freeze/unfreeze waveforms		Main menu
	NIBP start/stop key		Connector for satellite module rack
	Equipotential grounding		Video output
	USB connector		Network connector
	CIS connector		Auxiliary output
	Defibrillator		Zero key
	Check sensor		Calibrate key
	Measure/standby		Inserted direction
	Gas outlet		Serial number

	CIS connector		Manufacture date
	<p>The product bears CE mark indicating its conformity with the provisions of the Council Directive 93/42/EEC concerning medical devices and fulfils the essential requirements of Annex I of this directive.</p>		
	European community representative		
	ESD warning symbol for electrostatic sensitive devices.		
	Type CF applied part. Defibrillator-proof protection against electric shock.		
	Type BF applied part. Defibrillator-proof protection against electric shock.		
	<p>The following definition of the WEEE label applies to EU member states only.</p> <p>This symbol indicates that this product should not be treated as household waste. By ensuring that this product is disposed of correctly, you will help prevent bringing potential negative consequences to the environment and human health. For more detailed information with regard to returning and recycling this product, please consult the distributor from whom you purchased it.</p> <p>* For system products, this label may be attached to the main unit only.</p>		

2 The Basics

2.1 Monitor Description

2.1.1 Intended Use

This patient monitor is intended to be used for monitoring, displaying, reviewing, storing and transferring of multiple physiological parameters including ECG, heart rate (HR), respiration (Resp), temperature (Temp), SpO₂, pulse rate (PR), non-invasive blood pressure (NIBP), invasive blood pressure (IBP), cardiac output (C.O.), carbon dioxide (CO₂), oxygen (O₂), anesthetic gas (AG), impedance cardiograph (ICG), bispectral index (BIS), respiration mechanics (RM), continuous cardiac output (PiCCO), and central venous oxygen saturation (ScvO₂).

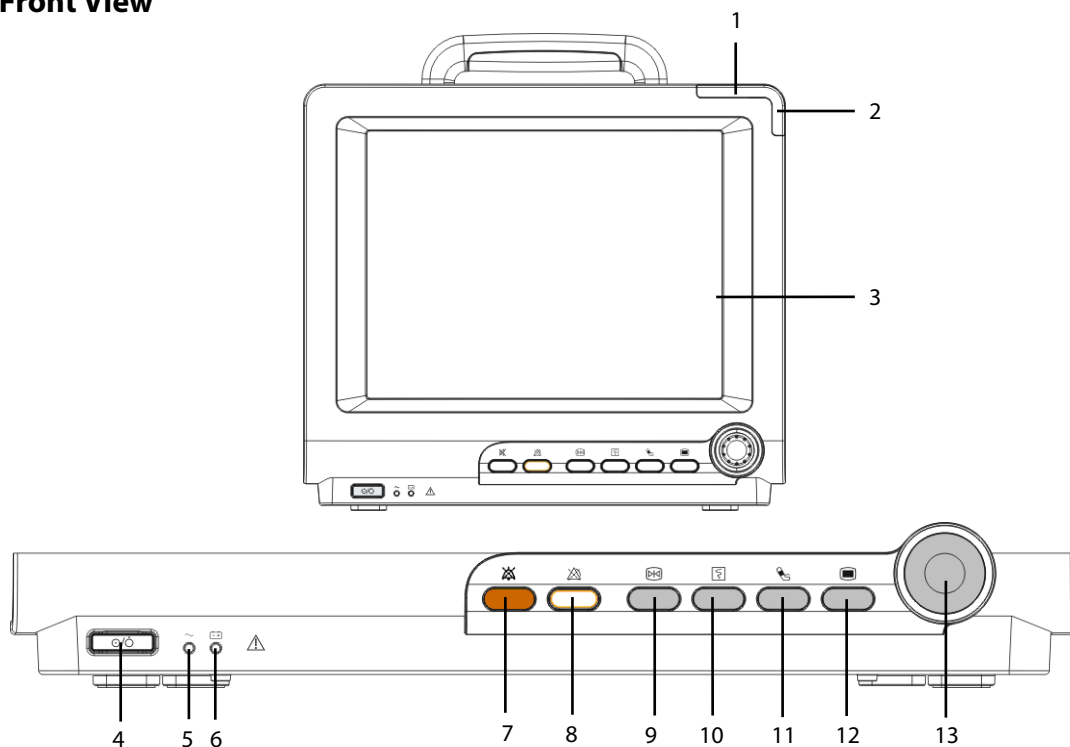
This monitor is to be used in healthcare facilities by clinical professionals or under their direction. It is not intended for helicopter transport, hospital ambulance, or home use.

WARNING

- **This patient monitor is intended for use only by clinical professionals or under their guidance. It must only be used by persons who have received adequate training in its use. Anyone unauthorized or untrained must not perform any operation on it.**
-

2.2 Main Unit

2.2.1 Front View



1. Alarm lamp

When a physiological alarm or technical alarm occurs, this lamp will flash as defined below.

- ◆ High level alarms: the lamp quickly flashes red.
- ◆ Medium level alarms: the lamp slowly flashes yellow.
- ◆ Low level physiological alarms: the lamp lights yellow without flashing.
- ◆ Low level technical alarms: the lamp does not light.

2. Technical alarm lamp

This lamp will light blue when a technical alarm occurs.

3. Display Screen

4. Power On/Off Switch


Press this switch to turn the patient monitor on. Press it again and hold for 2 seconds to turn the patient monitor off. An indicator is built in this switch. It turns on when the patient monitor is on and turns off when the patient monitor is off.


5. AC power LED

It turns on when AC power is connected.

6. Battery LED


- ◆ On: when at least a battery is installed to BeneView T5 monitor and the AC source is connected; when two batteries are installed to BeneView T8 monitor and the AC source is connected.
- ◆ Off: when no battery is installed, only one battery is installed to the BeneView T8 monitor, the installed battery is malfunction, or no AC source is connected when the patient monitor is power off.
- ◆ Flash: when the patient monitor operates on battery power.

7.  Press to silence all alarm sounds.

8.  Press to pause or restore alarms.

9.  Press to freeze or unfreeze waveforms.

10.  Press to start or stop recordings.

11.  Press to start or stop NIBP measurements.

12. 

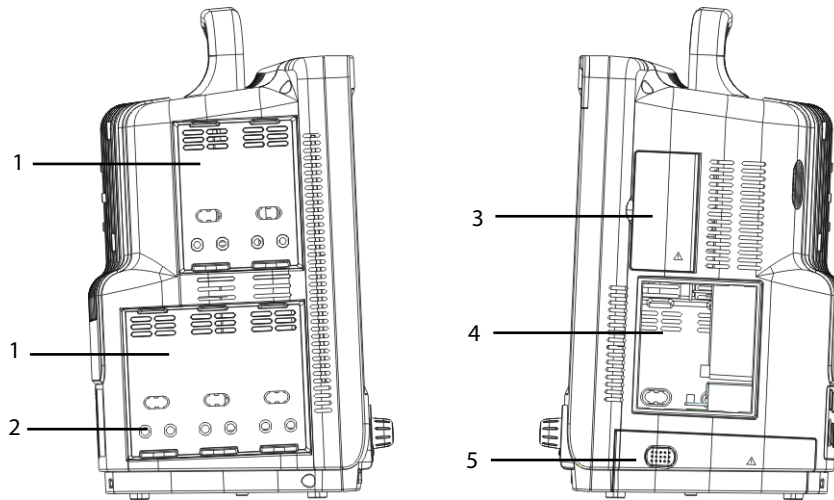
If no menu is displayed on the screen, pressing it will enter the main menu. If there is a menu displayed on the screen, pressing it will close that menu.

13. Knob

Rotate the Knob clockwise or anti-clockwise. With each click, the highlight jumps to the neighboring item. When you reach your desired item, press the Knob to select it.


2.2.2 Side View

BeneView T5



1. Integral Module Racks 2. Contact 3. Compartment for CF storage card slot
4. Recorder 5. Battery compartment

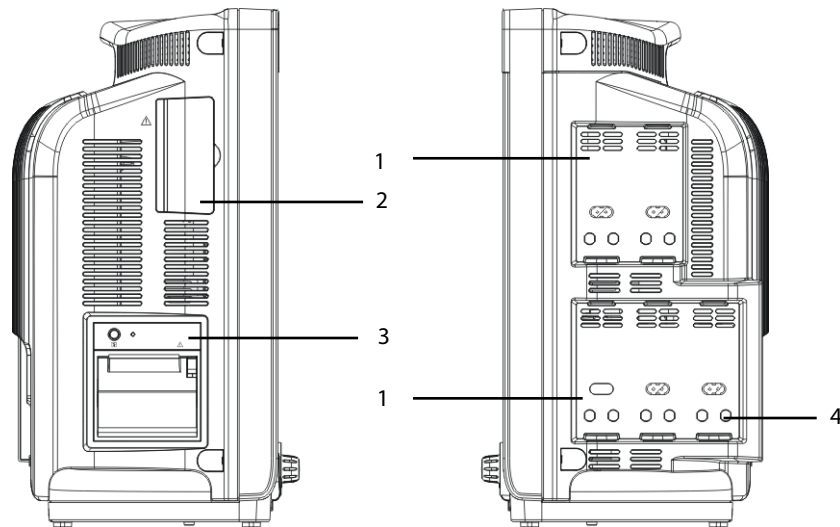BeneView T8




1. Integral Module Racks 2. Compartment for CF storage card slot
3. Recorder 4. Contact

NOTE

- To ensure a good contact, clean the contacts regularly, as dust and dirt may collect on them. When cleaning the contacts, wipe them with cotton, dampened with alcohol. (using forceps is recommended)


2.2.3 Rear View

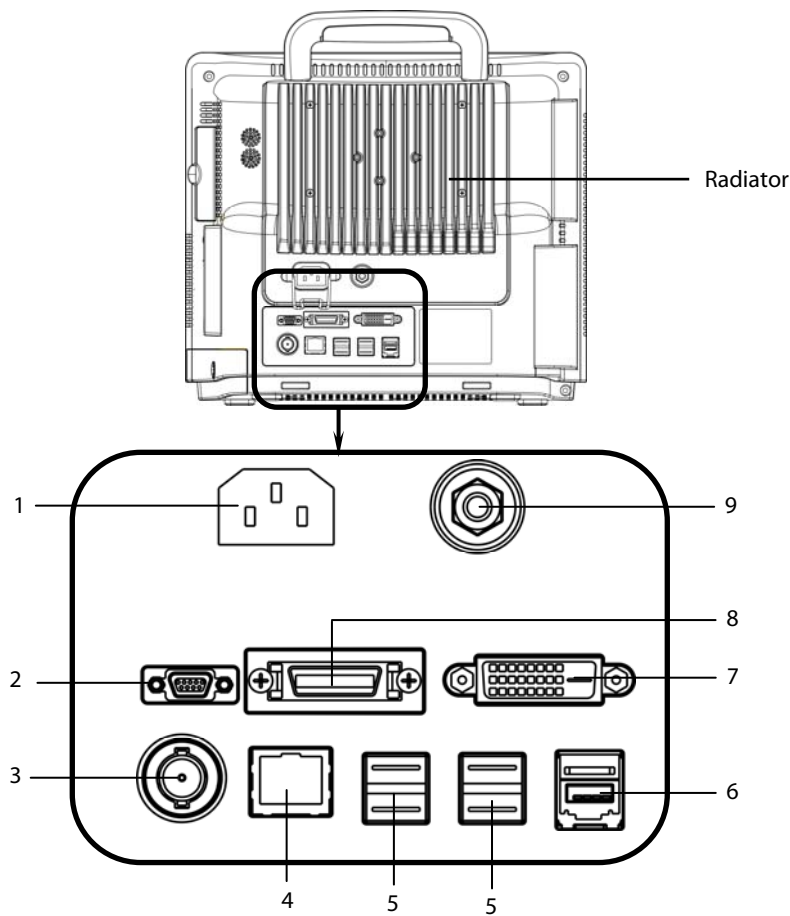BeneView T5




1. AC Power Input
2. Micro-D Connector: It outputs ECG, IBP and defibrillator synchronization signals simultaneously, among which the ECG signals supports pace pulses to be enhanced .
3. Nurse Call Connector: It connects the patient monitor to the hospital's nurse call system. Alarms indications are alerted to nurses through the nurse call system, if configured to do so.
4. Network Connector: It is a standard RJ45 connector that connects the patient monitor to the CMS.
5. USB Connectors: They connect such devices as the USB mice, USB keyboard, etc.
6. SMR Connector: It connects the satellite module rack (SMR).
7. Digital Video Interface (DVI): It connects a secondary display, which extends the display capability of your monitor. The contents displayed on the secondary display screen accords with those displayed on the monitor screen.
8. CIS Connector: It connects the patient monitor to the hospital's clinical information system (CIS) through an external CIS box.
9. Equipotential Grounding Terminal: When the patient monitor and other devices are to be used together, their equipotential grounding terminals should be connected together, eliminating the potential difference between them.

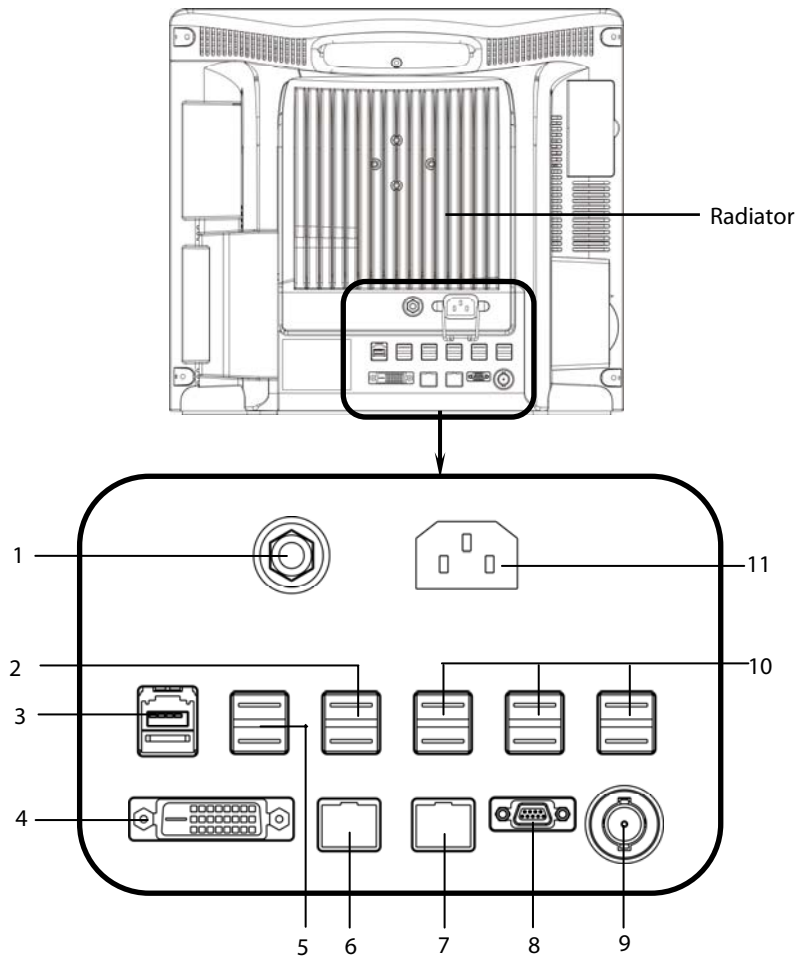BeneView T8




1. Equipotential Grounding Terminal: When the patient monitor and other devices are to be used together, their equipotential grounding terminals should be connected together, eliminating the potential difference between them.
2. USB Connectors: They are used for CIS maintenance.
3. SMR Connector: It connects the satellite module rack (SMR).
4. Digital Video Interface (DVI): It connects a secondary display, which extends the display capability of your monitor. The secondary display can be independently operated and controlled, and also display the contents different from the monitor screen.
5. USB Connectors: They connect the controlling devices (USB mouse and USB keyboard) of the secondary display.
6. CIS Connector: It is a standard RJ45 connector that connects the patient monitor to the hospital's clinical information system (CIS).
7. Network Connector: It is a standard RJ45 connector that connects the patient monitor to the CMS.
8. Micro-D Connector: It outputs ECG, IBP and defibrillator synchronization signals simultaneously, among which the ECG signals supports pace pulses to be enhanced.
9. Nurse Call Connector: It connects the patient monitor to the hospital's nurse call system. Alarms indications are alerted to nurses through the nurse call system, if configured to do so.
10. USB Connectors: They connect such devices as the USB mouse, USB keyboard, etc.
11. AC Power Input

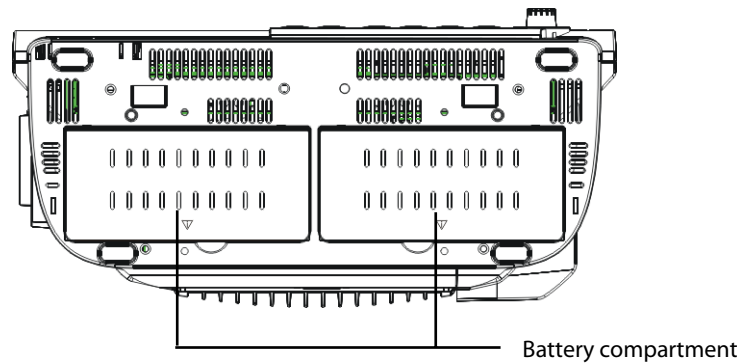2.2.4 Bottom View (BeneView T8)



2.3 Satellite Module Rack

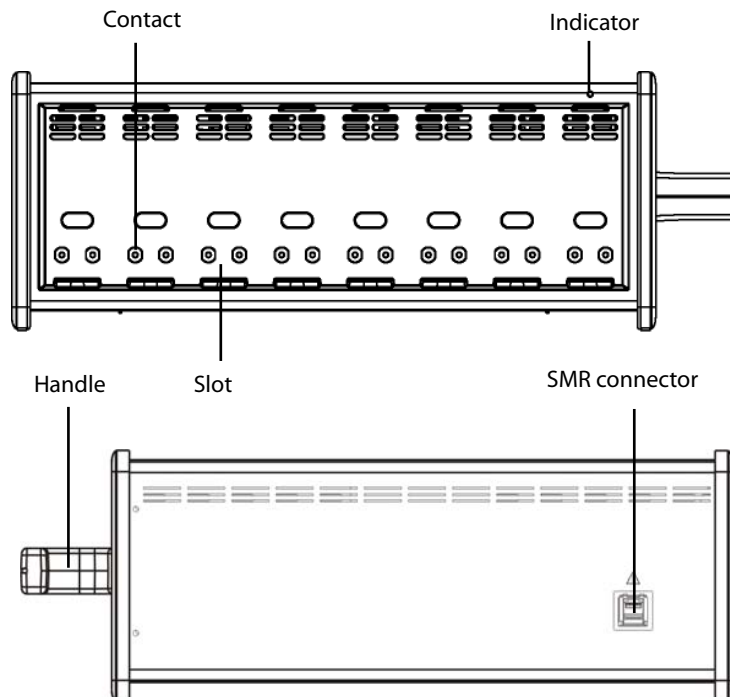The Satellite Module Rack (SMR) provides 8 slots for mounting measurement modules. The number of modules mounted in the SMR depends, as different modules may need different slots.



As shown in the figure above, there is an indicator telling the status of the SMR:

- On: when the SMR works normally.
- Off: when the SMR disconnects from the patient monitor, there is a problem with the power, or the patient monitor shuts down.

The SMR can be connected to the patient monitor through their SMR connectors via a SMR cable.

NOTE

- To ensure a good contact, clean the contacts regularly, as dust and dirt may collect on them. When cleaning the contacts, wipe them with cotton, dampened with alcohol. (using forceps is recommended)

2.4 Modules

As shown below, the patient monitor supports the following modules:



- SpO₂ module: Pulse oxygen saturation module.
- MPM: Multi-parameter module. It can simultaneously monitor ECG, respiration, SpO₂, temperature, NIBP and IBP.
- IBP module: Invasive blood pressure module.
- C.O. module: Cardiac output module.
- CCO/SvO₂ module: CCO/SvO₂ interface module, used to interface with Edwards Vigilance II® or Vigileo™ monitor.
- CO₂ module: Carbon dioxide module (including sidestream, microstream and mainstream).
- AG module: Anaesthesia gas module. The functions of the O₂ and BIS modules can be incorporated into it.
- BIS module: Bispectral index module.
- RM module: Respiration mechanics module.
- ICG module: Impedance cardiography module.
- PiCCO module : PiCCO module, used to measure cardiac output continuously.
- ScvO₂ module: Central venous oxygen saturation module.
- BeneLink module: BeneLink module is used for transmitting information from a connected external device to the BeneView patient monitor.
- BeneView T1: Used as a multi-measurement module for monitoring ECG, respiration, SpO₂, temperature, NIBP and IBP.

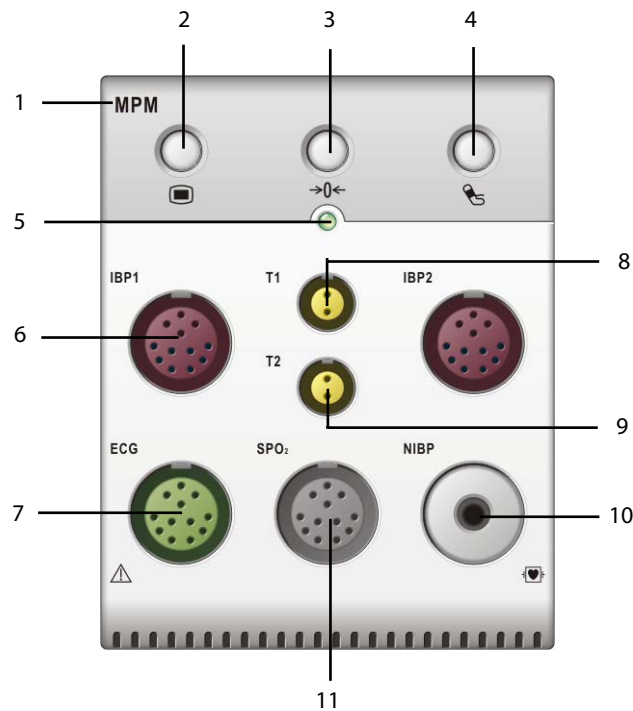
Under the maximum configuration, the patient monitor has one two-slot module rack, one three-slot module rack and one satellite module rack. The number of modules mounted in the patient monitor depends, as different modules may need different slots.

You can plug and unplug modules during patient monitoring. To plug a module, insert the module until the lever on the module clicks into place and then push the lock key at the bottom in position to lock the module. To remove a module, release the lock key, press the lever upwards and pull the module out.

Make sure that the indicator on the module lights on after the module is plugged in. Otherwise, re-plug the module until the indicator lights on.

2.4.1 Multi-Parameter Module


The multi-parameter module (MPM) incorporates multiple measurement modules. As shown below, the module name is located at the upper left corner, all hardkeys on the upper part, and all measurement connectors on the lower part. Other measurement modules look similar to the MPM.



1. Module name
2. Setup key: press to enter the **[MPM Setup]** menu.
3. Zero key: press to enter the **[Zero IBP]** menu.
4. NIBP start/stop key: press to start or stop NIBP measurements.
5. Indicator
 - ◆ On: when the patient monitor works correctly.
 - ◆ Flash: when the module is being initialized.
 - ◆ Off: when the module is either unconnected or broken.
6. Connector for IBP cable
7. Connector for ECG cable
8. Connector for Temp probe 1
9. Connector for Temp probe 2
10. Connector for NIBP Cuff
11. Connector for SpO₂ cable

2.4.2 BeneView T1

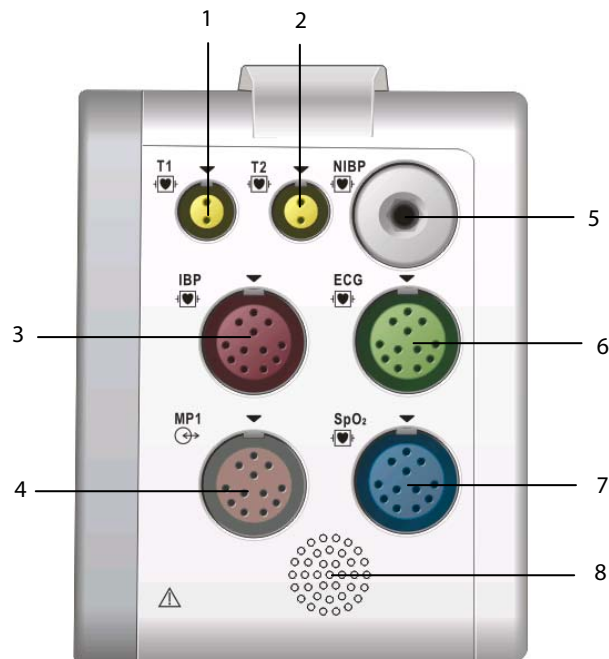BeneView T1 can be connected to a BeneView patient monitor as a multi-measurement module.



1. Connector for Temp probe 1
2. Connector for Temp probe 2
3. Connector for IBP cable
4. Multifunctional connector, connecting external parameter module and outputting analog and defib synchronization signal.
5. Connector for NIBP cuff
6. Connector for ECG cable
7. Connector for SpO₂ cable
8. Speaker

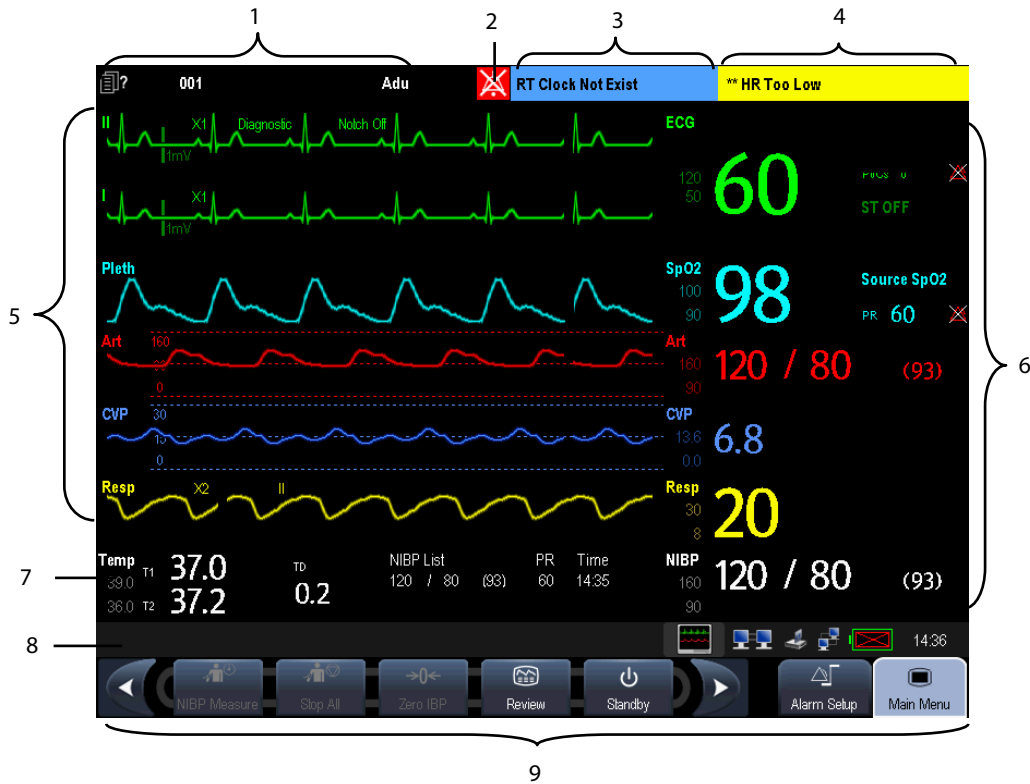
When the T1 is disconnected from BeneView T5 or T8, it can continue to monitor a patient as a stand-alone monitor running on battery power or external DC power supply. For details of using T1 as a stand-alone monitor, refer to ***BeneView T1 Operating Manual***.

NOTE

- **Micro-D connector is disabled when the multifunctional connector of T1 is in use.**
- **Please do not charge more than one BeneView T1 simultaneously with the module rack.**

2.5 Display Screen

This patient monitor adopts a high-resolution TFT LCD to display patient parameters and waveforms. A typical display screen is shown below.



1. Patient Information Area


This area shows the patient information such as department, bed number, patient name, patient category and paced status.

- ◆ : indicates that no patient is admitted or the patient information is incomplete.
- ◆ : indicates that the patient has a pacer.
- ◆ : indicates that the patient does not has a pacer or the patient's paced status is not selected.

If no patient is admitted, selecting this area will enter the **[Patient Setup]** menu. If a patient has been admitted, selecting this area will enter the **[Patient Demographics]** menu.

2. Alarm Symbols

- ◆ indicates alarms are paused.
- ◆ indicates alarm sounds are paused.
- ◆ indicates alarm sounds are turned off.

- ◆  indicates the system is in alarm off status.

3. Technical Alarm Area

This area shows technical alarm messages and prompt messages. When multiple messages come, they will be displayed circularly. Select this area and the technical alarm list will be displayed.

4. Physiological Alarm Area

This area shows physiological alarm messages. When multiple alarms occur, they will be displayed circularly. Select this area and the physiological alarm list will be displayed.

5. Waveform Area

This area shows measurement waveforms. The waveform name is displayed at the left upper corner of the waveform. Select this area and the corresponding measurement setup menu will be displayed.

6. Parameter Area A


This area shows measurement parameters. Each monitored parameter has a parameter window and the parameter name is displayed at the upper left corner. The corresponding waveform of each parameter is displayed in the same row in the waveform area. Select this area and the corresponding measurement setup menu will be displayed.

7. Parameter Area B


For the parameters displayed in this area, their corresponding waveform are not displayed.

8. Prompt Message Area

This area shows the prompt messages, network status icons, battery status icons, date and time, etc. For details about battery status symbols, refer to the chapter **32** Batteries.

- ◆  indicates patient monitor is connected to a wire network successfully.


- ◆  indicates the patient monitor has failed to connect a wire network.

- ◆  indicates the wireless function is working.

- ◆  indicates the wireless function is not working.

- ◆  indicates a CF storage card is inserted.

- ◆  indicates a USB disk is inserted.

- ◆  indicates a secondary display or remote display is connected.

- ◆  **[Screen Setup]** button

9. QuickKeys Area

This area contains QuickKeys that give you fast access to functions.

2.6 QuickKeys

A QuickKey is a configurable graphical key, located at the bottom of the main screen. They give you fast access to functions. Their availability and the order in which they appear on your screen, depend on how your patient monitor is configured.

The following QuickKeys can be displayed on the screen:

	Scroll left to display more QuickKeys.		Scroll right to display more QuickKeys.
	Enter the main menu		Enter standby mode
	Change alarm settings		Review the patient's data
	Enter the NIBP measurement menu		Stop all NIBP measurement
	Zero IBP		Silence all alarm sounds
	Pause or restore alarms		Change screen
	Enter the patient setup menu		Trigger a manual event
	Start the realtime print		Print Setup
	Have a split-screen view of minitrends		Enter the volume setup menu
	Load configurations		Have the CIS view
	Start cardiac output procedure (not available in USA)		View respiratory loops
	Perform calculations		Have a split-screen view of another patient's conditions
	Enter the full-screen 7-lead ECG screen		Have a split-screen view of oxyCRG trends
	Enter the [Parameters] menu		Enter the interpretation of resting 12-lead ECG screen (not available in USA)
	Start NIBP STAT measurement		Enter the [Unit Setup] menu
	Enter the PAWP measurement screen		Enter the CPB mode
	Enter the privacy mode		Enter the night mode

You can also select your desired QuickKeys to display on the screen.

1. Select [**Main Menu**]→[**Maintenance >>**]→[**Manage Configuration >>**]→enter the required password→[**Ok**].
2. In the [**Manage Configuration**] menu, select [**Edit Config.>>**].
3. In the pop-up menu, select the desired configuration and then select [**Edit**].
4. In the pop-up menu, select [**Screen Setup >>**].
5. In the [**Select QuickKeys**] screen, select your desired QuickKeys and the order of them.

FOR YOUR NOTES

3 Basic Operations

3.1 Installation

WARNING

- The equipment shall be installed by personnel authorized by us.
 - The software copyright of the equipment is solely owned by us. No organization or individual shall resort to juggling, copying, or exchanging it or to any other infringement on it in any form or by any means without due permission.
 - Devices connected to the equipment must meet the requirements of the applicable IEC standards (e.g. IEC 60950 safety standards for information technology equipment and IEC 60601-1 safety standards for medical electrical equipment). The system configuration must meet the requirements of the IEC 60601-1-1 medical electrical systems standard. Any personnel who connect devices to the equipment's signal input/output port is responsible for providing evidence that the safety certification of the devices has been performed in accordance to the IEC 60601-1-1. If you have any question, please contact us.
 - If it is not evident from the equipment specifications whether a particular combination with other devices is hazardous, for example, due to summation of leakage currents, please consult the manufacturers or else an expert in the field, to ensure the necessary safety of patients and all devices concerned will not be impaired by the proposed combination.
 - Not using screw and bracket specified by Mindray may cause the screw to touch the internal battery and lead to monitor damage.
-

3.1.1 Unpacking and Checking

Before unpacking, examine the packing case carefully for signs of damage. If any damage is detected, contact the carrier or us.

If the packing case is intact, open the package and remove the equipment and accessories carefully. Check all materials against the packing list and check for any mechanical damage. Contact us in case of any problem.

NOTE

- Save the packing case and packaging material as they can be used if the equipment must be reshipped.
-

WARNING

- When disposing of the packaging material, be sure to observe the applicable waste control regulations and keep it out of children's reach.
 - The equipment might be contaminated during storage and transport. Before use, please verify whether the packages are intact, especially the packages of single use accessories. In case of any damage, do not apply it to patients.
-
-

3.1.2 Environmental Requirements

The operating environment of the equipment must meet the requirements specified in this manual.

The environment where the equipment is used shall be reasonably free from noises, vibration, dust, corrosive, flammable and explosive substances. If the equipment is installed in a cabinet, sufficient space in front and behind shall be left for convenient operation, maintenance and repair. Moreover, to maintain good ventilation, the equipment shall be at least 2 inches (5cm) away from around the cabinet.

When the equipment is moved from one place to another, condensation may occur as a result of temperature or humidity difference. In this case, never start the system before the condensation disappears.



WARNING

- **Make sure that the operating environment of the equipment meets the specific requirements. Otherwise unexpected consequences, e.g. damage to the equipment, could result.**
-

3.2 Getting Started

3.2.1 Turning Power On

Once the patient monitor is installed, you can get ready for monitoring:

1. Before you start to make measurements, check the patient monitor, SMR and plug-in modules for any mechanical damage and make sure that all external cables, plug-ins and accessories are properly connected.
2. Plug the power cord into the AC power source. If you run the patient monitor on battery power, ensure that the battery is sufficiently charged.
3. Press the power on/off switch on the monitor's front. The start-up screens are displayed, and the technical alarm lamp and alarm lamp are lit in blue and yellow respectively. Then, the alarm lamp turns into red, and turns off together with the technical alarm lamp after the system gives a beep.
4. The monitor enters the main screen.



WARNING

- **Do not use the patient monitor for any monitoring procedure on a patient if you suspect it is not working properly, or if it is mechanically damaged. Contact your service personnel or us.**
-

3.2.2 Starting Monitoring

1. Decide which measurements you want to make.
 2. Connect the required modules, patient cables and sensors.
 3. Check that the patient cables and sensors are correctly connected.
 4. Check that the patient settings such as **[Patient Cat.]**, **[Paced]**, etc, are appropriate for your patient.
 5. Refer to the appropriate measurement section for details of how to perform the measurements you require.
-

3.3 Disconnecting from Power

To disconnect the patient monitor from the AC power source, follow this procedure:

1. Confirm that the patient monitoring is finished.
2. Disconnect patient cables and sensors from the patient.
3. Make sure to save or clear the patient monitoring data as required.
4. Press and hold the power on/off switch for above 2 seconds. The patient monitor shuts down and you can unplug the power cable.



CAUTION

- **Although not recommended, you can press and hold the power on/off switch for 10 seconds to forcibly shut down the monitor when it could not be shut down normally or under some special situations. This may cause loss of data of the patient monitor.**
-

3.4 Using a Mouse

You can use the USB mouse supplied with the equipment as a monitor input device. The USB mouse can be plugged and unplugged with the monitor on.

When you are using a mouse:

- By default, the left mouse-button is the primary button and the right one the secondary button.
- Clicking the primary button equals to pressing the knob or selecting the touchscreen.
- The secondary button is disabled.

You can also define the right mouse-button as the primary button by following this procedure:

1. Select [**Main Menu**]→[**Maintenance >>**]→[**User Maintenance >>**]→enter the required password.
2. Select [**Others >>**] to enter the [**Others**] menu.
3. Select [**Primary Button**] and then select [**Right**] from the popup list.




3.5 Using Keys

The monitor has three types of keys:

- Softkey: A softkey is a graphic key on the screen, giving you fast access to certain menus or functions. The monitor has three types of softkeys:
 - ◆ Parameter keys: Each parameter area or waveform area can be seen as a softkey. You can enter a parameter setup menu by selecting its corresponding parameter area or waveform area.
 - ◆ QuickKeys: QuickKeys are configurable graphical keys, located at the bottom of the main screen. For details, refer to the section **QuickKeys**.
- Hardkeys: A hardkey is a physical key on a monitoring device, such as the main menu hardkey on the monitor's front.
- Pop-Up Keys: Pop-up keys are task-related keys that appear automatically on the monitor screen when required. For example, the confirm pop-up key appears only when you need to confirm a change.


3.6 Using Keyboards

The on-screen keyboard enables you to enter information.


- Use the  key to delete the previously entered character.
- Use the  key to toggle between uppercase and lowercase letters.
- Select  to confirm what you have entered and close the on-screen keyboard.
- Press the switch button to print the special letter in the keyboards for Dutch, French, Spanish, and Portuguese. A special letter is composed of a normal letter and a symbol, for example, in French 'â' is a special letter composed of letter 'a' and symbol '^'. To print 'â', you should first press the switch button, and then the target special letter 'â'. The following table defines the switch buttons and the special letters corresponding with the keyboards for each language:


Language	Switch Button	Special Letter
Dutch	'	á, é, ú, í, ó, ý, ç, Á, É, Ú, Í, Ó, Ý, Ç
French	^	â, ê, û, î, ô
	¨	Ä, Ê, Ü, Î, Ö
Spanish	`	à, è, ù, ì, ò, À, È, Ù, Ì, Ò
	'	á, é, ú, í, ó, Á, É, Ú, Í, Ó
Portuguese	'	á, é, ú, í, ó, ý, Á, É, Ú, Í, Ó, Ý
	~	ã, õ, ñ, Ã, Õ, Ñ

3.7 Using the Touchscreen

Select screen items by pressing them directly on the patient monitor's screen. You can enable or disable touchscreen operation by pressing and holding the **[Main Menu]** QuickKey for 3 seconds. A padlock symbol  is displayed if touchscreen operation is disabled.

3.8 Setting the Screen

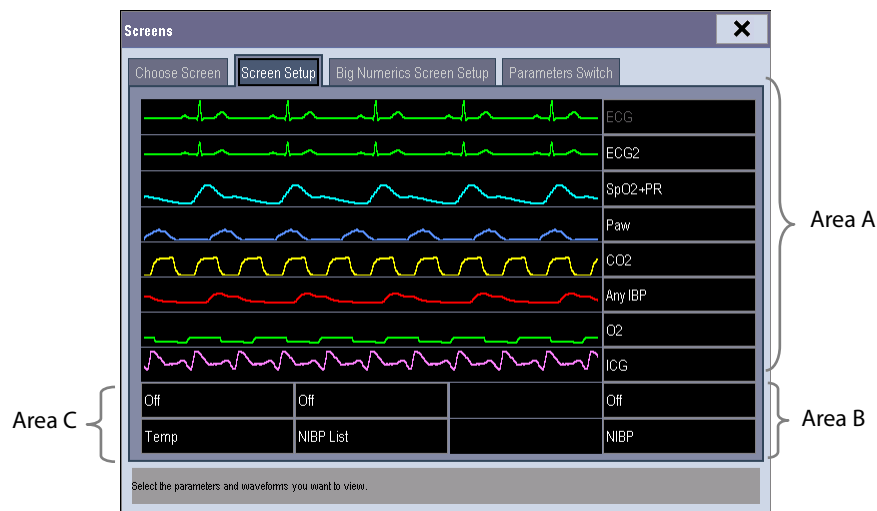You can enter the **[Screen Setup]** window as shown below by selecting the **[Screen Setup]** button  in the prompt message area. In this window, you can allocate the positions of the parameters and waveforms. The parameters or waveforms whose positions are not allocated will not be displayed.




The ECG parameter and the first ECG waveform always display in the first row. The configurable areas can be classified as Area A, Area B, and Area C.

- In Area A, you can choose to display the parameters (having waveforms) and their waveforms. Each parameter and the associated waveform are displayed in the same row.
- In Area B, you can choose to display the parameters and their waveforms. When there is no parameter displayed in area C, both the parameters and their waveforms will be displayed in area B. Otherwise, only the parameters will be displayed.
- In Area C, you can choose to display all the parameters whose associated waveforms will not be displayed.

The screen can automatically adjust to ensure the best view based on your screen setup.

If no corresponding parameter or waveform is displayed after the module is inserted, you should perform the following inspections:



- Check the connection between the module and lead, cable, sensor, or external device.
- Check whether there are the **[The display setup for XX is disabled]** message and the flashing **[Screen**


Setup] button  in the prompt message area. If yes, select this button to enter the **[Screen Setup]** window for the desired display configuration.

 **WARNING**

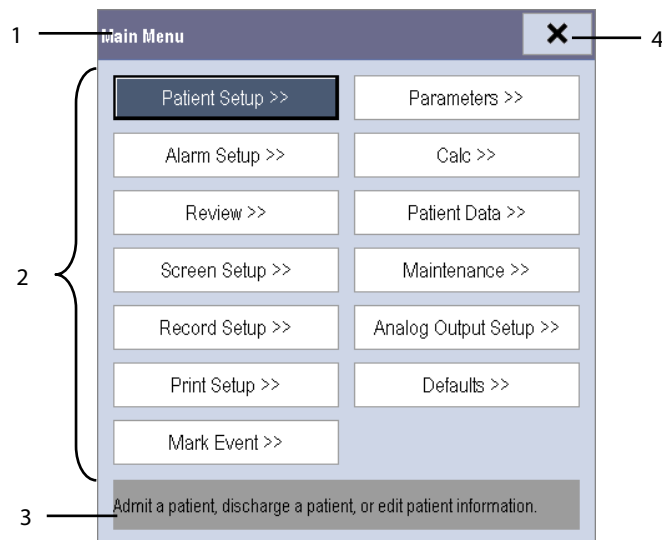
- **The parameters whose positions are not allocated in the [Screen Setup] window will not be displayed. However, the monitor can still give alarms of these parameters.**
-

3.9 Using the Main Menu

To enter the main menu, select the  on-screen QuickKey or the  hardkey on the monitor's front. Most of monitor operations and settings can be performed through the main menu.



Other menus are similar to the main menu and contain the following parts:

1. Heading: gives a sum-up for the current menu.
2. Main body: displays options, buttons, prompt messages, etc. The menu button with ">>" enlarges a secondary window to reveal more options or information.
3. Online help area: displays help information for the highlighted menu item.
4. : select to exit the current menu.

3.10 Setting Parameters

3.10.1 Switching the Parameters On/Off

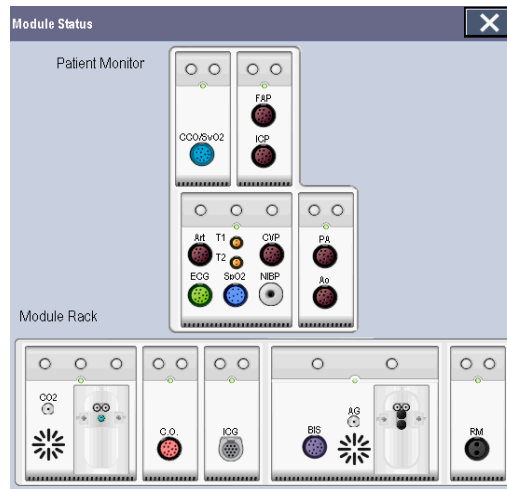
To switch the parameters on or off, select [Main Menu] → [Screen Setup >>] → [Screen Layout >>] → [Parameters Switch], or [Screen Layout] QuickKey → [Parameters Switch]. When a parameter is switched off, its corresponding parameter module stops working, and its parameter value and waveform are not shown on the monitor display.

NOTE

- **ECG is always selected, and you can not switch it off.**
-

3.10.2 Accessing the Parameters Menu

Select [**Parameters >>**] from the main menu or select the [**Parameters**] QuickKey at the bottom of the screen to enter the [**Parameters**] menu where you can get the access of each parameter's setup menu. You can further select [**Module Status >>**] to enter the menu as shown below. Your display may be configured to look slightly different depending on the modules mounted.



This menu displays the measurement modules mounted in the two-slot module rack, three-slot module rack and satellite module rack from top to bottom. Beside each measurement connector is the measurement label. The color in which a measurement connector appears matches the status as follows:



(colored) indicates that the module is turned on.



(grey) indicates that the module is turned off.



indicates a module name conflict.



indicates a module error.

3.10.3 Removing a Module Conflict

Besides three independent IBP modules and the IBP module on the MPM, the patient monitor supports only one more measurement module simultaneously. Otherwise, the message of module conflict will be prompted.

For example, if a CO₂ module (module A) is already loaded and then another CO₂ module (module B) is inserted, your patient monitor will then display module conflict. To use module A, just pull out module B. To use module B, pull both modules A and B out and then re-insert module B.

3.10.4 Removing a Label Conflict

Every label is unique and is assigned only once. The measurement label is stored inside the module. The system will prompt module name conflict when two measurement modules with the same name are used.

For example, an IBP module (module A) is already loaded and the Art label is used for module A. Then another IBP module (module B) is inserted and the Art label is also used for module B. In this case, your patient monitor will prompt the message of label conflict and display the **[Label]** menu.


- To use module A for Art measurement, just modify the label of module B on this channel in the **[Label]** menu. If the **[Label]** menu already exits inadvertently, you need to plug out and then plug in module B.
- To use module B for Art measurement, first exit the **[Label]** menu. Then select the Art parameter area on the screen and modify the label of module A on this channel in the popup menu. Finally, plug out and then plug in module B.

3.11 Using a CF Storage Card

A CF storage card is used to prevent data loss in case of a sudden power failure. The patient data such as trend data, waveform data, etc., will be automatically saved into the CF storage card during patient monitoring. In case of a sudden power failure, the patient data can be retrieved from the CF storage card after the patient monitor restarts.

To insert a CF storage card, open the compartment and then insert the card until the button flips out.

To remove the CF storage card, follow this procedure:

1. In the main menu, select **[Unload CF Storage Card]**, or **[Patient Data]**→**[Unload CF Storage Card]**. You can also click  icon in the lower right corner of the screen.
2. Select **[Ok]** from the popup menu to unload the CF storage card. A status message shown in the prompt message area will report completion of the unloading.
3. Press the button until the CF storage card flips out.

To browse the data saved in the CF storage card, follow this procedure:

1. Select **[Main Menu]**→**[Patient Data >>]**→**[History Data >>]**.
2. Select a patient whose data you want to view from the **[Patient Data List]** and then select **[Review]**.
3. Select **[Data Review]**.

As reviewing the history patient's data is just like reviewing the current patient's data, you can refer to the chapter **28 Review** for details.

NOTE

-
- **Data may be unable to be saved into the CF storage card when the patient monitor is just turned on.**
 - **If no CF storage card is used, all the data you have saved will get lost in case of monitor shut-down or sudden power interrupt.**
-



CAUTION

- **Unload the CF storage card before removing it from the patient monitor. Otherwise it may cause damage to the data in the card.**
 - **Use only the CF storage card specified by Mindray.**
 - **Please take measures against the static electricity such as Disposable Wrist Strap when you fetch the CF card.**
-

3.12 Changing General Settings

This chapter covers only general settings such as language, brightness, date and time, etc. Measurement settings and other settings can be referred to in respective sections.

3.12.1 Setting up a Monitor

In situations where you install a patient monitor or change the patient monitor's application site, you need to setup the patient monitor as follows:

1. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password.
2. In the **[User Maintenance]** menu, select, in turn, **[Monitor Name]**, **[Department]** and **[Bed No.]**, and then change their settings.

3.12.2 Changing Language

1. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password.
2. In the **[User Maintenance]** menu, select **[Language]** and then select the desired language.
3. Restart the patient monitor.

3.12.3 Adjusting the Screen Brightness

1. Select the **[Main Menu]**→**[Screen Setup >>]**→**[Brightness]**.
2. Select the appropriate setting for the screen brightness. 10 is the brightest, and 1 is the least bright.

If the patient monitor operates on battery power, you can set a less bright screen to prolong the operating time of the battery. When the patient monitor enters standby mode, the screen will change to the least brightness automatically.

3.12.4 Showing/Hiding the Help

The patient monitor provides online help information. The user can display or hide the help as required.

1. Select **[Main Menu]**→**[Screen Setup >>]**.
2. Select **[Help]** and toggle between **[On]** and **[Off]**.

3.12.5 Setting the Date and Time

1. Select [**Main Menu**] → [**Maintenance >>**] → [**System Time >>**].
2. Set the date and time.
3. Select [**Date Format**] and toggle between [**yyyy-mm-dd**], [**mm-dd-yyyy**] and [**dd-mm-yyyy**].
4. Select [**Time Format**] and toggle between [**24h**] and [**12h**].

If your patient monitor is connected to a central monitoring system (CMS), the date and time are automatically taken from that CMS. In that case, you cannot change the date and time settings on your patient monitor.



CAUTION

- **Changing date and time affects the storage of trends and events and may cause data missing.**
-

3.12.6 Adjusting Volume

Alarm Volume

1. Select the [**Volume Setup**] QuickKey, or [**Main Menu**] → [**Alarm Setup >>**] → [**Others**].
2. Select [**Alm Volume**] and then select the appropriate volume: X-10, in which X is the minimum volume, depending on the set minimum alarm volume (refer to the chapter Alarm), and 10 the maximum volume.

Key Volume

1. Select the [**Volume Setup**] QuickKey, or [**Main Menu**] → [**Screen Setup >>**].
2. Select [**Key Volume**] and then select the appropriate volume. 0 means off, and 10 the maximum volume.

QRS Volume

The QRS tone is derived from either the HR or PR, depending on which is currently selected as the alarm source in [**ECG Setup**] or [**SpO₂ Setup**]. When monitoring SpO₂, there is a variable pitch tone which changes as the patient's saturation level changes. The pitch of the tone rises as the saturation level increases and falls as the saturation level decreases. The volume of this tone is user adjustable.

1. Select the [**Volume Setup**] QuickKey, or the ECG parameter window → [**Others >>**], or the SpO₂ parameter window.
2. Select [**QRS Volume**] or [**Beat Vol**] and then select the appropriate volume. 0 means off, and 10 the maximum volume.

4 Managing Patients

4.1 Admitting a Patient


The patient monitor displays physiological data and stores them in the trends as soon as a patient is connected. This allows you to monitor a patient that is not admitted yet. However, it is recommended that you fully admit a patient so that you can clearly identify your patient, on recordings, reports and networking devices. To admit a patient:

1. Select the **[Patient Setup]** QuickKey, or **[Main Menu]→[Patient Setup >>]**.
2. Select **[Discharge Patient]** to clear any previous patient data. If you do not erase data from the previous patient, the new patient's data will be saved into the data of the previous patient. The monitor makes no distinction between the old and the new patient data.
3. If **[Discharge Patient]** button appears dimmed, directly select **[Admit Patient]** and then select:
 - ◆ **[Yes]** to apply the data saved in the patient monitor to the new patient, or
 - ◆ **[No]** to clear the data saved in the patient monitor.
4. In the **[Patient Demographics]** menu, enter the demographic details, of which:
 - ◆ **[Patient Cat.]** determines the way your patient monitor processes and calculates some measurements, and what safety and alarm limits are applied for your patient.
 - ◆ **[Paced]** determines whether to show pace pulse marks on the ECG waveform. When the **[Paced]** is set to **[No]**, pace pulse marks are not shown in the ECG waveform.
5. Select **[Ok]**.

WARNING

- **[Patient Cat.] and [Paced] will always contain a value, regardless of whether the patient is fully admitted or not. If you do not specify settings for these fields, the patient monitor uses the default settings from the current configuration, which might not be correct for your patient.**
 - **For paced patients, you must set [Paced] to [Yes]. If it is incorrectly set to [No], the patient monitor could mistake a pace pulse for a QRS and fail to alarm when the ECG signal is too weak.**
 - **For non-paced patients, you must set [Paced] to [No].**
-

4.2 Quick Admitting a Patient

Use **[Quick Admit]** only if you do not have the time or information to fully admit a patient. Complete the rest of the patient demographic details later. Otherwise, the  symbol will always be displayed in the patient information area.

1. Select the **[Patient Setup]** QuickKey, or **[Main Menu]→[Patient Setup >>]**.
2. Select **[Quick Admit]**. If a patient has been admitted at present, select **[OK]** to discharge the current patient. If no patient is admitted, you can choose either:
 - ◆ **[Yes]** to apply the data in your patient monitor to the new patient, or
 - ◆ **[No]** to clear any previous patient data.
3. Enter the patient category and paced status for the new patient, and then select **[Ok]**.

4.3 Querying and Obtaining Patient Information

The monitor can obtain patient information from HIS through eGateway. To query or obtain patient information from HIS,

1. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password→**[Gateway Communication Setting >>]**, and set**[IP Address]** and **[Port]**. Set **[ADT Query]** to **[On]**.
2. Click patient information area to enter the **[Patient Demographics]** menu.
3. Select **[Obtain Patient Info. >>]** to enter the **[Obtain Patient Information]** menu.
4. Input query condition and then select **[Query]**. The monitor will display the obtained patient information.
5. Select a patient and then click **[Import]**. Then the monitor will update the information of corresponding patient.
6. Select **X** to exit the **[Obtain Patient Information]** menu.

NOTE

- The option **[Obtain Patient Information]** is available in the **[Patient Setup]** menu only when **[ADT Query]** is set to **[On]**.
 - When obtaining patient information from HIS, the monitor only update patient information. The patient's monitoring data is not changed and the patient is not discharged.
-

4.4 Associating Patient Information

After associating patient information with HIS, the monitor will automatically update patient information if corresponding information in HIS has been is changed. The monitor can associate patient's MRN, first name, last name, date of birth, and gender with HIS.

NOTE

- A keyword takes effect only when being defined in eGateway. Refer to *eGateway Integration Manager Installation Guide* for details.
 - The monitor displays corresponding patient information only when all the keywords have been inputted.
-

4.5 Editing Patient Information

To edit the patient information after a patient has been admitted, or when the patient information is incomplete, or when you want to change the patient information:

1. Select the **[Patient Setup]** QuickKey, or **[Main Menu]**→**[Patient Setup >>]**.
2. Select **[Patient Demographics]** and then make the required changes.
3. Select **[Ok]**.

4.6 Discharging a Patient

To discharge a patient:

1. Select the [**Patient Setup**] QuickKey, or [**Main Menu**]→[**Patient Setup >>**].
2. Select [**Discharge Patient**]. In the popup menu, you can either:
 - ◆ Directly select [**Ok**] to discharge the current patient, or
 - ◆ Select [**Standby**] then [**Ok**]. The patient monitor enters the standby mode after discharging the current patient, or
 - ◆ Select [**Cancel**] to exit without discharging the patient.

NOTE

- **Discharging a patient clears all history data in the monitor.**
-

4.7 Transferring a Patient

You can transfer a patient with an MPM or BeneView T1 (referred to as T1 hereafter) to a new location without re-entering the patient demographic information or changing the settings. Transferring of patient data enables you to understand the patient's history condition. The patient data that can be transferred includes: patient demographics, trend data, alarm events and parameters alarm limits.

Select [**Others >>**] from [**User Maintenance**] menu. In the popup menu, you can set [**Transferred Data Length**]. The default is [**4 h**]. You can also set [**Data Transfer Method**]. The default is [**Off**].

WARNING

- **Do not discharge a patient before the patient is successfully transferred.**
 - **After a patient is successfully transferred, check if the patient settings (especially patient category, paced status and alarm limits settings, etc) on the monitor are appropriate for this patient.**
 - **Only when you open MPM transfer function and select [**Continue Module**], the IBP labels can be transferred along with the MPM.**
-
-

NOTE

- **The system automatically switches on the HR alarm and lethal arrhythmia alarm after transferring the patient data.**
-

4.7.1 Transferring Patient Data via MPM/T1

Familiarizing yourself with the data respectively stored in the patient monitor, T1 or MPM helps you understand the effects incurred by transferring patients with an MPM/T1.

Contents stored		In the patient monitor	In the MPM	In the T1
Data	Patient demographics (Name, Bed No., Gender, etc.)	Yes	Yes	Yes
	Trend data	Yes	Yes	Yes
	Calculation data (Dose calculations, oxygenation calculations, etc.)	Yes	No	No
	Event data (Marked events, alarm events, etc.)	Yes	No	Yes
Settings	Monitor settings (Alarm pause, alarm volume, etc.)	Yes	No	No
	Parameter settings (Alarm limits, etc.)	Yes	Yes	Yes

Before transferring a patient with an MPM/T1, set the destination monitor as follows:

1. Select [**Main Menu**]→[**Maintenance**]→[**User Maintenance >>**]→enter the required password.
2. Select [**Others >>**].
3. Set [**Data Transfer Method**] to [**Module**].
4. Set [**Apply Module Settings**] to [**On**]. If your patient monitor does not have this option, the system applies the MPM/T1's settings by default.

Then, follow this procedure to transfer the patient:

1. Disconnect MPM/T1 from the original monitor.
2. Connect MPM/T1 to the destination monitor.
3. If there is a mismatch between the MPM/T1 and monitor, the system will automatically display the [**Select Patient**] menu, from which you can choose the data set you want to continue using for this patient, either:
 - ◆ [**Continue Monitor**]:continue with the patient data and settings in the monitor, deleting all patient data and setting in MPM/T1 and copying all data in the monitor to MPM/T1.
 - ◆ [**Continue Module**]:continue with the patient data and settings in MPM/T1. Discharge the patient in the monitor. The monitor then automatically admits the patient and copies all data from MPM/T1.
 - ◆ [**New Patient**]: select this button if none of the information is correct. This deletes all data in the monitor and MPM/T1 and lets you admit a new patient on the monitor. In this case, you need to re-enter the patient demographics. The monitor will restore the settings according to the patient category.
 - ◆ [**Same Patient**]:select this button if the patient demographics are different, but it is the same patient. This merges the patient's trend data in the monitor and MPM/T1 and copies the settings in MPM/T1 to the monitor as well.
4. Select [**Yes**].

Operations	Examples of applications
Continue Monitor	1. Replace MPM/T1 during patient monitoring. 2. After the patient is admitted, connect the MPM/T1.
Continue Module	A patient is monitored using MPM/T1. You need to transfer the patient, e.g. from a ward (original monitor) to the operating room (destination monitor).
New Patient	Connect the MPM/T1 before admitting a new patient. However, the monitor and/or MPM/T1 store the previous patient's data and settings.
Same Patient	A patient is admitted by a monitor, to which MPM/T1 used in another monitor for monitoring this patient is connected.

4.7.2 Transferring Patient Data via Storage Medium

4.7.2.1 Transferring Data from the Monitor to Storage Medium

1. Select **[Main Menu]** → **[Patient Setup >>]**.
2. Select **[Transfer to Storage Medium]**. In the popup menu, you can:
 - ◆ Select **[Ok]** to transfer the patient data, or
 - ◆ Select **[Cancel]** to exit the menu.
3. Wait until the following message appears: **[Transfer to storage medium successful. Please remove the CF storage card.]** or **[Transfer to storage medium successful. Please remove the USB drive.]**.
4. Remove the CF storage card or USB drive from patient monitor.

4.7.2.2 Transferring Data from the Storage Medium to the Monitor

1. Connect the storage medium to the destination monitor.
2. In the popup menu, you can:
 - ◆ Select **[Transfer]** to transfer the patient data to the monitor, or
 - ◆ Select **[Cancel Transfer]** to cancel the operation of transferring patient data.
 - ◆ Select **[Unload CF Storage Card]** or **[Unload USB Drive]** to not transfer the patient data and to unload the card or USB drive.
3. After you select **[Transfer]**, in the popup menu you can further select the patient data contents that need to be transferred. **[Patient Demographics]** must be selected. After **[Ok]** is selected, the monitor compares the patient information stored in both the storage medium and monitor and deals with the patient data based on the following.
 - ◆ Different Patients: The monitor erases all the current patient data, transfers the patient data from the storage medium, and loads the configuration according to the patient category.
 - ◆ Same Patient: In the popup dialog box, you can:
 - ◆ Select **[Yes]** to merge the patient data in the monitor and storage medium.
 - ◆ Select **[No]** to erase all the current patient data in the monitor and to transfer the patient data from the storage medium.

Wait until the following message appears: **[Transfer from storage medium successful.]**.

 **WARNING**

- **The USB drive you use may have write-protect function. In this case, please make sure the USB drive for data transfer is in read/write mode.**
 - **Do not remove the storage medium during data transfer process. Otherwise, data files may be damaged.**
-
-

4.8 Connecting to a Central Monitoring System

If your patient monitor is connected to a central monitoring system (CMS):

- All patient information, measurement data and settings on the patient monitor can be transferred to the CMS.
- All patient information, measurement data and settings can be displayed simultaneously on the patient monitor and CMS. For some functions such as editing patient information, admitting a patient, discharging a patient, starting/stopping NIBP measurements, etc., bi-directional control can be achieved between your patient monitor and the CMS.

For details, refer to the CMS's instructions for use.

5 Managing Configurations

5.1 Introduction

When performing continuous monitoring on a patient, the clinical professional often needs to adjust the monitor's settings according to the patient's condition. The collection of all these settings is called a configuration. Allowing you to configure the monitor more efficiently, the monitor offers different sets of configuration to suit different patient categories and departments. You can change some settings from a certain set of configuration and then save the changed configuration as a user configuration.

The default configurations provided for your monitor are department-oriented. You can choose either from:

- General
- OR
- ICU
- NICU
- CCU

Each department has three different sets of configurations tailored for adult, pediatric and neonatal patients.

WARNING

- **The configuration management function is password protected. The configuration management tasks must be performed by clinical professionals.**
-

The system configuration items can be classified as:

- Parameter configuration items

These items relates to parameters, e.g., waveform gain, alarm switch, alarm limits..

- Conventional configuration items


These items define how the monitor works, e.g., screen layout, record, print and alarm settings.

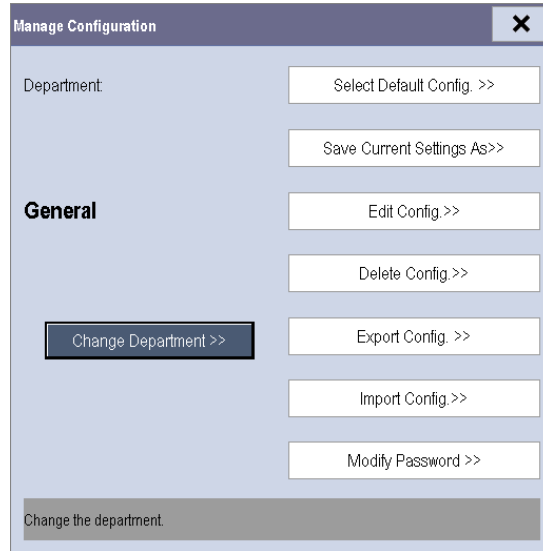
- User maintenance items

These items relates to user maintenance settings, e.g., unit setup, time format and data format.

For the important configuration items and their default values and user maintenance items, see appendix **Configuration Default Information**.


5.2 Entering the [Manage Configuration] Menu

1. Press the  hardkey on the monitor's front to enter the main menu.
2. Select [**Maintenance >>**]→[**Manage Configuration >>**]. Enter the required password and then select [**Ok**].



5.3 Changing Department

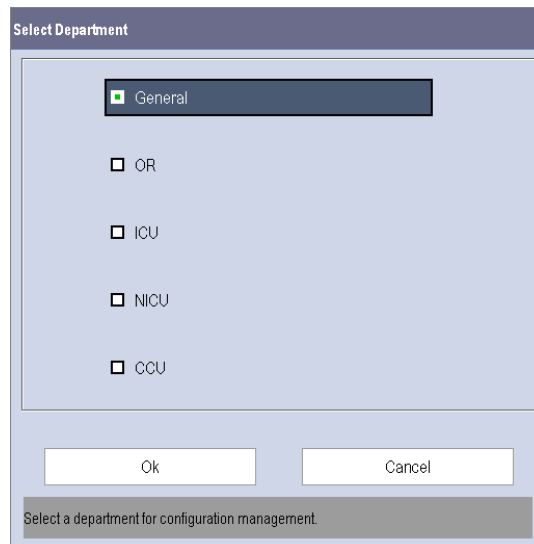If the current department configuration is not the one you want to view, you can select [**Change Department >>**] in the [**Manage Configuration**] menu and then choose the one you want for viewing as shown below.



NOTE

- **Changing the department will delete all current user configurations. Please act with caution.**

5.4 Setting Default Configuration

The monitor will load the pre-set default configuration in the following cases.

- The patient monitor restarts after quitting over 120 seconds.
- A patient is admitted.
- A patient is discharged.
- Patient data is cleared.
- Patient category is changed.

To set default configuration:

1. Select [**Select Default Config. >>**] in the [Manage Configuration] menu.
2. In the [**Select Default Config.**] menu, select [**Load the Latest Config.**] or [**Load Specified Config.**].

When you select [**Load Specified Config.**], the configuration (adult, pediatric or neonate) to be restored is subject to the patient category. This configuration can be either factory configuration or saved user configuration. Take adult as an example, select [**Default Adu Config.**] and toggle between [**Defaults**] or user configuration(s).

NOTE

- **To know what configuration is restored when the patient monitor starts, enter the main screen to check the prompt information at the lower part of the screen (displayed for about 10 seconds).**
-

5.5 Saving Current Settings

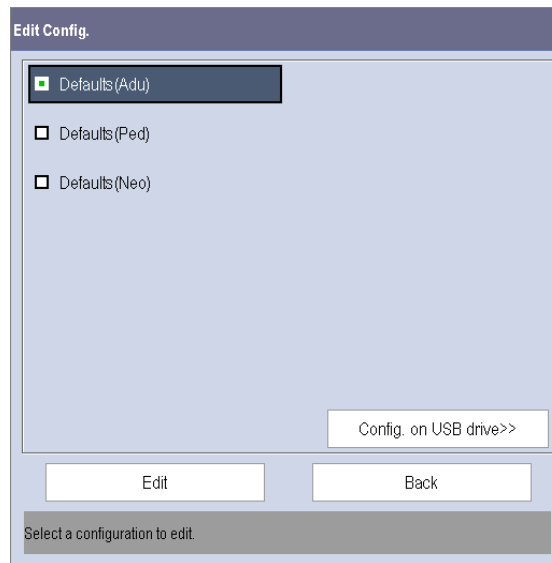
Current settings can be saved as user configuration. Up to 10 user configurations can be saved.

To save current settings:


1. Select [**Save Current Settings As >>**] in the [**Manage Configuration**] menu.
2. In the popup dialog box, enter the configuration name and then select [**Ok**].

5.6 Editing Configuration

1. Select [**Edit Config. >>**] in the [**Manage Configuration**] menu. The following menu appears.



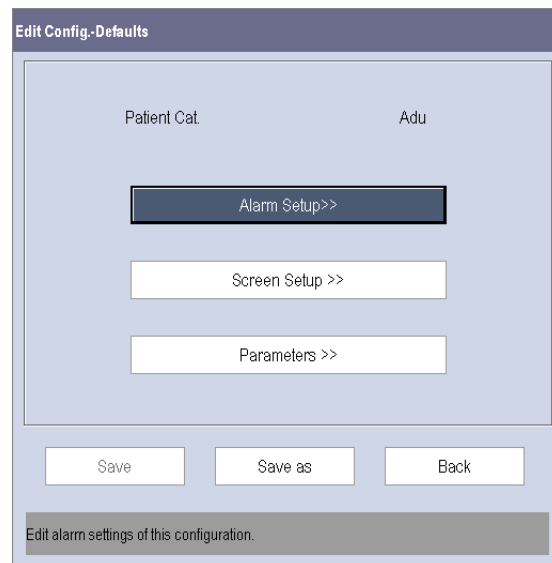2. The popup menu shows the existing configurations on the monitor. Selecting [**Config. on USB drive >>**] will show the existing configurations on the USB drive. Select the desired configuration and then select the [**Edit**] button. The following menu appears.



3. Select [**Alarm Setup >>**], [**Screen Setup >>**] or [**Parameter >>**] to enter the corresponding menu in which settings can be changed. The changed items of alarm setup will be marked in red.
4. You can select [**Save**] or [**Save as**] to save the changed configuration. Select [**Save**] to overwrite the original configuration. Select [**Save as**] to save the changed configuration in another name.

5.7 Deleting a Configuration

1. Select [**Delete Config. >>**] in the [**Manage Configuration**] menu.
2. The popup menu shows the existing user configurations on the monitor. Selecting [**Config. on USB drive >>**] will show the existing user configurations on the USB drive. Select the user configurations you want to delete and then select [**Delete**].
3. Select [**Yes**] in the popup.

5.8 Transferring a Configuration

When installing several monitors with identical user configuration it is not necessary to set each unit separately. An USB drive may be used to transfer the configuration from monitor to monitor.

To export the current monitor's configuration:

1. Connect the USB drive to the monitor's USB port.
2. Select [**Export Config. >>**] in the [**Manage Configuration**] menu.
3. In the [**Export Config.**] menu, select the configurations and [**User Maintenance Settings**] to export. Then select the [**Export**] button. A status message will report completion of the transfer.

To import the configuration on the USB drive to the monitor:

1. Connect the USB drive to the monitor's USB port.
2. Select [**Import Config. >>**] in the [**Manage Configuration**] menu.
3. In the [**Import Config.**] menu, select the configurations and [**User Maintenance Settings**] to import. Then select the [**Import**] button. A status message will report completion of the transfer.

5.9 Loading a Configuration

You may make changes to some settings during operation. However, these changes or the pre-selected configuration may not be appropriate for the newly admitted patient. Therefore, the monitor allows you to load a desired configuration so as to ensure that all the settings are appropriate for your patient.

To load a configuration,

1. Select [**Load Configuration >>**] from the main menu.
2. The popup menu shows the existing configurations on the monitor. Selecting [**Config. on USB drive >>**] will show the existing configurations on the USB drive.
3. Select a desired configuration.
4. Select [View] to view the configuration details. In the popup menu, you can select [**Alarm Setup >>**], [**Screen Setup >>**] or [**Parameter >>**] to view the corresponding contents. The alarm setup items which are different than those currently used are marked in red.
5. Select [**Load**] to load this configuration.

NOTE

- **The monitor may configure some settings by default when you load a configuration of different version with current configuration.**
-

5.10 Restoring the Latest Configuration Automatically

During operation, you may make changes to some settings. However, these changes may not be saved as user configuration. To prevent the changes from losing in case of a sudden power failure, the patient monitor stores the configuration in real time. The saved configuration is the latest configuration.

The monitor restore the latest configuration if restarts within 60 seconds after the power failure. And it will restore the default configuration rather than the latest configuration if restarts 120 seconds later after the power failure. The monitor may load either the latest configuration or the default configuration if restarts from 60-120 seconds after the power failure.

5.11 Modifying Password

To modify the password for accessing the **[Manage Configuration]** menu,

1. Select **[Modify Password >>]** in the **[Manage Configuration]** menu.
2. Input a new password in the popup menu.
3. Select **[Ok]**.

6 User Screens

6.1 Tailoring Your Screens

You can tailor your patient monitor's screens by setting:

- Waveform sweep mode
- Wave line size
- The color in which each measurement's numerics and waveform are displayed
- The parameter to be monitored.

Changing some settings may be hazardous. Therefore, those settings are password-protected and can be modified by authorized personnel only. Once change is made, those who use the patient monitor should be notified.

6.1.1 Setting the Waveform Sweep Mode

1. Select **[Main Menu]**→**[Screen Setup >>]**.
2. Select **[Sweep Mode]** and toggle between **[Refresh]** and **[Scroll]**.
 - ◆ **[Refresh]**: The waveforms keep stationary, being refreshed from left to right by a moving "erase bar".
 - ◆ **[Scroll]**: The waveforms move from the right to the left with time passing by.

6.1.2 Changing the Wave Line Size

1. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password.
2. Select **[Others >>]**.
3. Select **[Wave Line]** and toggle between **[Thick]**, **[Mediate]** and **[Thin]**.

6.1.3 Changing Measurement Colors

1. Select **[Main Menu]**→**[Screen Setup >>]**→**[Measurement Color Setup >>]**.
2. Select the color box next to your desired measurement and then select a color from the popup menu.

6.1.4 Changing Screen Layout

Select the **[Screens]** QuickKey, or **[Main Menu]**→**[Screen Setup >>]**→**[Screen Layout >>]** to enter the **[Screens]** menu.



- You can choose the desired screen type in the **[Choose Screen]** window.
- You can select the parameters and waveforms you want to view in the **[Screen Setup]** window. For details, please refer to the section **Setting the Screen**.
- You can select the parameters you want to view on big numerics screen in the **[Big Numerics Screen Setup]** window.
- You can switch on or off the connected parameter modules in the **[Parameters Switch]** window. If a parameter module is switched off, parameter values and waveforms will not display on the screen.

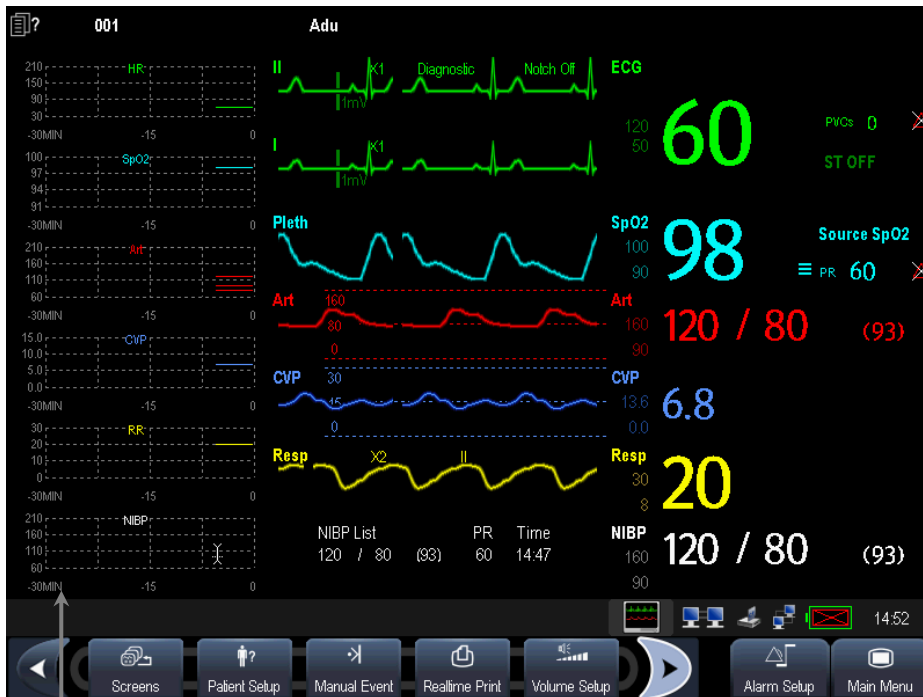
6.2 Viewing Minitrends

6.2.1 Having a Split-Screen View of Minitrends

You can split the normal screen so that one part of the screen, on the left hand side, continuously shows graphic minitrends beside waveforms as shown in the figure below.


To have a split-screen view of minitrends, you can:

- Select [**Minitrends**] QuickKey, or
- Select [**Screens**] QuickKey → [**Minitrends Screen**] → , or
- Select [**Main Menu**] → [**Screen Setup >>**] → [**Screen Layout >>**] → [**Minitrends Screen**] → .



Minitrend View

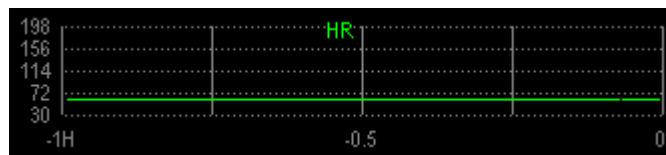The split-screen view provides minitrends for multiple parameters. In each field, the label, scale and time are respectively displayed at the top, left, and bottom as shown below.



6.2.2 Setting Minitrends

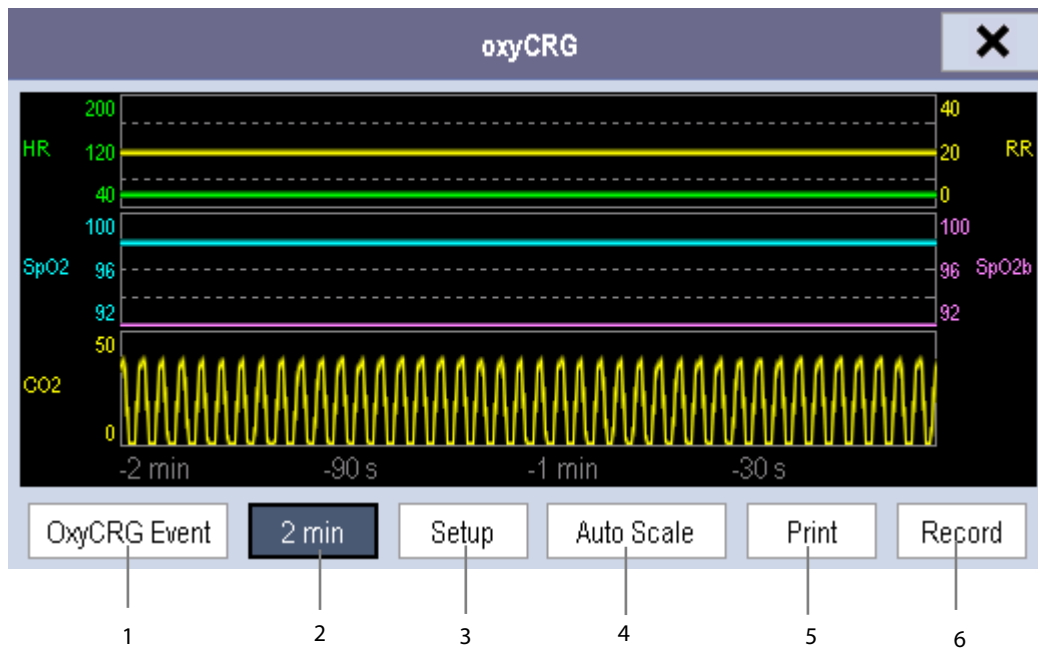
Select the minitrends area. From the pop-up [**Minitrend Setup**] menu, you can:

- Select the parameters to be displayed, or
- Select [**Minitrend Length**] and then select the appropriate setting.

6.3 Viewing OxyCRG

To have a split screen view of OxyCRG, you can:

- Select [**OxyCRG**] QuickKey, or
- Select [**Screens**] QuickKey→[**OxyCRG Screen**]→**X**, or.
- Select [**Main Menu**]→[**Screen Setup >>**]→[**Screen Layout >>**]→[**OxyCRG Screen**]→**X**.



The split-screen view covers the lower part of the waveform area and shows HR trend, SpO₂ trend, SpO₂b trend, RR trend and a compressed wave (CO₂ wave or Resp wave). At the bottom, there are controls:

1. OxyCRG Event
You can enter the [**Review**] menu by selecting the [**OxyCRG Event**] button.
2. Trend length list box
In the trend length list box, you can select [**1 min**], [**2 min**], [**4 min**], or [**8 min**].
3. Setup
Select [**Setup**] button to enter [**Setup**] menu, in which you can select the parameters for display, the time length to be saved before and after an event, and the scale of the graphic trends and waveform.
4. Auto Scale
Select [**Auto Scale**] button, and the system automatically adjusts the scaling.
5. Print
Select [**Print**] to print out the realtime OxyCRG.
6. Record
Through this button, you can print out the currently displayed OxyCRG trends by the recorder.

6.4 Viewing Other Patients

6.4.1 Care Group

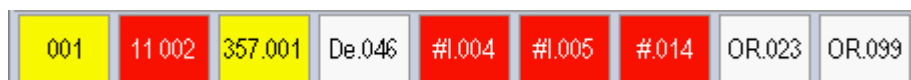
You can select other patient monitors (including telemetry) connected to the same LAN into a Care Group. This lets you:

- View information on the monitor screen from another bed in the same Care Group.
- Be notified of physiological and technical alarm conditions at the other beds in the same Care Group.

You can select up to 10 patient monitors for BeneView T5 and 16 for BeneView T8 in a Care Group. To have a Care Group:

1. Open the **[View Other Patient]** window by:
 - ◆ Selecting **[Others]** QuickKey, or
 - ◆ Selecting **[Screens]** QuickKey→**[View Others Screen]**→**X**, or
 - ◆ Selecting **[Main Menu]**→**[Screen Setup >>]**→**[Screen Layout >>]**→**[View Others Screen]**→**X**.
2. Select **[Setup]** in the **[View Other Patient]** window.
3. Select the desired patient monitors from the **[Connected Monitor List]**, and then select the **X** button. The selected patient monitors constitute a Care Group.

6.4.2 Viewing the Care Group Overview Bar



The Care Group overview bar locates at the bottom of the **[View Other Patient]** window. In the overview bar, the department and bed label for any Care Group beds are displayed. For telemetry, # is displayed before the department label. The color in which a Care Group bed appears matches its status:

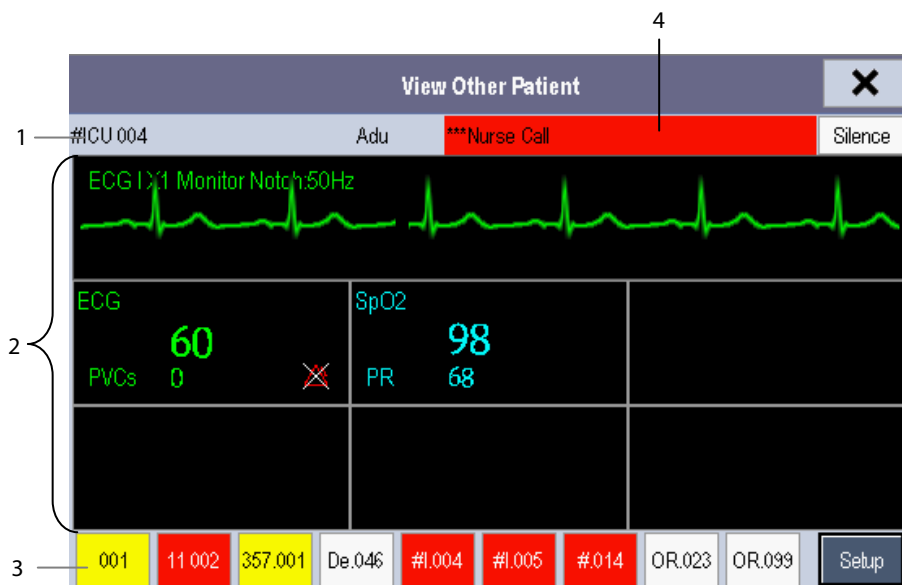
- Red: indicates the bed is giving high-level physiological alarms or the telemetry is giving alarm, such as nurse call or event.
- Yellow: indicates the bed is giving medium- or low-level physiological alarms, or medium-level technical alarms.
- Blue: indicates the bed is giving low-level technical alarms.
- Grey: indicates the bed fails to be networked or stays in the standby mode.

You can view a Care Group bed's alarms by selecting it from the care group, and as well you can select the **[View This Patient]** button to view this bed in the **[View Other Patient]** window.

For more details about Care Group alarms, refer to the **Alarms** chapter.

6.4.3 Understanding the View Other Patient Window

When you first open the [View Other Patient] window, the patient monitor automatically selects a monitor from the network to display in the [View Other Patient] window.



The [View Other Patient] window covers the lower part of the waveform area and consists of:

1. Information Area: shows the patient information (including department, bed number, patient name, etc.), and network status symbol.
2. View Area: shows physiological waveforms and parameters. You can switch a waveform area to a parameter area by selecting your desired waveform area and then selecting [Switch to Parameter Area], or switch a parameter area to a waveform area by selecting your desired parameter area and then selecting [Switch to Waveform Area].
3. Care Group Overview Bar.
4. Message Area: shows physiological, technical and prompt messages from the currently viewed patient monitor. It also shows the alarm given by the telemetry such as nurse call or event. By selecting this area, you can enter the [Alarm Information List] to view all physiological, technical and prompt messages coming from the currently viewed patient.

Additionally, you can change a waveform or parameter for viewing

- To change a waveform for viewing, select the waveform segment where you want a new waveform to appear and then select the waveform you want from the popup menu.
- To change a parameter for viewing, select the parameter window where you want a new parameter to appear and then select the parameter you want from the popup menu.

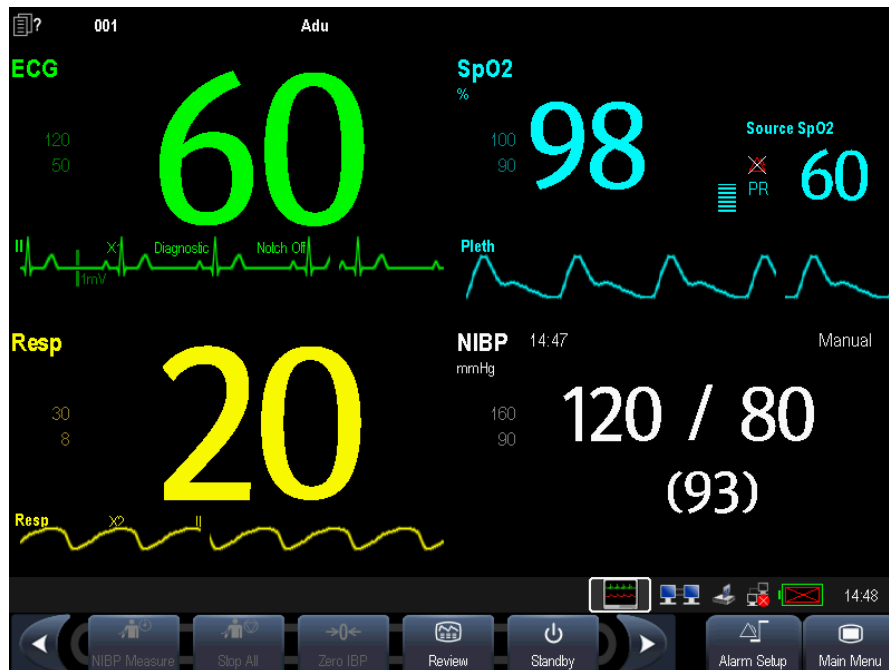
WARNING

- The data presented in the [View Other Patient] window have delay. Do not rely on this window for realtime data.
-

6.5 Understanding the Big Numerics Screen

To enter the big numerics screen:

1. Select the **[Screens]** QuickKey, or **[Main Menu]**→**[Screen Setup >>]**→**[Screen Layout >>]**.
2. Select **[Big Numerics]**→**X**.



You can select your desired parameters to display in this screen: select the **[Screens]** QuickKey→**[Big Numerics Screen Setup]** and then select the parameters you want. For parameters having a waveform, the waveform will also be displayed.

7 Alarms

Alarms, triggered by a vital sign that appears abnormal or by technical problems of the patient monitor, are indicated to the user by visual and audible alarm indications.

WARNING

- **A potential hazard can exist if different alarm presets are used for the same or similar equipment in any single area, e.g. an intensive care unit or cardiac operating room.**
 - **If your patient monitor is connected to a CMS, remote suspension, inhibition, silence and reset of monitor alarms via the CMS may cause a potential hazard.**
-
-

7.1 Alarm Categories

By nature, the patient monitor's alarms can be classified into three categories: physiological alarms, technical alarms and prompt messages.

1. Physiological alarms

Physiological alarms, also called patient status alarms, are triggered by a monitored parameter value that violates set alarm limits or an abnormal patient condition. Physiological alarm messages are displayed in the physiological alarm area.

2. Technical alarms

Technical alarms, also called system status alarms, are triggered by a device malfunction or a patient data distortion due to improper operation or mechanical problems. Technical alarm messages are displayed in the technical alarm area.

Apart from the physiological and technical alarm messages, the patient monitor shows some messages telling the system status or patient status. Messages of this kind are included into the prompt message category and usually displayed in the prompt information area. Some prompt messages that indicate the arrhythmia events are displayed in the physiological alarm area. For some measurements, their related prompt messages are displayed in their respective parameter windows.

7.2 Alarm Levels

By severity, the patient monitor's alarms can be classified into three categories: high level, medium level and low level..

	Physiological alarms	Technical alarms
High level	Indicate that your patient is in a life threatening situation, such as Asystole, Vfib/Vtac and so forth, and an emergency treatment is demanded.	Indicate a severe device malfunction or an improper operation, which could make it possible that the monitor cannot detect critical patient status and thus threaten the patient's life, such as low battery.
Medium level	Indicate that your patient's vital signs appear abnormal and an immediate treatment is required.	Indicate a device malfunction or an improper operation, which may not threaten the patient's life but may compromise the monitoring of vital physiological parameters.
Low level	Indicate that you patient's vital signs appear abnormal and an immediate treatment may be required.	Indicate a device malfunction or an improper operation, which may compromise a certain monitoring function but will not threaten the patient's life.

7.3 Alarm Indicators

When an alarm occurs, the patient monitor will indicate it to the user through visual or audible alarm indications.

- Alarm lamp
- Alarm message
- Flashing numeric
- Audible alarm tones

7.3.1 Alarm Lamp

If a technical alarm occurs, the technical alarm lamp will turn blue. If a technical alarm or physiological alarm occurs, the alarm lamp will flash. The flashing color and frequency match the alarm level as follows:

- High level alarms: the lamp quickly flashes red.
- Medium level alarms: the lamp slowly flashes yellow.
- Low level physiological alarms: the lamp turns yellow without flashing.
- Low level technical alarms: the lamp does not light.

7.3.2 Alarm Message

When an alarm occurs, an alarm message will appear in the technical or physiological alarm area. For physiological alarms, the asterisk symbols (*) before the alarm message match the alarm level as follows:

- High level alarms: ***
- Medium level alarms: **
- Low level alarms: *

Additionally, the alarm message uses different background color to match the alarm level:

- High level alarms: red
- Medium level alarms: yellow
- Low level physiological alarms: yellow
- Low level technical alarms: blue

You can view the alarm messages by selecting the physiological or technical alarm area.

7.3.3 Flashing Numeric

If an alarm triggered by an alarm limit violation occurs, the numeric of the measurement in alarm will flash every second, and the corresponding alarm limit will also flash at the same frequency indicating the high or low alarm limit is violated.

7.3.4 Audible Alarm Tones

This monitor has three choices of alarm tones and patterns: ISO, Mode 1 and Mode 2. For each pattern, the alarm tones identify the alarm levels as follows:





- ISO pattern:
 - ◆ High level alarms: triple+double+triple+double beep.
 - ◆ Medium level alarms: triple beep.
 - ◆ Low level alarms: single beep.
- Mode 1:
 - ◆ High level alarms: high-pitched single beep.
 - ◆ Medium level alarms: double beep.
 - ◆ Low level alarms: low-pitched single beep.
- Mode 2:
 - ◆ High level alarms: high-pitched triple beep.
 - ◆ Medium level alarms: double beep.
 - ◆ Low level alarms: low-pitched single beep.

NOTE

-
- **When multiple alarms of different levels occur simultaneously, the patient monitor will select the alarm of the highest level and give visual and audible alarm indications accordingly.**
-

7.3.5 Alarm Status Symbols

Apart from the aforementioned alarm indicators, the patient monitor still uses the following symbols telling the alarm status:

-  indicates alarms are paused.
-  indicates alarm sound is silenced.
-  indicates the alarm sound is turned off.
-  indicates individual measurement alarms are turned off or the system is in alarm off status.

7.4 Alarm Tone Configuration


7.4.1 Setting the Minimum Alarm Volume

1. Select [**Main Menu**]→[**Maintenance >>**]→[**User Maintenance >>**]→enter the required password.
2. Select [**Alarm Setup >>**] to enter the [**Alarm Setup**] menu.
3. Select [**Minimum Alarm Volume**] and toggle between 0 and 10.

The minimum alarm volume refers to the minimum value you can set for the alarm volume, which is not affected by user or factory default configurations. The setting of minimum alarm volume remains unchanged when the patient monitor shuts down and restarts.

7.4.2 Changing the Alarm Volume

1. Select the [**Volume Setup**] QuickKey or the [**Alarm Setup**] QuickKey→[**Others**], or [**Main Menu**]→[**Alarm Setup >>**]→[**Others**].
2. Select the appropriate volume from [**Alm Volume**]: X-10, in which X is the minimum volume, depending on the set minimum alarm volume, and 10 the maximum volume.
3. Select [**High Alarm Volume**] to set the volume of the high priority alarm as [**Alm Volume+0**], [**Alm Volume+1**] or [**Alm Volume+2**].
4. Select [**Reminder Vol**] to set the volume of the reminder tone as [**High**], [**Med**] or [**Low**].

When alarm volume is set to 0, the alarm sound is turned off and a  symbol appears on the screen.

7.4.3 Setting the Interval between Alarm Sounds

You cannot change the interval between alarm tones if you choose mode 1 or 2 as your desired alarm tone pattern. For these two patterns, the interval between alarm tones identifies the alarm levels as follows:

- Mode 1:
 - ◆ Interval between high level alarm tones: continuously.
 - ◆ Interval between medium level alarm tones: 5 s.
 - ◆ Interval between low level alarm tones: 20 s.
- Mode 2:
 - ◆ Interval between high level alarm tones: 1 s.
 - ◆ Interval between medium level alarm tones: 5 s.
 - ◆ Interval between low level alarm tones: 20 s.

If you choose the ISO pattern, you can change the interval between alarm tones. To change the interval between alarm tones:

1. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password.
2. Select **[Alarm Setup >>]** to enter the **[Alarm Setup]** menu.
3. Select **[High Alarm Interval (s)]**, **[Med Alarm Interval (s)]** and **[Low Alarm Interval (s)]** in turn and then select the appropriate settings.



WARNING

- **When the alarm sound is switched off, the patient monitor will give no audible alarm tones even if a new alarm occurs. Therefore the user should be very careful about whether to switch off the alarm sound or not.**
 - **Do not rely exclusively on the audible alarm system for patient monitoring. Adjustment of alarm volume to a low level may result in a hazard to the patient. Always keep the patient under close surveillance.**
-

7.4.4 Changing the Alarm Tone Pattern

To change the alarm tone pattern:

1. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password.
2. Select **[Alarm Setup >>]** to enter the **[Alarm Setup]** menu.
3. Select **[Alarm Sound]** and toggle between **[ISO]**, **[Mode 1]** and **[Mode 2]**.

User or factory default configurations exert no impact on the setup of alarm tone pattern. The alarm tone pattern remains unchanged after the monitor restarts.

7.4.5 Setting the Reminder Tones

When the alarm volume is set to zero, or the alarm tone is silenced or turned off, the patient monitor issues a periodical reminder tone.

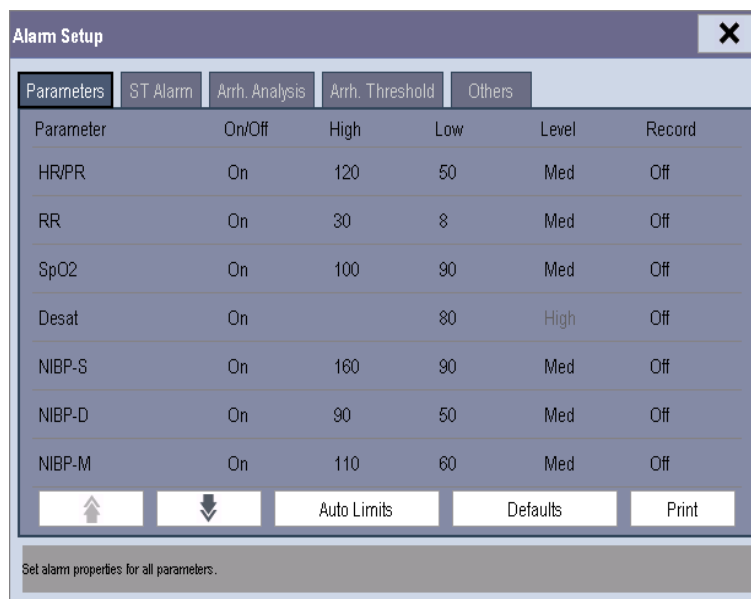
1. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password.
2. Select **[Alarm Setup >>]** to enter the **[Alarm Setup]** menu.
 - ◆ To switch the reminder tones on or off, select **[Reminder Tones]** and toggle between **[On]** and **[Off]**.
 - ◆ To set the interval between reminder tones, select **[Reminder Interval]** and toggle between **[1min]**, **[2min]** and **[3min]**.

In addition, you can set the volume of alarm reminder tones. To set the volume of alarm reminder tones, select **[Main Menu]**→**[Alarm Setup >>]**→**[Others]** or the **[Alarm Setup]** QuickKey→**[Others]**. Then, select **[Reminder Vol]** and toggle between **[High]**, **[Medium]** and **[Low]**.

7.5 Understanding the Alarm Setup Menu

Select the **[Alarm Setup]** QuickKey or **[Main Menu]**→**[Alarm Setup >>]** to enter the **[Alarm Setup]**, where you can:

- Set alarm properties for all parameters.
- Change ST alarm settings.
- Change arrhythmia alarm settings.
- Set the threshold for some arrhythmia alarms.
- Change other settings.



The screenshot shows the 'Alarm Setup' window with a close button (X) in the top right corner. Below the title bar are five tabs: 'Parameters', 'ST Alarm', 'Arrh. Analysis', 'Arrh. Threshold', and 'Others'. The 'Parameters' tab is active, displaying a table with the following data:

Parameter	On/Off	High	Low	Level	Record
HR/PR	On	120	50	Med	Off
RR	On	30	8	Med	Off
SpO2	On	100	90	Med	Off
Desat	On		80	High	Off
NIBP-S	On	160	90	Med	Off
NIBP-D	On	90	50	Med	Off
NIBP-M	On	110	60	Med	Off

At the bottom of the table are five buttons: an up arrow, a down arrow, 'Auto Limits', 'Defaults', and 'Print'. Below the buttons is a footer text: 'Set alarm properties for all parameters.'

Please refer to the **ECG** section for how to change ST alarm settings, how to change arrhythmia alarm settings and how to set the threshold for some arrhythmia alarms.

7.5.1 Setting Alarm Properties for All Parameters

In the main menu, select **[Alarm Setup >>]** → **[Parameters]**. You can review and set alarm limits, alarm switches, alarm level and alarm recordings for all parameters.

When a measurement alarm occurs, automatic recording of all the measurement numerics and related waveforms is possible when the measurement's **[On/Off]** and **[Record]** are set on.

WARNING

- **Make sure that the alarm limits settings are appropriate for your patient before monitoring.**
- **Setting alarm limits to extreme values may cause the alarm system to become ineffective. For example, High oxygen levels may predispose a premature infant to retrolental fibroplasia. If this is a consideration do NOT set the high alarm limit to 100%, which is equivalent to switching the alarm off.**

7.5.2 Adjusting Alarm Limits Automatically

The monitor can automatically adjust alarm limits according to the measured vital signs, using the auto limits function. When auto limits are selected, the monitor calculates safe auto limits based on the latest measured values.

To get accurate auto alarm limits, you need to collect a set of measured vital signs as a baseline. Then, in the main menu, select **[Alarm Setup >>]** → **[Parameters]** → **[Auto Limits]** → **[Ok]**. The monitor will create new alarm limits based on the measured values.

Before applying these automatically created alarm limits, confirm if they are appropriate for your patient in the mass alarm setup menu. If not, you can adjust them manually. These alarm limits will remain unchanged until you select auto limits again or adjust them manually.

The monitor calculates the auto limits based on the following rules.

Module	Parameter	Low alarm limit		High alarm limit		Auto alarm limits range
		Adult/ pediatric	Neonate	Adult/ pediatric	Neonate	
ECG	HR/PR	HR × 0.8 or 40bpm (whichever is greater)	(HR – 30) or 90bpm (whichever is greater)	HR × 1.25 or 240bpm (whichever is smaller)	(HR + 40) or 200bpm (whichever is smaller)	Adult/pediatric: 35 to 240 Neonate: 55 to 225
Resp	RR	RR × 0.5 or 6 rpm (whichever is greater)	(RR – 10) or 30 rpm (whichever is greater)	RR × 1.5 or 30 rpm (whichever is smaller)	(RR + 25) or 85 rpm (whichever is smaller)	Adult/pediatric: 6 to 55 Neonate: 10 to 90
SpO ₂	SpO ₂	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the measurement range

Module	Parameter	Low alarm limit		High alarm limit		Auto alarm limits range
		Adult/ pediatric	Neonate	Adult/ pediatric	Neonate	
NIBP	NIBP-S	$(SYS \times 0.68 + 10)$ mmHg	$(SYS - 15)$ or 45mmHg (whichever is greater)	$(SYS \times 0.86 + 38)$ mmHg	$(SYS + 15)$ or 105mmHg (whichever is smaller)	Adult: 45 to 270 Pediatric: 45 to 185 Neonate: 35 to 115
	NIBP-D	$(Dia \times 0.68 + 6)$ mmHg	$(Dia - 15)$ or 20mmHg (whichever is greater)	$(Dia \times 0.86 + 32)$ mmHg	$(Dia + 15)$ or 80mmHg (whichever is smaller)	Adult: 25 to 225 Pediatric: 25 to 150 Neonate: 20 to 90
	NIBP-M	$(Mean \times 0.68 + 8)$ mmHg	$(Mean - 15)$ or 35mmHg (whichever is greater)	$(Mean \times 0.86 + 35)$ mmHg	$(Mean + 15)$ or 95mmHg (whichever is smaller)	Adult: 30 to 245 Pediatric: 30 to 180 Neonate: 25 to 105
Temp	T1	$(T1 - 0.5)^\circ C$	$(T1 - 0.5)^\circ C$	$(T1 + 0.5)^\circ C$	$(T1 + 0.5)^\circ C$	1 to 49 °C
	T2	$(T2 - 0.5)^\circ C$	$(T2 - 0.5)^\circ C$	$(T2 + 0.5)^\circ C$	$(T2 + 0.5)^\circ C$	1 to 49 °C
	TD	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the measurement range
IBP: ART/ Ao/ UAP/ BAP/ FAP/ LV/ P1-P4 (Arterial pressure)	IBP-S	$(SYS \times 0.68 + 10)$ mmHg	$(SYS - 15)$ or 45mmHg (whichever is greater)	$(SYS \times 0.86 + 38)$ mmHg	$(SYS + 15)$ or 105mmHg (whichever is smaller)	Adult: 45 to 270 Pediatric: 45 to 185 Neonate: 35 to 115
	IBP-D	$(Dia \times 0.68 + 6)$ mmHg	$(Dia - 15)$ or 20mmHg (whichever is greater)	$(Dia \times 0.86 + 32)$ mmHg	$(Dia + 15)$ or 80mmHg (whichever is smaller)	Adult: 25 to 225 Pediatric: 25 to 150 Neonate: 20 to 90
	IBP-M	$(Mean \times 0.68 + 8)$ mmHg	$(Mean - 15)$ or 35mmHg (whichever is greater)	$(Mean \times 0.86 + 35)$ mmHg	$(Mean + 15)$ or 95mmHg (whichever is smaller)	Adult: 30 to 245 Pediatric: 30 to 180 Neonate: 25 to 105
PA	IBP-S	$SYS \times 0.75$	$SYS \times 0.75$	$SYS \times 1.25$	$SYS \times 1.25$	3 to 120mmHg
	IBP-D	$Dia \times 0.75$	$Dia \times 0.75$	$Dia \times 1.25$	$Dia \times 1.25$	
	IBP-M	$Mean \times 0.75$	$Mean \times 0.75$	$Mean \times 1.25$	$Mean \times 1.25$	

Module	Parameter	Low alarm limit		High alarm limit		Auto alarm limits range
		Adult/ pediatric	Neonate	Adult/ pediatric	Neonate	
IBP: CVP/ ICP/ LAP/ RAP/ UVP/ P1-P4 (Venous pressure)	IBP-M	Mean × 0.75	Mean × 0.75	Mean × 1.25	Mean × 1.25	3 to 40mmHg
CO ₂	EtCO ₂	0 to 32mmHg: remains the same	0 to 32mmHg: remains the same	0 to 32mmHg: remains the same	0 to 32mmHg: remains the same	Same as the measurement range
		32 to 35mmHg: 29mmHg	32 to 35mmHg: 29mmHg	32 to 35mmHg: 41mmHg	32 to 35mmHg: 41mmHg	
		35 to 45mmHg: (etCO ₂ -6) mmHg	35 to 45mmHg: (etCO ₂ -6) mmHg	35 to 45mmHg: (etCO ₂ +6) mmHg	35 to 45mmHg: (etCO ₂ +6) mmHg	
		45 to 48mmHg:39 mmHg	45 to 48mmHg:39 mmHg	45 to 48mmHg:51 mmHg	45 to 48mmHg:51 mmHg	
		>48mmHg: remains the same	>48mmHg: remains the same	>48mmHg: remains the same	>48mmHg: remains the same	
	FICO ₂	N/A	N/A	Same as the default alarm limit	Same as the default alarm limit	Same as the measurement range
	awRR	awRR × 0.5 or 6 rpm (whichever is greater)	(awRR – 10) or 30 rpm (whichever is greater)	awRR × 1.5 or 30 rpm (whichever is smaller)	(awRR+25) or 85 rpm (whichever is smaller)	Adult/pediatric: 6 to 55 Neonate: 10 to 90
AG	EtCO ₂ (AG)	Same as CO ₂ module				
	FICO ₂ (AG)					

Module	Parameter	Low alarm limit		High alarm limit		Auto alarm limits range
		Adult/ pediatric	Neonate	Adult/ pediatric	Neonate	
	awRR	awRR × 0.5 or 6 rpm (whichever is greater)	awRR – 10 or 30 rpm (whichever is greater)	awRR × 1.5 or 30 rpm (whichever is smaller)	awRR+25 or 85 rpm (whichever is smaller)	Adult/pediatric: 6 to 55 Neonate: 10 to 90
	FiAA/ EtAA	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the measurement range
	FiO ₂ / EtCO ₂	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the measurement range
	FiN ₂ O/ EtN ₂ O	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the default alarm limit	Same as the measurement range
C.O.	BT	Adult: (BT – 1)°C	N/A	Adult: (BT – 1)°C	N/A	Same as the measurement range
ICG	C.I.	N/A				
	TFC					
RM	RR(RM)	awRR × 0.5 or 6 rpm (whichever is greater)	N/A	awRR × 1.5 or 30 rpm (whichever is smaller)	N/A	Adult/pediatric: 6 to 55 Neonate: 10 to 90
	PEEP	(PEEP – 5) cmH ₂ O	N/A	(PEEP + 5) cmH ₂ O	N/A	Same as the measurement range
	PIP	(PIP – 10) cmH ₂ O	N/A	(PIP + 10) cmH ₂ O	N/A	Same as the measurement range
	MVe	(MVe – 2) L/min	N/A	(MVe + 2) L/min	N/A	Same as the measurement range
BIS	BIS	N/A				
CCO	CCO/ CCI, EDV/ EDVI, SVR/ SVRI, SV/SVI, RVEF	N/A				
SvO ₂	SvO ₂	(SvO ₂ – 5)%	N/A	(SvO ₂ + 5)%	N/A	Same as the measurement range
	ScvO ₂	(ScvO ₂ – 5)%	N/A	(ScvO ₂ + 5)%	N/A	Same as the measurement range

7.5.3 Setting Alarm Delay Time

You can set the alarm delay time for over-limit alarms of continuously measured parameters. If the alarm-triggered condition disappears within the delay time, the patient monitor will not give the alarm. You can set the **[Alarm Delay]**, in the **[Others]** window of **[Alarm Setup]** menu. Alarm delay is not applied to the following physiological alarms:

- Apnea
- ST alarms
- Arrhythmia alarms
- ECG Weak Signal
- Resp Artifact
- SpO₂ Desat
- No Pulse
- Nellcor SpO₂ over alarm limits
- FiO₂ Shortage
- Measurements of noncontinuous parameters over alarm limits
- HR over alarm limits
- Anesthetic Mixture's MAC>3

You can set **[Apnea Delay]** and **[ST Alarm Delay]** separately in the **[Others]** window of **[Alarm Setup]** menu.


7.5.4 Setting Recording Length

You can change the length of the recorded waveforms. In the **[Others]** window of the **[Alarm Setup]** menu, select **[Recording Length]** and toggle between **[8 s]**, **[16 s]** and **[32 s]**:

- **[8 s]**: 4 seconds respectively before and after the alarm or manual event trigger moment.
- **[16 s]**: 8 seconds respectively before and after the alarm or manual event trigger moment.
- **[32 s]**: 16 seconds respectively before and after the alarm or manual event trigger moment.

7.5.5 Entering CPB Mode

When performing Cardiopulmonary bypass (CPB), you can set the patient monitor to enter CPB mode in order to reduce unnecessary alarms. The CPB mode is activated only if you select **[OR]**. To select **[OR]**,

1. Press the  hardkey on the monitor's front panel to enter **[Main Menu]**.
2. Select **[Maintenance >>]**→**[Manage Configuration >>]**. Enter the required password and then select **[OK]**.
3. Select **[Change Department >>]**→**[OR]**.

In the CPB mode, all the physiological alarms are switched off except for the following alarms.


- BIS-related alarms
- FiCO₂/EtCO₂ too high (for CO₂ module and AG module)
- FiO₂/EtO₂ too high or too low
- FiAA/EtAA too high (AA represents the anaesthetic gas)
- FiN₂O/EtN₂O too high


In CPB mode, [**CPB Mode**] is displayed in the physiological alarm area with red background color.

To enter CPB mode:

Select the [**CPB Mode**] Quickkey or select [**Enter CPB Mode**] in the [**Others**] window of the [**Alarm Setup**] menu. Then select [**Ok**] in the popup dialog box.


7.6 Pausing Alarms

If you want to temporarily prevent alarms from sounding, you can pause alarms by pressing the  hardkey on the monitor's front. When alarms are paused:

- No alarm lamps flash and no alarms are sounded.
- No numeric and alarm limit flash.
- No alarm messages are shown.
- The remaining pause time is displayed in the physiological alarm area.
- The  alarms paused symbol is displayed in the sound symbol area.

The patient monitor enters into the alarm paused status as soon as it is turned on. The alarm pause time is fixed to be 2 minutes.


When the alarm pause time expires, the alarm paused status is automatically cancelled and the alarm tone will sound.


You can also cancel the alarm paused status by pressing the  hardkey.


You can set the alarm pause time as desired. The default alarm pause time is 2 minutes.

1. Select [**Main Menu**] → [**Maintenance >>**] → [**User Maintenance >>**] → enter the required password.
2. Select [**Alarm Setup >>**] → [**Alarm Pause Time**] and then select the appropriate setting from the popup list.

7.7 Switching Off All Alarms

If [Alarm Pause Time] is set to [Permanent], the patient monitor will enter into the alarm off status after the  hardkey is pressed. During the alarm off status,




- As for physiological alarms, no alarm lamps flash and no alarms are sounded.
- As for physiological alarms, no numeric and alarm limit flash.
- No physiological alarm messages are shown.
- [Alarm Off] is displayed in the physiological alarm area with red background.
- As for technical alarms, no alarms are sounded.
- The  alarm off symbol is displayed in the sound symbol area.

You can cancel the alarm off status by pressing the  hardkey.

WARNING

- **Pausing or switching off alarms may result in a hazard to the patient. Please be very careful.**
-

7.8 Silencing the Alarm Sound

You can silence all alarm sounds by pressing the  hardkey on the monitor's front. In that case, the alarm lamp flashing and alarm tones are cleared and  appears in the sound symbol area. After the physiological alarm is silenced,  appears before the alarm message and the numeric and alarm limit still flash. For the performance after the technical alarm is silenced, please refer to the **Silencing Technical Alarms** section.

The alarm silenced status will be automatically cancelled if you switch the patient monitor to other alarm statuses or when a new physiological or technical alarm occurs.

7.9 Latching Alarms

The alarm latching setting for your patient monitor defines how the alarm indicators behave when you do not acknowledge them. When alarms are set to non-latching, their alarm indications end when the alarm condition ends. If you switch alarm latching on, all visual and audible alarm indications last until you acknowledge the alarms, except that the measurement numeric and violated alarm limit stop flashing as soon as the initial alarm condition goes away.


To set alarms to latching or non-latching:


1. Select [**Main Menu**]→[**Maintenance >>**]→[**User Maintenance >>**]→enter the required password.
2. Select [**Alarm Setup >>**].
3. Select [**Latching Alarms**] and toggle between [**High only**], [**Hi&Med**], [**All**] and [**Off**]. If you select [**High only**], only high priority alarm are latched; if you select [**Hi&Med**], both high priority alarms and mediate priority alarms are latched; if you select [**All**], all alarms are latched; if you select [**Off**], the alarm latching is turned off.



NOTE

- **Changing of alarm priority may affect the latching status of corresponding alarm. Please determine if you need to reset the latching status for the specific alarm when you have changed its alarm priority.**
-

7.10 Silencing Technical Alarms

For some technical alarms, their alarm lamp flashing and alarm tones are cleared and the alarm messages change to prompt messages after the  hardkey is pressed. After the patient monitor restores the normal status, the patient monitor can give alarm indications correctly when these alarms are triggered again.

For some technical alarms, all their alarm indications are cleared after the  hardkey is pressed. After the patient monitor restores the normal status, the patient monitor can give alarm indications correctly when these alarms are triggered again.

For some other technical alarms, their alarm lamp flashing and alarm tones are cleared and  appears before the alarm message after the  hardkey is pressed. After the patient monitor restores the normal status, the patient monitor can give alarm indications correctly when these alarms are triggered again.

7.11 Testing Alarms

When the monitor starts up, a selftest is performed. In the meantime, the start-up screens are displayed, and the technical alarm lamp and alarm lamp are lit in blue and yellow respectively. Then, the alarm lamp turns into red, and turns off together with the technical alarm lamp after the system gives a beep. This indicates that the visible and audible alarm indicators are functioning correctly.

For further testing of individual measurement alarms, perform the measurement on yourself (for example SpO₂ or CO₂) or use a simulator. Adjust alarm limits and check that appropriate alarm behaviour is observed.

7.12 When an Alarm Occurs

When an alarm occurs, observe the following steps and take proper actions:

1. Check the patient's condition.
2. Confirm the alarming parameter or alarm category.
3. Identify the source of the alarm.
4. Take proper action to eliminate the alarm condition.
5. Make sure the alarm condition is corrected.

For troubleshooting specific alarms, see appendix **Alarm Messages**.

7.13 Using Care Group Alarms

7.13.1 Care Group Auto Alarms

When a Care Group is set up on your monitor, a flashing symbol will appear beside the QuickKeys area if any monitor in your Care Group, which is not currently viewed by your monitor, is alarming. The alarm symbol is shown as below.



When a patient monitor in the Care Group is disconnected, the flashing symbol is shown as below.



The department and bed label of the alarming monitor appear on the symbols. You can enter the view other patient window by pressing the symbol.

7.13.2 Setting Care Group Alert Tone

When a monitor in the Care Group issue an alarm, your patient monitor prompts you by giving an alert tone. To set the alert tone,

1. In the main menu, select [**Screen Setup >>**]→[**Screen Layout >>**].
2. In the [**Choose Screen**] window, select [**View Others Screen**].
2. In the view other patient window, select [**Setup>>**], and set [**Alert Tone**] to [**Repeat**], [**Once**] or [**Off**].

7.13.3 Silencing Care Group Alarms

You can silence the alarm sound of the currently viewed bed in the view other patient window. This function can be set in the [**Alarm Setup**] menu from the [**User Maintenance**] menu only.

When the alarm silence function for other patients is active and the currently viewed bed is in normal alarm status or alarm sound off status, press the [**Silence**] button in the view other patient window. The currently viewed bed will then enter into the alarm silenced status.

Note that this button is disabled when the currently viewed bed is in alarms off or paused status.

 **WARNING**

-
- **Silencing care group alarms may cause a potential hazard. Please act with caution.**
-

8 Monitoring ECG

8.1 Introduction

The electrocardiogram (ECG) measures the electrical activity of the heart and displays it on the patient monitor as a waveform and a numeric. This patient monitor measures ECG using the MPM module or the BeneView T1. ECG monitoring provides two algorithms:

1. Mindray algorithm

The Mindray algorithm enables 3-, 5- and 12-lead ECG monitoring, ST-segment analysis, arrhythmia analysis and interpretation of resting 12-lead ECG.

2. Mortara algorithm

The Mortara algorithm enables 3-, 5- and 12-lead ECG monitoring, ST-segment analysis and arrhythmia analysis.

You can select either algorithm as required. For the patient monitor incorporating the Mortara algorithm, interpretation of resting 12-lead ECG in the Mindray algorithm is optional. The MPM module or the BeneView T1 incorporating Mortara algorithm is labelled with the logo of Mortara.

8.2 Safety



WARNING

- **Use only ECG electrodes and cables specified by the manufacturer.**
 - **When connecting electrodes and/or patient cables, make sure that the connectors never come into contact with other conductive parts, or with earth. In particular, make sure that all of the ECG electrodes are attached to the patient, to prevent them from contacting conductive parts or earth.**
 - **Periodically inspect the electrode application site to ensure skin quality. If the skin quality changes, replace the electrodes or change the application site.**
 - **Use defibrillator-proof ECG cables during defibrillation.**
 - **Do not touch the patient, or table, or instruments during defibrillation.**
-
-

NOTE

- **After defibrillation, the screen display recovers within 10 seconds if the correct electrodes are used and applied in accordance with the manufacturer's instructions for use.**
 - **Interference from a non-grounded instrument near the patient and electrosurgery interference can cause problems with the waveform.**
-

8.3 Preparing to Monitor ECG

8.3.1 Preparing the Patient and Placing the Electrodes

1. Prepare the patient's skin. Proper skin preparation is necessary for good signal quality at the electrode, as the skin is a poor conductor of electricity. To properly prepare the skin, choose flat areas and then follow this procedure:
 - ◆ Shave hair from skin at chosen sites.
 - ◆ Gently rub skin surface at sites to remove dead skin cells.
 - ◆ Thoroughly cleanse the site with a mild soap and water solution. We do not recommend using ether or pure alcohol, because this dries the skin and increases the resistance.
 - ◆ Dry the skin completely before applying the electrodes.
2. Attach the clips or snaps to the electrodes before placing them.
3. Place the electrodes on the patient.
4. Attach the electrode cable to the patient cable and then plug the patient cable into the ECG connector on the MPM or the BeneView T1.

8.3.2 Choosing AHA or IEC Lead Placement

1. Select the ECG parameter window or waveform area to enter the **[ECG Setup]** menu.
2. Select **[Others >>]**→**[Lead Set]** and then select **[3-lead]**, **[5-lead]**, **[12-lead]** or **[Auto]** according to the applied electrodes.
3. Select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password
4. Select **[Others >>]**→**[ECG Standard]** and then select **[AHA]** or **[IEC]** according to the standard that is applied for your hospital.

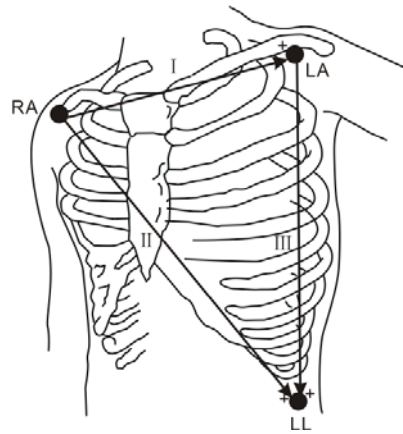
8.3.3 ECG Lead Placements

The electrode placement illustrations in this chapter adopt the AHA standard.

3-Leadwire Electrode Placement

Following is an electrode configuration when using 3 leadwires:


- RA placement: directly below the clavicle and near the right shoulder.
- LA placement: directly below the clavicle and near the left shoulder.
- LL placement: on the left lower abdomen.



5-Leadwire Electrode Placement

Following is an electrode configuration when using 5 leadwires:

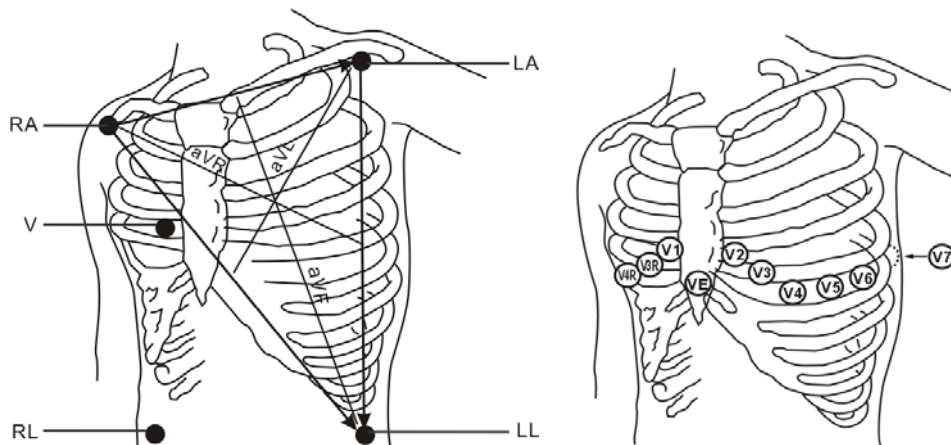- RA placement: directly below the clavicle and near the right shoulder.
- LA placement: directly below the clavicle and near the left shoulder.
- RL placement: on the right lower abdomen.
- LL placement: on the left lower abdomen.
- V placement: on the chest.




The chest (V) electrode can be placed on one of the following positions:

- V1 placement: on the fourth intercostal space at the right sternal border.
- V2 placement: on the fourth intercostal space at the left sternal border.
- V3 placement: midway between the V2 and V4 electrode positions.
- V4 placement: on the fifth intercostal space at the left midclavicular line.
- V5 placement: on the left anterior axillary line, horizontal with the V4 electrode position.
- V6 placement: on the left midaxillary line, horizontal with the V4 electrode position.
- V3R-V6R placement: on the right side of the chest in positions corresponding to those on the left.
- VE placement: over the xiphoid process.
- V7 placement: on posterior chest at the left posterior axillary line in the fifth intercostal space.
- V7R placement: on posterior chest at the right posterior axillary line in the fifth intercostal space.

12-Leadwire Electrode Placement

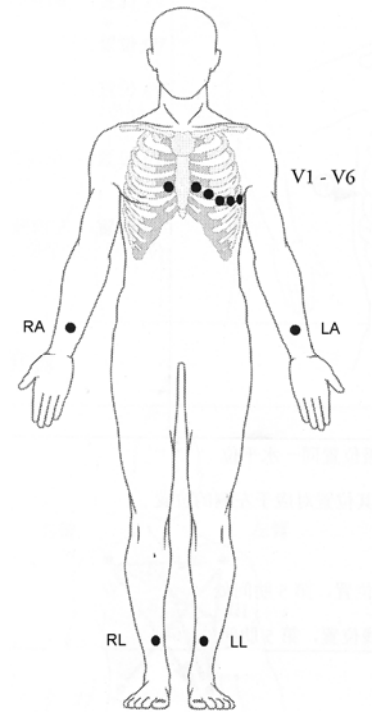12-lead ECG uses 10 electrodes, which are placed on the patient's four limbs and chest. The limb electrodes should be placed on the soft skin and the chest electrodes placed according to the physician's preference.



Lead Placement for Surgical Patients



The surgical site should be taken into consideration when placing electrodes on a surgical patient. e.g. for open-chest surgery, the chest electrodes can be placed on the lateral chest or back. To reduce artifacts and interference from electrosurgical units, you can place the limb electrodes close to the shoulders and lower abdomen and the chest electrodes on the left side of the mid-chest. Do not place the electrodes on the upper arm. Otherwise, the ECG waveform will be very small.



WARNING

- When using electrosurgical units (ESU), patient leads should be placed in a position that is equal distance from the Electrosurgery electrotome and the grounding plate to avoid burns to the patient. Never entangle the ESU cable and the ECG cable together.
- When using electrosurgical units (ESU), never place ECG electrodes near to the grounding plate of the ESU, as this can cause a lot of interference on the ECG signal.

8.3.4 Checking Paced Status

It is important to set the paced status correctly when you start monitoring ECG. The paced symbol  is displayed when the [Paced] status is set to [Yes]. The pace pulse markers “|” are shown on the ECG wave when the patient has a paced signal. If [Paced] is set to [No] or the patient's paced status is not selected, the symbol  will be shown in the patient information area.

To change the paced status, you can select either:

- the patient information area, or
- [Main Menu]→[Patient Setup]→[Patient Demographics], or,
- the ECG parameter window or waveform area→[Others >>],

and then, select [Paced] from the popup menu and toggle between [Yes] and [No].


If you do not set the paced status, the patient monitor issues a prompt tone when pace pulse is detected. At the same time, the paced symbol flashes and the message “Please confirm the pace of patient” appears in the ECG waveform area. Then, please check and set the paced status of the patient.

 **WARNING**

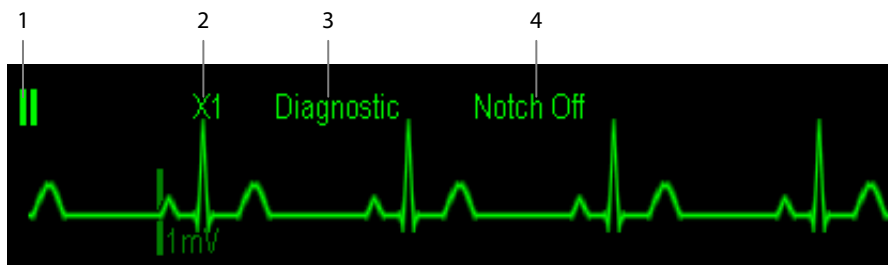- For paced patients, you must set [Paced] to [Yes]. If it is incorrectly set to [No], the patient monitor could mistake a pace pulse for a QRS and fail to alarm when the ECG signal is too weak. Do not rely entirely on rate meter alarms when monitoring patients with pacemakers. Always keep these patients under close surveillance.
 - For non-paced patients, you must set [Paced] to [No].
 - The auto pacer recognition function is not applicable to neonatal patients.
-

8.4 Understanding the ECG Display

Your display may be configured to look slightly different.



1. Lead label of the displayed wave
2. ECG gain
3. ECG filter label
4. Notch filter status

Besides, when a paced signal has been detected, the pace pulse marks “|” are shown on the ECG wave if the **[Paced]** has been set to **[Yes]**.



1. Current heart rate alarm limits
2. Current heart rate
3. Heart beat symbol

For 12-lead ECG display screen, refer to the section **12-Lead ECG Monitoring**.

8.5 Changing ECG Settings

8.5.1 Accessing ECG Menus

By selecting the ECG parameter window or waveform area, you can access the **[ECG Setup]** menu.

8.5.2 Setting Pacemaker Rate (For Mortara only)

Some pacemaker pulses can be difficult to reject. When this happens, the pulses are counted as a QRS complex and could result in an incorrect HR and failure to detect some arrhythmias. You can set **[Pacemaker Rate]** to the pacemaker's rate in the **[ECG Setup]** menu. In this way, the patient monitor can calculate HR and detect arrhythmias more accurately. When **[Paced]** is set to **[No]**, the pacemaker rate cannot be set.

8.5.3 Choosing the Alarm Source

In most cases the HR and PR numerics are identical. In order to avoid simultaneous alarms on HR and PR, the monitor uses either HR or PR as its active alarm source. To change the alarm source, select **[Alm Source]** in the **[ECG Setup]** menu and then select either:

- **[HR]**: if you want the HR to be the alarm source for HR/PR.
- **[PR]**: if you want the PR to be the alarm source for HR/PR.
- **[Auto]**: If the **[Alm Source]** is set to **[Auto]**, the patient monitor will use the heart rate from the ECG measurements as the alarm source whenever a valid heart rate is available. If the heart rate becomes unavailable, for example the ECG module is turned off or becomes disconnected, the patient monitor will automatically switch to PR as the alarm source.

8.5.4 Setting the ECG Lead Set

You can set the **[Lead Set]** by selecting **[ECG Setup]**→**[Others>>]**. You can set the **[Lead Set]** as **[Auto]** if the auto lead detection function is available.

8.5.5 Choosing an ECG Display Screen

When monitoring with a 5-lead or 12-lead set, you can select the **[Screens]** Quickkey. In the **[Choose Screen]** window, choose the screen type as:

- **[Normal Screen]**: The ECG waveform area shows 2 ECG waveforms.
- **[ECG 7-Lead Full-Screen]**: The whole waveform area shows 7 ECG waveforms only.
- **[ECG 7-Lead Half-Screen]**: The upper half part of the whole waveform area displays 7 ECG waveforms.

When monitoring with a 12-lead set, you can also choose the screen type as **[ECG 12-Lead Full-Screen]**.

When the screen type is set to **[Normal Screen]** and **[Sweep Mode]** is set to **[Refresh]**, cascaded ECG waveforms can be displayed. To cascade ECG waveforms:

1. Select the **[Screens]** Quickkey→**[Screen Setup]**.
2. Select **[ECG1 Casc.]** in the second row. A cascaded waveform is displayed in two waveform positions.

8.5.6 Changing the ECG Filter Settings

The ECG filter setting defines how ECG waves are smoothed. To change the filter setting, select **[Filter]** from **[ECG Setup]** and then select the appropriate setting.

- **[Monitor]**: Use under normal measurement conditions.
- **[Diagnostic]**: Use when diagnostic quality is required. The unfiltered ECG wave is displayed so that changes such as R-wave notching or discrete elevation or depression of the ST segment are visible.
- **[Surgery]**: Use when the signal is distorted by high frequency or low frequency interference. High frequency interference usually results in large amplitude spikes making the ECG signal look irregular. Low frequency interference usually leads to wandering or rough baseline. In the operating room, the surgery filter reduces artifacts and interference from electrosurgical units. Under normal measurement conditions, selecting **[Surgery]** may suppress the QRS complexes too much and then interfere with ECG analysis.
- **[ST]**: Use when ST monitoring is applied.

 **WARNING**

- **The [Diagnostic] filter is recommended when monitoring a patient in an environment with slight interference only.**
-

8.5.7 Setting the Notch Filter

The notch filter removes the line frequency interference. Only when **[Filter]** is set to **[Diagnostic]**, the **[Notch Filter]** is adjustable.

1. Select the ECG parameter window or waveform area to enter its setup menu. Then select **[Others >>]**.
2. Set **[Notch Filter]** to
 - **[Strong]** when there is strong interference (such as spikes) with the waveform.
 - **[Weak]** when there is weak interference with the waveform.
 - **[Off]** to turn the notch filter off.
3. When **[Notch Filter]** is set on, select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password.
4. Select **[Others >>]**→**[Notch Freq.]** and then select **[50Hz]** or **[60Hz]** according to the power line frequency.

8.5.8 Changing the Pacer Reject Settings

Select **[ECG Setup]**→**[Others>>]**→**[Pacer Reject]**, and toggle between **[On]** and **[Off]**.

When **[Paced]** is set to **[Yes]**:

- When **[Pacer Reject]** is switched on, the pace pulses are not counted as extra QRS complexes.
- The pace pulse marks “|” are shown on the ECG wave when pace pulses are detected.

When **[Paced]** is set to **[No]**, the pace markers are not shown on the ECG wave, and the options of **[Pacer Reject]** are invalid.

8.5.9 About the Defibrillator Synchronization

If a defibrillator is connected, a defibrillator synchronization pulse (100 ms, +5V) is outputted through the Defib. Sync Connector every time when the patient monitor detects an R-wave. The defibrillator synchronization function is always enabled.

WARNING

- **Improper use of a defibrillator may cause injury to the patient. The user should determine whether to perform defibrillation or not according to the patient's condition.**
 - **Before defibrillation, the user must ensure both defibrillator and monitor has passed the system test and can be safely used jointly.**
 - **Before defibrillation, make sure that the [Filter] is set to [Diagnostic].**
 - **After defibrillation is finished, select the filter mode as required.**
-

8.5.10 Changing ECG Wave Settings

In the [ECG Setup] menu:

- You can select [ECG], [ECG1], or [ECG2] to select a lead to view. The waveform of selected lead should have the following characteristics:
 - ◆ The QRS should be either completely above or below the baseline and it should not be biphasic.
 - ◆ The QRS should be tall and narrow.
 - ◆ The P-waves and T-waves should be less than 0.2mV.
- If the wave is too small or clipped, you can change its size by selecting an appropriate [Gain] setting. If you select [Auto] from [Gain], the patient monitor will automatically adjust the size of the ECG waves. In normal screen, only the selected ECG wave's size is adjusted. In other screens, all ECG waves' size is adjusted simultaneously.
- You can change the wave sweep speed by selecting [Sweep] and then selecting the appropriate setting.

8.5.11 Enabling Smart Lead Off

When the smart lead off function is set on and there is a "lead off" in the lead that has an ECG waveform in filter mode and notch status, if another lead is available, this available lead automatically becomes that lead. The system will re-calculate HR and analyze and detect arrhythmia. When the "lead off" condition is corrected, the leads are automatically switched back.

To switch on/off the smart lead off function, select [Others >>] from the [ECG Setup] menu; select [Smart Lead Off] and toggle between [On] and [Off] from the popup menu.

8.5.12 Setting the Alarm Level for ECG Lead Off Alarms

Select [Alarm Setup >>] from the [User Maintenance] menu. You can set [ECGLeadOff Lev.] from the popup menu.

8.5.13 Adjusting QRS Volume

QRS sounds are produced based on the alarm source. To adjust the QRS volume, select [Others >>] from the [ECG Setup] menu; select [QRS Volume] from the popup menu and select the appropriate setting. When valid SpO₂ measured value is available, the system adjusts the pitch tone of QRS sound based on the SpO₂ value.

8.6 About ST Monitoring

- Mortara ST segment analysis is not intended for neonatal patients.
- ST segment analysis calculates ST segment elevations and depressions for individual leads and then displays them as numerics in the ST1 and ST2 areas.
- A positive value indicates ST segment elevation; a negative value indicates ST segment depression.
- Measurement unit of the ST segment: mV or mm. You can set the unit in the [Unit Setup] menu from the [User Maintenance] menu.
- Measurement range of the ST segment: -2.0 mV to +2.0 mV.



WARNING

- **The ST algorithm has been tested for accuracy of the ST segment data. The significance of the ST segment changes need to be determined by a clinician.**
-

8.6.1 Switching ST On and Off

To switch ST monitoring on or off:

1. In the [ECG Setup] menu, select [ST Analysis >>].
2. Select [ST Analysis] to toggle between [On] and [Off].

Reliable ST monitoring can hardly be ensured if:

- You are unable to get a lead that is not noisy.
- Arrhythmias such as atrial fib/flutter cause irregular baseline.
- The patient is continuously ventricularly paced.
- The patient has left bundle branch block.

In these cases, you may consider switching ST monitoring off.

8.6.2 Changing ST Filter Settings

ST-segment analysis can be carried out only when the filter mode is set to [Diagnostic] or [ST]. When ST-segment analysis is switched on, [Filter] will automatically switch to [ST] if it is not [Diagnostic] or [ST]. When ST-segment analysis is switched off, the filter mode automatically switches to previous manual setting.

However, if you switch [Filter] to [Monitor] or [Surgery], ST-segment analysis will turn off automatically. In case that you change [Monitor] or [Surgery] to [Diagnostic] or [ST], ST-segment analysis keeps off, you can turn it on manually.

8.6.3 Understanding the ST Display

8.6.3.1 ST Numerics

This example shows ST numerics with 5-lead ECG. Your monitor screen may look slightly different from the illustration.

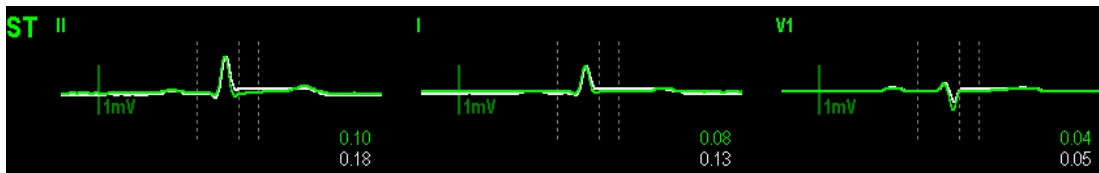


8.6.3.2 ST Segment


ST segment shows a QRS complex segment for each measured ST lead. The current ST segment is drawn in the same color as the ECG wave, usually green, superimposed over the stored reference segment, drawn in a different color. The information is updated once every ten seconds.

To display the ST segment on normal screen:

1. Enter the **[ST Analysis]** menu. Set **[ST Analysis]** to **[On]**.
2. Enter the **[Screen Setup]** window of **[Screens]** menu. Set **[ST Segment]** to be displayed.



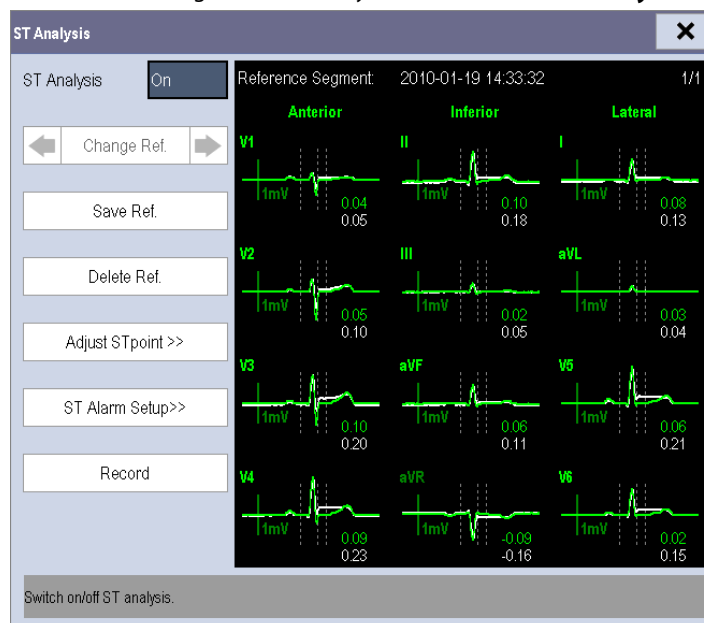Select the ST parameter window or ST segment area and you can enter the **[ST Analysis]** menu.



8.6.4 Saving the Current ST Segment as Reference

Select **[Save Ref.]** in the **[ST Analysis]** menu to save the current segment as reference. Up to 20 reference segment groups can be saved.

NOTE

- If the memory is full and you do not delete a group before saving a new one, the oldest saved group is deleted automatically.

8.6.5 Changing the Reference Segment

Select the ◀ and ▶ arrow keys beside the [Change Ref.] to switch between different reference segment groups.

8.6.6 Deleting a Reference Segment

To delete the current ST reference segment, select [Delete Ref.] in the [ST Analysis] menu and then select [Ok] in the popup.

8.6.7 Recording the ST Segment

To record the current ST segment and reference segment, select [Record] in the [ST Analysis] menu.

8.6.8 Changing the ST Alarm Limits

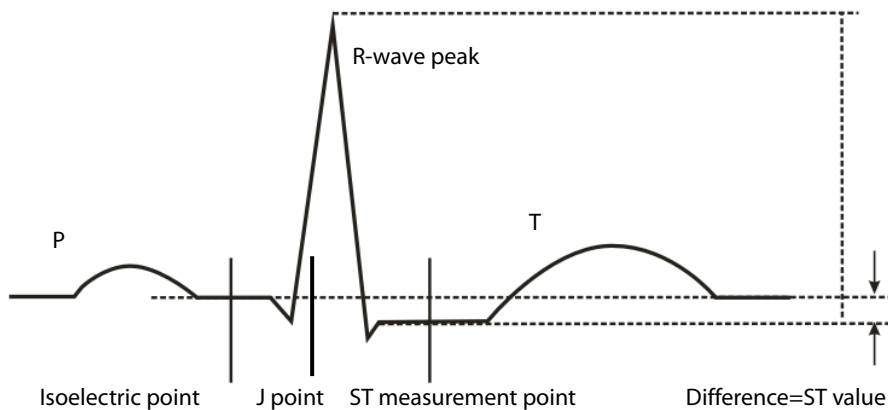
High and low ST alarm limits can be set individually for each ECG lead. Alarm limits can also be set separately for single-lead and multi-lead ST monitoring. You can select [ST Alarm Setup >>] from [ST Analysis] menu and then change ST alarm settings for each lead.

8.6.9 Setting the ST Alarm Delay Time

You can set the ST alarm delay time from the [Others] window of [Alarm Setup] menu.

8.6.10 Adjusting ST Measurement Points

As shown in the figure below, the ST measured for each beat complex is the vertical difference between two measurement points with the R-wave peak as the baseline for the measurement.



The ISO and ST points need to be adjusted when you start monitoring and if the patient's heart rate or ECG morphology changes significantly. Exceptional QRS complexes are not considered for ST-segment analysis.



WARNING

- Always make sure that the positions of ST measurement points are appropriate for your patient.

To adjust the ST measurement points:

1. In the [ST Analysis] menu, select [Adjust ST Point >>]. In the [Adjust ST Point] window, three vertical lines represent the ISO, J and ST point positions respectively.
2. Select [View Leads] and use the Knob to select an ECG lead with obvious J point and R wave.
3. Select [ISO], [J] or [ST Point] and then use the Knob to adjust the position of each point.
 - ◆ The ISO-point (isoelectric) position is given relative to the R-wave peak. Position the ISO-point in the middle of the flattest part of the baseline (between the P and Q waves of in front of the P wave).
 - ◆ The J-point position is given relative to the R-wave peak and helps locating the ST-point. Position the J-point at the end of the QRS complex and the beginning of the ST segment.
 - ◆ The ST-point is positioned a fixed distance from the J-point. Move the J-point to position the ST-point at the midpoint of the ST segment. Position the ST-point relative to the J-point at either [J+60/80ms], [J+40ms], [J+60ms] or [J+80ms]. When [J+60/80ms] is selected, the ST-point will be positioned 80 ms (heart rate 120 bpm or less) or 60 ms (heart rate more than 120 bpm) from the J-point.

8.7 About Arrhythmia Monitoring

Arrhythmia analysis provides information about your patient's condition, including heart rate, PVC rate, rhythm and ectopics.



WARNING

- **Arrhythmia analysis program is intended to detect ventricular arrhythmias. It is not designed to detect atrial or supraventricular arrhythmias. It may incorrectly identify the presence or absence of an arrhythmia. Therefore, a physician must analyze the arrhythmia information with other clinical findings.**
- **Mortara arrhythmia algorithm is not intended for neonatal patients.**

8.7.1 Understanding the Arrhythmia Events

Mindray algorithm

Arrhythmia message	Description	Category
Asystole	No QRS detected within the set time threshold in absence of ventricular fibrillation or chaotic signal.	Lethal arrhythmia
Vfib/Vtac	A fibrillatory wave for 6 consecutive seconds. A dominant rhythm of adjacent Vs and a HR > the V-Tac HR limit.	
Vtac	The consecutive PVCs > Vtac PVCs limit, and the HR > the Vtac HR limit.	
Vent. Brady	The consecutive PVCs ≥ the Vbrd threshold and the ventricular HR < the Vbrd Rate threshold.	
Extreme Tachy	The heart rate is greater than the extreme tachycardia limit.	
Extreme Brady	The heart rate is less than the extreme bradycardia limit.	
PVCs	PVCs/min exceeds high limit	Nonlethal arrhythmia
PNP	No pace pulse detected for 1.75 x average R-to-R intervals following a QRS complex (for paced patients only).	

Arrhythmia message	Description	Category
PNC	No QRS complex detected for 300 milliseconds following a pace pulse (for paced patients only).	
PVC	One PVC detected in normal heartbeats.	
Couplet	Paired PVCs detected in normal heartbeats.	
VT > 2	More than 2 consecutive PVCs within the last minute.	
Bigeminy	A dominant rhythm of N, V, N, V, N, V.	
Trigeminy	A dominant rhythm of N, N, V, N, V, N, V.	
R on T	R on T detected in normal heartbeats.	
Missed Beats	No beat detected for 1.75 x average R-R interval for HR <120, or No beat for 1 second with HR > 120 (for non-paced patients only), or No beat detected for more than the set pause threshold.	
Brady	The average heart rate is less than the bradycardia limit.	
Tachy	The average heart rate is greater than the tachycardia limit.	
Vent. Rhythm	The consecutive PVCs > the Vbrd PVCs limit, and the HR is between Vbrd Rate limit and the Vtac Rate limit.	
Multif. PVC	Multiform PVCs detected in Multif. PVC's Window (which is adjustable).	
Nonsus. Vtac	The consecutive PVCs < the Vtac PVCs limit but > 2, and HR > the Vtac Rate limit.	
Pause	No QRS detected within the set time threshold of pause.	
Irr. Rhythm	Consistently irregular rhythm.	

Mortara algorithm

Arrhythmia Message	Description	Category
Asystole	No QRS complex detected within the set time threshold (in absence of ventricular fibrillation or chaotic signals).	Lethal arrhythmia
Vfib	Ventricular fibrillation occurs and persists for 6 seconds.	
Vtac	Ventricular HR is greater or equal to the preset threshold and the number of consecutive PVCs is greater than the preset threshold.	
PVCs	PVCs/min exceeds high limit	Nonlethal arrhythmia
PNP	No pace pulse detected for (60*1000/pace rate +90) milliseconds following a QRS complex or a pacer pulse (for paced patients only).	
PNC	No QRS complex detected for 300 milliseconds following a pace pulse (for paced patients only).	
Multif. PVC	More than 2 PVCs of different forms occur in the predefined search window (3-31).	
Couplet	Paired PVCs are detected.	
VT > 2	Ventricular HR is greater than or equal to the preset threshold and the number of PVCs is greater than or equal to 3 but less than the preset threshold.	
Vent. Rhythm	Ventricular HR is less than the preset threshold and the number of PVCs is greater than or equal to 3.	
Bigeminy	A dominant rhythm of N, V, N, V, N, V.	
Trigeminy	A dominant rhythm of N, N, V, N, V, N, V.	
R on T	R on T is detected.	
Irr. Rhythm	Consistently irregular rhythm	

Arrhythmia Message	Description	Category
Missed Beats	No beat detected for 1.75x average R-R interval for HR <120, or No beat for 1 second with HR >120 (for non-paced patients only), or No beat detected for more than the set pause threshold.	
Brady	The HR is less than the set bradycardia low limit.	
Tachy	The HR is greater than the set tachycardia high limit.	

8.7.2 Changing Arrhythmia Alarm Settings

To change arrhythmia alarm settings, select the ECG parameter area or waveform area → [ECG Setup] → [Arrh. Analysis >>]. In the pop-up menu, you can set the [Alm Lev] to [High], [Med], [Low] or [Message], or switch on lethal arrhythmia analysis alarms only or switch on/off all arrhythmia analysis alarms. In the [Alarm Setup] menu from the [User Maintenance] menu, you can enable/disable turning off lethal arrhythmia analysis alarms.



WARNING

- If you switch off all arrhythmia analysis alarms, the monitor cannot give any arrhythmia analysis alarm. Always keep the patient under close surveillance.

8.7.3 Changing Arrhythmia Threshold Settings

Select the ECG parameter window or waveform area → [Arrh. Analysis >>] → [Arrh. Threshold], and you can then change threshold settings for some arrhythmia alarms. In case an arrhythmia violates its threshold, an alarm will be triggered. The asystole delay time relates to ECG relearning. When HR is less than 30 bpm, it is recommended to set the asystole delay time to 10 seconds.

Mindray algorithm

Arrh. event	Range	Default	Step	Unit
PVCs High	1 to 100	10	1	/min
Asys. Delay	3 to 10	5	1	s
Tachy High	60 to 300	Adult: 120 Pediatric: 160 Neonate: 180	5	bpm
Brady Low	15 to 120	Adult: 50 Pediatric: 75 Neonate: 90	5	bpm
Extreme Tachy	120 to 300	Adult: 160 Pediatric: 180 Neonate: 200	5	bpm
Extreme Brady	15 to 60	Adult: 35 Pediatric: 50 Neonate: 60	5	bpm
Multif. PVC's Window	3 to 31	15	1	/min
Vtac Rate	100 to 200	Adult, pediatric: 130 Neonate: 160	5	bpm

Arrh. event	Range	Default	Step	Unit
Vtac PVCs	3 to 99	6	1	/min
Pause Time	1.5, 2.0, 2.5	2	/	s
Vbrd PVCs	3 to 99	5	1	/min
Vbrd Rate	15 to 60	40	5	bpm

Mortara algorithm

Arrh. event	Range	Default	Step	Unit
PVCs High	1 to 100	10	1	/min
Asys. Delay	2 to 10	5	1	s
Vtac Rate	100 to 200	130	5	bpm
Vtac PVC	3 to 12	6	1	beats
Multif. PVC	3 to 31	15	1	beats
Tachy High	Adult: 100 to 300 Pediatric: 160 to 300	Adult: 100 Pediatric: 160	5	bpm
Brady Low	Adult: 15 to 60 Pediatric: 15 to 80	Adult: 60 Pediatric: 80	5	bpm

8.7.4 Setting the Extended Arrh. (For Mindray Algorithm Only)

The following arrhythmia events are defined as extended arrhythmia:

- Extreme Tachy
- Extreme Brady
- Vent. Brady
- Nonsus. Vtac
- Multif. PVC
- Irr. Rhythm
- Pause

You can select **[Main Menu]**→**[Maintenance >>]**→**[User Maintenance >>]**→enter the required password→select **[Alarm Setup >>]**, and set **[Extended Arrh.]** to **[Enable]** or **[Disable]**. When **[Extended Arrh.]** is set to **[Disable]**, the patient monitor does not analysis the extended arrhythmia events and corresponding alarms are not given.



- Set **[Extended Arrh.]** to **[Disable]** when the patient monitor is connected to the Central Monitoring System of version prior to 06.01.00. Failure to do so may cause the Central Monitoring System unable to display extended arrhythmia related alarms normally when extended arrhythmia occurs.
-

8.7.5 Reviewing Arrhythmia Events

Please refer to the **Review** chapter.

8.8 ECG Relearning

8.8.1 Initiating an ECG Relearning Manually

During ECG monitoring, you may need to initiate an ECG relearning when the patient's ECG template changes dramatically. A change in the ECG template could result in:

- incorrect arrhythmia alarms
- loss of ST measurement, and/or
- inaccurate heart rate

ECG relearning allows the monitor to learn the new ECG template so as to correct arrhythmia alarms and HR value, and restore ST measurements. To initiate relearning manually, select the ECG parameter window or waveform area → **[Relearn]**. When the patient monitor is learning, the message **[ECG Learning]** is displayed in the technical alarm area.

CAUTION

- **Take care to initiate ECG relearning only during periods of normal rhythm and when the ECG signal is relatively noise-free. If ECG learning takes place during ventricular rhythm, the ectopics may be incorrectly learned as the normal QRS complex. This may result in missed detection of subsequent events of V-Tach and V-Fib.**
-

8.8.2 Automatic ECG Relearning

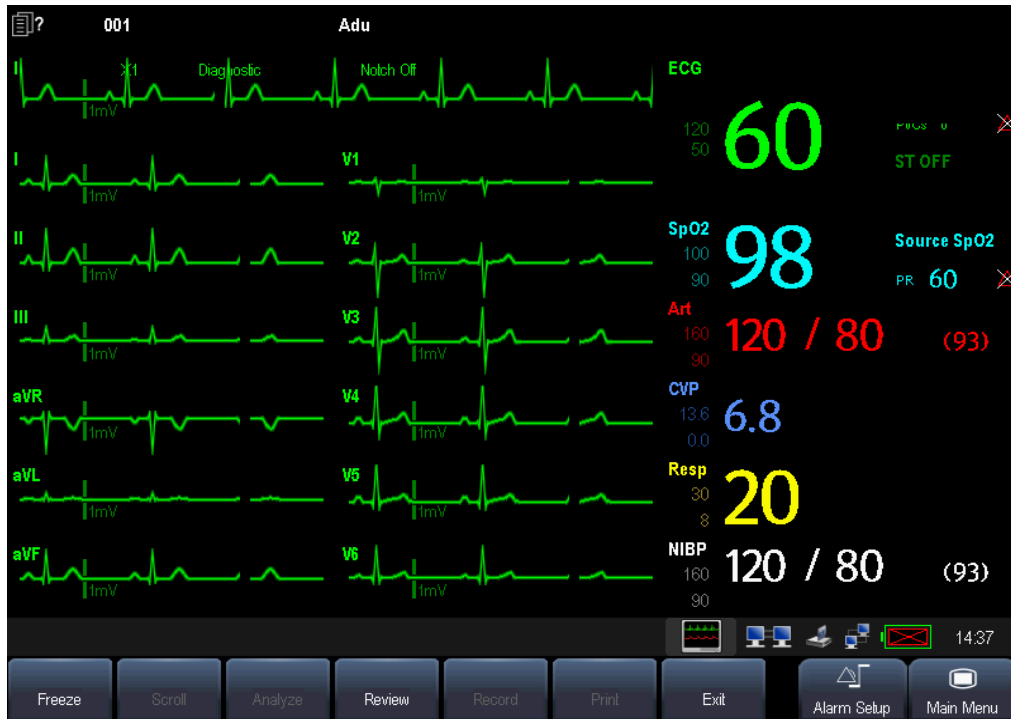
ECG relearning is initiated automatically whenever:

- The ECG lead or lead label is changed
- The ECG lead is re-connected
- A new patient is admitted
- After the calibration is completed, select **[Stop Calibrating ECG]**
- A switch happens between the options of screen type during 5/12-lead ECG monitoring.
- The paced status of the patient is changed.

8.9 12-Lead ECG Monitoring

8.9.1 Entering the 12-lead ECG Monitoring Screen

1. Refer to the section **8.3.3 ECG Lead Placements** for placing the electrodes.
2. In the **[ECG Setup]** menu, select **[Lead Set]**→**[12-Lead]**. Select **[Screens]** Quickkey→**[ECG 12-Lead Full-Screen]**.




There are totally 12 ECG waves and 1 rhythm wave displayed on the screen. The rhythm lead is ECG I before entering the 12-lead ECG monitoring screen. The ST numerics are displayed in three groups:

- ST Ant (anterior): V1, V2, V3, V4
- ST Inf (inferior): II, III, aVF, (aVR)
- ST Lat (lateral): I, aVL, V5, V6

Although aVR is displayed in the ST Inf group, it is not an inferior lead.

Additionally, the 12-lead ECG monitoring has the following features:

- The **[Filter]** mode is automatically switched to **[Diagnostic]** when the patient monitor accesses the 12-lead full-screen; the **[Filter]** mode resumes to the configuration before accessing the 12-lead full screen when the patient monitor exit the 12-lead full screen.
- In the adult mode, the  hardkey on the monitor's front is disabled.

8.9.2 Interpretation of resting 12-lead ECG



WARNING

- Interpretation of resting 12-lead ECG is restricted to adult patients only.

You can only start a interpretation of resting 12-lead ECG 11 seconds after entering the 12-lead ECG monitoring screen. Otherwise, the prompt message **[Not enough data. Cannot analyze.]** will be displayed. To start a interpretation of resting 12-lead ECG, select **[Freeze]** and then **[Analyze]**. The following screen will be displayed. In this screen, you can:

- Select **[Record Result]** to print out the interpretation of resting 12-lead ECG results by the recorder.
- Select **[Record Wave]** to print out the interpretation of resting 12-lead ECG results and waves by the recorder.
- Select **[Print Report]** to print out the interpretation of resting 12-lead ECG report by the printer.

Resting 12-Lead ECG [X]

Analyzing Time	2009-08-07	14:18:06
Heart Rate	60	QT/QTc 346/ 346
PR Interval	176	P/QRS/T Axis 53/ 44/ 54
QRS Duration	72	RV5/SV1 1.09/ -0.58
		RV5+SV1 1.67

CODE	Diagnosis
210000	Sinus Rhythm

Record Result Record Wave Print Report

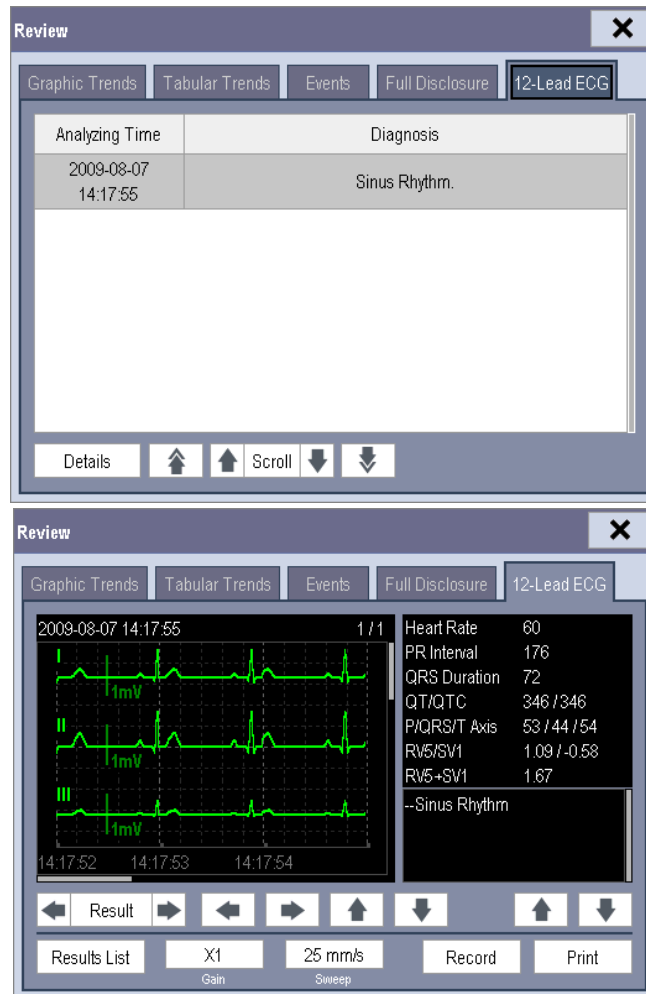
Print the interpretation of resting 12-lead ECG by recorder.

Besides, after selecting **[Freeze]**, you can:

- Browse the frozen ECG waves by selecting **[Scroll]** and rotating the Knob, or selecting the ◀ or ▶ button beside **[Scroll]**.
- Print out the currently frozen waves by selecting **[Record]**.

8.9.3 Reviewing Interpretation of resting 12-lead ECG Results

In the 12-lead ECG monitoring screen, you can review previous 12-lead ECG analyses by selecting **[Review]**.



In this review window, you can: switch between details and results list by selecting **[Details]** or **[Results List]**. When viewing the details, you can:

- Select **◀** or **▶** beside **[Result]** to switch between results.
- Adjust **[Gain]** and **[Sweep]**.
- Select **[Record]** to print out the currently displayed interpretation of resting 12-lead ECG results by the recorder.
- Select **[Print]** to print out the currently displayed interpretation of resting 12-lead ECG results by the printer.

FOR YOUR NOTES

9 Monitoring Respiration (Resp)

9.1 Introduction

Impedance respiration is measured across the thorax. When the patient is breathing or ventilated, the volume of air changes in the lungs, resulting in impedance changes between the electrodes. Respiration rate (RR) is calculated from these impedance changes, and a respiration waveform appears on the patient monitor screen.

9.2 Safety Information

WARNING

- When monitoring the patient's respiration, do not use ESU-proof ECG cables.
 - If you do not set the detection level for the respiration correctly in manual detection mode, it may not be possible for the monitor to detect apnea. If you set the detection level too low, the monitor is more likely to detect cardiac activity, and to falsely interpret cardiac activity as respiratory activity in the case of apnea.
 - The respiration measurement does not recognize the cause of apneas. It only indicates an alarm if no breath is detected when a pre-adjusted time has elapsed since the last detected breath. Therefore, it cannot be used for diagnostic purpose.
 - If operating under conditions according to the EMC Standard EN 60601-1-2 (Radiated Immunity 3V/m), field strengths above 1V/m may cause erroneous measurements at various frequencies. Therefore it is recommended to avoid the use of electrically radiating equipment in close proximity to the respiration measurement unit.
-

9.3 Understanding the Resp Display



By selecting the waveform area or parameter area, you can enter the **[Resp Waveform]** menu. By selecting the Resp parameter window, you can enter the **[Resp Setup]** menu.

NOTE

- Respiration monitoring is not for use on the patients who are very active, as this will cause false alarms.
-

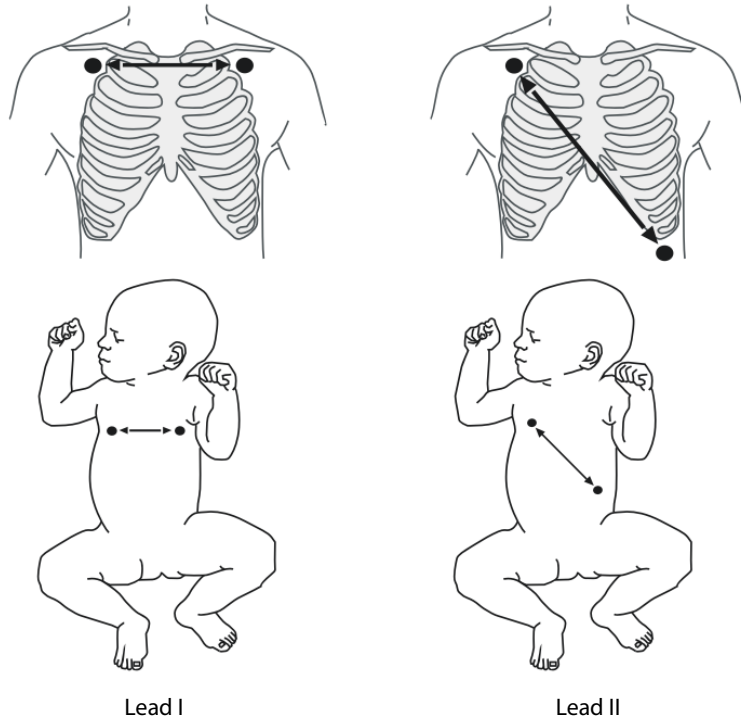
9.4 Placing Resp Electrodes

As the skin is a poor conductor of electricity, preparing the skin is necessary for a good Respiration signal. You can refer to the ECG section for how to prepare the skin.

As the Respiration measurement adopts the standard ECG electrode placement, you can use different ECG cables (3-lead, 5-lead or 12-lead). Since the respiration signal is measured between two ECG electrodes, if a standard ECG electrode placement is applied, the two electrodes should be RA and LA of ECG Lead I, or RA and LL of ECG Lead II.

NOTE

- **To optimize the respiration waveform, place the RA and LA electrodes horizontally when monitoring respiration with ECG Lead I; place the RA and LL electrodes diagonally when monitoring respiration with ECG Lead II.**



9.4.1 Optimizing Lead Placement for Resp

If you want to measure Resp and you are already measuring ECG, you may need to optimize the placement of the two electrodes between which Resp will be measured. Repositioning ECG electrodes from standard positions results in changes in the ECG waveform and may influence ST and arrhythmia interpretation.

9.4.2 Cardiac Overlay

Cardiac activity that affects the Resp waveform is called cardiac overlay. It happens when the Resp electrodes pick up impedance changes caused by the rhythmic blood flow. Correct electrodes placement can help to reduce cardiac overlay: avoid the liver area and the ventricles of the heart in the line between the respiratory electrodes. This is particularly important for neonates.

9.4.3 Abdominal Breathing

Some patients with restricted movement breathe mainly abdominally. In these cases, you may need to place the left leg electrode on the left abdomen at the point of maximum abdominal expansion to optimise the respiratory wave.

9.4.4 Lateral Chest Expansion

In clinical applications, some patients (especially neonates) expand their chests laterally, causing a negative intrathoracic pressure. In these cases, it is better to place the two respiration electrodes in the right midaxillary and the left lateral chest areas at the patient's maximum point of the breathing movement to optimise the respiratory waveform.

9.5 Choosing the Respiration Lead

In the **[Resp Setup]** menu, set **[Resp Lead]** to **[I]**, **[II]** or **[Auto]**.

9.6 Changing the Apnea Alarm Delay



The apnea alarm is a high-level alarm used to detect apneas. You can set the apnea alarm delay time after which the patient monitor alarms if the patient stops breathing. In the **[Resp Setup]** menu, select **[Apnea Delay]** and then select the appropriate setting. The **[Apnea Delay]** of Resp, CO₂, AG, and RM module keeps consistent with each other.

9.7 Changing Resp Detection Mode

In the **[Resp Setup]** menu, select **[Detection Mode]** and toggle between **[Auto]** and **[Manual]**.

- In auto detection mode, the patient monitor adjusts the detection level automatically, depending on the wave height and the presence of cardiac artifact. Note that in auto detection mode, the detection level (a dotted line) is not displayed on the waveform.

Use auto detection mode for situations where:

- ◆ The respiration rate is not close to the heart rate.
- ◆ Breathing is spontaneous, with or without continuous positive airway pressure (CPAP).
- ◆ Patients are ventilated, except patients with intermittent mandatory ventilation (IMV).
- In manual detection mode, you adjust the dotted detection level line to the desired level by selecting **[Upper Line]** or **[Lower Line]** and then selecting  or  beside them. Once set, the detection level will not adapt automatically to different respiration depths. It is important to remember that if the depth of breathing changes, you may need to change the detection level.

Use manual detection mode for situations where:

- ◆ The respiration rate and the heart rate are close.
- ◆ Patients have intermittent mandatory ventilation.
- ◆ Respiration is weak. Try repositioning the electrodes to improve the signal.

In Auto Detection Mode, if you are monitoring Resp and ECG is switched off, the monitor cannot compare the ECG and Resp rates to detect cardiac overlay. The respiration detection level is automatically set higher to prevent the detection of cardiac overlay as respiration.

In Manual Detection Mode, cardiac overlay can in certain situations trigger the respiration counter. This may lead to a false indication of a high respiration or an undetected apnea condition. If you suspect that cardiac overlay is being registered as breathing activity, raise the detection level above the zone of cardiac overlay. If the Resp wave is so small that raising the detection level is not possible, you may need to optimize the electrode placement as described in the section "Lateral Chest Expansion".

9.8 Changing Resp Wave Settings

WARNING

- **When monitoring in manual detection mode, make sure to check the respiration detection level after you have increased or decreased the size of the respiration wave.**
-


In the [**Resp Setup**] menu, you can:

- Select [**Gain**] and then select an appropriate setting. The bigger the gain is, the larger the wave amplitude is.
- Select [**Sweep**] and then select an appropriate setting. The faster the wave sweeps, the wider the wave is.

9.9 Setting RR Source

To set RR source:

1. Enter the [**Resp Setup**] menu.
2. Select [**RR Source**] and then select a source or [**Auto**] from the dropdown list.

The dropdown list displays the currently available RR source. When you select [**Auto**], the system will automatically select the RR source according to the priority. When the current RR source does not have valid measurement, the system will automatically switch the [**RR Source**] to [**Auto**]. RR source switches back to impedance respiration if you press the  silence hardkey on the monitor's front during an apnea alarm.

The priority of RR source is (from high to low): CO₂ measurement, RM measurement and impedance respiration measurement.

The **[RR Source]** settings of Resp, CO₂, AG and RM module are linked.

The RR source options and description are shown in the table below.

Option	Description
Auto	RR source is automatically selected according to the priority.
CO ₂	RR source is from CO ₂ measurement.
RM	RR source is from RM measurement.
ECG	RR source is from impedance respiration measurement.

9.10 Setting alarm properties

Select **[Alarm Setup >>]** from the **[Resp Setup]** menu. In the popup menu, you can set alarm properties for this parameter.

FOR YOUR NOTES

10 Monitoring PR

10.1 Introduction

The pulse numeric counts the arterial pulsations that result from the mechanical activity of the heart. You can display a pulse from any measured SpO₂ or any arterial pressure (see the IBP section). The displayed pulse numeric is color-coded to match its source.



1. PR: detected beats per minute.
2. PR Source

10.2 Setting the PR Source

The current pulse source is displayed in the PR parameter area. The pulse rate chosen as pulse source:

- is monitored as system pulse and generates alarms when you select PR as the active alarm source;
- is stored in the monitor's database and reviewed in the graphic/tabular trends; in trend graphs, as the PR curve is in the same color with the PR source, it is unlikely to distinguish the PR source;
- is sent via the network to the central monitoring system, if available.

To set which pulse rate as PR source:

1. Enter the **[SpO₂ Setup]** menu.
2. Select **[PR Source]** and then select a label or **[Auto]** from the popup menu.

The popup menu displays the currently available PR sources from top to bottom by priority. When you select **[Auto]**, the system will automatically select the first option as the PR source from the popup menu. When the current PR source is unavailable, the system will automatically switch **[PR Source]** to **[Auto]**. When you select **[IBP]**, the system will automatically select the first pressure label as the PR source from the popup menu.

10.3 Selecting the Active Alarm Source

In most cases the HR and pulse numerics are identical. In order to avoid simultaneous alarms on HR and Pulse, the monitor uses either HR or Pulse as its active alarm source. To change the alarm source, select **[Alm Source]** in the **[ECG Setup]** or **[SpO₂ Setup]** menu and then select either:

- **[HR]**: The monitor will use the HR as the alarm source for HR/pulse.
- **[PR]**: The monitor will use the PR as the alarm source for HR/pulse.
- **[Auto]**: If the **[Alm Source]** is set to **[Auto]**, the monitor will use the heart rate from the ECG measurement as the alarm source whenever the ECG measurement is switched on and a valid heart rate is available. If the heart rate becomes unavailable, for example if leads becomes disconnected, and a pulse source is switch on and available, the monitor will automatically switch to Pulse as the alarm source. When the Leads Off condition is corrected, the monitor will automatically switch back to the heart rate as the alarm source.

10.4 QRS Tone

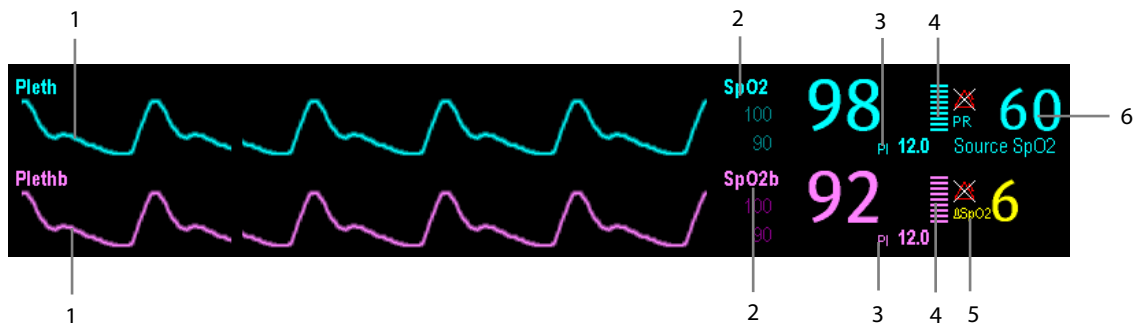
When PR is used as the alarm source, the PR source will be used as a source for the QRS tone. You can change the QRS volume by adjusting **[Beat Vol]** in the **[SpO₂ Setup]** menu. When a valid SpO₂ value exists, the system will adjust the pitch tone of QRS volume according to the SpO₂ value.

11 Monitoring SpO₂

11.1 Introduction

SpO₂ monitoring is a non-invasive technique used to measure the amount of oxygenated haemoglobin and pulse rate by measuring the absorption of selected wavelengths of light. The light generated in the probe passes through the tissue and is converted into electrical signals by the photodetector in the probe. The SpO₂ module processes the electrical signal and displays a waveform and digital values for SpO₂ and pulse rate.

This device is calibrated to display functional oxygen saturation. It provides four measurements:



1. Pleth waveform (Pleth/Plethb): visual indication of patient's pulse. The waveform is not normalized.
2. Oxygen saturation of arterial blood (SpO₂/SpO_{2b}): percentage of oxygenated hemoglobin in relation to the sum of oxyhemoglobin and deoxyhemoglobin. SpO₂ measurement is obtained through the MPM module, and SpO_{2b} measurement is obtained through the SpO₂ module.
3. Perfusion index (PI): gives the numerical value for the pulsatile portion of the measured signal caused by arterial pulsation. PI is an indicator of the pulsatile strength. You can also use it to assess the quality of SpO₂ measurement. Above 1 is optimal, between 0.3 and 1 is acceptable. Below 0.3 indicates low perfusion; reposition the SpO₂ sensor or find a better site. If low perfusion persists, choose another method to measure oxygen saturation if possible.
PI is available for Mindray SpO₂ module or Masimo SpO₂ module.
4. Perfusion indicator: the pulsatile portion of the measured signal caused by arterial pulsation.
5. SpO₂ difference (Δ SpO₂): Δ SpO₂ = | SpO₂ - SpO_{2b} | .
6. Pulse rate (derived from pleth wave): detected pulsations per minute.

In the case that you need to measure SpO₂ and SpO_{2b}, select the same type of modules. Otherwise, the SpO₂ module for SpO_{2b} is closed automatically. For example, if MPM module configuring Mindray SpO₂ and Masimo SpO₂ module are applied simultaneously, Masimo SpO₂ module is closed automatically.

11.2 Safety



WARNING

- **Use only SpO₂ sensors specified in this manual. Follow the SpO₂ sensor's instructions for use and adhere to all warnings and cautions.**
 - **When a trend toward patient deoxygenation is indicated, blood samples should be analyzed by a laboratory co-oximeter to completely understand the patient's condition.**
 - **Do not use SpO₂ sensors during magnetic resonance imaging (MRI). Induced current could potentially cause burns. The sensor may affect the MRI image, and the MRI unit may affect the accuracy of the oximetry measurements.**
 - **Prolonged continuous monitoring may increase the risk of undesirable changes in skin characteristics, such as irritation, reddening, blistering or burns. Inspect the sensor site every two hours and move the sensor if the skin quality changes. Change the application site every four hours. For neonates, or patients with poor peripheral blood circulation or sensitive skin, inspect the sensor site more frequently.**
-

11.3 Identifying SpO₂ Connectors

To identify which SpO₂ connector is incorporated into your MPM, BeneView T1, or SpO₂ module, see the company logo located at the right upper corner. The color of the cable connector matches the company as shown below:

- Mindray SpO₂ connector: a blue connector without logo.
- Masimo SpO₂ connector: a purple connector with a logo of Masimo SET.
- Nellcor SpO₂ connector: a grey connector with a logo of Nellcor.

The connectors for these three SpO₂ sensors are mutually exclusive.

11.4 Applying the Sensor

1. Select an appropriate sensor according to the module type, patient category and weight.
2. Remove colored nail polish from the application site.
3. Apply the sensor to the patient.
4. Select an appropriate adapter cable according to the connector type and plug this cable into the SpO₂ connector.
5. Connect the sensor cable to the adapter cable.

11.5 Changing SpO₂ Settings

11.5.1 Accessing SpO₂ Menus

By selecting the SpO₂ parameter window or waveform area, you can access the **[SpO₂ Setup]** or **[SpO₂b Setup]** menu.

11.5.2 Adjusting the Desat Alarm

The desat alarm is a high level alarm notifying you of potentially life threatening drops in oxygen saturation. Select **[Alarm Setup >>]** from the **[SpO₂ Setup]** or **[SpO₂b Setup]** menu. From the popup menu, you can set low alarm limit, alarm switch, and alarm recording for **[Desat]** or **[Desatb]**. When the SpO₂ value is below the desat alarm limit and desat alarm switch is set on, the message **[SpO₂ Desat]** or **[SpO₂b Desat]** is displayed.

11.5.3 Setting SpO₂ Sensitivity

For Masimo SpO₂ module, you can set **[Sensitivity]** to **[Normal]** or **[Maximum]** in the **[SpO₂ Setup]** or **[SpO₂b Setup]** menu. When the **[Sensitivity]** is set to **[Maximum]**, the patient monitor is more sensitive to minor signals. When monitoring critically ill patients whose pulsations are very weak, it is strongly recommended that the sensitivity is set to **[Maximum]**. When monitoring neonatal or non-critically ill patients who tend to move a lot, noise or invalid signals may be caused. In this case, it is recommended that the sensitivity is set to **[Normal]** so that the interference caused by motion can be filtered and therefore the measurement stability can be ensured. The settings of sensitivity in the **[SpO₂ Setup]** and **[SpO₂b Setup]** menus are linked.

11.5.4 Changing Averaging Time

The SpO₂ value displayed on the monitor screen is the average of data collected within a specific time. The shorter the averaging time is, the quicker the patient monitor responds to changes in the patient's oxygen saturation level. Contrarily, the longer the averaging time is, the slower the patient monitor responds to changes in the patient's oxygen saturation level, but the measurement accuracy will be improved. For critically ill patients, selecting shorter averaging time will help understanding the patient's state.

To set the averaging time:

- For Mindray SpO₂ module, select **[Sensitivity]** in the **[SpO₂ Setup]** or **[SpO₂b Setup]** menu and then toggle between **[High]**, **[Med]** and **[Low]**, which respectively correspond to 7 s, 9 s and 11 s.
- For Masimo SpO₂ module, select **[Averaging]** in the **[SpO₂ Setup]** or **[SpO₂b Setup]** menu and then toggle between **[2-4 s]**, **[4-6 s]**, **[8 s]**, **[10 s]**, **[12 s]**, **[14 s]** and **[16 s]**.

11.5.5 Monitoring SpO₂ and NIBP Simultaneously

When monitoring SpO₂ and NIBP on the same limb simultaneously, you can switch **[NIBP Simul]** on in the **[SpO₂ Setup]** or **[SpO₂b Setup]** menu to lock the SpO₂ alarm status until the NIBP measurement ends. If you switch **[NIBP Simul]** off, low perfusion caused by NIBP measurement may lead to inaccurate SpO₂ readings and therefore cause false physiological alarms.

11.5.6 Sat-Seconds Alarm Management

With traditional alarm management, high and low alarm limits are set for monitoring oxygen saturation. During monitoring, as soon as an alarm limit is violated, an audible alarm immediately sounds. When the patient % SpO₂ fluctuates near an alarm limit, the alarm sounds each time the limit is violated. Such frequent alarm can be distracting. Nellcor’s Sat-Seconds alarm management technique is used to reduce these nuisance alarms.

The Sat-Seconds feature is available with the Nellcor SpO₂ module to decrease the likelihood of false alarms caused by motion artifacts. To set the Sat-Seconds limit, select [**Sat-Seconds**] in the [**SpO₂ Setup**] menu and then select the appropriate setting.


With Sat-Seconds alarm management, high and low alarm limits are set in the same way as traditional alarm management. A Sat-Seconds limit is also set. The Sat-Seconds limit controls the amount of time that SpO₂ saturation may be outside the set limits before an alarm sounds. The method of calculation is as follows: the number of percentage points that the SpO₂ saturation falls outside the alarm limit is multiplied by the number of seconds that it remains outside the limit. This can be stated as the equation:

$$\text{Sat-Seconds} = \text{Points} \times \text{Seconds}$$

Only when the Sat-Seconds limit is reached, the monitor gives a Sat-Seconds alarm. For example, the figure below demonstrates the alarm response time with a Sat-Seconds limit set at 50 and a low SpO₂ limit set at 90%. In this example, the patient % SpO₂ drops to 88% (2 points) and remains there for 2 seconds. Then it drops to 86% (4 points) for 3 seconds, and then to 84% (6 points) for 6 seconds. The resulting Sat-Seconds are:

% SpO ₂	Seconds	Sat-Seconds
2×	2=	4
4×	3=	12
6×	6=	36
Total Sat-Seconds=		52

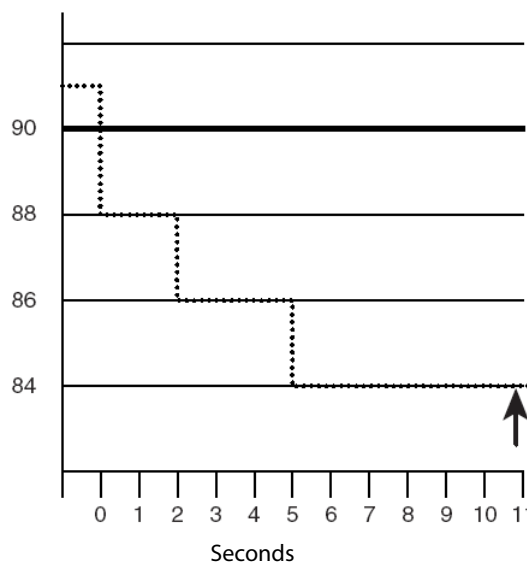After approximately 10.9 seconds, a Sat-Second alarm would sound, because the limit of 50 Sat-Seconds would have been exceeded.



Saturation levels may fluctuate rather than remaining steady for a period of several seconds. Often, the patient % SpO₂ may fluctuate above and below an alarm limit, re-entering the non-alarm range several times. During such fluctuation, the monitor integrates the number of %SpO₂ points, both positive and negative, until either the Sat-Seconds limit is reached, or the patient%SpO₂ re-enters the non-alarm range and remains there.

11.5.7 Changing the Speed of the Pleth/Plethb Wave

In the [**SpO₂ Setup**] or [**SpO₂b Setup**] menu, select [**Sweep**] and then select the appropriate setting. The faster the waveform sweeps, the wider the waveform is.

11.5.8 Setting the Alarm Level for SpO₂ Sensor Off Alarm

Select [**Alarm Setup >>**] from the [**User Maintenance**] menu. You can set the [**SpO₂ SensorOff Lev.**] in the popup menu.

11.5.9 Setting the SpO₂ Tone Mode

Select [**Others >>**] from the [**User Maintenance**] menu. In the popup menu, you can set [**SpO₂ Tone**] as [**Mode 1**] or [**Mode 2**].



- **The same SpO₂ tone mode shall be used for the same patient monitors in a single area.**
-

11.6 Measurement Limitations

If you doubt the measured SpO₂, check patient vital signs first. Then check the patient monitor and SpO₂ sensor. The following factors may influence the accuracy of measurement:

- Ambient light
- Physical movement (patient and imposed motion)
- Diagnostic testing
- Low perfusion
- Electromagnetic interference, such as MRI environment
- Electrosurgical units
- Dysfunctional haemoglobin, such as carboxyhemoglobin (COHb) and methemoglobin (MetHb)
- Presence of certain dyes, such as methylene and indigo carmine
- Inappropriate positioning of the SpO₂ sensor, or use of incorrect SpO₂
- Drop of arterial blood flow to immeasurable level caused by shock, anemia, low temperature or vasoconstrictor.

11.7 Masimo Information



- Masimo Patents

This device is covered under one or more the following U.S.A. patents: 5,758,644, 6,011,986, 6,699,194, 7,215,986, 7,254,433, 7,530,955 and other applicable patents listed at: www.masimo.com/patents.htm.

- No Implied License

Possession or purchase of this device does not convey any express or implied license to use the device with unauthorized sensors or cables which would, alone, or in combination with this device, fall within the scope of one or more of the patents relating to this device.

11.8 Nellcor Information



- Nellcor Patents

This device may be covered by one or more of the following US patents and foreign equivalents: 5,485,847, 5,676,141, 5,743,263, 6,035,223, 6,226,539, 6,411,833, 6,463,310, 6,591,123, 6,708,049, 7,016,715, 7,039,538, 7,120,479, 7,120,480, 7,142,142, 7,162,288, 7,190,985, 7,194,293, 7,209,774, 7,212,847, 7,400,919.

- No Implied License

Possession or purchase of this device does not convey any express or implied license to use the device with unauthorized replacement parts which would, alone, or in combination with this device, fall within the scope of one or more of the patents relating to this device.

12 Monitoring NIBP

12.1 Introduction

The MPM and BeneView T1 uses the oscillometric method for measuring the non-invasive blood pressure (NIBP). This measurement can be used for adults, pediatrics and neonates.

Automatic non-invasive blood pressure monitoring uses the oscillometric method of measurement. To understand how this method works, we'll compare it to the auscultative method. With auscultation, the clinician listens to the blood pressure and determines the systolic and diastolic pressures. The mean pressure can then be calculated with reference to these pressures as long as the arterial pressure curve is normal.

Since the monitor cannot hear the blood pressure, it measures cuff pressure oscillation amplitudes. Oscillations are caused by blood pressure pulses against the cuff. The oscillation with the greatest amplitude is the mean pressure. This is the most accurate parameter measured by the oscillometric method. Once the mean pressure is determined, the systolic and diastolic pressures are calculated with reference to the mean.

Simply stated, auscultation measures systolic and diastolic pressures and the mean pressure is calculated. The oscillometric method measures the mean pressure and determines the systolic and diastolic pressures.

As specified by IEC 60601-2-30/EN60601-2-30, NIBP measurement can be performed during electro-surgery and discharge of defibrillator.

NIBP diagnostic significance must be decided by the doctor who performs the measurement.

NOTE

-
- **Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultatory method or an intra-arterial blood pressure measurement device, within the limits prescribed by the American National Standard, Manual, electronic, or automated sphygmomanometers.**
-

12.2 Safety



WARNING

- **Be sure to select the correct patient category setting for your patient before measurement. Do not apply the higher adult settings for pediatric or neonatal patients. Otherwise it may present a safety hazard.**
 - **Do not measure NIBP on patients with sickle-cell disease or any condition where skin damage has occurred or is expected.**
 - **Use clinical judgement to determine whether to perform frequent unattended blood pressure measurements on patients with severe blood clotting disorders because of the risk of hematoma in the limb fitted with the cuff.**
 - **Do not use the NIBP cuff on a limb with an intravenous infusion or arterial catheter in place. This could cause tissue damage around the catheter when the infusion is slowed or blocked during cuff inflation.**
 - **If you doubt the NIBP readings, determines the patient's vital signs by alternative means and then verify that the monitor is working correctly.**
-

12.3 Measurement Limitations

Measurements are impossible with heart rate extremes of less than 40bpm or greater than 240bpm, or if the patient is on a heart-lung machine.

The measurement may be inaccurate or impossible:

- If a regular arterial pressure pulse is hard to detect
- With excessive and continuous patient movement such as shivering or convulsions
- With cardiac arrhythmias
- Rapid blood pressure changes
- Severe shock or hypothermia that reduces blood flow to the peripheries
- Obesity, where a thick layer of fat surrounding a limb dampens the oscillations coming from the artery

12.4 Measurement Methods

There are three methods of measuring NIBP:

- Manual: measurement on demand.
- Auto: continually repeated measurements at set intervals.
- STAT: continually rapid series of measurements over a five minute period, then return to the previous mode.

12.5 Setting Up the NIBP Measurement


12.5.1 Preparing to Measure NIBP

1. Power on the monitor.
2. Verify that the patient category is correct. Change it if necessary.
3. Plug the air tubing into the NIBP connector on the MPM module or BeneView T1.
4. Select a correct sized cuff and then apply it as follows:
 - ◆ Determine the patient's limb circumference.
 - ◆ Select an appropriate cuff by referring to the limb circumference marked on the cuff. The width of the cuff should be 40% (50% for neonates) of the limb circumference, or 2/3 of the upper arm's length. The inflatable part of the cuff should be long enough to encircle at least 50% to 80% of the limb.
 - ◆ Apply the cuff to an upper arm or leg of the patient and make sure the Φ marking on the cuff matches the artery location. Do not wrap the cuff too tightly around the limb. It may cause discoloration, and ischemia of the extremities. Make sure that the cuff edge falls within the marked range. If it does not, use a larger or smaller cuff that will fit better.
5. Connect the cuff to the air tubing and make sure that the bladder inside the cover is not folded and twisted.

NOTE

-
- The use of the equipment is restricted to one patient at a time;
-

12.5.2 Starting and Stopping Measurements

Select the [**NIBP Measure**] QuickKey and you can start the desired measurement from the popup menu. You can select [**Stop All**] QuickKey to stop all NIBP measurements. You can also start and stop measurements by using the  hardkey on either the monitor's front panel or the MPM module.

12.5.3 Correcting the Measurement if Limb is not at Heart Level

The cuff should be applied to a limb at the same level as the patient's heart. If the limb is not at the heart level, to the displayed value:

- Add 0.75 mmHg (0.10 kPa) for each centimetre higher, or
- Deduct 0.75 mmHg (0.10 kPa) for each centimeter lower.

12.5.4 Enabling NIBP Auto Cycling and Setting the Interval

1. Select the NIBP parameter window to enter the [**NIBP Setup**] menu.
2. Select [**Interval**] and then select a desired time interval. Selecting [**Manual**] switches to manual mode.
3. Start a measurement manually. The monitor will then automatically repeat NIBP measurements at the set time interval.

12.5.5 Starting a STAT Measurement

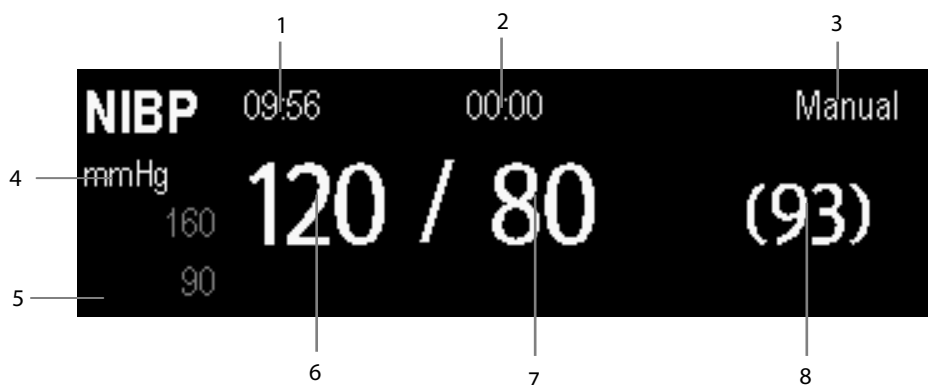
1. Select the NIBP parameter window to enter the [NIBP Setup] menu.
2. Select [NIBP STAT]. The STAT mode initiates 5 minutes of continuous, sequential, automatic NIBP measurements.

WARNING

- **Continuous non-invasive blood pressure measurements may cause purpura, ischemia and neuropathy in the limb with the cuff. Inspect the application site regularly to ensure skin quality and inspect the extremity of the cuffed limb for normal color, warmth and sensitivity. If any abnormality occurs, move the cuff to another site or stop the blood pressure measurements immediately.**
-

12.6 Understanding the NIBP Numerics

The NIBP display shows numerics only as below. Your display may be configured to look slightly different.



1. Time of last measurement
2. Time remaining to next measurement
3. Measurement mode
4. Unit of pressure: mmHg or kPa
5. Prompt message area: shows NIBP-related prompt messages
6. Systolic pressure
7. Diastolic pressure
8. Mean pressure obtained after the measurement and cuff pressure obtained during the measurement

12.7 Changing NIBP Settings

By selecting the NIBP parameter window, you can enter the [NIBP Setup] menu.

12.7.1 Setting the Initial Cuff Inflation Pressure

You can set the initial cuff inflation pressure manually. In the [NIBP Setup] menu, select [Initial Pressure] and then select the appropriate setting.

12.7.2 Setting NIBP Alarm Properties

Select [**Alarm Setup >>**] from the [**NIBP Setup**] menu. You can set the alarm properties for this parameter in the popup menu.

12.7.3 Switching On NIBP End Tone

The monitor can issue a reminder tone at the completion of NIBP measurement. The NIBP end tone is off by default. You can switch it on by accessing the [**NIBP Setup**] menu.

12.7.4 Displaying NIBP List

Select [**Screens**] QuickKey→[**Screen Setup**]. You can set [**NIBP List**] to be displayed at the bottom area of the screen. Then, multiple sets of most recent NIBP measurements will be displayed. And PR displayed is derived from NIBP.

NIBP List	PR	Time
120 / 80 (93)	60	17:15
120 / 80 (93)	60	16:23
120 / 80 (93)	60	16:09
120 / 80 (93)	60	14:24


You can not display NIBP list in some screens such as the big numerics screen and the interpretation of resting 12-lead ECG screen.

12.7.5 Setting the Pressure Unit

Select [**Unit Setup >>**] from the [**User Maintenance**] menu. In the popup menu, select [**Press. Unit**] and toggle between [**mmHg**] and [**kPa**].

12.8 Assisting Venous Puncture

You can use the NIBP cuff to cause sub-diastolic pressure to block the venous blood vessel and therefore help venous puncture.

1. Select [**VeniPuncture >>**] from the [**NIBP Setup**] menu. In the popup menu, verify that the [**Cuff Press.**] value is appropriate. Change it if necessary.
2. Select [**VeniPuncture**].
3. Puncture vein and draw blood sample.
4. Select the  hardkey on the monitor's front, or the [**Stop All**] QuickKey to deflate the cuff. The cuff deflates automatically after a set time if you do not deflate it.

During measurement, the NIBP display shows the inflation pressure of the cuff and the remaining time in venous puncture mode.

FOR YOUR NOTES

13 Monitoring Temp

13.1 Introduction

You can simultaneously monitor two temperature sites using the MPM or the BeneView T1.

13.2 Safety

WARNING

- **Verify that the probe detection program works correctly before monitoring. Plug out the temperature probe cable from the T1 or T2 connector, and the monitor can display the message [T1 Sensor Off] or [T2 Sensor Off] and give alarm tones correctly.**
-

13.3 Making a Temp Measurement

1. Select an appropriate probe for your patient.
2. If you are using a disposable probe, connect the probe to the temperature cable.
3. Plug the probe or temperature cable to the temperature connector.
4. Attach the probe to the patient correctly.
5. Check that the alarm settings are appropriate for this patient.

13.4 Understanding the Temp Display

The temperature monitoring is displayed on the monitor as three numerics: T1, T2 and TD. By selecting this area, you can enter the **[Alarm Setup]** menu.



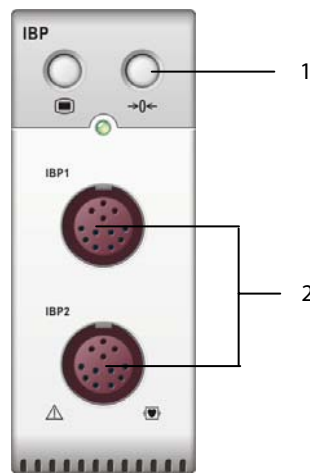
13.5 Setting the Temperature Unit

Select [**Unit Setup >>**] from the [**User Maintenance**] menu. In the popup menu, select [**Temp Unit**] and toggle between [**°C**] and [**°F**].

14 Monitoring IBP

14.1 Introduction

You can measure invasive blood pressure using the MPM, BeneView T1, or the pressure plug-in module. The monitor can monitor up to 8 invasive blood pressures and displays the systolic, diastolic and mean pressures and a waveform for each pressure.



1. Zero key

2. Connector for IBP cable

14.2 Safety

WARNING

- Use only pressure transducers specified in this manual. Never reuse disposable pressure transducers.
 - Make sure that the applied parts never contact other conductive parts.
 - To reduce the hazard of burns during high-frequency surgical procedure, ensure that the monitor's cables and transducers never come into contact with the high-frequency surgical units.
 - When using accessories, their operating temperature should be taken into consideration. For details, refer to instructions for use of accessories.
-
-

14.3 Setting Up the Pressure Measurement

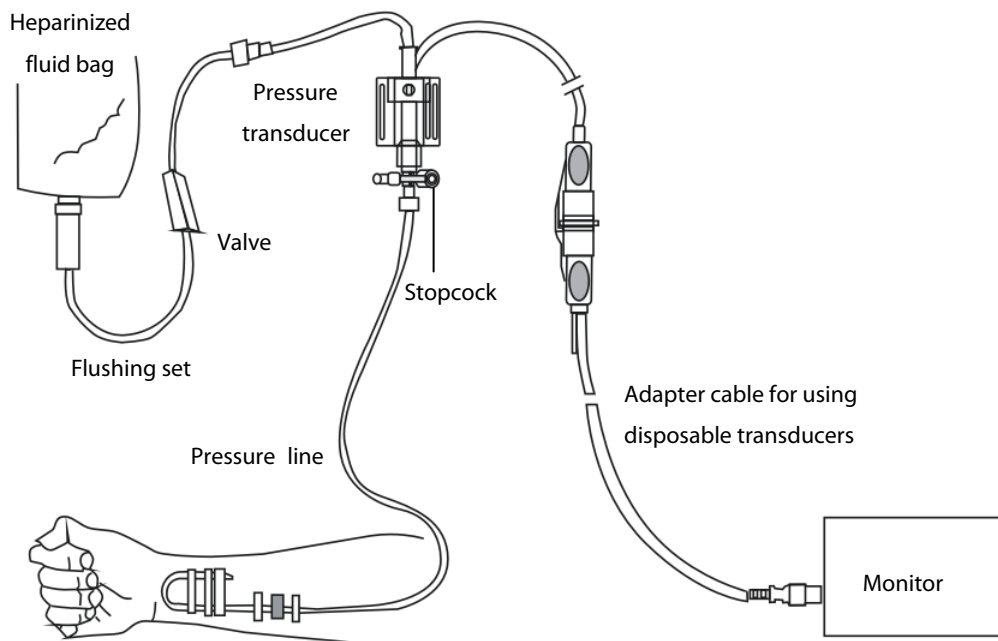
1. Plug the pressure cable into the IBP connector.
2. Prepare the flush solution.
3. Flush the system to exhaust all air from the tubing. Ensure that the transducer and stopcocks are free of air bubbles.



WARNING

- **If air bubbles appear in the tubing system, flush the system with the infusion solution again. Air bubble may lead to wrong pressure reading.**

4. Connect the pressure line to the patient catheter.
5. Position the transducer so that it is level with the heart, approximately at the level of the midaxillary line.
6. Select the appropriate label.
7. Zero the transducer. After a successful zeroing, turn off the stopcock to the atmosphere and turn on the stopcock to the patient.




WARNING

- **If measuring intracranial pressure (ICP) with a sitting patient, level the transducer with the top of the patient's ear. Incorrect leveling may give incorrect values.**

14.4 Understanding the IBP Display

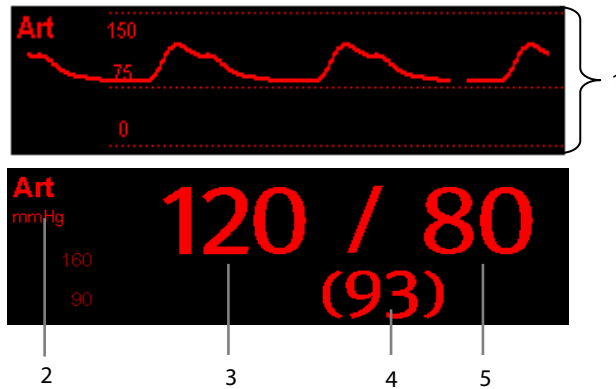The IBP measurement is displayed on the monitor as a waveform and numeric pressures. The figure below shows the waveform and numerics for the Art pressure. For different pressures, this display may be slightly different.




1. Waveform
2. Pressure unit
3. Systolic pressure
4. Mean pressure
5. Diastolic pressure

For some pressures, the parameter window may show the mean pressure only. For different pressures, their default unit may be different. If the Art and ICP pressures are measured simultaneously, the ICP parameter area will display numeric CPP, which is obtained by subtracting ICP from the Art mean.

14.5 Changing IBP Settings

14.5.1 Changing a Pressure for Monitoring

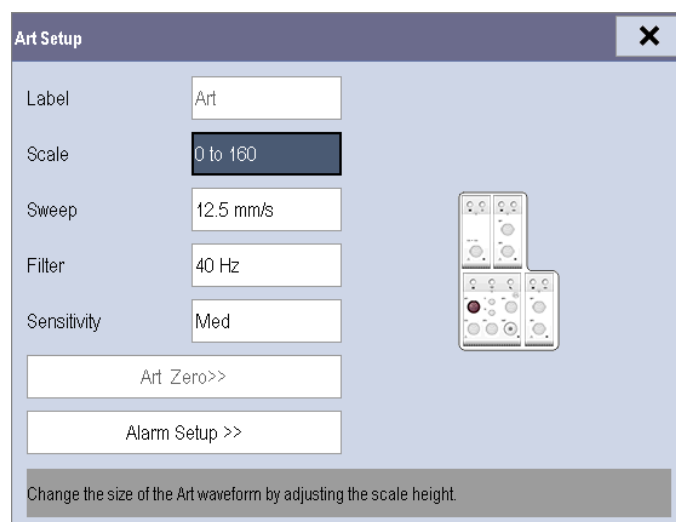1. Select the pressure you want to change to enter its setup menu. In the menu, there is a figure showing the current IBP measurement connector.



2. Select [**Label**] and then select your desired label from the list. The already displayed labels cannot be selected.

Label	Description	Label	Description
PA	Pulmonary artery pressure	CVP	Central venous pressure
Ao	Aortic pressure	LAP	Left atrial pressure
UAP	Umbilical arterial pressure	RAP	Right atrial pressure
BAP	Brachial arterial pressure	ICP	Intracranial pressure
FAP	Femoral arterial pressure	UVP	Umbilical venous pressure
Art	Arterial blood pressure	LV	Left ventricular pressure
P1 to P4	Non-specific pressure label		

14.5.2 Setting Alarm Properties

Select [**Alarm Setup >>**] from the parameter setup menu. You can set alarm properties for this parameter in the popup menu.

14.5.3 Changing Averaging Time

The IBP value displayed on the monitor screen is the average of data collected within a specific time. The shorter the averaging time is, the quicker the patient monitor responds to changes in the patient's blood pressure. Contrarily, the longer the averaging time is, the slower the patient monitor responds to changes in the patient's blood pressure, but the measurement accuracy will be improved. For critically ill patients, selecting shorter averaging time will help understanding the patient's state.

To set the averaging time, in the parameter setup menu, select [**Sensitivity**] and toggle between [**High**], [**Med**] and [**Low**], the corresponding averaging time is about 1 s, 8 s and 12 s respectively.

14.5.4 Setting the Pressure Unit

Select [**Unit Setup >>**] from the [**User Maintenance**] menu. In the popup menu, select [**Press. Unit**] and toggle between [**mmHg**] and [**kPa**]. Select [**CVP Unit**] and toggle between [**mmHg**], [**cmH₂O**] and [**kPa**].

14.5.5 Setting Up the IBP Wave

In the setup menu for the pressure, you can:

- Select [**Sweep**] and then select the appropriate setting. The faster the wave sweeps, the wider the wave is.
- Select [**Scale**] and then select the appropriate setting. If [**Auto**] is selected, the size of the pressure's waveform will be adjusted automatically.
- Select [**Filter**] and then select the desired option.

14.6 Measuring PAWP

Pulmonary Artery Wedge Pressure (PAWP) values, used to assess cardiac function, are affected by fluid status, myocardial contractility, and valve and pulmonary circulation integrity.

Obtain the measurement by introducing a balloon-tipped pulmonary artery flotation catheter into the pulmonary artery. When the catheter is in one of the smaller pulmonary arteries, the inflated balloon occludes the artery allowing the monitor to record changes in the intrathoracic pressures that occur throughout the respiration cycle.

The pulmonary wedge pressure is the left ventricular end diastolic pressure when the airway pressure and valve function are normal. The most accurate PAWP values are obtained at the end of the respiration cycle when the intrathoracic pressure is fairly constant and the artifact caused by respiration is minimal.

14.6.1 Preparing to Measure PAWP

1. Prepare the same accessories as in the C.O. measurement. Connect the parts such as catheter, syringe, etc. following the C.O. measurement steps and use the balloon inflation port.
2. Connect the PAWP cable into the IBP connector on the monitor. Since PAWP is measured on PA, selecting **[PA]** as the IBP label is recommended.
3. Select the PA parameter window or waveform area to enter its setup menu. Then, select **[PAWP]** to enter the PAWP measurement window. You can also enter the PAWP measurement window from the P1-P4 parameter window.

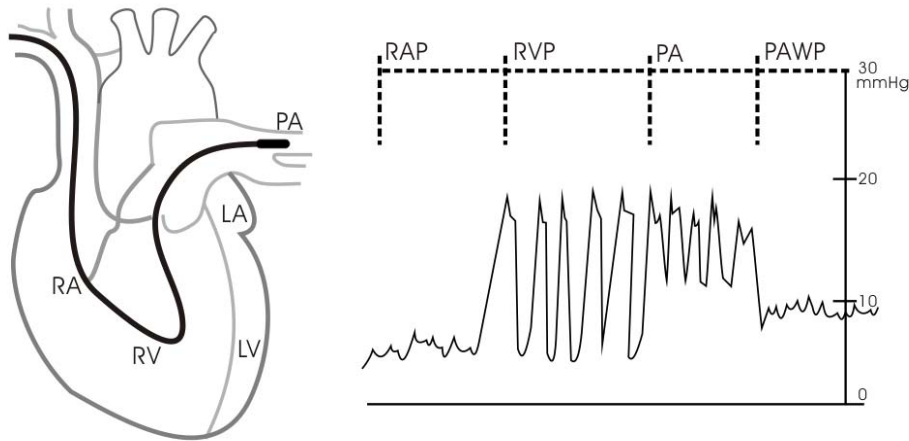


NOTE

- After entering the PAWP measurement window, the monitor will turn off the PA alarm automatically.

14.6.2 Setting Up the PAWP Measurement

1. Wedge the flotation catheter into the pulmonary artery. Then inflate the balloon and pay attention to PA waveform changes on the screen.



2. After obtaining a stable PAWP waveform, press the **[Freeze]** key to freeze the waveform and deflate the balloon.
3. You can adjust the PAWP scale to an appropriate position by adjusting ▲ or ▼ beside the **[Adjust]** button. Press the **[Confirm]** key to save one PAWP measurement.
4. If you need to start a new measurement, select **[Next Measure]**.

WARNING

- **Prolonged inflation can cause pulmonary hemorrhage, infarction or both. Inflate the balloon for the minimum time necessary to get an accurate measurement.**
 - **If the PAWP is greater than the PA (systolic), deflate the balloon and report the incident in accordance with hospital policy. Because the pulmonary artery could be accidentally ruptured, and the PAWP value derived will not reflect the patient's hemodynamic state, but will merely reflect the pressure in the catheter or balloon.**
-

14.6.3 Understanding the PAWP Setup Menu

Select **[Setup]** to enter the **[PAWP Setup]** menu. In this menu, you can:

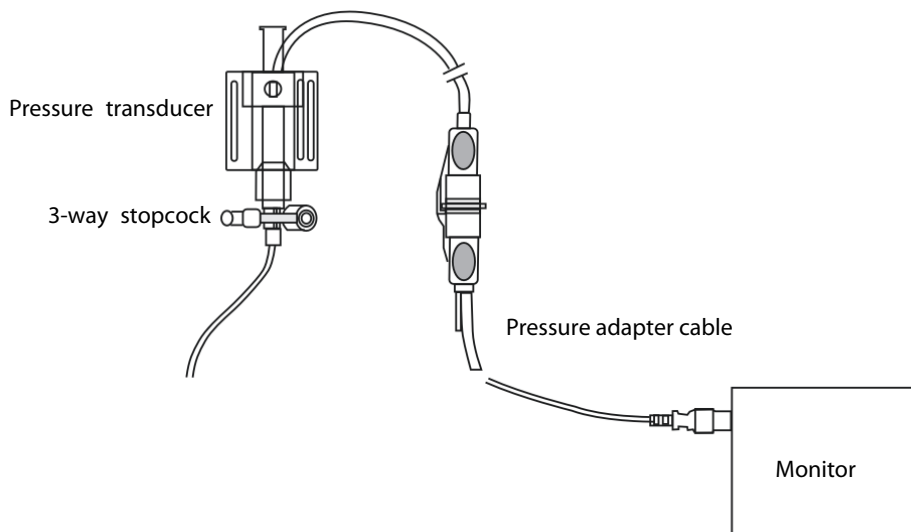
- Select a ECG lead wave as the first reference wave.
- Select a respiration wave as the second reference wave.
- Select a sweep speed for the displayed waveform.
- Change the size of the PA waveform by adjusting the scale height.

14.7 Zeroing the Transducer

To avoid inaccurate pressure readings, the monitor requires a valid zero. Zero the transducer in accordance with your hospital policy (at least once per day). Zero whenever:

- A new transducer or adapter cable is used.
- You reconnect the transducer cable to the monitor.
- The monitor restarts.
- You doubt the readings.

1. Turn off the stopcock to the patient.



2. Vent the transducer to the atmospheric pressure by turning on the stopcock to the air.
3. Press the **→0←** hardkey on the module, or, in the setup menu for the pressure (e.g. Art), select [**Art Zero** >>]→[**Zero**]. During zero calibration, the [**Zero**] button appears dimmed. It recovers after the zero calibration is completed. To zero all IBP channels, select [**Zero IBP**] hotkey, and then select [**Zero All Channels**] in the popup menu.
4. After the zero calibration is completed, close the stopcock to the air and open the stopcock to the patient.

NOTE

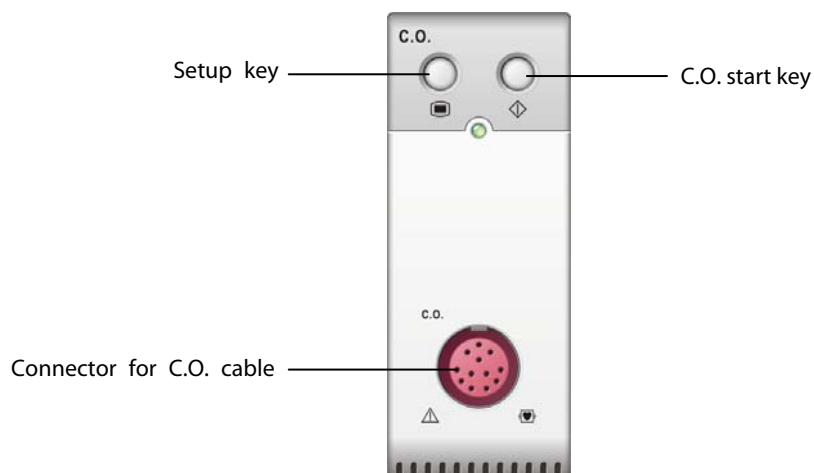
-
- **Your hospital policy may recommend that the ICP transducer is zeroed less frequently than other transducers.**
-

FOR YOUR NOTES

15 Monitoring Cardiac Output

15.1 Introduction

The cardiac output (C.O.) measurement invasively measures cardiac output and other hemodynamic parameters using the right heart (atria) thermodilution method. A cold solution of known volume and temperature is injected into the right atrium through the proximal port of a pulmonary artery (PA) catheter. The cold solution mixes with the blood in the right ventricle and the change in blood temperature is measured with a thermistor at the distal end of the catheter in the pulmonary artery. The temperature change is displayed as a curve in the C.O. split screen, and the monitor calculates the C.O. value from this curve. The C.O. value is inversely proportional to the area under the curve. As cardiac output varies continuously, a series of measurements must be carried out to achieve a reliable C.O. average value. Always use the average of multiple thermodilution measurements for therapy decisions. The monitor is capable of storing 6 measurements.



NOTE

- **C.O. monitoring is restricted to adult patients only.**
-

15.2 Understanding the C.O. Display

The C.O. measurement is displayed on the monitor as numeric C.O., C.I. and TB in the C.O. parameter window as shown below. To enter the [C.O. Setup] menu, select the C.O. parameter window.



1. Cardiac output
2. Time at which the C.O. average is calculated
3. Cardiac index
4. Blood temperature

15.3 Influencing Factors

The factors that affect cardiac output are:

- temperature of injectate solution,
- volume of injectate solution,
- patient's baseline blood temperature,
- patient's inspiratory/expiratory cycle,
- placement of catheter with relation to proximity of lung field,
- the catheter itself,
- the patient rhythm and hemodynamic status, and
- any other rapid IV solutions which are infused while the C.O. measurement is being performed.

Followings are some technique suggestions to obtain accurate C.O.:

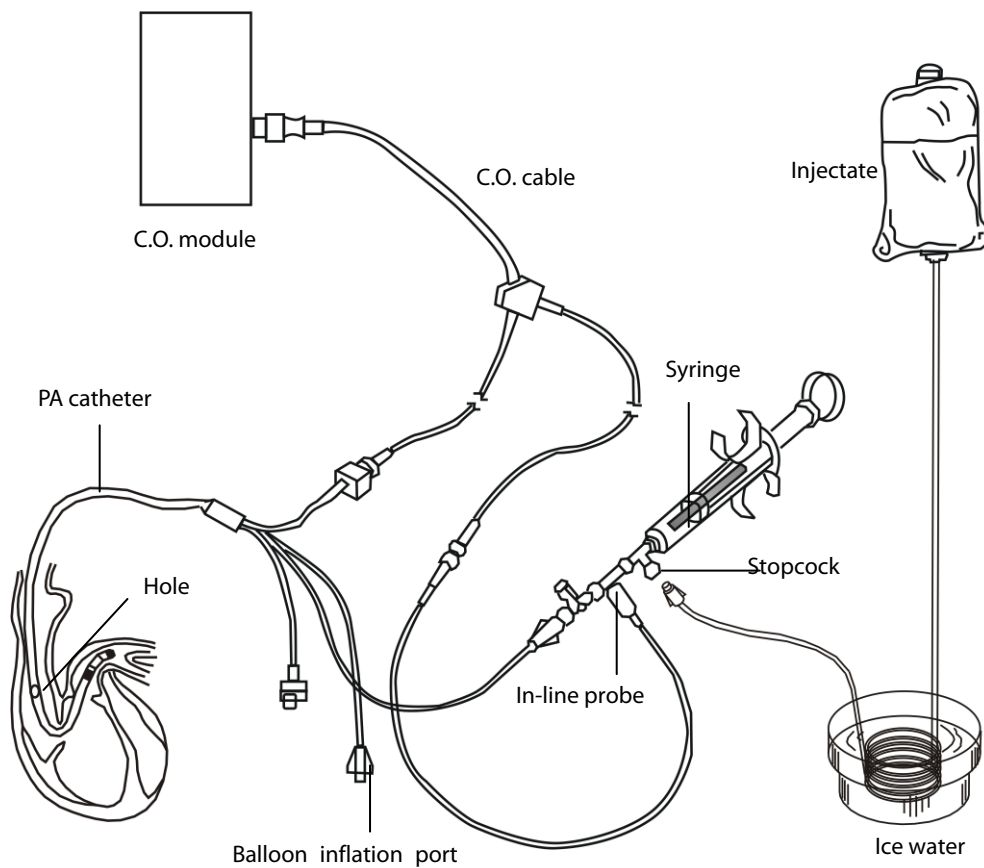
- Injectate solution must be cooler than the patient's blood.
- Inject solution rapidly and smoothly.
- Inject at end expiration.

15.4 Setting Up the C.O. Measurement

WARNING

- Use only accessories specified in this manual. Make sure that the accessories never come into contact with conductive parts.
-


1. Connect the C.O. cable to the C.O. connector.
2. Interconnect the C.O. module, catheter and syringe as shown below. Make sure that:
 - ◆ The module is securely inserted.
 - ◆ The PA catheter is in place in the patient.
 - ◆ The C.O. cable is properly connected to the module.



NOTE

- The above picture is connecting illustration when TI sensor PN 6000-10-02079 is used. The connection may be different if other TI sensors are used.
-

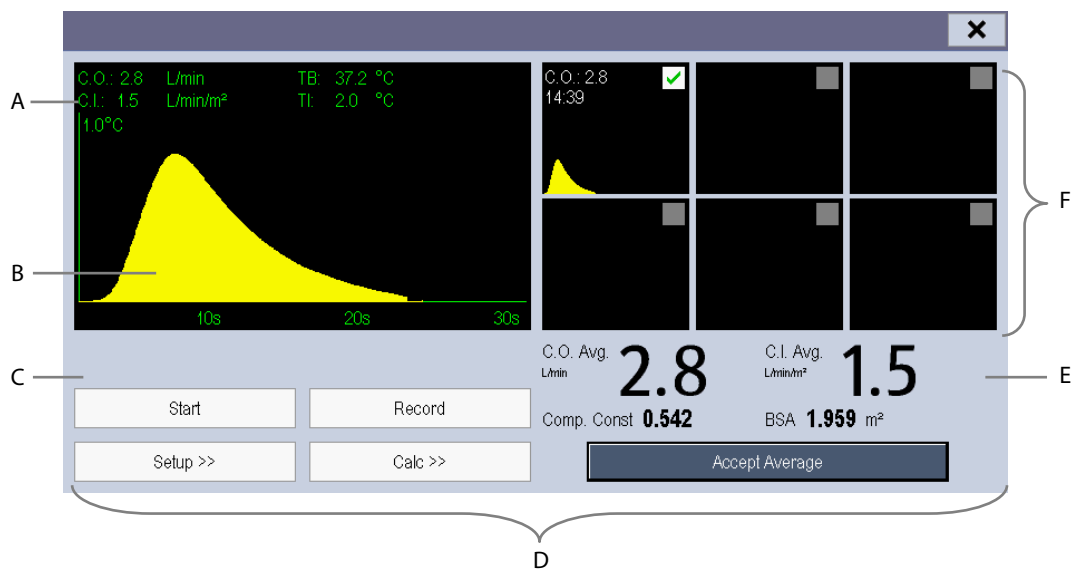3. Select the C.O. parameter window to enter the **[C.O. Setup]** menu. Check if the height and weight are appropriate for your patient. Change if necessary.
4. In the **[C.O. Setup]** menu:
 - ◆ Check that the correct computation constant is entered. Refer to the Instruction for Use of pulmonary artery catheter to determine the **[Comp. Const]** according to the entered injectate volume and temperature. To change the computation constant, select **[Comp. Const]** and then enter the correct value. When a new catheter is used, the computation constant should be adjusted in accordance with the manufacturer's instructions for use.
 - ◆ Set the **[Auto TI]** to **[Manual]** or **[Auto]**. If you select **[Auto]**, the system automatically detects the injectate temperature, and the **[Manual TI]** is disabled. If you select **[Manual]**, you need to enter the injectate temperature at **[Manual TI]** manually.
 - ◆ Set the **[Measuring Mode]** to **[Auto]** or **[Manual]**. In **[Auto]** mode, the monitor automatically takes the C.O. measurement after establishing a baseline blood temperature. In this mode, it is not necessary to select the **[Start]** button in the C.O. measurement window. In **[Manual]** mode, the monitor takes the C.O. measurement after **[Start]** button is selected.
5. Select **[Enter C.O. Screen]** to enter the C.O. measurements window.



- A. Currently measured numeric
- B. Currently measured C.O. curve
- C. Prompt message area
- D. Buttons
- E. Averaged values
- F. Measurement windows

6. Proceed as follows.

- ◆ In [**Manual**] measure mode, when you see the message [**Ready for new set of measurement**], select the [**Start**] button and then inject the solution quickly. As shown in the figure above, during the measurement, the currently measured thermodilution curve is displayed. At the end of the measurement, the thermodilution curve is transferred to one of the 6 measurement windows and the monitor prompts you to wait for a certain period of time before starting a new measurement.
- ◆ In [**Auto**] measure mode, the C.O. measurements can be performed consecutively, without the need for pressing the [**Start**] button between measurements. A new thermodilution measurement is possible as soon as [**Inject now!**] is displayed on the screen. The patient monitor automatically detects further thermodilution measurements.

7. Consecutively take 3 to 5 single measurements as instructed by Step 6.

A maximum of 6 measurements can be stored. If you perform more than six measurements without rejecting any, the oldest will automatically be deleted when a seventh curve is stored. Select from the 6 measurement curves and the system will automatically calculate and display the averaged C.O. and C.I. values. Then select the [**Accept Average**] button to accept and store the averaged values.

When injecting, the stopcock to the PA catheter is open and the stopcock to the injectate solution is closed. After the measurement is completed, turn off the stopcock to the PA catheter and turn on the stopcock to the injectate solution, and then draw the injectate solution into the injectate syringe.

In the buttons area, you can:

- Select [**Start**] to start a C.O. measurement.
- Select [**Stop**] to stop the current measurement.
- Select [**Cancel**] during a measurement to cancel the measurement.
- Select [**Record**] to print out the curves selected for average calculation, numerics and averaged values by the recorder.
- Select [**Setup >>**] to access the [**C.O. Setup**] menu.
- Select [**Calc >>**]→[**Hemodynamic >>**] to access the [**Hemodynamic Calculation**] menu.

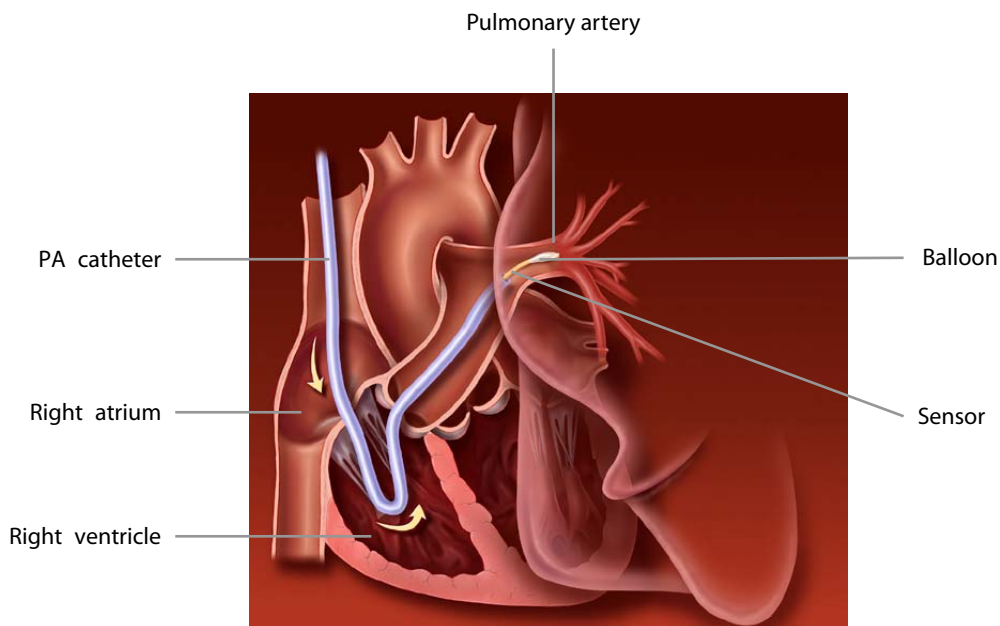
The system can automatically adjust the X-axis scale range to 30 s or 60 s and Y-axis scale range to 0.5°C, 1.0°C, or 2.0°C.

NOTE

- Starting measuring without blood temperature being stable yet may cause measuring failure.
 - During the cardiac output measurement, blood temperature alarms are inactive.
 - Please refer to the Instructions for Use of pulmonary artery catheter delivered with the patient monitor to determine the [Comp. Const] and the volume of injectate.
-

15.5 Measuring the Blood Temperature

As shown below, the blood temperature is measured with a temperature sensor at the distal end of the catheter in the pulmonary artery. During C.O. measurements, blood temperature alarms are suppressed to avoid false alarms. They will automatically recover as soon as the C.O. measurements are completed.



15.6 Changing C.O. Settings

15.6.1 Setting the Temperature Unit

Select [**Unit Setup >>**] from the [**User Maintenance**] menu. In the popup menu, select [**Temp Unit**] to toggle between [**°C**] and [**°F**].

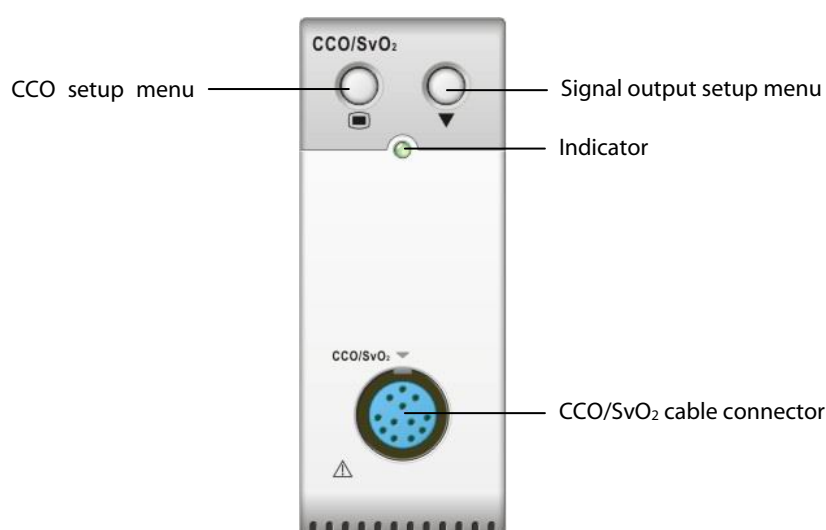
15.6.2 Setting Alarm Properties

Select [**Alarm Setup >>**] from the [**C.O. Setup**] menu. You can set alarm properties for this parameter in the popup menu.

16 Monitoring CCO/SvO₂

16.1 Introduction

Edwards Vigilance II® monitor / Vigileo™ monitor measures continuous cardiac output (CCO), mixed venous oxygen saturation (SvO₂), central venous oxygen saturation (ScvO₂) etc. It also calculates hemodynamic and oxygenation parameters. This patient monitor can be connected to the Vigilance II® monitor / Vigileo™ monitor and can display, store, and review the measured and calculated parameter values from the Vigilance II® monitor / Vigileo™ monitor. This patient monitor can also give alarms of these measured parameters. You must set alarm on/off, alarm limits, alarm level, and alarm record separately on this monitor. The alarm is Off by default.



16.2 Safety

WARNING

- **The Vigilance II® monitor and Vigileo™ monitor are manufactured by Edwards Lifesciences. This company provides the technology of measuring and calculating the relevant parameters. We only provide the connection between this patient monitor and Vigilance II® monitor/ Vigileo™ monitor.**
 - **If you have any doubts about the operation and maintenance of the Vigilance II® monitor/ Vigileo™ monitor, please read the Vigilance II® monitor/ Vigileo™ monitor Operator's Manual or contact Edwards Lifesciences (www.edwards.com) directly.**
 - **Fully observe the Vigilance II® monitor/ Vigileo™ monitor Operator's Manual to make settings and to connect the monitor with the patient.**
-

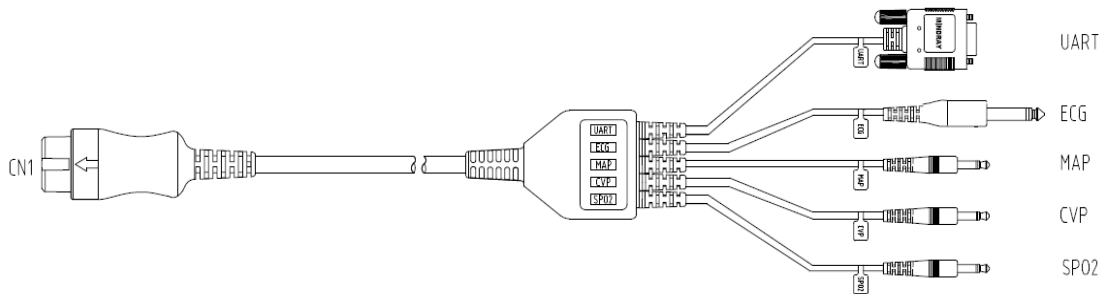
16.3 Automatic Communication Detection

The relevant parameter window is not displayed on the screen if this patient monitor detects communication failure between the CCO/SvO₂ module and Vigilance II® monitor / Vigileo™ monitor automatically.

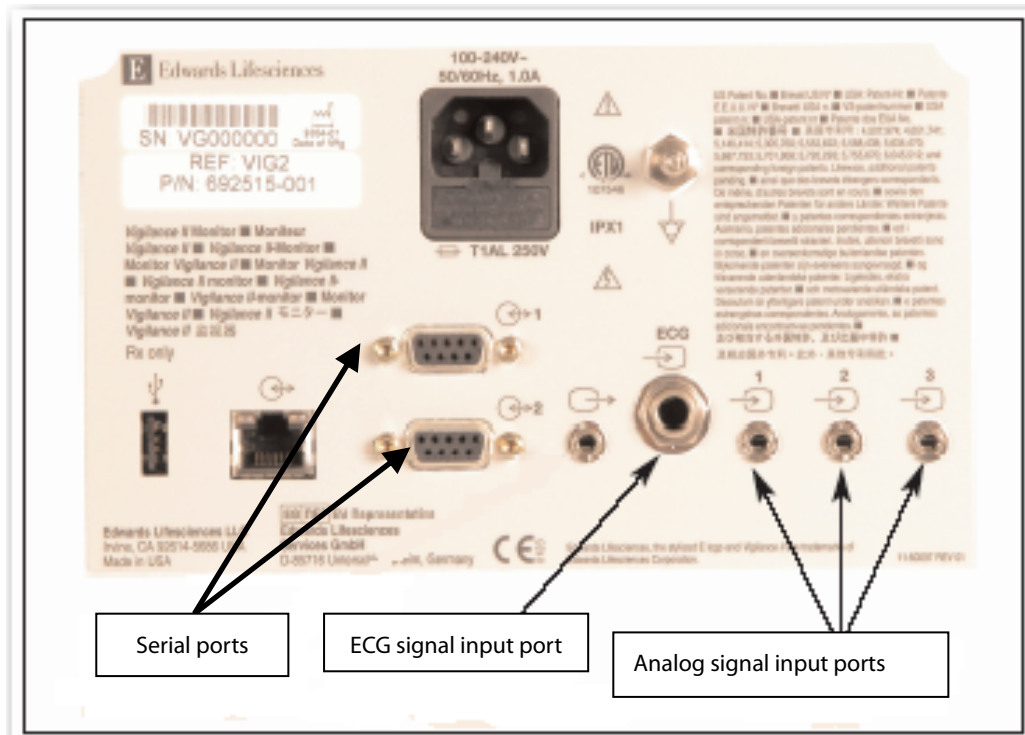
16.4 Connecting the Device

16.4.1 Connecting the Vigilance II® Monitor

The following figure shows how to connect this patient monitor to the Vigilance II® monitor through cables.

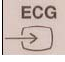


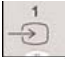
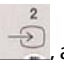

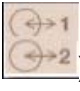
The following figure shows the rear housing of the Vigilance II® monitor.



To connect the Vigilance II® monitor,

1. Connect CN1 with the CCO/SvO₂ connector on the patient monitor.

2. Insert the ECG signal end into the ECG signal input port marked  on the rear housing of the Vigilance II® monitor.

3. Insert the MAP signal end into the analog signal input port 1 marked , the CVP signal end into port 2 marked , and SPO₂ signal end into port 3 marked  respectively on the rear housing of the Vigilance II® monitor.
4. Insert UART into either of the serial ports (marked ) on the rear housing of the Vigilance II® monitor
5. Set the Vigilance II® monitor as follows:
 - Access the **[Serial Port Setup]** menu.
 - ◆ Set **[Device]** to **[IFMout]**, **[Baud Rate]** to **[19200]**, **[Parity]** to **[None]**, **[Stop Bits]** to **[1]**, **[Data Bits]** to **[8]**, and **[Flow Control]** to **[2 s]**.
 - Access the **[Analog Input Setup]** menu.
 - ◆ For port 1, set **[Parameter]** to **[MAP]**, **[Voltage Range]** to **[0-5 v]**, **[Full Scale Range]** to **500 mmHg** (66.7 kPa), **[Simulated High Value]** to **500 mmHg** (66.7 kPa), and **[Simulated Low Value]** to **0 mmHg** (0.0 kPa).
 - ◆ For port 2, set **[Parameter]** to **[CVP]**, **[Voltage Range]** to **[0-5 v]**, **[Full Scale Range]** to **100 mmHg** (13.3 kPa), **[Simulated High Value]** to **100 mmHg** (13.3 kPa), and **[Simulated Low Value]** to **0 mmHg** (0.0 kPa).
 - ◆ For port 3, set **[Parameter]** to **[SaO₂]**, **[Voltage Range]** to **[0-10 v]**, **[Full Scale Range]** to **[100%]**, **[Simulated High Value]** to **[100%]**, and **[Simulated Low Value]** to **[0%]**.

Refer to the Vigilance II® Operator's Manual for the operation of the monitor.

WARNING

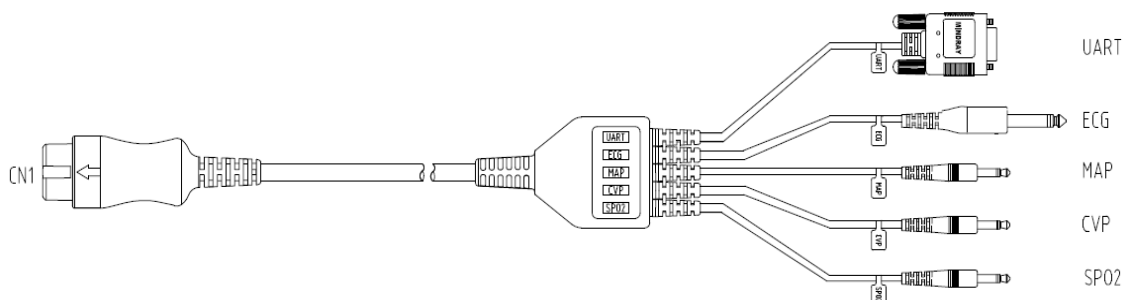
- **Calibrate the Vigilance II® monitor before monitoring. Refer to the Vigilance II® Operator's Manual for the calibration of the monitor.**

Notes


- **For the Vigilance II® monitor, [Flow Control] must be set to 2 seconds.**

16.4.2 Connecting the Vigileo™ Monitor

The following figure shows how to connect this patient monitor to the Vigileo™ monitor through cables.



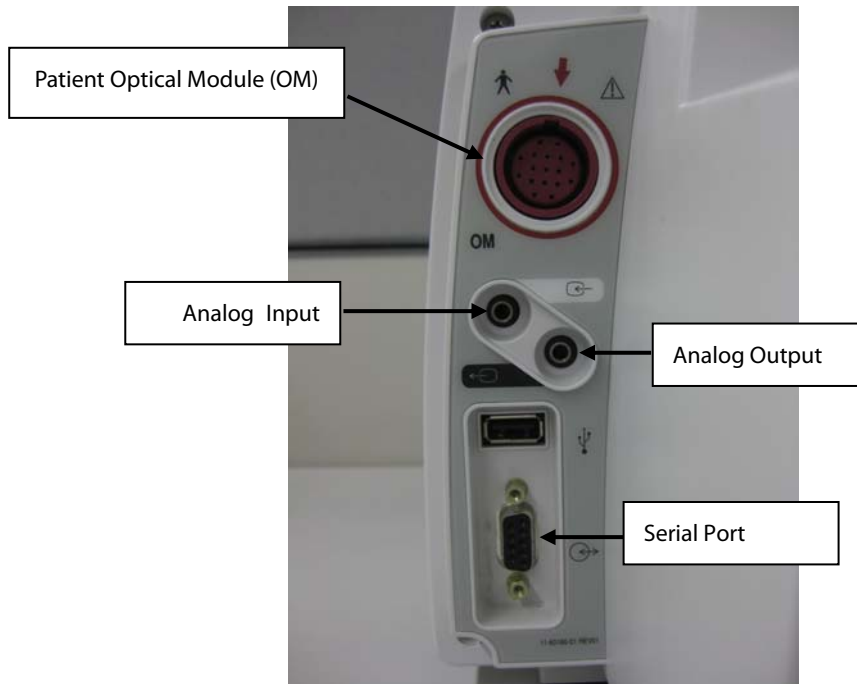The following figure shows the rear housing of the Vigileo™ monitor.



To connect the Vigileo™ monitor,

1. Connect CN1 with the CCO/SvO₂ connector on the patient monitor.
2. Insert the CVP signal end into the analog signal input port on the rear housing of the Vigileo™ monitor.
3. Insert UART into the serial port on the rear housing of the Vigileo™ monitor.
4. Set the Vigileo™ monitor as follows:
 - Access the [**Serial Port Setup**] menu.
 - ◆ Set [**Device**] to [**IFMout**], [**Baud Rate**] to [**19200**], [**Parity**] to [**None**], [**Stop Bits**] to [**1**], [**Data Bits**] to [**8**], and [**Flow Control**] to [**2 seconds**].
 - Access the [**Analog Input Port Setup**] menu.
 - ◆ Set [**Parameter**] to [**CVP**], [**Voltage Range**] to [**0-5 v**], [**Full Scale Range**] to **100 mmHg** (13.3 kPa), [**Simulated High Value**] to **100 mmHg** (13.3 kPa), and [**Simulated Low Value**] to **0 mmHg** (0.0 kPa).

Refer to the Vigileo™ Operator's Manual for the operation of the monitor.

 **WARNING**

- **Calibrate the Vigileo™ monitor before monitoring. Refer to the Vigileo™ Operator's Manual for the calibration of the monitor.**
-

Notes

- **For the Vigileo™ monitor, [Flow Control] must be set to 2 seconds.**
-

16.5 Understanding CCO Parameters

When the patient monitor connects Vigilance II® monitor, by selecting the CCO parameter window→[**Hemodynamic Parameters >>**], you can view the hemodynamic parameters for evaluation of the patient's hemodynamic status.

Abbreviation	Unit	Full spelling
CCO	L/min	continuous cardiac output
CCI	L/min/m ²	continuous cardiac index
C.O.	L/min	cardiac output
C.I.	L/min/m ²	cardiac index
EDV	ml	end diastolic volume
EDVI	ml/m ²	end diastolic volume index
SV	ml/b	stroke volume
SVI	ml/b/m ²	stroke volume index
SVR	DS/cm ⁵ or kPa-s/l	systemic vascular resistance
SVRI	DS-m ² /cm ⁵ or kPa-s-m ² /l	systemic vascular resistance index
RVEF	%	right ventricular ejection fraction
BT	°C or °F	blood temperature
ESV	ml	end systolic volume
ESVI	ml/m ²	end systolic volume index
CVP	cmH ₂ O, kPa or mmHg	central venous pressure
MAP	mmHg or kPa	mean arterial pressure
HR	rpm	heart rate

When the patient monitor connects Vigileo™ monitor, by selecting the CCO parameter window→[**Hemodynamic Parameters >>**], you can view the hemodynamic parameters for evaluation of the patient's hemodynamic status.

Abbreviation	Unit	Full spelling
CCO	L/min	continuous cardiac output
CCI	L/min/m ²	continuous cardiac index
SV	ml/b	stroke volume
SVI	ml/b/m ²	stroke volume index
SVV	%	stroke volume variation
SVR	DS/cm ⁵ or kPa-s/l	systemic vascular resistance
SVRI	DS-m ² /cm ⁵ or kPa-s-m ² /l	systemic vascular resistance index
CVP	cmH ₂ O, kPa or mmHg	central venous pressure
MAP	mmHg or kPa	mean arterial pressure
HR	rpm	heart rate

16.6 Understanding the CCO Display

- When the patient monitor connects Vigilance II® monitor:

In the continuous measurement mode, the CCO parameter window displays the values of one primary parameter and up to three secondary parameters. You can select the desired parameters to be displayed through the menu. The default secondary parameters displayed are SVR, EDV and SV.



In the intermittent measurement mode, the CCO parameter window displays the values of two primary parameters and two secondary parameters. You can also select the desired parameter to be displayed through the menu.

- When the patient monitor connects Vigileo™ monitor:

The CCO parameter window displays the values of one primary parameter and up to three secondary parameters. You can select the desired parameters to be displayed through the menu. The default secondary parameters displayed are SVR, SVV and SV.

16.7 Changing CCO Settings

16.7.1 Selecting Vascular Resistance Unit

1. Access the [**CCO Setup**] menu.
2. Select [**SVR Unit**] and toggle between [**DS/cm5**] and [**kPa-s/l**].

16.7.2 Selecting the Displayed Parameters

1. Access the [**CCO Setup**] menu.
2. Select [**Select Parameters >>**].
3. Select the parameters to be displayed from the pop-up menu.

16.7.3 Checking the C.O. Measurements

When the patient monitor connects Vigilance II® monitor, you can check the C.O. measurements in the intermittent measurement mode.

1. Access the [**CCO Setup**] menu.
2. Select [**C.O. Measurements >>**].

16.7.4 Setting Signal Output

- When the patient monitor connects Vigilance II® monitor:

This patient monitor outputs analog signals for the Vigilance II® monitor. You can select [**Signal Output Setup >>**] from the [**CCO Setup**] menu to set the source of MAP signals. You can also select [**Simulated High Value**] or [**Simulated Low Value**] to provide simulated high value or low value signals for calibrating the Vigilance II® monitor. Refer to the Vigilance II® Operator's Manual for the calibration of the monitor.

- When the patient monitor connects Vigileo™ monitor:

Select [**Signal Output Setup >>**] from the [**CCO Setup**] menu. In the popup menu, you can select [**Simulated High Value**] or [**Simulated Low Value**] to provide simulated high value or low value signals for calibrating the Vigileo™ monitor. Refer to the Vigileo™ Operator's Manual for the calibration of the monitor.

16.7.5 Selecting Alarm Properties

You can select [**Alarm Setup >>**] from the [**CCO Setup**] menu to set the alarm properties for the relevant parameters.



- **Because the alarm limits of the relevant measured parameters can be set on this patient monitor, the alarms of these parameters on this patient monitor may be different from those on the Vigilance II® / Vigileo™ monitor. Please pay special attention to the alarms on the Vigilance II® / Vigileo™ monitor.**
 - **The alarm of the relevant measured parameters on this patient monitor is Off by default. Please pay special attention to the alarms on the Vigilance II® / Vigileo™ monitor.**
-

16.8 Understanding SvO₂ Parameters

When the patient monitor connects Vigilance II® monitor, by selecting the SvO₂ parameter window → [**SvO₂ Setup**] → [**Oxygenation Parameters >>**], you can view all the oxygenation parameters.

Abbreviation	Unit	Full spelling
SvO ₂	%	mixed venous oxygen saturation
ScvO ₂	%	central venous oxygen saturation
SaO ₂	%	arterial oxygen saturation
DO ₂	ml/min	oxygen delivery
VO ₂	ml/min	oxygen consumption
O ₂ El	%	oxygen extraction index

16.9 Understanding the SvO₂ Display

The parameter window displays the primary parameter, secondary parameter and SQI bargraph.



16.10 Changing SvO₂ Settings

16.10.1 Setting Signal Output

This patient monitor outputs analog signals for the Vigilance II® monitor. You can select [**Signal Output Setup >>**] from the [**SvO₂ Setup**] menu to set the source of MAP signals. You can also select [**Simulated High Value**] or [**Simulated Low Value**] to provide simulated high value or low value signals for the Vigilance II® monitor. Refer to the Vigilance II® Operator's Manual for the calibration of the monitor.

16.10.2 Selecting Alarm Properties

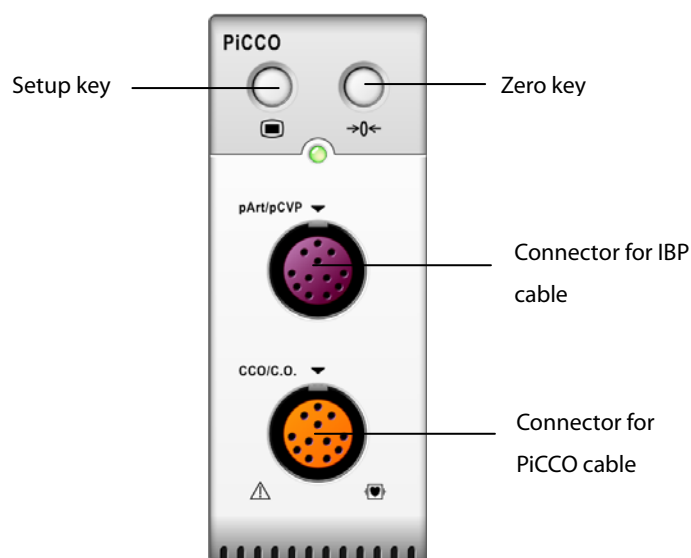
When the patient monitor connects Vigilance II® monitor, you can select [**Alarm Setup >>**] from the [**SvO₂ Setup**] menu to set the alarm properties for the relevant parameters.

When the patient monitor connects Vigileo™ monitor, select SvO₂ or ScvO₂ parameter area. You can set the alarm properties for the relevant parameters in the popup menu.

17 Monitoring PiCCO

17.1 Introduction

The PiCCO method combines transpulmonary thermodilution and pulse contour analysis on the blood pressure waveform. A cold bolus (e.g. normal saline 0.9%) with a known volume and temperature is injected into the right atrium through a central venous catheter. The cold bolus mixes with the blood in the heart and the change in blood temperature is measured with a thermistor at the distal end of the arterial thermodilution catheter placed in one of the bigger systemic arteries, for example, the femoral artery. The monitor uses the transpulmonary thermodilution method to measure C.O., GEDV (Global End Diastolic Volume) and EVLW (Extra Vascular Lung Water). With the C.O. value measured with the transpulmonary thermodilution method and the result of the pulse contour analysis, a patient-specific calibration factor is calculated. The monitor uses this value to compute CCO and the other continuous hemodynamic parameters.



17.2 Safety Information

WARNING

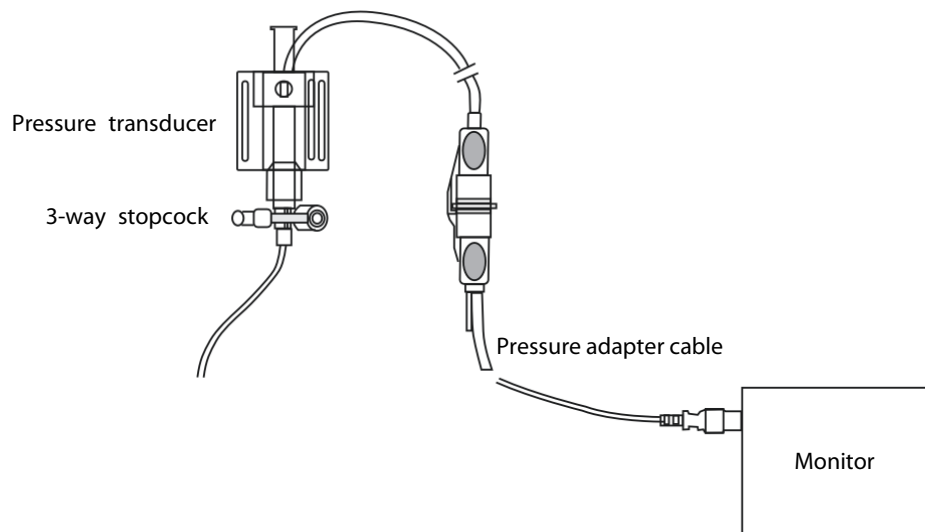
- **PiCCO monitoring is restricted to adult and pediatric patients.**
 - **Use only pressure transducers specified in this manual. Never reuse disposable pressure transducers.**
 - **Make sure that the applied parts never contact other conductive parts.**
 - **To reduce the hazard of burns during high-frequency surgical procedure, ensure that the monitor's cables and transducers never come into contact with the high-frequency surgical units.**
 - **When using accessories, their operating temperature should be taken into consideration. For details, refer to instructions for use of accessories.**
-

17.3 Zeroing the Transducer

To avoid inaccurate pressure readings, the monitor requires a valid zeroing. Zero the transducer in accordance with your hospital policy (at least once per shift). Zero whenever:

- A new transducer or adapter cable is used.
- You reconnect the transducer cable to the monitor.
- The monitor restarts.
- You doubt the readings.


1. Turn off the stopcock to the patient.



2. Vent the transducer to the atmospheric pressure by turning on the stopcock to the air.
3. Press the **→0←** hardkey on the module, or, in the setup menu for the pressure (e.g. pArt), select **[pArt Zero >>]→[Zero]**. During zero calibration, the **[Zero]** button appears dimmed. It recovers after the zero calibration is completed. To zero all IBP channels, select **[Zero IBP]** hotkey, and then select **[Zero All Channels]** in the popup menu.
4. After the zero calibration is completed, close the stopcock to the air and open the stopcock to the patient.

17.4 Setting up the PiCCO Measurements

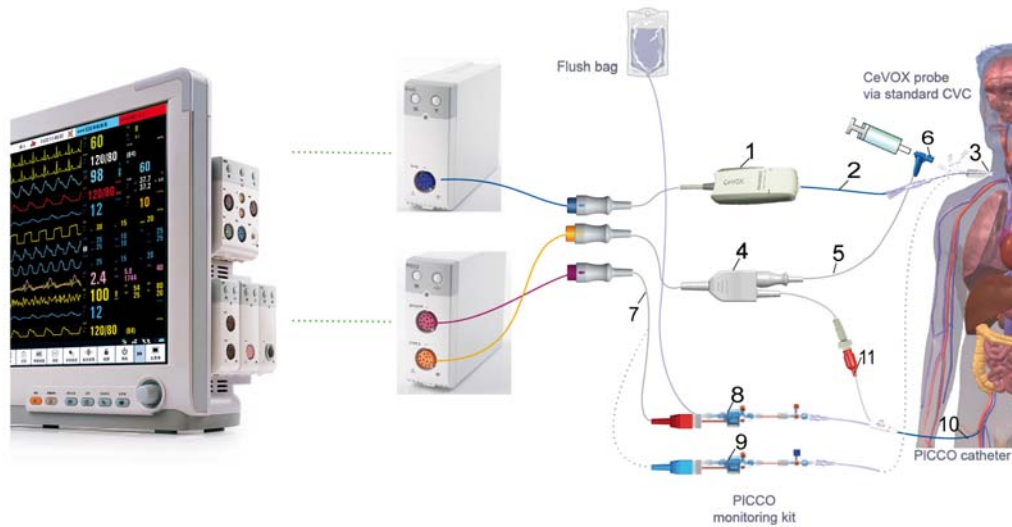Please refer to the following figure and procedure to set up the PiCCO measurements:



1. Optical module
2. CeVOX fiberoptic probe
3. Central venous catheter
4. PiCCO cable
5. Injectate temperature sensor cable
6. Injectate temperature sensor
7. IBP cable
8. Arterial pressure transducer
9. CVP transducer
10. Arterial thermodilution catheter
11. Blood temperature sensor

17.5 Preparation for PiCCO Measurements

1. Place the arterial thermodilution catheter.

WARNING

- **The arterial thermodilution catheter must be placed in one of the bigger systemic arteries, for example, the femoral, the brachial or the axillary artery.**
 - **You must use the approved catheters and puncture locations.**
-

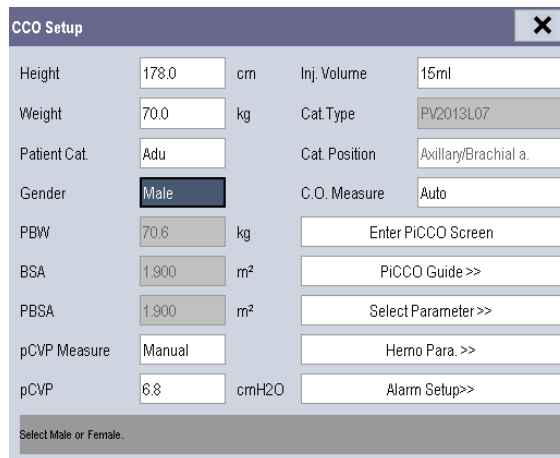
2. Place the central venous catheter.
3. Connect the injectate temperature sensor to the central venous catheter.

4. Plug the PiCCO cable into the CCO/C.O. connector on the PiCCO module, and connect the following devices to the PiCCO cable:
 - ◆ Injectate temperature sensor probe
 - ◆ Blood temperature sensor connector.
5. Connect one end of the arterial pressure transducer to the arterial thermodilution catheter and the other end to the IBP cable marked with pArt.

WARNING

- **If air bubbles appear in the tubing system, flush the system with the infusion solution again. Air bubble may lead to wrong pressure reading.**

6. Connect one end of the CVP transducer to the central venous catheter and the other end to the IBP cable marked with pCVP (neglect this procedure if CVP measurement is not performed). Then plug the IBP cable to the pArt/pCVP connector on the PiCCO module.
7. Access the **[CCO Setup]** menu by selecting **[PiCCO Measurement]** → **[Setup>>]** → **[CCO Setup]**. You can also select **[Main Menu]** → **[Parameters]** → **[CCO Setup>>]** to access the **[CCO Setup]** menu.



The screenshot shows the 'CCO Setup' menu with the following fields and values:

Height	178.0	cm	Inj. Volume	15ml
Weight	70.0	kg	Cat. Type	PV2013L07
Patient Cat	Adu		Cat. Position	Axillary/Brachial a.
Gender	Male		C.O. Measure	Auto
PBW	70.6	kg	Enter PiCCO Screen	
BSA	1.900	m ²	PiCCO Guide >>	
PBSA	1.900	m ²	Select Parameter >>	
pCVP Measure	Manual		Hemo Para. >>	
pCVP	6.8	cmH2O	Alarm Setup >>	

At the bottom, there is a note: "Select Male or Female."

8. Check that the correct arterial catheter constant is displayed at **[Cat.Type]** in **[CCO Setup]** menu. The monitor can recognize the arterial catheter automatically when the PiCCO cable is connected to the CCO/C.O. connector.

NOTE

- **If the catheter constant is not recognized, enter the correct value for the catheter in the [Cat.Type] edit box. The catheter constant is usually written either on the catheter or on the catheter packaging.**

9. Set up the patient information in **[CCO Setup]** menu.

NOTE

- **Correct input of height, weight, category and gender is mandatory for the accuracy of the displayed parameters as well as for the correct indexing of some parameters.**
- **Input a proper pCVP value in the [CCO Setup] menu if CVP is not measured. The system adopts 5mmHg by default if the pCVP value is neither measured nor input manually.**

10. Enter the **[CCO Setup]** menu to select the injectate volume. If the injectate volume is not selected, the system sets the volume by default, which is 15ml for adult and 10 ml for pediatric. The following table displays the recommended injectate volume depending on body weight and ELWI (Extravascular Lung Water Index):

	ELWI < 10	ELWI > 10	ELWI < 10
Patient Weight (kg)	Iced Injectate	Iced Injectate	Room Temperature Injectate
<3	2ml	2ml	3ml
<10	2ml	3ml	3ml
<25	3ml	5ml	5ml
<50	5ml	10ml	10ml
<100	10ml	15ml	15ml
≥100	15ml	20ml	20ml

11. Set up the C.O. measure mode by selecting **[C.O. Measure]** from the **[CCO Setup]** menu, and toggling between **[Auto]** and **[Manual]**.
- If you select **[Manual]**, you should start each measurement manually by pressing the **[Start]** key in the **[PiCCO Measurement]**.
 - If you select **[Auto]**, the C.O. measurements can be performed consecutively, without the need for pressing the **[Start]** key.

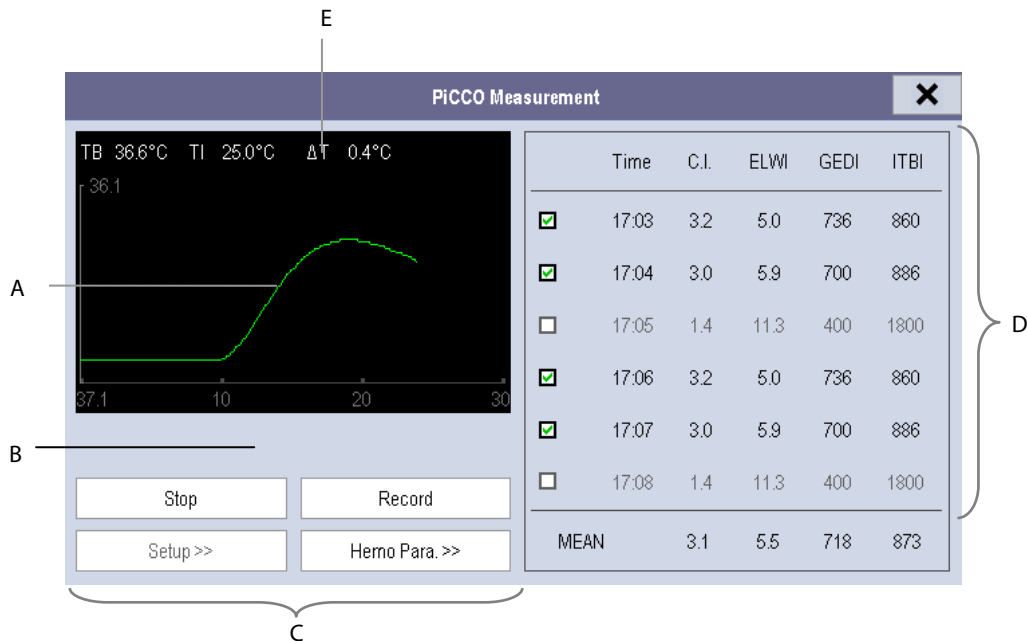
NOTE

- **Steps 8 to 10 can also be conducted with the [C.O. Measure (Transpulmonary) Setup Guide] menu, which can be accessed by selecting [PiCCO Guide>>] in [CCO Setup] menu. In order to ensure correct PiCCO calibration, please be sure the information you have entered is correct.**

17.6 Performing PiCCO Measurements and CCO Calibration

Please perform the PiCCO measurements according to the following procedure:

Open the [PiCCO Measurement] menu.



- A. Thermodilution curve
 - B. Prompt message area
 - C. Buttons
 - D. History window
 - E. Measurement quality: ΔT
2. Select the [Start] button and inject the bolus rapidly (<7sec) and smoothly as soon as the message [Inject xx ml!] and prompt tone appear. As shown in the figure above, during the measurement, the currently measured thermodilution curve is displayed. At the end of the measurement, the measured values are displayed in the history window and the monitor prompts you to wait for a certain period of time before starting a new measurement. The ΔT value should be greater than 0.15°C to ensure high accuracy. A low ΔT can be caused by a very high ELWI or an extreme low CI. If ΔT is too low, you can try to increase it by
 - Injecting more volume (remember to reenter the injectate volume in [CCO Setup] menu before injecting).
 - Injecting colder bolus.
 - Injecting the bolus in a shorter time.
 3. Perform 3 to 5 single measurements direct after each other within a maximum of 10 minutes as described in Step 2. A new measurement is available when you see the blood temperature is steady in the [PiCCO Measurement] window.
 - If you've selected [Manual] measure in the [CCO Setup] menu, you should repeat Step 2 manually.
 - If you've selected [Auto] measure in the [CCO Setup] menu, the C.O. measurements can be performed consecutively, without the need for pressing the [Start] button between measurements. A new thermodilution

measurement is possible as soon as [**Inject xx ml!**] is displayed on the screen. The patient monitor automatically detects further thermodilution measurements.

4. A maximum of 6 measurements can be stored. If you perform more than six measurements without rejecting any, the oldest will be automatically deleted when a seventh curve is stored. Select the measurement values and the system will automatically perform calibration and calculate the averaged CCO and CCI values.

In the buttons area, you can:

- Select [**Stop**] during a measurement to stop the measurement.
- Select [**Record**] to print out the curves selected for average calculation, numerics and averaged values by the recorder.
- Select [**Setup >>**] to access the [**C.O. Setup**] menu.
- Select [**Hemo Para.>>**] to access the [**Hemodynamic Parameters**] menu.

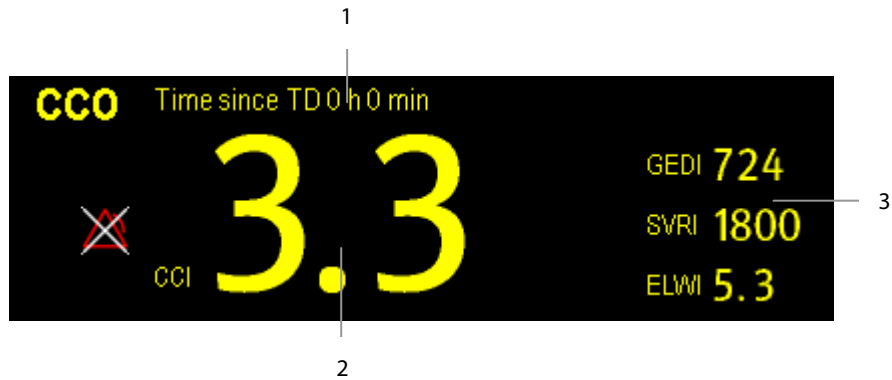


CAUTION

- **Three to five single thermodilution measurements within 10 minutes are recommended. For a stable patient it is recommended to perform a thermodilution measurement every 8 hours. For an unstable patient it may be necessary to perform thermodilution measurements more frequently in order to determine the patient's volume status and to recalibrate the continuous determination of C.O..**
 - **As the pulse contour cardiac output of children has not been sufficiently validated thus far, the C.O. should be checked by thermodilution before therapeutic interventions.**
 - **If the system can not get a reliable pArt value during a C.O. measure, the corresponding C.O. value is invalid for PiCCO calibration.**
 - **Recalibration is recommended with significant changes in hemodynamic conditions, such as volume shifts or changes to medication.**
 - **If the option of the continuous CVP measurement is not used, CVP should be updated as soon as a new value is obtained to accurately calculate SVR and CCO.**
 - **If the displayed continuous parameters are not plausible, they should be checked by a thermodilution measurement. The PiCCO measurement will be recalibrated automatically.**
 - **Faulty measurements can be caused by incorrectly placed catheters, interfering signal transmission e.g. of arterial pressure, defective connections or sensors, or by electromagnetic interference (e.g. electric blankets, electric coagulation).**
 - **Aortic aneurysms may cause the displayed blood volume (GEDV/ITBV) derived by thermodilution measurement to be erroneously high if the arterial thermodilution catheter is placed in the femoral artery.**
-

17.7 Understanding the Displayed PiCCO Parameters


17.7.1 Understanding the CCO Display



1. Prompt message: the time since previous TD measurement
2. Label and value for main parameter
3. Labels and values for secondary parameters

17.7.2 Understanding the pArt Display

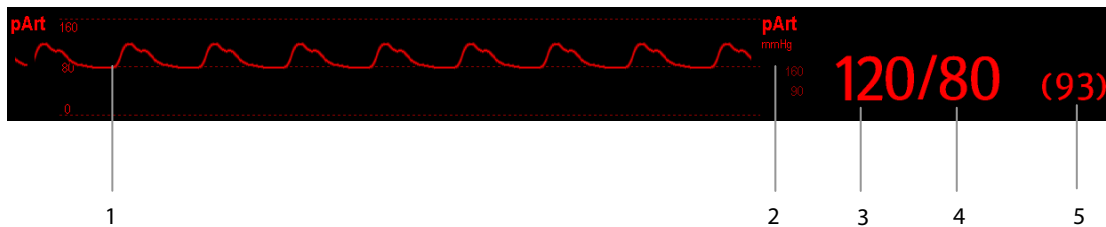
The artery pressure is displayed on the monitor as a waveform and numeric pressures. The figure below shows the pArt waveform and numerics.



1. Waveform
2. Pressure unit
3. Systolic pressure
4. Diastolic pressure
5. Mean pressure

17.7.3 Understanding the pCVP Display

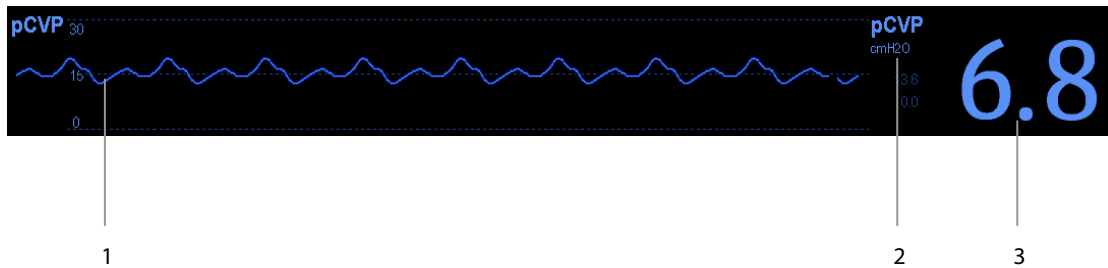The central venous pressure is displayed on the monitor as a waveform and numeric pressures. The figure below shows the pCVP waveform and numerics.



1. Waveform
2. Pressure unit
3. Central venous pressure

17.8 Understanding PiCCO Parameters

You can enter the **[Hemodynamic Parameters]** menu either by:


- Accessing the **[CCO Setup]** menu and selecting **[Hemo Para.>>]**, or
- Accessing the **[PiCCO Measurement]** menu and selecting **[Hemo Para.>>]**.

17.8.1 Spider Vision


17.8.1.1 Spider Vision Diagram

The spider vision diagram shows all continuous parameters in dynamic conjunction.

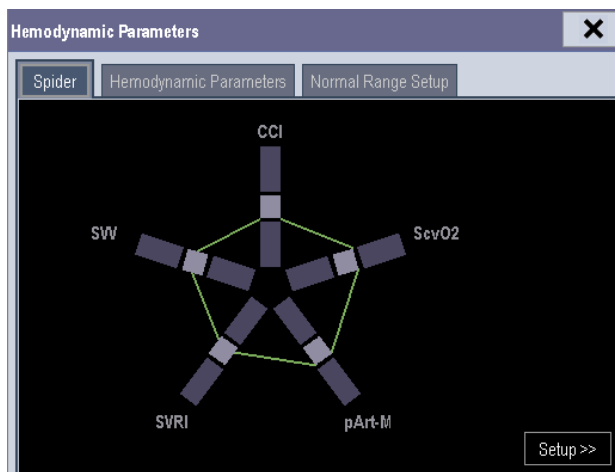
Each spider leg is divided into 3 segments indicating different value ranges for the respective parameters. The segment in the middle indicates the normal range for the respective parameter. The outer segment will be highlighted when corresponding parameter value exceeds the upper limit. The inner segment will be highlighted when its corresponding parameter value exceeds the lower limit.



The diagram is displayed GREEN when all displayed parameters are within the normal range.

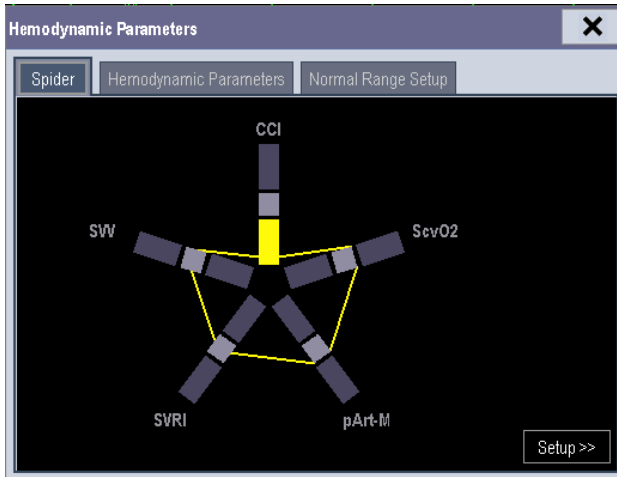

The diagram is displayed YELLOW immediately when one of the displayed parameters goes outside the normal range.

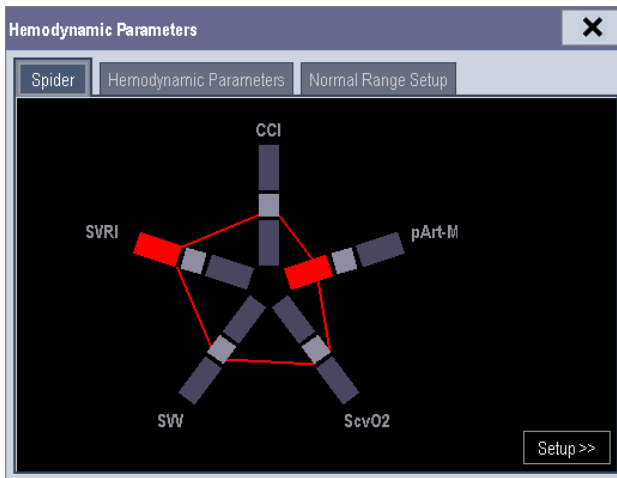


The diagram appears RED when two or more displayed parameters are outside the normal range.

The parameter whose default normal range is changed will be marked with the symbol



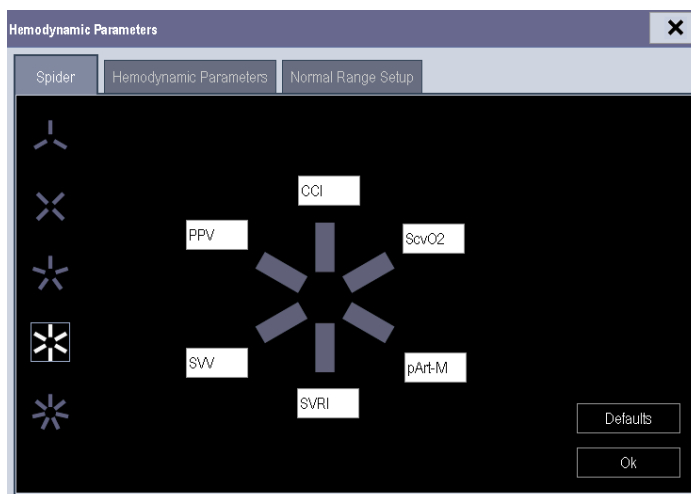17.8.1.2 Spider Configuration



The spider vision diagram can be configured individually. You can select [**Setup>>**] in the spider vision screen and set the diagram by the following procedure:

1. Select the number of spider legs (3to7).
2. Select the parameter to be displayed.

17.8.2 Hemodynamic Parameters

Select [**Hemodynamic Parameters**] tab from the [**Hemodynamic Parameters**] menu to view the patient's hemodynamic parameters. In the [**Hemodynamic Parameters**] menu, you can select [**Range**] to view the referential normal range of each parameter. If a parameter value exceeds its normal range, the system will add a " ↑ " or " ↓ " to the right of the parameter.

The screenshot shows a window titled "Hemodynamic Parameters" with three tabs: "Spider", "Hemodynamic Parameters", and "Normal Range Setup". The "Hemodynamic Parameters" tab is active, displaying a table of parameters. The parameters are grouped into four columns: Output, Preload Volume, Afterload, and Oxygenation Parameters. The values are as follows:

Output	Preload Volume	Afterload	Oxygenation Parameters
CCO 5.9 L/min	GEDV 1347 ml	pArt-M 90 mmHg	ScvO2 80 %
CCI 3.1 L/min/m ²	GEDI 718 ml/m ²	pArt-S 120 mmHg	Hb 15.0 g/dl
SV 95 ml	ITBV 1659 ml	pArt-D 75 mmHg	DO2 855 ml/min
SVI 50 ml/m ²	ITBI 873 ml/m ²	Organ Function	DO2I 450 ml/min/m ²
HR 60 bpm	SW 6 %	EVLW 392 ml	VO2 247 ml/min
Contractility	PPV 6 %	ELWI 5.5 ml/kg	VO2I 130 ml/min/m ²
GEF 35 ↑ %	Afterload	CPO 1.1 W	SaO2 95.0 %
CFI 5.2 1/min	SVR 1000 DS/cm ⁵	CPI 0.6 W/m ²	
dPmx 243 mmHg/s	SVRI 1900 DS·m ² /cm ⁵	PVPI 1.5	
		TB 37.0 °C	

At the bottom of the window, there are two buttons: "Range" and "Record".

	Abbreviation	Full Spelling	Unit	Default Normal Range
Output	CCO	Continuous Cardiac Output	L/min	/
	CCI	Continuous Cardiac Index	L/min/m ²	3.0-5.0
	SV	Stroke Volume	ml/b	/
	SVI	Stroke Volume Index	ml/b/m ²	40-60
	HR	Heart Rate	bpm	60-80
Contractility	GEF	Global Ejection Fraction	%	25-35
	CFI	Cardiac Function Index	1/min	4.5-6.5
	dPmx	Left Ventricular Contractility	mmHg/s	/
Preload Volume	GEDV	Global End Diastolic Volume	ml	/
	GEDI	Global End Diastolic Volume Index	ml/m ²	680-800
	ITBV	Intrathoracic Blood Volume	ml	/
	ITBI	Intrathoracic Blood Volume Index	ml/m ²	850-1000
	SVV	Stroke Volume Variation	%	0-10
	PPV	Pulse Pressure Variation	%	0-10
Afterload	SVR	Systemic Vascular Resistance	DS/cm ⁵ or kPa·s/l	/
	SVRI	Systemic Vascular Resistance Index	DS·m ² /cm ⁵ or kPa·s·m ² /l	1700-2400
	pArt-M	Mean Artery Pressure	mmHg/kPa or cmH ₂ O	70-90
	pArt-D	Diastolic Artery Pressure	mmHg/kPa or cmH ₂ O	60-80
	pArt-S	Systolic Artery Pressure	mmHg/kPa or cmH ₂ O	100-140

	Abbreviation	Full Spelling	Unit	Default Normal Range
Organ Function	EVLW	Extravascular Lung Water	ml	/
	ELWI	Extravascular Lung Water Index	ml/kg	3.0-7.0
	CPO	Cardiac Power Output	W	/
	CPI	Cardiac Power Index	W/ m ²	0.5-0.7
	PVPI	Pulmonary Vascular Permeability Index	no unit	1.0-3.0
	TB	Blood Temperature	°C	/
Oxygenation Parameters	ScvO ₂	Central Venous Oxygen Saturation	%	70-80
	Hb	Hemoglobin	g/dl	/
	DO ₂	Oxygen Delivery	ml/min	/
	DO ₂ I	Oxygen Delivery Index	ml/min/m ²	400-650
	VO ₂	Oxygen Consumption	ml/min	/
	VO ₂ I	Oxygen Consumption Index	ml/min/m ²	125-175
	SaO ₂	Arterial Oxygen Saturation	%	90-100

17.8.3 Normal Range Setup

You can select [**Normal Range Setup**] tab from the [**Hemodynamic Parameters**] menu to set up the normal ranges for 20 parameters. The system adopts the default normal ranges for the parameters if the ranges are not set up manually. Please refer to the above table for the hemodynamic parameters to see the default normal ranges of the hemodynamic parameters.

NOTE

- The normal ranges are based upon clinical experience and can vary from patient to patient. The stated values are therefore offered without guarantee. Indexed parameters are related to body surface area, predicted body weight or predicted body surface area and can also be displayed as absolute values.
- The values listed are not recommended for use on a specific patient. The treating physician is in any case responsible for determining and utilizing the appropriate diagnostic and therapeutic measures for each individual patient.

17.9 Changing PiCCO Settings

17.9.1 Selecting the Displayed Parameters

Select [**Select Parameter>>**] from the [**CCO Setup**] menu. In the pop-up menu, select the parameters to be displayed.

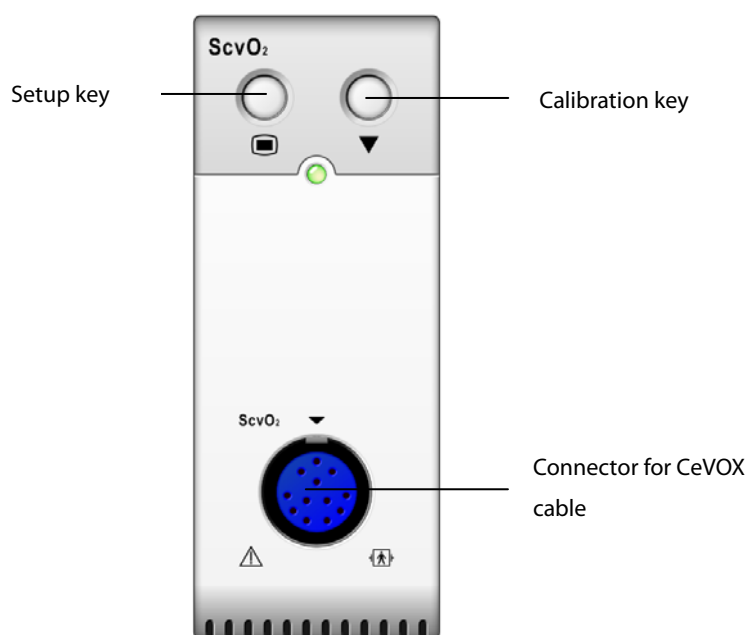
17.9.2 Selecting Alarm Properties

Select [**Alarm Setup >>**] from the [**CCO Setup**] menu to set the alarm properties for the relevant parameters.

18 Monitoring ScvO₂

18.1 Introduction

Central venous oxygen saturation (ScvO₂) is measured across spectrophotometry. Spectrophotometry involves the use of light emitting diodes (LED) that produce light of various wavelengths in red and infrared spectra. The light is transmitted to the blood through a fiberoptic in the probe, reflected off the red blood cells and transmitted back through a separate fiberoptic to an optical module. The central venous oxygen saturation is calculated through the analysis of the reflected spectra.



18.2 Safety Information

WARNING

- ScvO₂ monitoring is restricted to adult and pediatric patients.
-

18.3 Performing ScvO₂ Measurements

Please refer to the following procedure to perform the ScvO₂ measurements:

1. Apply the central venous catheter.
2. Place one end of the fiberoptic probe into the central venous catheter through the distal lumina, and connect the other end to the CeVOX optical module. Then plug the CeVOX cable into the ScvO₂ module.
3. If you see the message [**Calibration Required**], calibrate the ScvO₂ before performing the measurements. For detailed information on ScvO₂ calibration, please see **18.4 ScvO₂ Calibration**.
4. Check the reading in the ScvO₂ parameter window.



WARNING

- **To avoid installation failure, ensure that proper fiberoptic probe is selected.**
 - **Incorrect placement of the fiberoptic probe can lead to vessel perforation. Therefore check the correct position of the probe as indicated in the probe's instructions for use.**
-

18.4 ScvO₂ Calibration

Regular in vivo calibration is required using blood gas analysis of a central venous blood sample to ensure accurate measurement of continuous ScvO₂. For optimal accuracy, it is recommended that an in vivo calibration be performed at least every 24 hours or if hemoglobin is changing (for more details, check the notes below). Please refer to the following procedure to perform calibration:

1. Check central venous catheter and CeVOX probe for proper placement.
2. Check the quality of the signal. The Signal Quality Indicator (SQI) is used for assessing the quality of fiberoptical signals during probe placement, calibration and measurement. The signal quality is indicated by bars of different height levels. Generally, the higher the level, the better the signal.
3. Withdraw a sufficient amount of central venous blood from the side port of the CeVOX probe to avoid intermixture of infusion/injection with the withdrawn blood.
4. Slowly withdraw 2ml blood from the side port of the CeVOX probe. Do not pull too strongly in order to avoid a hemolysis.
5. Immediately confirm by pressing the [**Sample drawn**] button.
6. If necessary put blood sample on ice and perform an analysis by a blood gas analysis device or a laboratory oximeter.
7. Input lab values for Hb/Hct and ScvO₂ and press [**Calibrate**] to confirm.

NOTE

-
- **The SQI signal can be affected by the presence of electrosurgical units. Keep electrocautery equipment and cables away from the monitor and use separate power socket if possible.**
 - **To achieve optimal accuracy, it is recommended that the entered hemoglobin and hematocrit values are updated when there is a change of 6 % or more in hematocrit, or of 1.8 g/dl (1.1 mmol/l) or more in hemoglobin. A change in hemoglobin may also affect SQI.**
 - **Dye (e.g. Indocyanine Green) or other substances, containing dyes which usually modify the light absorption capacities, can lead to faulty measurement values of the oxygen saturation.**
-

18.5 Understanding the ScvO₂ Display



18.6 Understanding ScvO₂ Parameters

Apart from ScvO₂, the patient monitor can also monitor DO₂, VO₂, DO₂I, and VO₂I. You can access the [**ScvO₂ Calibration**] menu from the [**ScvO₂ Setup**] menu and input a SaO₂ value in [**SaO₂**] edit box. The patient monitor will calculate the values for oxygenation parameters automatically, and displays these parameters at [**Oxygenation Parameters**] in the [**ScvO₂ Setup**] menu. If a parameter value exceeds its normal range, the system will add a " ↑ " or " ↓ " to the right of the parameter.

WARNING

- **The patient monitor may only be regarded as a device providing early warning. If there is an indication of a trend towards de-oxygenation of the patient, blood samples must be taken and tested on a laboratory oximeter in order to arrive at a decision concerning the condition of the patient.**
-

18.7 Changing ScvO₂ Settings

18.7.1 Selecting Hb/Hct

1. Open the [**ScvO₂ Setup**] menu.
2. Select [**Hb/Hct**] and toggle between [**Hb**] and [**Hct**].

18.7.2 Selecting Alarm Properties

Select [**Alarm Setup >>**] from the [**ScvO₂ Setup**] menu to set the alarm properties for the relevant parameters.

19 Monitoring Carbon Dioxide

19.1 Introduction

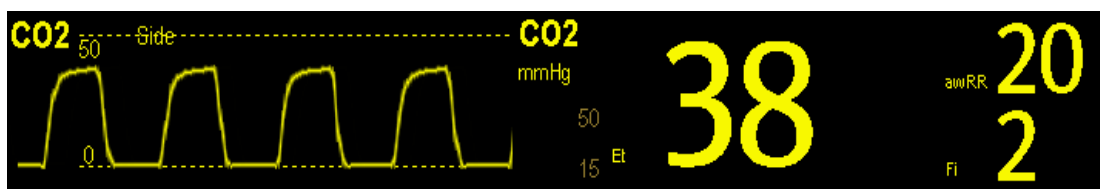
CO₂ monitoring is a continuous, non-invasive technique for determining the concentration of CO₂ in the patient's airway by measuring the absorption of infrared (IR) light of specific wavelengths. The CO₂ has its own absorption characteristic and the amount of light passing the gas probe depends on the concentration of the measured CO₂. When a specific band of IR light is passed through respiratory gas samples, some of IR light will be absorbed by the CO₂ molecules. The amount of IR light transmitted after it has been passed through the respiratory gas sample is measured with a photodetector. From the amount of IR light measured, the concentration of CO₂ is calculated.

There are two methods for measuring CO₂ in the patient's airway:

1. Mainstream measurement uses a CO₂ sensor attached to an airway adapter directly inserted into the patient's breathing system.
2. Sidestream/Microstream measurement samples expired patient gas at a constant sample flow from the patient's airway and analyzes it with a CO₂ sensor built into the CO₂ module.


The measurement provides:

1. A CO₂ waveform
2. End tidal CO₂ value (EtCO₂): the CO₂ value measured at the end of the expiration phase.
3. Fraction of inspired CO₂ (FiCO₂): the smallest CO₂ value measured during inspiration.
4. Airway respiration rate (awRR): the number of breaths per minute, calculated from the CO₂ waveform.



19.2 Identifying CO₂ Modules

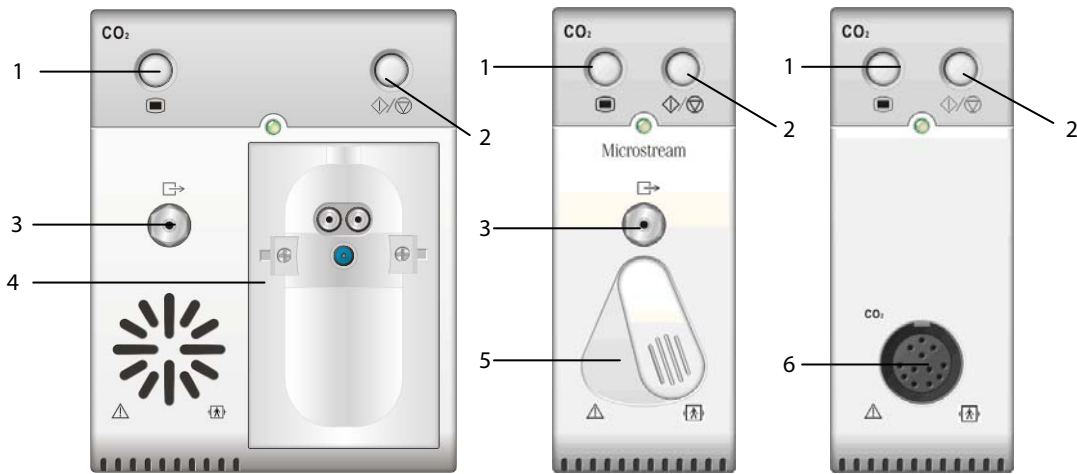From left to right are sidestream CO₂ module, microstream CO₂ module and mainstream CO₂.




1. Setup key to enter the CO₂ setup menu
2. Measure/standby
3. Gas outlet
4. Slot for CO₂ watertrap
5. Connector for sampling line
6. Connector for CO₂ transducer

If you measure CO₂ using the AG module, see the section **Monitoring AG**.

19.3 Preparing to Measure CO₂

19.3.1 Using a Sidestream CO₂ Module

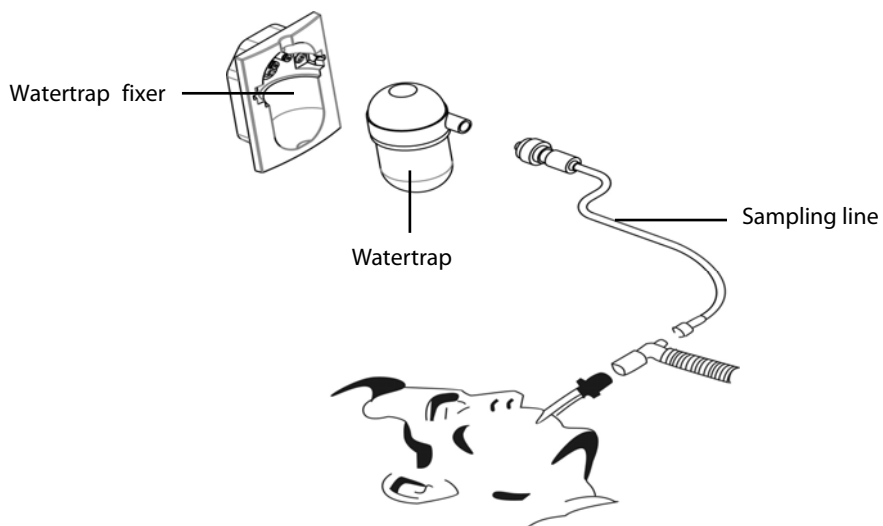1. Attach the watertrap to the module and then connect the CO₂ components as shown below.



2. By default, the sidestream CO₂ module is in measure mode. The **[CO₂ Startup]** message appears on the screen when the CO₂ module is plugged.
3. After start-up is finished, the CO₂ module needs time to warm up to reach the operating temperature. The message **[CO₂ Sensor Warmup]** is displayed. If you perform CO₂ measurements during warm-up, the measurement accuracy may be compromised.
4. After warm-up is finished, you can perform CO₂ measurements.

NOTE

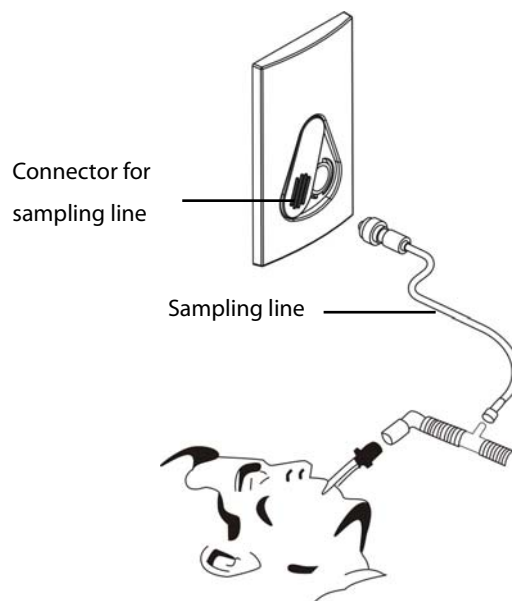
-
- **Do not apply adult watertrap to the neonate patient. Otherwise, patient injury could result.**
 - **To extend the lifetime of the watertrap and module, disconnect the watertrap and set the operating mode to standby mode when CO₂ monitoring is not required.**
-

CAUTION

-
- **The watertrap collects water drops condensed in the sampling line and therefore prevents them from entering the module. If the collected water reaches a certain amount, you should drain it to avoid blocking the airway.**
 - **The watertrap has a filter preventing bacterium, water and secretions from entering the module. After a long-term use, dust or other substances may compromise the performance of the filter or even block the airway. In this case, replace the watertrap. It is recommended to replace the watertrap once every two months, or when the watertrap is found leaky, damaged or contaminated.**
-

19.3.2 Using a Microstream CO₂ Module


1. Connect the sampling line to the module and then connect the CO₂ components as shown below.



2. By default, the microstream CO₂ module is in measure mode. The message **[CO₂ Sensor Warmup]** appears on the screen when the CO₂ module is plugged.
3. After warm-up, you can perform CO₂ measurements.

19.3.3 Using a Mainstream CO₂ Module

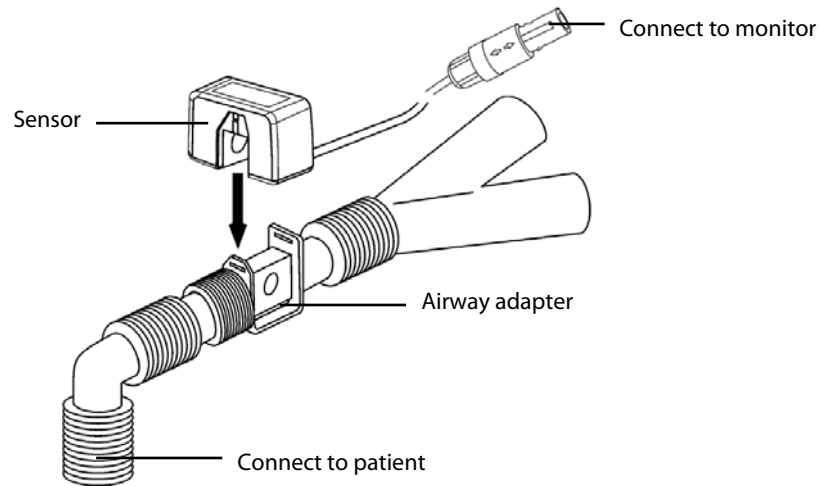1. Connect the sensor to the module.
2. By default, the mainstream CO₂ module is in measure mode. The message [**CO₂ Sensor Warmup**] appears on the screen when the CO₂ module is plugged.
3. After warm-up is finished, connect the transducer to the airway adapter.
4. Perform a zero calibration per the **Zeroing the Sensor** section.
5. After the zero calibration is finished, connect the airway as shown below.



6. Make sure there are no leakages in the airway and then start a measurement.

NOTE

-
- **Always position the sensor with the adapter in an upright position to avoid collection of fluids on the windows of the adapter. Large concentrations of fluids at this point will obstruct gas analysis.**
-

19.4 Changing CO₂ Settings

19.4.1 Accessing CO₂ Menus



By selecting the CO₂ parameter window or waveform, you can access the [**CO₂ Setup**] menu.

19.4.2 Entering the Standby Mode

The standby mode of the CO₂ module relates to the standby mode of the monitor as follows:

- If the monitor enters the standby mode, the CO₂ module will also enter the standby mode.
- If the monitor exits the standby mode, the CO₂ module will also exit the standby mode.
- If the CO₂ module enters or exits the standby mode, it will not affect the monitor.

To enter or exit the standby mode manually,

- select the   hardkey on the module, or
- select [**Operating Mode**] in the [**CO₂ Setup**] menu and then toggle between [**Standby**] and [**Measure**].

When you set the sidestream CO₂ module to the standby mode, the CO₂ gas sample intake pump automatically sets the sample flow rate to zero. When exiting the standby mode, the CO₂ module continues to work at the preset sample flow rate with no need to warm up again. After nearly 1 minute, the module enters the full accuracy mode.

For the sidestream CO₂ module, you can set the delay time. After the delay time the CO₂ module enters the standby mode if no breath is detected.

For the microstream CO₂ module, you can also set a period of time after which the CO₂ module enters the standby mode if no breath is detected since the CO₂ module is powered on or the CO₂ module switches to the measuring mode or the automatic standby time is re-set. To set the standby time, in the **[CO₂ Setup]** menu, select **[Auto Standby]** and then select the appropriate setting.

19.4.3 Setting the CO₂ Unit

Select **[Unit Setup >>]** from the **[User Maintenance]** menu. In the popup menu, select **[CO₂. Unit]** and toggle between **[mmHg]**, **[%]** and **[kPa]**.

19.4.4 Setting up Gas Compensations



- **Make sure that the appropriate compensations are used. Inappropriate compensations may cause inaccurate measurement values and result in misdiagnosis.**
-

For the sidestream CO₂ module:

1. Select **[CO₂ Setup]**.
2. According to the actual condition, set the concentration required for the following compensations:
 - ◆ **[O₂ Compen]**
 - ◆ **[N₂O Compen]**
 - ◆ **[Des Compen]**

For the microstream CO₂ module, gas compensations are not required.

For the mainstream CO₂ module, in the **[CO₂ Setup]** menu, respectively select:

- **[Balance Gas]** and toggle between **[Room Air]** and **[N₂O]**. Select **[Room Air]** when air predominates in the ventilation gas mixture and select **[N₂O]** when N₂O predominates in the ventilation gas mixture and select **[He]** when He predominates in the ventilation gas mixture.
- **[O₂ Compen]** and then select **[Off]** or an appropriate setting according to the amount of O₂ in the ventilation gas mixture. When the amount of O₂ is less than 30%, you'd better switch this compensation off.
- **[AG Compen]** and enter the concentration of anesthetic gas present in the ventilation gas mixture. This could compensate for the effect of AG on the readings.

19.4.5 Setting up Humidity Compensation

Sidestream and microstream CO₂ modules are configured to compensate CO₂ readings for either Body Temperature and Pressure, Saturated Gas (BTPS), to account for humidity in the patient's breath, or Ambient Temperature and Pressure, Dry Gas (ATPD).

1. ATPD: $P_{CO_2} (mmHg) = CO_2 (vol\%) \times P_{amb} / 100$
2. BTPS: $P_{CO_2} (mmHg) = CO_2 (vol\%) \times (P_{amb} - 47) / 100$

Where, P_{CO_2} = partial pressure, $vol\%$ = CO₂ concentration, P_{amb} = ambient pressure, and unit is mmHg.

As the mainstream CO₂ module has a built-in heating component to prevent water vapour from condensing, setting humidity compensation is not needed. For the sidestream and microstream CO₂ module, you can set the humidity compensation on or off according to the actual condition. To set the humidity compensation:

1. In the [**CO₂ Setup**] menu, select [**BTPS Compen**].
2. Select either [**On**] for BTPS or [**Off**] for ATPD, depending on which compensation applies.

19.4.6 Setting the Apnea Alarm Delay

In the [**CO₂ Setup**] menu, select [**Apnea Delay**] and then select the appropriate setting. The monitor will alarm if the patient has stopped breathing for longer than the preset apnea time. The [**Apnea Delay**] of Resp, CO₂, AG, and RM module keeps consistent with each other.

WARNING

- **The respiration measurement does not recognize the cause of apneas. It only indicates an alarm if no breath is detected when a pre-adjusted time has elapsed since the last detected breath. Therefore, it cannot be used for diagnostic purpose.**
-

19.4.7 Choosing a Time Interval for Peak-Picking

For microstream and mainstream CO₂ modules, you can select a time interval for picking the highest CO₂ as the EtCO₂ and the lowest as the FiCO₂.

In the [**CO₂ Setup**] menu, select [**Max Hold**] and toggle between [**Single Breath**], [**10 s**], [**20 s**] and [**30 s**].

- [**Single Breath**]: EtCO₂ and FiCO₂ are calculated for every breath.
- [**10 s**] or [**20 s**]: EtCO₂ and FiCO₂ are calculated using 10, 20 or 30 seconds of data.

19.4.8 Setting the Flow Rate

For the sidestream CO₂ module, you can change the sampling rate of respiratory gas in the patient's airway by setting the flow rate. To set the flow rate, enter the [**CO₂ Setup**] menu and select an appropriate setting from [**Flow Rate**].

 **WARNING**

- **Please consider the patient's actual bearing capability and select the appropriate flow rate when setting the flow rate.**
-

19.4.9 Setting up the CO₂ Wave

In the [CO₂ Setup] menu, you can:

- Select [Wave Type] and toggle between [Draw] and [Fill]:
 - ◆ [Draw]: The CO₂ wave is displayed as a curved line.
 - ◆ [Fill]: The CO₂ wave is displayed as a filled area.
- Select [Sweep] and then select the appropriate setting. The faster the wave sweeps, the wider the wave is.
- Change the size of the CO₂ waveform by adjusting the wave [Scale].

19.5 Setting RR Source

To set RR source:

1. Enter the [CO₂ Setup] menu.
2. Select [RR Source] and then select a source or [Auto] from the dropdown list.

The [RR Source] settings of Resp, CO₂, AG and RM module are linked. For details, please refer to the section **Setting RR Source** of chapter *Resp*.

19.6 Setting Barometric Pressure Compensation

Both sidestream and microstream CO₂ modules have the function of automatic barometric pressure compensation (the system automatically measures the barometric pressure which the patient monitor is exposed to). However, the mainstream CO₂ module does not have such function. For the mainstream CO₂ module, the default barometric pressure is 760 mmHg. You must modify the barometric pressure based on the actual situation as follows:

1. Select [Main Menu]→[Maintenance >>]→[User Maintenance >>]→enter the required password→[Maintain CO₂ >>]→[Calibrate CO₂ >>].
2. Select [Barometric Pressure] and then enter the value of barometric pressure to which the patient monitor is exposed to.

 **WARNING**

- **Be sure to set the barometric pressure properly before using the mainstream CO₂ module. Improper settings will result in erroneous CO₂ reading.**
-

19.7 Measurement Limitations

The following factors may influence the accuracy of measurement:

- Leaks or internal venting of sampled gas
- Mechanical shock
- Cyclic pressure up to 10 kPa (100 cmH₂O)
- Other sources of interference, if any

19.8 Leakage test

When the modules need maintenance, the monitor will prompt on the CO₂ parameter window: [**Need maintenance. Enter CO₂ setup menu.**] Then, you can access [**CO₂ Setup**]→ [**Maintain CO₂**], and perform leakage test according to the prompt messages on the menu.

19.9 Troubleshooting the Sidestream CO₂ Sampling System

When the sampling system of the sidestream CO₂ module works incorrectly, check if the sampling line is kinked. If not, remove it from the watertrap. If the monitor gives a message indicating the airway still works incorrectly, it indicates that the watertrap must have been blocked, and you should replace with a new one. Otherwise, you can determine that the sampling line must have been blocked. Replace with a new sampling line.

19.10 Removing Exhaust Gases from the System



- **Anesthetics: When using the Sidestream or Microstream CO₂ measurement on patients who are receiving or have recently received anesthetics, connect the outlet to a scavenging system, or to the anesthesia machine/ventilator, to avoid exposing medical staff to anesthetics.**
-

To remove the sample gas to a scavenging system, connect an exhaust tube to the gas outlet connector of the module.

19.11 Zeroing the Sensor

The zero calibration eliminates the effect of baseline drift during CO₂ measurement exerted on the readings and therefore maintains the accuracy of the CO₂ measurements.

19.11.1 For Sidestream and Microstream CO₂ Modules

For sidestream and microstream CO₂ modules, a zero calibration is carried out automatically when necessary. You can also start a manual zero calibration if necessary. To manually start a zero calibration, select [**Maintain CO₂ >>**] from the [**User Maintenance**] menu. Then select [**Calibrate CO₂ >>**] → [**Start Zero Cal.**]. Disconnecting the patient airway is not required when performing a zero calibration.

19.11.2 For Mainstream CO₂ Modules

For mainstream CO₂ modules, zero the sensor whenever:

- A new adapter is used;
- You reconnect the sensor to the module;
- You see the message [**CO₂ Zero Required**]. In this case, check the airway adapter for any blockage, e.g. mucus, etc. If a blockage is detected, clear or replace the adapter.

To zero the sensor, follow this procedure:

1. Connect the sensor to the module.
2. In the [**CO₂ Setup**] menu, set the [**Operating Mode**] to [**Measure**]. The message [**CO₂ Sensor Warmup**] is displayed.
3. After warm-up is finished, connect the sensor to a clean, dry airway adapter. The adapter should be vented to the air and isolated from CO₂ sources, such as ventilator, the patient's breathing, your own breathing, etc.
4. Select [**Start Zero Cal.**] in the [**CO₂ Setup**] menu. The message [**CO₂ Zero Running**] is displayed.
5. It takes about 15 to 20 seconds. The message disappears when the zero calibration is completed.

WARNING

- **When perform a zero calibration during the measurement, disconnect the transducer from the patient's airway first.**
 - **Please do not rely on the readings during zeroing.**
-

19.12 Calibrating the Sensor

For sidestream or microstream CO₂ modules, a calibration should be performed once every year or when the readings go far beyond the range. For mainstream CO₂ modules, no calibration is required. For details, refer to the chapter **0** Maintenance.

19.13 Oridion Information

Microstream

This trademark is registered in Israel, Japan, German and America.

Oridion Patents

The capnography component of this product is covered by one or more of the following US patents: 6,428,483; 6,997,880; 6,437,316; 7,488,229; 7,726,954 and their foreign equivalents. Additional patent applications pending.

No Implied License

Possession or purchase of this device does not convey any express or implied license to use the device with unauthorized CO₂ sampling consumables which would, alone, or in combination with this device, fall within the scope of one or more of the patents relating to this device and/or CO₂ sampling consumable.

20 Monitoring tcGas

20.1 Introduction

This patient monitor can connect a Radiometer TCM monitor for continuous transcutaneous blood gas monitoring.

This patient monitor can display, store and review measurements from TCM monitor, as well as present related alarms. On this patient monitor, you can separately set the level of tcGas related alarms and switch on or off alarm recording; you can also view TCM monitor settings of alarm limits and alarm switch.

This patient monitor can integrate the following TCM monitors:

- TCM CombiM
- TCM TOSCA

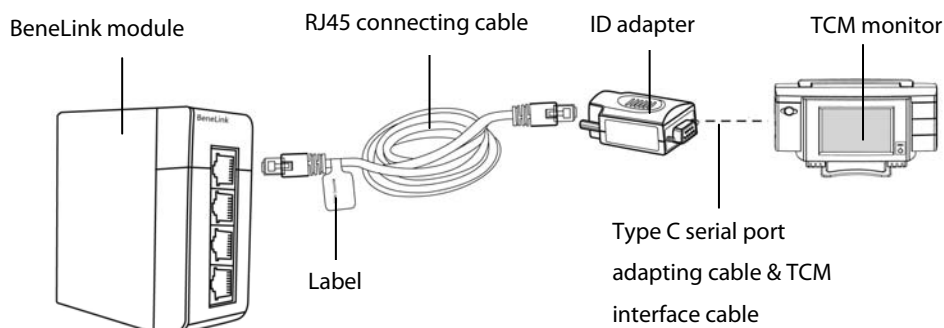
20.2 Safety

WARNING

- **TCM CombiM monitor and TCM TOSCA monitor are manufactured by Radiometer Medical ApS. This company provides the technology for measuring tcGas parameters. We only provide the connection between this patient monitor and TCM monitor.**
 - **If you have any doubts about the operation and maintenance of the TCM monitor, please refer to TCM monitor operator's manual or directly contact Radiometer Medical (www.radiometer.com).**
 - **Fully observe TCM monitor operator's manual to make settings and to connect the monitor with a patient.**
-

20.3 Connecting a TCM monitor

The TCM monitor connects with BeneLink module through an ID adapter, see the picture below.



Please refer to the following procedure to connect the TCM monitor:

1. Insert a BeneLink module into a BeneView patient monitor module rack.
2. Connect the ID adapter that matches the TCM monitor to the BeneLink module with an RJ45 connecting cable.
3. Connect the ID adapter to the serial port (COM port) of the TCM monitor with Mindray type C serial port adapting cable (PN: 009-001769-00) and an interface cable provided with the TCM monitor.
4. Stick a label indicating device name to the RJ45 connecting cable at the end nearby the BeneLink module. When the BeneLink module is connected to several external devices, you can tell the devices easily with these labels.
5. Turn on both monitors.

20.4 tcGas Parameters

TCM CombiM monitor provides the following measurements:

- tcpCO₂
- tcpO₂
- Power
- Tsensor


In which, tcpCO₂ and tcpO₂ are primary parameters, and the others are secondary parameters.

TCM TOSCA provides the following measurements:

- tcpCO₂
- SpO₂
- PR
- Power
- Tsensor

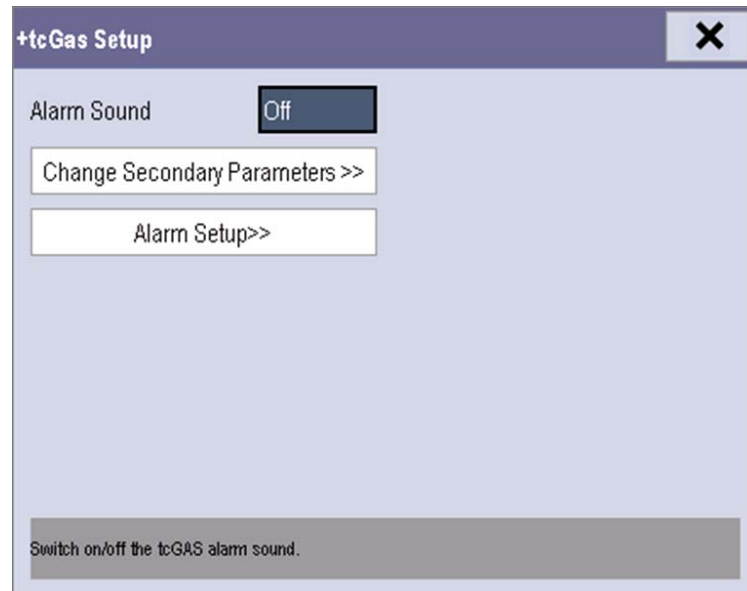
In which, tcpCO₂ is primary parameter, and the others are secondary parameters.

20.5 Displaying tcGas Parameters

To display tcGas parameters on this patient monitor, select the **[Screen Setup]** button  to enter the **[Screens]** window, and then select **[Screen Setup]**. You can choose where to display the tcGas parameters on the screen.

20.6 Enter the tcGas Setup menu

You can access the **[+tcGas Setup]** menu by selecting the tcGas area or selecting **[Main Menu]** → **[Parameters >>]** → **[+tcGas Setup >>]**.



In the **[+tcGas Setup]** menu, you can

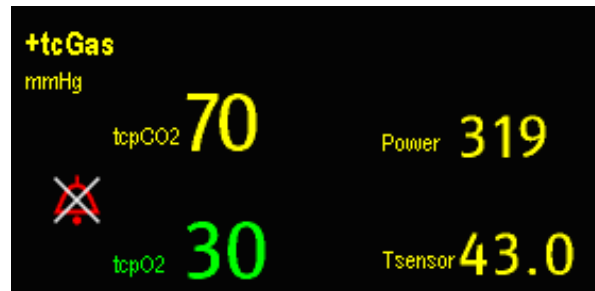
- Toggle **[Alarm Sound]** between **[On]** and **[Off]** to switch on or off tcGas alarms on this patient monitor.
- Choose the secondary parameters to be displayed. The tcGas area can display maximum three secondary parameters.
For TCM CombiM monitor, only two secondary parameters, Power and Tsensor, are measured, so in **[+tcGas Setup]** menu the option **[Change Secondary Parameters >>]** is not available.
- Set alarm level for tcGas parameters, switch on or off alarm record.

20.7 Setting tcpCO₂/tcpO₂ Unit


You can enter the [User Maintenance] menu to [Unit Setup >>] to set [tcpCO₂/tcpO₂ Unit] to [mmHg] or [kPa].

20.8 tcGas Display

If TCM CombiM monitor is connected, the tcGas area is shown as follows:



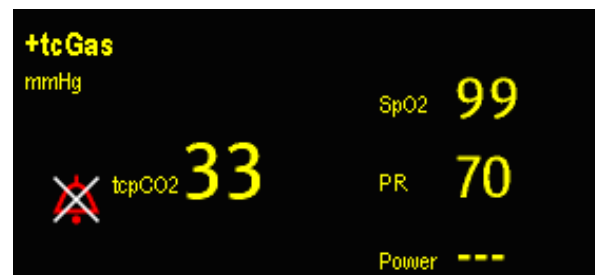If TCM TOSCA monitor is connected, the tcGas area is shown as follows:



21 Monitoring AG

21.1 Introduction

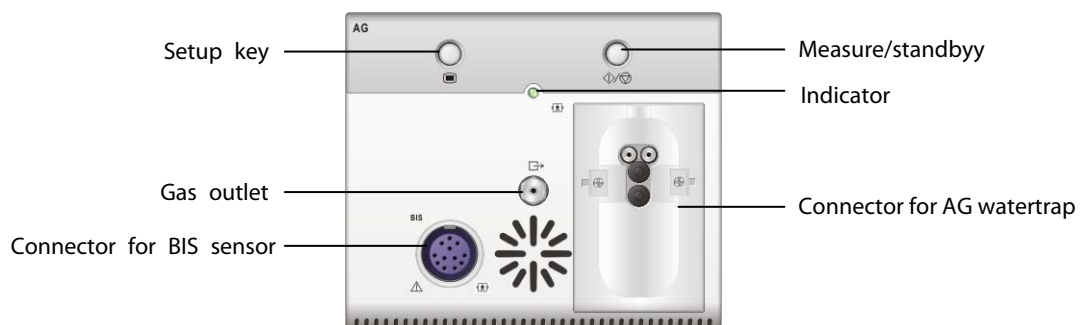
The anaesthetic gas (AG) module measures the patient's anesthetic and respiratory gases, and incorporates the features of the O₂ module and BIS module as well.

The AG module determines the concentration of certain gases using the infrared (IR) light absorption measurement. The gases that can be measured by the AG module absorb IR light. Each gas has its own absorption characteristic. The gas is transported into a sample cell, and an optical IR filter selects a specific band of IR light to pass through the gas. For multiple gas measurement, there are multiple IR filters. The higher the concentration of gas in a given volume the more IR light is absorbed. This means that higher concentration of IR absorbing gas cause a lower transmission of IR light. The amount of IR light transmitted after it has been passed though an IR absorbing gas is measured. From the amount of IR light measured, the concentration of gas present can be calculated.

Oxygen does not absorb IR light as other breathing gases and is therefore measured relying on its paramagnetic properties. Inside the O₂ sensor are two nitrogen-filled glass spheres mounted on a strong rare metal taut-band suspension. This assembly is suspended in a symmetrical non-uniform magnetic field. In the presence of paramagnetic oxygen, the glass spheres are pushed further away from the strongest part of the magnetic field. The strength of the torque acting on the suspension is proportional to the oxygen concentration. From the strength of the torque, the concentration of oxygen is calculated.

21.2 Identifying AG Modules

AG module can identify two anesthetic gases in a mixture automatically and distinguish between them according to their contributions to the MAC value for display as the primary and secondary anesthetic agent.




For details on BIS, refer to the chapter **23 Monitoring BIS**.

NOTE

- **The AG module is configured with automatic barometric pressure compensation function.**

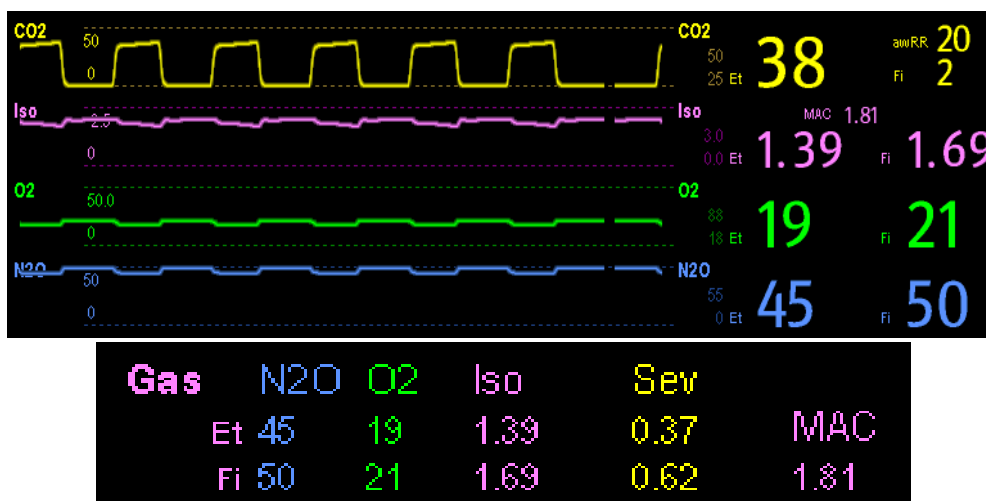21.3 Understanding the AG Display



The AG module can send waves and numerics for all measured anesthetic gases for display on the monitor, including:

- CO₂, O₂, N₂O and AA waves
- awRR: airway respiratory rate
- MAC: minimum alveolar concentration
- End tidal (Et) and fraction of inspired (Fi) numerics for CO₂, O₂, N₂O and AA

Where AA represents Des (desflurane), Iso (isoflurane), Enf (enflurane), Sev (sevoflurane), or Hal (halothane).

The AA waveform area displays the primary anesthetic gas's waveform. When O₂ module does not exist, no O₂ waveform will be displayed. When O₂ module exists, the O₂ waveform will be displayed only when the O₂ waveform is currently switched on.

WARNING

- **To avoid explosion hazard, do not use flammable anesthetic agent such as ether and cyclopropane for this equipment.**
-

21.4 MAC Values

Minimum alveolar concentration (MAC) is the minimum concentration of the agent in the alveoli. It is a basic index to indicate the depth of anesthesia. The standard ISO 21647 defines MAC as this: alveolar concentration of an inhaled anesthetic agent that, in the absence of other anesthetic agents and at equilibrium, prevents 50% of patients from moving in response to a standard surgical stimulus.

Minimum alveolar concentration (MAC) values are listed below:

Agent	Des	Iso	Enf	Sev	Hal	N2O
1 MAC	7.3%*	1.15%	1.7%	2.1%	0.77%	105%**

* The data is taken from a patient of 25 years old.

** indicates 1 MAC nitrous oxide can only be reached in hyperbaric chamber.

NOTE

- The MAC values shown in the table above are those published by the U.S. Food and Drug Administration for a healthy 40-year-old adult male patient.
- In actual applications, the MAC value may be affected by age, weight and other factors.

The formula to calculate the MAC value is as follows:

$$MAC = \sum_{i=0}^{N-1} \frac{EtAgent_i}{AgentVol_i}$$

Where N is the number of all agents (including N₂O) that the AG module can measure, EtAgent_i is the concentration of each agent and AgentVol_i is the concentration of each agent at 1 MAC.

For example, the AG module measures there are 4% of Des, 0.5% of Hal and 50% of N₂O in the patient's end-tidal gas:

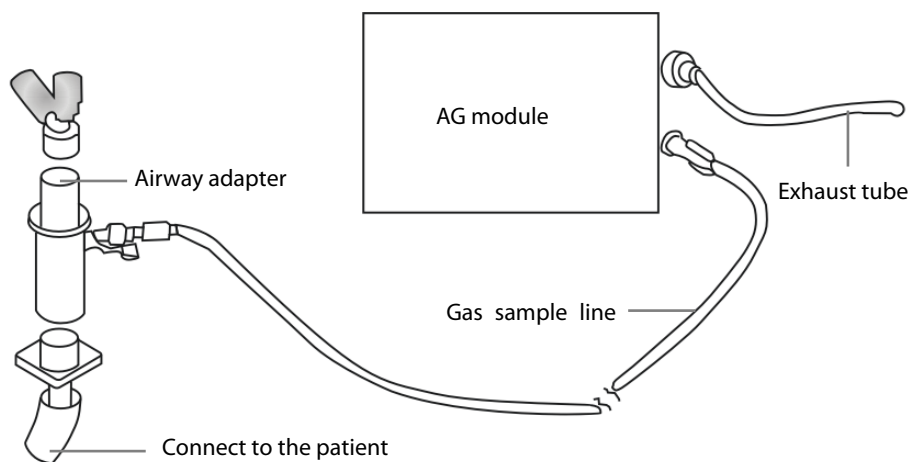
$$MAC = \frac{4.0\%}{7.3\%} + \frac{0.5\%}{0.77\%} + \frac{50\%}{105\%} = 1.67$$

NOTE

- The formula mentioned above is intended for adult patients only.

21.5 Preparing to Measure AG

1. Select an appropriate watertrap according to patient category and attach it to the module.
2. Connect the gas sample line to the connector of the watertrap.
3. Connect the other end of the gas sampling line to the patient via the airway adapter.
4. Connect the gas outlet to a scavenging system using an exhaust tube.



5. Insert the AG module into the SMR or the patient monitor and the patient monitor will prompt **[AG Startup]**. Within 10 minutes after startup is finished, the AG module enters the iso accuracy mode. After that, the module enters the full accuracy mode.

 **CAUTION**

- Position the airway adapter so that the part connecting to the gas sample line is pointing upwards. This prevents condensed water from passing into the gas sample line and causing an occlusion.
 - The watertrap collects water drops condensed in the sampling line and therefore prevents them from entering the module. If the collected water reaches to a certain amount, you should drain it to avoid blocking the airway.
 - The watertrap has a filter preventing bacterium, water and secretions from entering the module. After a long-term use, dust or other substances may compromise the performance of the filter or even block the airway. In this case, replace the watertrap. Replacing the watertrap once a month is recommended.
-

 **WARNING**

- Make sure that the connections are tight. Any leak in the system can result in erroneous readings due to ambient air mixing with patient gases.
 - Do not apply adult watertrap to the neonate patient. Otherwise, patient injury could result.
 - Using high-frequency electrosurgical units may increase the risk of skin burn. In this case, do not use antistatic or conductive respiratory tubing.
-

21.6 Changing AG Settings

21.6.1 Setting Gas Unit

For N₂O and AA, the unit of the measured gas is fixed to “%”.

Select [**Unit Setup >>**] from the [**User Maintenance**] menu. In the popup menu, you can select [**CO₂ Unit**] or [**O₂ Unit**] and toggle between [**mmHg**], [%] and [**kPa**].

21.6.2 Setting the Apnea Alarm Delay

In the [**AG Setup**] menu, select [**Apnea Delay**] and select the appropriate setting. The monitor will alarm if the patient has stopped breathing for longer than the preset apnea time.

The [**Apnea Delay**] of Resp, CO₂, AG, and RM module keeps consistent with each other.

 **WARNING**

- The respiration measurement does not recognize the cause of apneas. It only indicates an alarm if no breath is detected when a pre-adjusted time has elapsed since the last detected breath. Therefore, it cannot be used for diagnostic purpose.
-

21.6.3 Changing the Sample Flow Rate

In the setup menu for any gas, select [**Flow Rate**] and then choose either:

- [**High**]: 200 ml/min for adult and pediatric patients, and 120 ml/min for neonatal patients.
 - [**Med**]: 150 ml/min for adult and pediatric patients, and 90 ml/min for neonatal patients.
 - [**Low**]: 120 ml/min for adult and pediatric patients, and 70 ml/min for neonatal patients.
-

21.6.4 Setting up the O₂ Compensation

If the AG module does not incorporate the O₂ module, you need to manually select [**O₂ Compen**] and then select [**Off**] or an appropriate setting according to the amount of O₂ in the ventilation gas mixture. When the amount of O₂ is less than 30%, you'd better switch this compensation off.

If the AG module incorporates the O₂ module, the system will directly use the O₂ concentration detected by the O₂ module to make compensation. At this time, the [**O₂ Compen**] in the setup menu for any gas is fixed to [**Off**].

21.6.5 Entering the Standby Mode

For the AG module, the default operating mode is measure. When you set the AG module to the standby mode, the AG gas sample intake pump automatically sets the sample flow rate to zero. When exiting the standby mode, the AG module continues to work at preset sample flow rate with no need to warm up again. After nearly 1 minute, the module enters the full accuracy mode. The standby mode of the AG module relates to the standby mode of the monitor as follows:

- If the monitor enters the standby mode, the AG module will also enter the standby mode.
- If the monitor exits the standby mode, the AG module will also exit the standby mode.
- If the AG module enters or exits the standby mode, it will not affect the monitor.

To enter or exit the standby mode manually, in the agent's setup menu, select [**Operating Mode**] and then toggle between [**Standby**] and [**Measure**]. You can also set a period of time after which the AG module enters the standby mode automatically if no breath is detected since the last detected breath. To set the standby time, in the agent's setup menu, select [**Auto Standby (min)**] and then select the appropriate setting.

21.6.6 Setting up the AG Wave

In the [**AG Setup**] menu, you can:

- Select [**CO₂ Wave Type**] and toggle between [**Draw**] and [**Fill**]:
 - ◆ [**Draw**]: The CO₂ wave is displayed as a curved line.
 - ◆ [**Fill**]: The CO₂ wave is displayed as a filled area.
- Select [**Sweep**] and then select the appropriate setting. The faster the wave sweeps, the wider the wave is.
- Change the size of the waveform by adjusting the scale.

21.6.7 Setting RR Source

To set RR source:

1. Enter the [**AG Setup**] menu.
2. Select [**RR Source**] and then select a source or [**Auto**] from the dropdown list.

The [**RR Source**] settings of Resp, CO₂, AG and RM module are linked. For details, please refer to the section **Setting RR Source** of chapter *Resp*.

21.7 Changing the Anesthetic Agent

When the anesthetic agent used on the patient is changed, the AG module can detect the mixed anesthetic gas during the transition of two anesthetic agents. The time required for completing the replacement of anesthetic agent depends on anesthesia type (low flow or high flow) and the characteristics of anesthetic agents (pharmacokinetics). During the transition of two anesthetic agents, the patient monitor gives no prompt messages and the MAC value displayed may be inaccurate.

The AG module can identify two anesthetic agents automatically. When the proportion of the primary and secondary anesthetic agents in the mixture changes, the AG module can distinguish between them according to their contributions to the MAC value. Then the primary and secondary anesthetic agents will be exchanged for display.

21.8 Measurement Limitations

The following factors may influence the accuracy of measurement:

- Leaks or internal venting of sampled gas
- Mechanical shock
- Cyclic pressure up to 10 kPa (100 cmH₂O)
- Other sources of interference, if any

21.9 Troubleshooting

21.9.1 When the Gas Inlet is Blocked

If the gas inlet (including watertrap, sampling line and airway adapter) is occluded by condensed water, the message [**AG Airway Occluded**] will appear.

To remove the occlusion:

- Check the airway adapter for an occlusion and replace if necessary.
- Check the sampling line for an occlusion or kinking and replace if necessary.
- Check the watertrap for a build up of water. Empty the watertrap. If the problem persists, replace the watertrap.

21.9.2 When an Internal Occlusion Occurs

Condensed water may enter the module and cause contamination and/or internal occlusions. In this case, the message [**AG Airway Occluded**] will be displayed.

To remove the occlusion:

- Check for any occlusion in the gas inlet and/or outlet system.
- If the problem persists, internal occlusions may exist. Contact your service personnel.

21.10 Removing Exhaust Gases from the System

WARNING

- **Anesthetics: When using the AG measurement on patients who are receiving or have recently received anesthetics, connect the outlet to a scavenging system, or to the anesthesia machine/ventilator, to avoid exposing medical staff to anesthetics.**
-
-

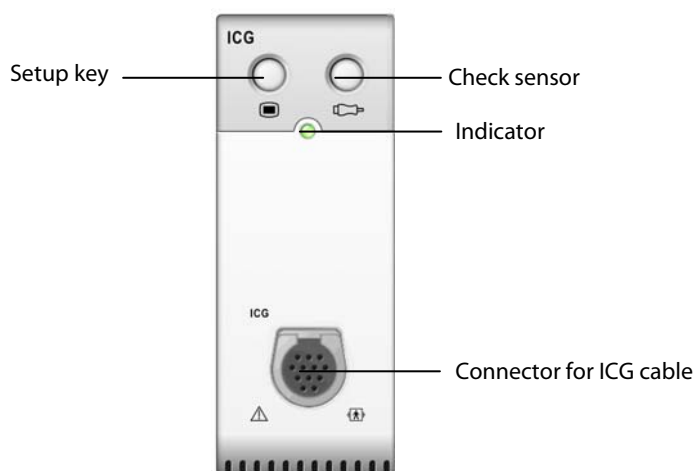
To remove the sample gas to a scavenging system, connect an exhaust tube to the gas outlet connector of the module.

FOR YOUR NOTES

22 Monitoring ICG

22.1 Introduction

Impedance cardiography (ICG) measures a patient's hemodynamic status using a safe, non-invasive method based on thoracic electrical bioimpedance (TEB) technology. ICG uses four pairs of sensors to transmit a small electrical signal through the thorax. As velocity and volume of blood in the aorta change, the ICG measures the changes in impedance from systole to diastole to calculate hemodynamic parameters.



22.2 Safety

WARNING

- **Apply ICG monitoring to adult patients in height of 122 to 229 cm, weight of 30 to 159 kg (67 to 341 pounds) only.**
 - **ICG monitoring should not be used concurrently on patients with minute ventilation pacemakers when the MV sensor function is activated.**
 - **During ICG monitoring, make sure that the conductive paste on the ICG sensors never come into contact with other conductive parts.**
 - **ICG sensors are for single patient use only.**
-

22.3 Understanding ICG Parameters

By selecting the ICG parameter window → [ICG Setup] → [Hemodynamic Parameters >>], you can view the hemodynamic parameters for evaluation of the patient's hemodynamic status.

22.3.1 Measured Parameters

Abbreviation	Unit	Full spelling
ACI	/100s ²	acceleration index
VI	/1000s	velocity index
PEP	ms	Pre-ejection period
LVET	ms	Left ventricular ejection time
TFI	Ω	Thoracic fluid index
TFC	/kΩ	thoracic fluid content
HR*	bpm	heart rate


*The HR value is directly derived from the ICG module.

22.3.2 Calculated Parameters

Abbreviation	Unit	Full spelling
BSA	m ²	Body surface area
C.O.	L/min	Cardiac output
C.I.	L/min/m ²	Cardiac index
SV	ml	Stroke volume
SVI	ml/m ²	Stroke volume index
SVR	DS/cm ⁵	Systemic vascular resistance
SVRI	DS·m ² /cm ⁵	Systemic vascular resistance index
PVR	DS/cm ⁵	Pulmonary vascular resistance
PVRI	DS·m ² /cm ⁵	Pulmonary vascular resistance index
LCW	kg·m	Left cardiac work
LCWI	kg·m/m ²	Left cardiac work index
LVSW	g·m	Left ventricular stroke work
LVSWI	g·m/m ²	Left ventricular stroke work index
STR	none	Systolic time ratio
VEPT	ml	Volume of electrically participating tissue

22.4 Understanding the ICG Display

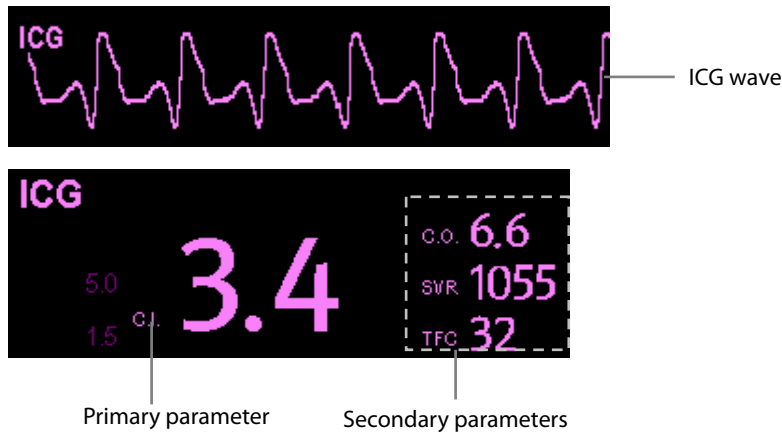The ICG monitoring provides a continuous display of the impedance waveform and four numerics. Of four numerics, one is the primary parameter C.I. and the other three are secondary parameters. The secondary parameters are user-selectable, and C.O., SVR and TFC are the defaults.



By selecting the ICG waveform area or ICG parameter window, you can access the **[ICG Setup]** menu.

22.5 ICG Limitations

The measurement accuracy may be compromised when patients present with the following conditions or anomalies:

- Septic shock.
- Aortic valve regurgitation.
- Severe hypertension (Art mean > 130 mmHg).
- The patient's weight and height are out of range.
- Connection to an intra-aortic balloon pump.
- With excessive and continuous patient movements such as shivering.
- Signal interference from cable connections and/or power cords.
- Open-chest surgeries that could result in changes in the normal pattern of blood flow and/or the electrical current through the chest cavity.

22.6 Preparing to Monitor ICG

1. Insert the ICG module into the monitor.
2. Connect the patient cable to the ICG module.
3. Prepare the patient's skin and place ICG sensors on the patient.
4. Connect the ICG sensor connector end to the patient cable lead wires.
5. Enter the patient information.

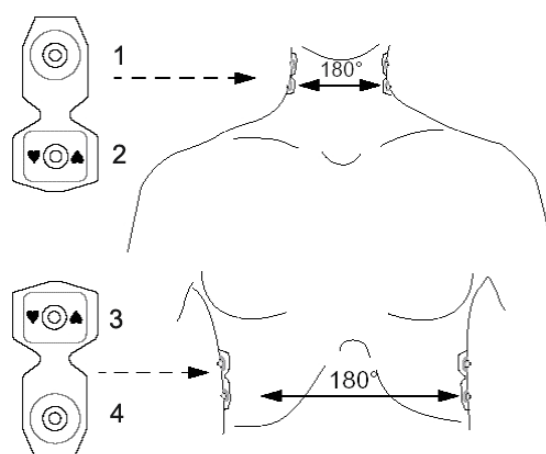
22.6.1 Preparing the Patient

Proper skin preparation is necessary for good signal quality at the sensor, as the skin is a poor conductor of electricity. To properly prepare the skin, choose flat, non-muscular areas and then follow this procedure:

1. Shave hair from skin at chosen sites.
2. Gently rub skin surface with a gauze pad at sites to remove dead skin cells.
3. Thoroughly cleanse the site with a mild soap and water solution. Be sure to remove all oily residue, dead skin cells, and abrasives. Leftover abrasion particles can be a source of noise.
4. Dry the skin completely before applying the sensor.

22.6.2 Placing ICG Sensors

Appropriate sensor placement is important for good signal quality and accurate measurements. Attach ICG sensors to the patient as shown below:



1. Two neck sensors placed on each side of the neck, with the rectangular shaped end of the sensor (2) placed at the base (or root) of the neck and the circular shaped end of the sensor (1) placed directly superior and in line with the earlobe.
2. Two thorax sensors placed on each side of the thorax, with the rectangular shaped end of the sensor (3) at the level with the xyphoid process and the circular shaped end of the sensor (4) directly inferior and in line with the midaxillary line. Each pair of sensors should be opposite directly to each other (180°) as shown in the figure above.

22.6.3 Setting up the Patient Information

1. Enter the [**ICG Setup**] menu.
2. Select [**Patient Demographics >>**].
3. Select [**Height**] and [**Weight**] and then select the appropriate settings. The patient's height and weight are important for ICG monitoring. The system will automatically check them when an ICG module is connected. If no values are entered or the entered values do not meet the requirements, corresponding prompt messages will be given in the ICG parameter window.
4. If the mean arterial blood pressure (Art mean) is not obtained automatically from either IBP or NIBP module, then enter Art mean. Enter CVP and PAWP, obtained from invasive catheters or enter an assumed value. (Note: CVP and PAWP are used only in the calculation of SVR, SVRI, LCW, LCWI, LVSW, and LVSWI and the value of CVP and PAWP do not normally have a significant effect on the calculated parameters.)

22.7 Changing ICG Settings

22.7.1 ICG Averaging

The ICG value is the average of multiple measurements. You can select an interval (heart beats) for averaging ICG, ranging from 5-60 beats.

1. Enter the **[ICG Setup]** menu.
2. Select **[Averaging]** and then select the appropriate setting. The greater the averaging interval is, the less the ICG value is affected by human interference and vice versa.

22.7.2 Selecting Secondary Parameters


C.O., SVR and TFC are the default three secondary parameters. You can also select your desired secondary parameter for display.

1. Enter the **[ICG Setup]** menu.
2. Select **[Change Secondary Parameters >>]**.
3. Select three parameters from the popup menu.

22.7.3 Checking Sensors

During ICG monitoring, the ICG sensors should be checked regularly to ensure that no sensor becomes disconnected. During sensor checking, the ICG waveform is displayed as a straight line and the message **[ICG Sensor Check]** is displayed. Once a disconnected sensor is detected, a prompt message in which the sensor's application site is indicated will be displayed.

To initiate a sensor check:

- In the **[ICG Setup]** menu, select **[Check Sensor]**.
- Press the  hardkey.

22.7.4 Changing the ICG Wave Speed

1. Enter the **[ICG Setup]** menu.
2. Select **[Sweep]** and then select the appropriate setting. The faster the wave sweeps, the wider the wave is.

FOR YOUR NOTES

23 Monitoring BIS

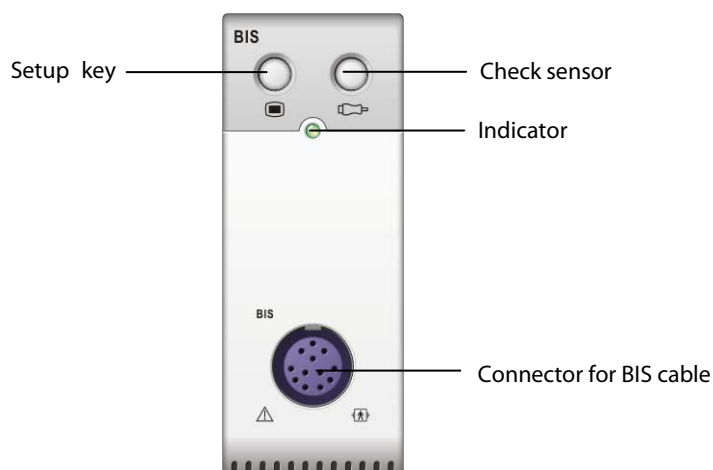
23.1 Introduction

Bispectral index (BIS) monitoring is for use on adult and pediatric patients within a hospital or medial facility providing patient care to monitor the state of the brain by data acquisition of EEG signals.

The BIS, a processed EEG variable, may be used as an aid in monitoring the effects of certain anesthetic agents. Use of BIS monitoring to help guide anesthetic administration may be associated with the reduction of the incidence of awareness with recall during general anesthesia or sedation.

BISx is for brain's single side BIS monitoring. BISx4 is for brain's single side or both sides BIS monitoring. BISx4 can be used for brain's both sides BIS monitoring only when BIS Bilateral Sensor is connected.

The BISx or BISx4 equipment must be used under the direct supervision of a licensed healthcare practitioner or by personnel trained in its proper use.



23.2 Safety Information

For patients with neurological disorders, patients taking psychoactive medication, and children below the age of 1 year, BIS values should be interpreted cautiously.

WARNING

- The conductive parts of sensors and connectors should not come into contact with other conductive parts, including earth.
 - To reduce the hazard of burns in the high-frequency surgical neutral electrode connection, the BIS sensor should not be located between the surgical site and the electro-surgical unit return electrode.
 - The BIS sensor must not be located between defibrillator pads when a defibrillator is used on a patient connected to the patient monitor.
 - The BIS component using on our monitor is purchased from Aspect Medical System. It is important to recognize this index is derived using solely that company's proprietary technology. Therefore, it is recommended that clinicians have reviewed applicable information on its utility and/or risks in published articles and literature/web site information from Aspect Medical Systems, Inc. or contact that company itself at www.aspectmedical.com, if you have clinical-based BIS questions relating to this module portion of the patient monitor. Failure to do so could potentially result in the incorrect administration of anesthetic agents and/or other potential complications of anesthesia or sedation. We recommend that clinicians also review the following practice advisory (that includes a section on BIS monitoring): The American Society of Anesthesiologists, Practice Advisory for Intraoperative Awareness and Brain Function Monitoring (Anesthesiology 2006;104:847-64). Clinicians are also recommended to maintain current knowledge of FDA or other federal-based regulatory, practice or research information on BIS and related topics.
 - The Bispectral Index is a complex technology, intended for use only as an adjunct to clinical judgment and training.
 - The clinical utility, risk/benefit and application of the BIS component have not undergone full evaluation in the pediatric population.
-

23.3 Understanding the BIS Display

23.3.1 BIS Parameter Area

For brain's single side BIS monitoring, the BIS parameter area displays the following parameters:



1. Bispectral Index (BIS)

The BIS numeric reflects the patient's level of consciousness. It ranges from 100 for wide awake to 0 in the absence of brain activity.

BIS numeric	Description
100	The patient is widely awake.
70	The patient is underdosed but still unlikely to become aware.
60	The patient is under general anesthesia and loses consciousness.
40	The patient is overdosed and in deep hypnosis.
0	The EEG waveform is displayed as a flat line, and the patient has no electrical brain activity.

2. Electromyograph (EMG)

EMG bar graph reflects the electrical power of muscle activity and high frequency artifacts. The power range is 30-55 dB. When the EMG indicator is low, it indicates that EMG activity is low. BIS monitoring conditions are optimal when the bar is empty.

1 bar represents power in the 31-35 range.

2 bars represent power in the 36-40 range.

3 bars represent power in the 41-45 range.

4 bars represent power in the 46-50 range.

5 bars represent power greater than 51.

- EMG > 55 dB: this is an unacceptable EMG.
- EMG < 55 dB: this is an acceptable EMG.
- EMG ≤ 30 dB: this is an optimal EMG.

3. Suppression Ratio (SR)

SR numeric is the percentage of time over the last 63-second period during which the EEG is considered to be in a suppressed state.

4. Spectral Edge Frequency (SEF)

The SEF is a frequency below which 95% of the total power is measured.

5. Signal Quality Index (SQI)

The SQI numeric reflects signal quality and provides information about the reliability of the BIS, SEF, TP, and SR numerics during the last minute. Signal quality is optimal when all five bars of the SQI icon are filled with color. SQI ranges from 0-100%.

1 bar represents SQI in the 1%-20% range.

2 bars represent SQI in the 21%-40% range.

3 bars represent SQI in the 41%-60% range.

4 bars represent SQI in the 61%-80% range.

5 bars represent SQI in the 81%-100% range.

- ◆ 0 to 15%: the numerics cannot be derived.
- ◆ 15% to 50%: the numerics cannot be reliably derived.
- ◆ 50% to 100%: the numerics are reliable.

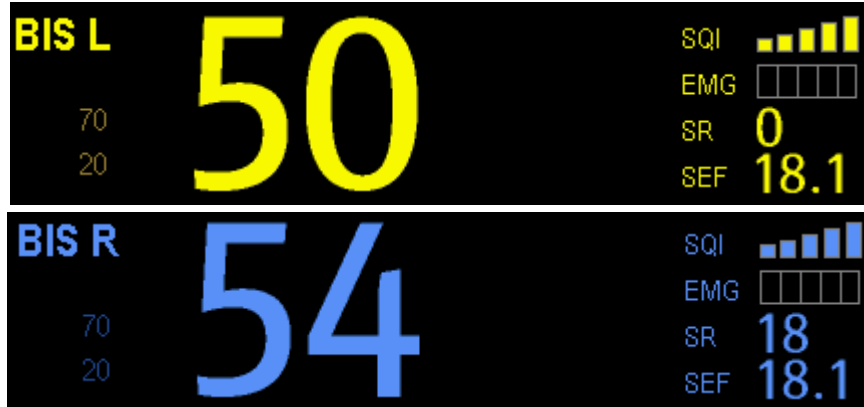
6. Total Power (TP)

TP numeric which only monitors the state of the brain indicates the power in the frequency band 0.5-30Hz. The useful range is 40-100db.

7. Burst Count (BC)

A burst means a period (at least 0.5 second) of EEG activity followed and preceded by inactivity. The BC numeric helps you quantify suppression by measuring the number of EEG bursts per minute. This parameter is intended for the BIS module with the Extend Sensor or Bilateral Sensor only. BC numeric is valid only when $SQI \geq 15\%$ and $SR \geq 5\%$.

For brain's both sides BIS monitoring, the BIS parameter area displays the following parameters (L: Left brain hemisphere; R: Right brain hemisphere):



1. BIS L BIS R
2. EMG L EMG R
3. SR L SR R
4. SEF L SEF R
5. SQI L SQI R
6. TP L TP R
7. BC L BC R
8. sBIS L sBIS R

sBIS (BIS Variability Index)

This numeric represents the standard deviation of the BIS variable over the last three minutes.

9. sEMG L sEMG R

sEMG (EMG Variability Index)

This numeric represents the standard deviation of the EMG value over the last three minutes.

10. ASYM

Asymmetry (ASYM) is a processed variable indicating the percentage of EEG power present in left or right hemispheres with respect to total (left and right) EEG power.

Designation 'L' of the asymmetry data indicates asymmetry to the left side.

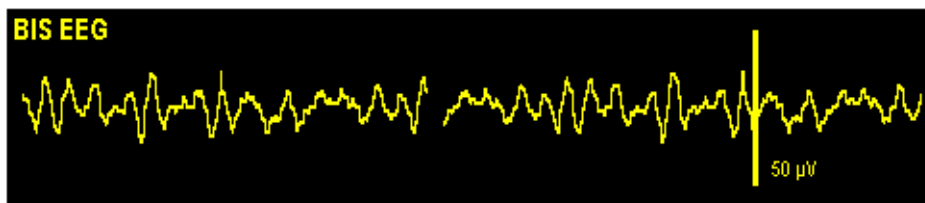
Designation 'R' of the asymmetry data indicates asymmetry to the right side.

23.3.2 BIS Waveform Area


The BIS waveform area allows you to view either EEG waveform or BIS trend. A secondary parameter's trend line can also be displayed together with BIS trend line.

1. Enter the **[BIS Setup]** menu.
2. Select **[Display]** and then select the desired option..


◆ **[EEG]**



◆ **[BIS Trend]**

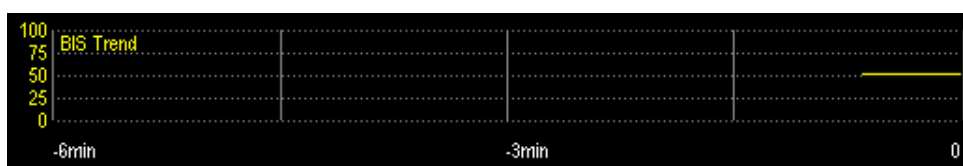

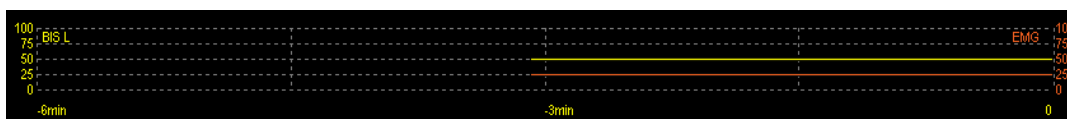- ◆ The available options for BIS trend superimpose display include: **[BIS+EMG Trend]**, **[BIS+SQI Trend]**, **[BIS+SR Trend]**, **[BIS L+R Trend]**, **[BIS+BC Trend]**, **[BIS+sBIS Trend]** or **[BIS+sEMG Trend]**, depending on the sensor type.




23.3.3 BIS Expand View

When BIS Bilateral Sensor is used for bilateral monitoring, BIS expand view can be displayed.


1. Enter the **[BIS Setup]** menu.
2. Select **[BIS Expand View >>]**.
3. Select **[Display]** and then toggle between **[EEG]**, **[BIS Trend]** and **[DSA]**.

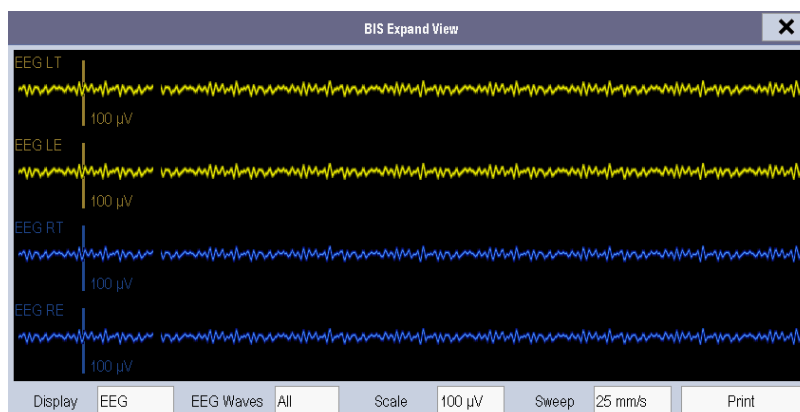23.3.3.1 Displaying EEG Waveforms




You can select the EEG waveforms to be displayed. You can also select the desired scale and sweep speed.

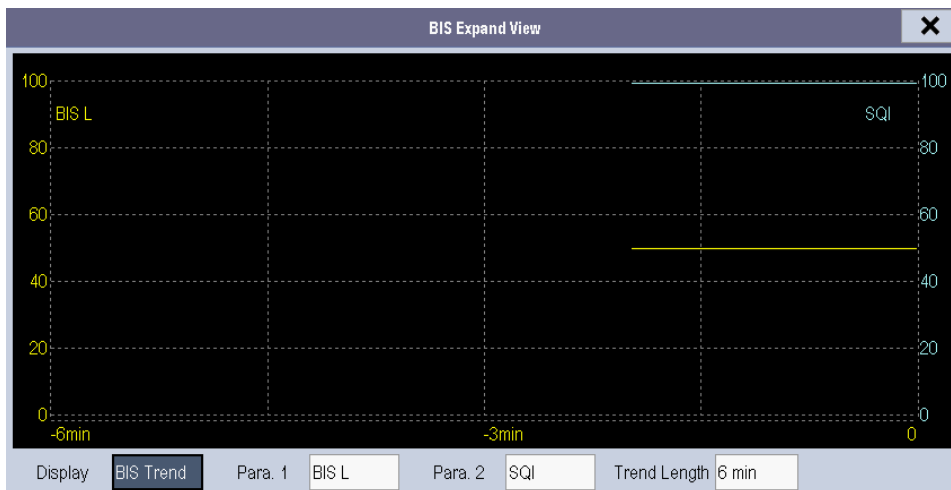23.3.3.2 Displaying BIS Trend



You can the desired trend lines to be displayed and set the time scale. The artifact mark is displayed at the bottom to indicate SQI value. When $SQI < 15\%$, the artifact mark is yellow and the corresponding trend lines of BIS, SR, BC and sBIS are not displayed. When $15\% \leq SQI < 50\%$, the artifact mark is brown.

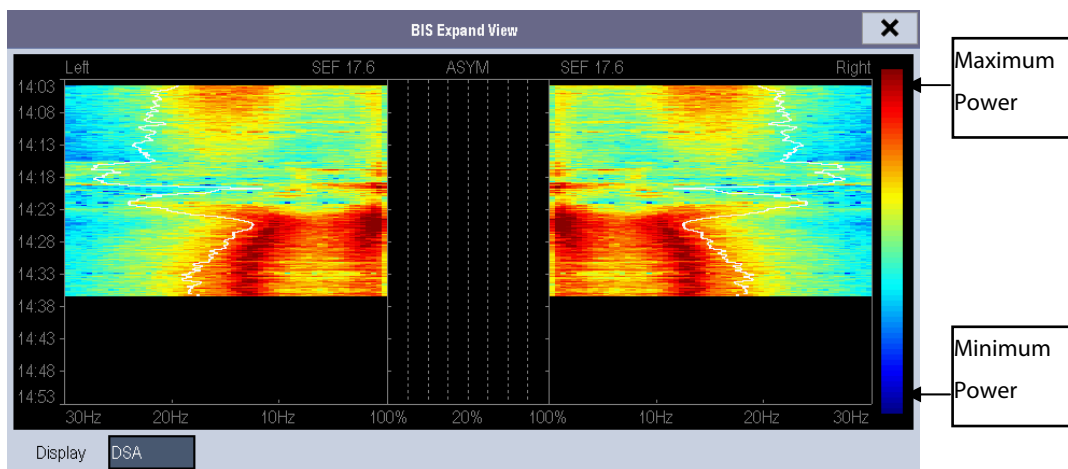23.3.3.3 Displaying DSA




The Density Spectral Array (DSA) shows changes in the power spectrum distribution over a certain time period. The DSA represents the power spectra ranging from 49-94 dB. The color bar to the right of the time scale shows the range of colors used to indicate minimum and maximum power. The frequency scale is shown on the horizontal axis with a range from 0-30 Hz.

A white Spectral Edge line is superimposed on the graph where 95% of the total power lies on one side of the line (toward the inside of the graph) and 5% lies on the other. The Spectral Edge Frequency value (SEF) displays above the graph.

The ASYM graph in the center of the screen shows the degree of asymmetry in EEG power between the left and right hemispheres. The ASYM scale begins at 20% at the center line and runs left or right to 100%. Asymmetry data less than 20% are not displayed on the graph, but are available in the tabular trends.

23.4 Setting up the BIS Measurement

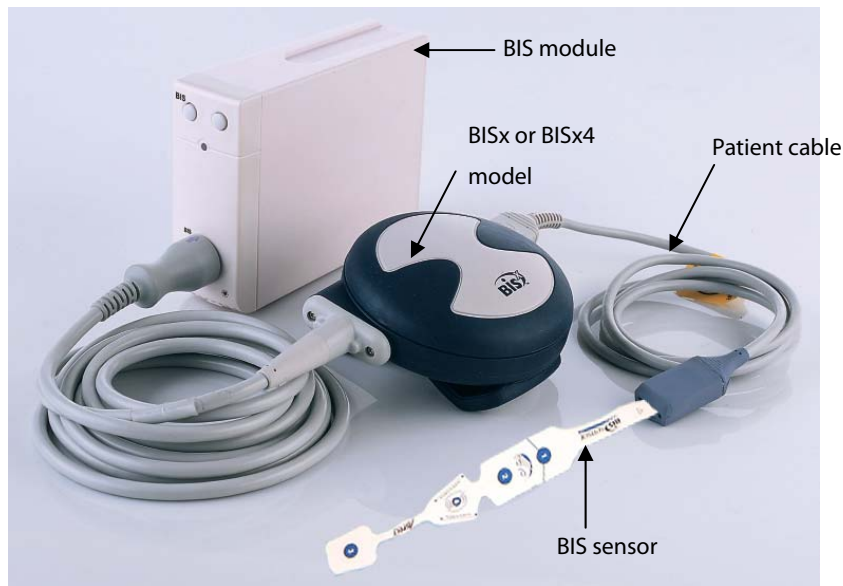1. Connect the BISx or BISx4 model to the BIS module.



2. Use the attachment clip to secure the BISx or BISx4 model near, but not above the level of the patient's head.
3. Connect the BISx or BISx4 model to the patient cable.
4. Attach the BIS sensor to the patient following the instructions supplied with sensor.

NOTE

-
- **Make sure the patient's skin is dry. A wet sensor or a salt bridge could result in erroneous BIS and impedance values.**
-

5. Connect the BIS sensor to the patient interface cable.

CAUTION

-
- **Do not attach the BISx or BISx4 model to the patient's skin for long time. Otherwise, the BISx or BISx4 heats while on the patient and may cause discomfort.**
-

23.5 Auto Impedance Check

By default, this check is switched on. It checks:

- The combined impedance of the signal electrodes plus the reference electrode. This is done automatically and continuously and does not affect the EEG wave. As long as the impedances are within the valid range, there is no prompt message of this check or its results.
- The impedance of the ground electrode. This is done every ten minutes and takes approximately four seconds. It causes an artifact in the EEG wave, and the message [**BIS Ground Checking**] is displayed on the monitor during the check. If the ground electrode does not pass this check, another check is initiated. This continues until the ground electrode passes the check.

If the auto impedance check interferes with other measurements, it can be switched off. To do this:



1. Select [**Sensor Check**] in the [**BIS Setup**] menu to open the sensor check window.
2. Set [**Automatic Check**] to [**Off**].

CAUTION

- **Switching the auto impedance check off will disable automatic prompt to the user of impedance value changes, which may lead to incorrect BIS values. Therefore, this should only be done if the check interferes with or disturbs other measurements.**
-

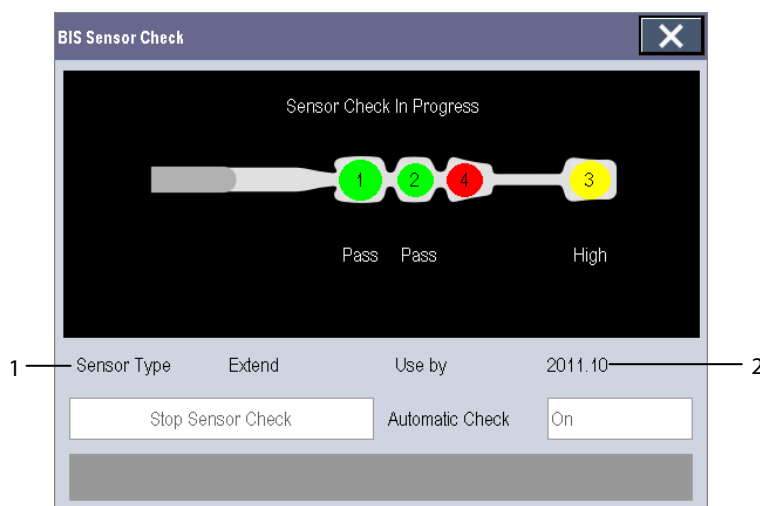
23.6 Sensor Check

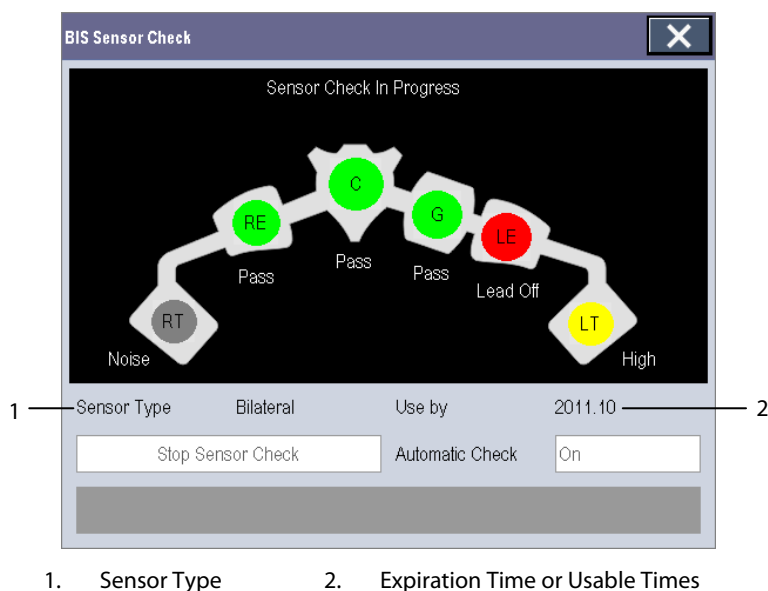
This measures the exact impedance of each individual electrode. It causes a disturbed EEG wave, and a prompt message is displayed on the monitor

- The sensor check is automatically initiated when a sensor is connected. To manually start a sensor check, you can either:
 - ◆ Press the  hardkey on the BIS module.
 - ◆ Select [**Sensor Check**] in the [**BIS Setup**] menu.
 - ◆ Select [**Start Sensor Check**] in the BIS sensor window.
- The sensor check stops automatically if the impedances of all electrodes are within the valid range. To manually stop a sensor check, you can either:
 - ◆ Press the  hardkey on the BIS module.
 - ◆ Select [**Stop Sensor Check**] in the sensor check window.

23.7 BIS Sensor Check Window

To open the sensor check window, select [**Sensor Check**] in the [**BIS Setup**] menu. The graphic in the BIS sensor check window automatically adapts to show the type of sensor you are using, show each electrode as required. Each symbol in the graphic represents an electrode and illustrates the most recently-measured impedance status of the electrodes.





Different colors indicate different statuses. The electrode status is displayed below each electrode:

Color	Status	Description	Action
Red	[Lead Off]	Electrode falls off and has no skin contact	Reconnect electrode, or check the sensor-to-skin contact. If necessary, clean and dry skin.
Grey	[Noise]	The EEG signal is too noisy. Impedance cannot be measured	Check the sensor-to-skin contact. If necessary, clean and dry skin.
Yellow	[High]	The impedance is above the limit	
Green	[Pass]	The impedance is within valid range	No action necessary.

Although BIS may still be measured when the electrode status is **[Noise]** or **[High]**, for best performance, all electrodes should be in **[Pass]** status.

23.8 Choosing the BIS Smoothing Rate

To change the smoothing rate:

1. Select the BIS parameter window to enter the **[BIS Setup]** menu.
2. Select **[Smoothing Rate]** and then toggle between **[10 s]**, **[15 s]** and **[30 s]**

The smoothing rate defines how the monitor averages the BIS value. With the smoothing rate becoming smaller, the monitor provides increased response to changes in the patient's state. Contrarily, the monitor provides a smoother BIS trend with decreased variability and sensitivity to artifacts.

NOTE

- When **[Smoothing Rate]** is set as **[10 s]** or **[30 s]**, **sBIS** and **sEMG** are displayed as invalid values.

23.9 Changing the Secondary Parameters

You can choose the desired secondary parameters for display on the screen.

1. Enter the [**BIS Setup**] menu.
2. Select [**Change Secondary Parameter**>>] and then select at most 2 desired parameters from the popup menu.

23.10 Changing the EEG Wave Size

1. Enter the [**BIS Setup**] menu.
2. Select [**EEG**] from [**Display**].
3. Select [**Scale**] and then select the appropriate setting.

23.11 Changing the Speed of the EEG Wave

1. Enter the [**BIS Setup**] menu.
2. Select [**EEG**] from [**Display**].
3. Select [**Sweep**] and then select the appropriate setting. The faster the wave sweeps, the wider the wave is.

23.12 Setting the Trend Length

1. Enter the [**BIS Setup**] menu.
2. Select a BIS trend option from [**Display**].
3. Select [**Trend Length**] and then select the appropriate BIS time length setting.

23.13 Switching the Filter On or Off

1. Enter the [**BIS Setup**] menu.
2. Select [**Filter**] and then toggle between [**On**] and [**Off**]. The default is [**On**].

The filter screens out undesirable interference from the raw EEG wave display. The notch filter includes filters for both 50 and 60 Hz. Filter settings do not affect processing of the trend variables (i.e., BIS, EMG, and SR).

24 Monitoring NMT

24.1 Introduction

This patient monitor can connect a Organon TOF-Watch® SX monitor for NMT(neuromuscular transmission) monitoring . This patient monitor can display, store and review measurements from TOF-Watch® SX monitor, as well as present related alarms. On this patient monitor, you can separately set the level of NMT related alarms and switch on or off alarm recording; you can also view TOF-Watch® SX monitor settings of alarm limits and alarm switch.

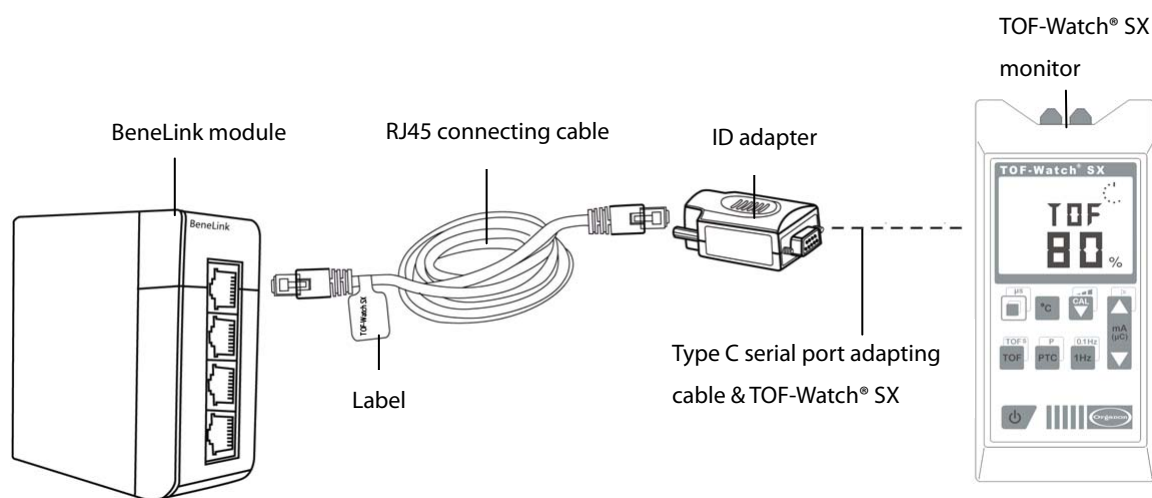
24.2 Safty

WARNING

- **TOF-Watch® SX monitor is manufacutred by Organon. This company provides the technology for measuring NMT parameters. We only provide the connection between this patient monitor and TOF-Watch® SX monitor.**
 - **If you have any doubts about the operation and maintenance of the TOF-Watch® SX monitor, please refer to TOF-Watch® SX monitor operator’s manual or directly contact Organon .**
 - **Fully observe TOF-Watch® SX monitor operator’s manual to make settings and to connect the monitor with a patient.**
-

24.3 Connecting a TOF-Watch® SX monitor

The TOF-Watch® SX monitor connects with BeneLink module through an ID adapter, see the picture below.



Please refer to the following procedure to connect the TOF-Watch® SX monitor:

1. Insert a BeneLink module into a BeneView patient monitor module rack.
2. Connect the ID adapter that matches the TOF-Watch® SX monitor to the BeneLink module with an RJ45 connecting cable.
3. Connect the ID adapter to the TOF-Watch® SX interface with Mindray type C serial port adapting cable (PN: 009-001769-00).
4. Connect the TOF-Watch® SX interface to the TOF-Watch® SX monitor.
5. Stick a label indicating device name to the RJ45 connecting cable at the end nearby the BeneLink module. When the BeneLink module is connected to several external devices, you can tell the devices easily with these labels.
6. Turn on both monitors.

24.4 NMT Parameters

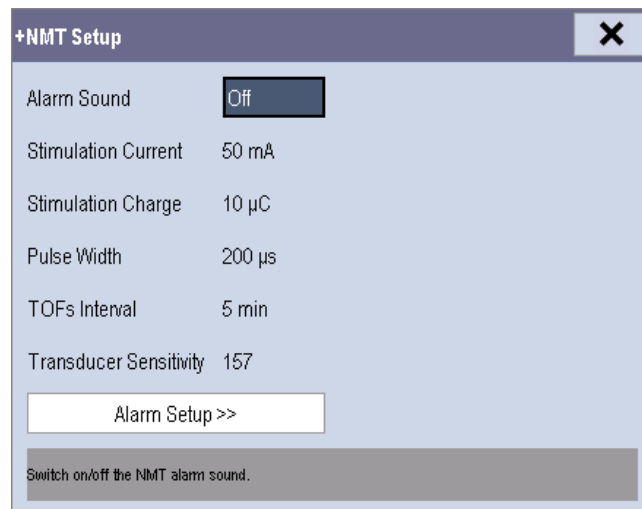
TOF-Watch® SX monitor provides the following measurements:

- TOF-Ratio
- TOF-Count
- PTC
- Single
- Tskin

24.5 Accessing the NMT Setup menu

You can access the [+NMT Setup] menu by selecting the NMT area or selecting [Main Menu] → [Parameters


>>] → [+NMT Setup>>].



In the [+NMT Setup] menu, you can

- Toggle [**Alarm Sound**] between [**On**] and [**Off**] to switch on or off NMT alarms on this patient monitor.
- View the setup as follows:
 - ◆ Stimulation Current
 - ◆ Stimulation Charge
 - ◆ Pulse Width
 - ◆ TOFs Interval
 - ◆ Transducer Sensitivity
- Set alarm level for TOF-Ratio and TOF-Count, switch on or off alarm record.

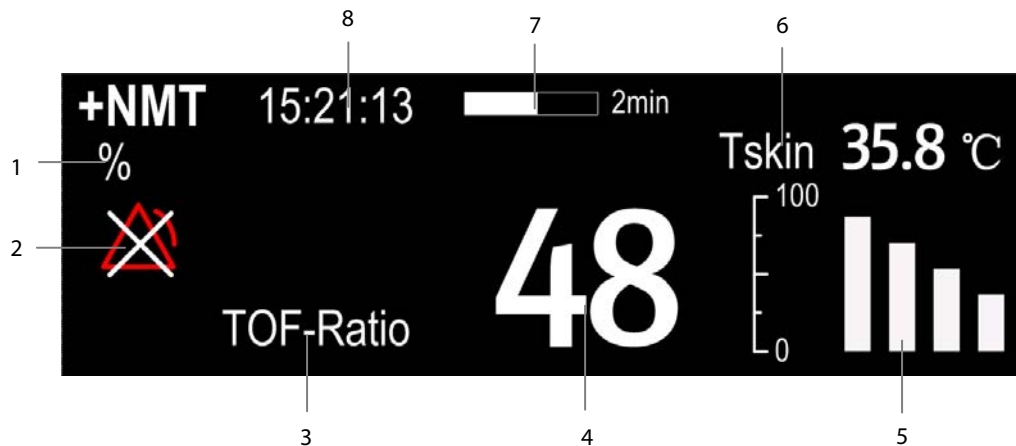24.6 NMT Display



1. Parameter unit
2. Alarm status
3. Parameter label
4. Parameter measurement
5. Response amplitude of stimulation
6. Skin temperature
7. Measurement countdown
8. Time of last measurement

In the case that you take a measurement in TET50Hz mode, TET100Hz mode, DBS3.3 mode or DBS3.2 mode, only mode label is displayed in the NMT parameter area, which is shown as follows:



Measuring mode


25 Monitoring RM

25.1 Introduction

 **WARNING**

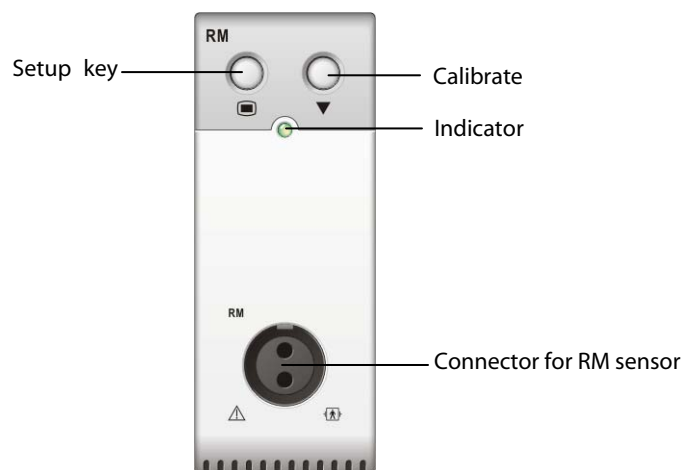- **RM monitoring is not intended for neonatal patients.**
-
-

In the respiratory mechanics measurement, the airway pressures are measured, from the part between the patient circuit and intubation tube, using a flow sensor between the Y-piece of patient circuit and the patient connection. The pressure is transferred to the monitor through the tube and measured by a pressure transducer in the RM module. The pressure difference together with the gas concentration information is used to calculate flow. The volume information is obtained by integrating the flow signal. From these three parameter, other parameter such as RR, I:E, Compl, etc. are derived.



The RM monitoring enables clinicians to understand the ventilator operation and patient respiratory status.

RM monitoring displays the following waveforms and loops:

- Flow waveform
- Paw waveform
- Vol waveform
- FV (flow-volume) loop
- PV (paw-volume) loop

RM monitoring provides values for 15 parameters. The 15 parameters can be classified into 4 categories:

1. Paw parameters

- ◆ PIP: peak inspiratory pressure (unit: cmH₂O)
- ◆ Pplat: pressure (unit:cmH₂O)
- ◆ PEEP: positive end expiratory pressure (unit: cmH₂O)
- ◆ Pmean: mean pressure (unit: cmH₂O)

2. Flow parameters

- ◆ PIF: peak inspiratory flow (unit: L/min)
- ◆ PEF: peak expiratory flow (unit: L/min)

3. Vol parameters

- ◆ TVi: inspiratory tidal volume (unit: ml)
- ◆ TVe: expiratory tidal volume (unit: ml)
- ◆ MVi: inspirator minute volume (L)
- ◆ MVe: expiratory minute volume (L)

4. Other parameters

- ◆ RR: respiratory rate (unit: rpm)
- ◆ I: E: ratio of the inspiratory and expiratory time
- ◆ Compl: compliance (unit: ml/cmH₂O)
- ◆ FEV1.0: first second forced expiratory volume ratio (unit: %)
- ◆ RSBI: rapid shallow breathing index (unit: rpm/L)

25.2 Safety Information

WARNING

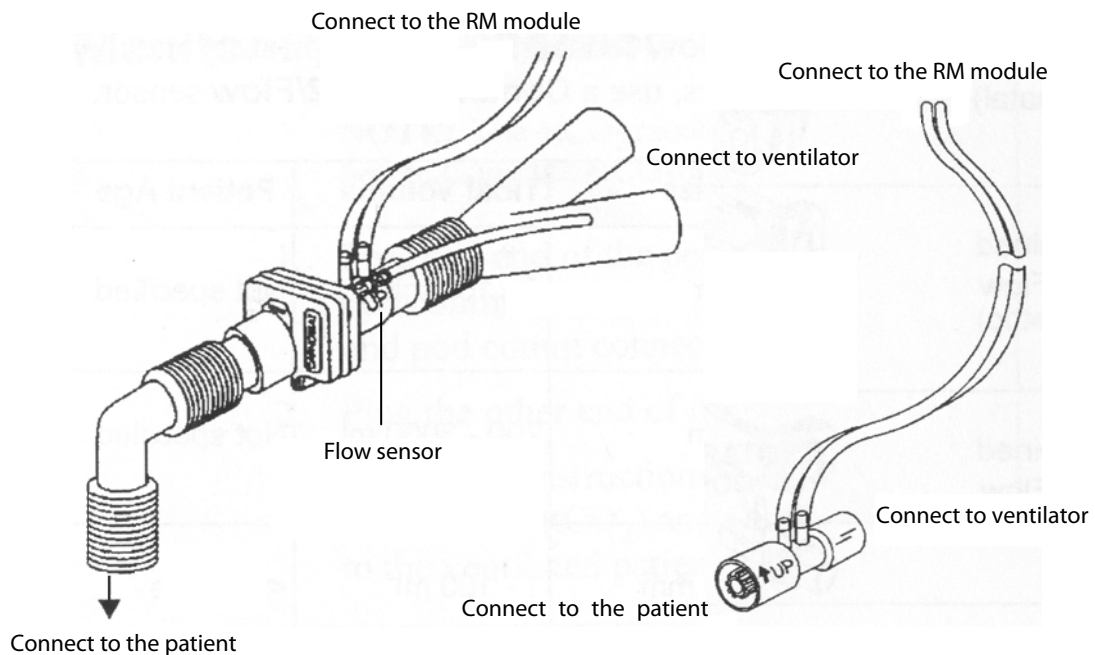
- **Check for leaks in the breathing circuit system, as they may significantly affect respiratory mechanics readings.**
 - **Match the airway adapter you select to the appropriate patient category. Improper sensor selection may produce excessive ventilation resistance or introduce excessive airway deadspace, as well as inaccurate scales and alarm limits.**
 - **Periodically check the flow sensor and tubing for excessive moisture or secretion build-up and purge if necessary.**
-

NOTE

- To avoid the affects of excessive moisture in the measurement circuit, insert the flow sensor airway adapter in the breathing circuit with the tubes upright.
- Do not place the airway adapter between the endotracheal tube and an elbow as this may allow patient secretions to block the adapter windows.
- Measurement values provided by a ventilator may differ significantly from the values provided by the RM module, due to different locations of the flow sensor.

25.3 Preparing to Monitor RM


1. Select an appropriate flow sensor in accordance with the patient category.
2. Connect the small tubes of the flow sensor to the RM connector of the module using a color-coded adapter.
3. Insert the flow sensor between the Y-piece of the patient circuit and the patient connection.



4. Calibrate the flow sensor: select the RM parameter window to open the [**Calibrate RM**] menu when you see the prompt message [**Calibration Required**] appears on the RM parameter window. Then calibrate the flow sensor according to the procedure described in **25.8 Calibrating the Flow Sensor**.
5. Select [**Sensor Type**] in the [**Calibrate RM**] menu and then choose [**Infant One-time**], [**Disposable**] or [**Reusable**] according to the selected sensor.

25.4 Understanding the RM Display

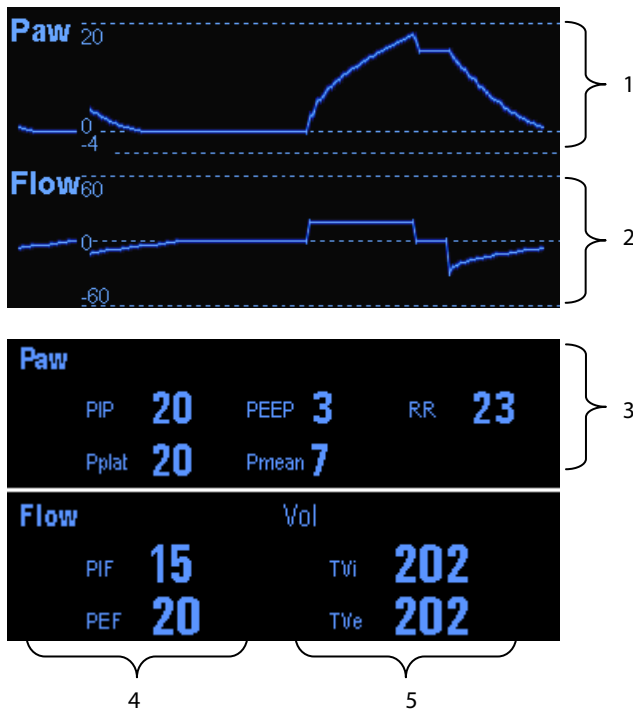The RM display shows either the Paw and Flow waveforms, or the Paw and Vol waveforms in the waveform area.



1. Paw waveform
2. Flow waveform
3. Paw parameter window
4. Flow parameter window
5. Vol parameter window

25.5 Changing RM Settings

25.5.1 Accessing RM Menus

- By selecting the RM parameter window or waveform area, you can access the **[RM Setup]** menu.
- By selecting the Paw wave, you can access the **[Paw Waveform]** menu.
- By selecting the Flow wave, you can access the **[Flow Waveform]** menu.
- By selecting the Vol wave, you can access the **[Vol Waveform]** menu.

25.5.2 Setting the Apnea Alarm Delay

In the **[RM Setup]** menu, select **[Apnea Delay]** and then select the appropriate setting. The monitor will alarm if the patient has stopped breathing for longer than the preset apnea time. The **[Apnea Delay]** of Resp, CO₂, AG, and RM module keeps consistent with each other.

WARNING

- **The respiration measurement does not recognize the cause of apneas. It only indicates an alarm if no breath is detected when a pre-adjusted time has elapsed since the last detected breath. Therefore, it cannot be used for diagnostic purpose.**
-

25.5.3 Selecting TV or MV for Display

To select tidal volume (TV) or minute volume (MV) for display in the Vol parameter window, in the **[RM Setup]** menu, select **[TV/MV]** and toggle between **[TV]** and **[MV]**. By default, the Vol parameter window displays TV values.

25.5.4 Selecting Flow or Vol Waveform for Display

To select Flow or Vol waveform for display:

1. Enter the **[RM Setup]** menu.
2. Select **[Flow/Vol]** and toggle between **[Flow]** and **[Vol]**.

25.5.5 Changing the Wave Sweep Speed

1. Enter the **[RM Setup]** menu.
2. Select **[Sweep]** and select the appropriate setting. The faster the wave sweeps, the wider the wave is.

25.5.6 Changing the Wave Scale

1. Select **[Wave Scale >>]** from the **[RM Setup]** menu..
2. Select the appropriate settings in the popup menu.

25.5.7 Setting RR Source

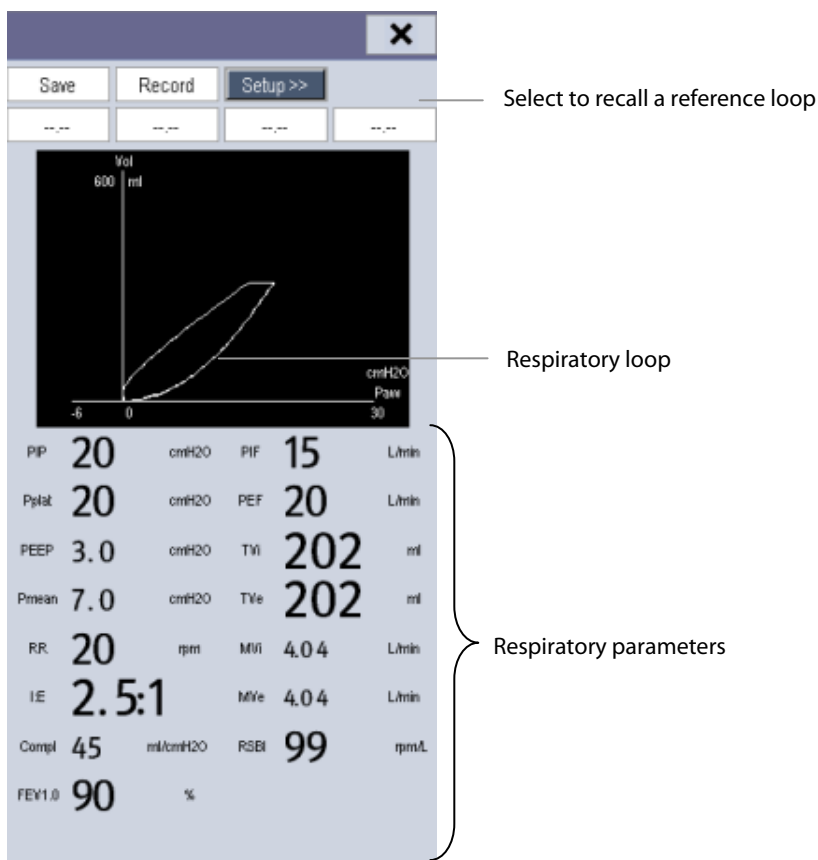
To set RR source:

1. Enter the **[RM Setup]** menu.
2. Select **[RR Source]** and then select a source or **[Auto]** from the dropdown list.

The **[RR Source]** settings of Resp, CO₂, AG and RM module are linked. For details, please refer to the section **Setting RR Source** of chapter *Resp*.

25.6 Understanding the Respiratory Loops

Select **[Respiratory Loop]** in the **[RM Setup]** menu. The following window will be displayed.



In this window, you can:

- Select **[Save]** to save the respiratory loops in the current respiratory cycle as the reference loops. Up to 4 groups of respiratory loops can be saved, and the saving time is displayed above the respiratory loops.
- Change the respiratory loops displayed on the screen: select **[Setup >>]** → **[Display Loop]** and then toggle between **[PV Loop]** and **[FV Loop]**.
- Turn on/off reference loop: select **[Setup >>]** → **[Reference Loop]**, and then toggle between **[On]** and **[Off]**.
- Change the size of the PV and FV loops: select **[Setup >>]**, and then adjust the **[Paw Scale]**, **[Vol Scale]** or **[Flow Scale]**.
- Select parameters for display: select **[Setup >>]** → **[Select RM Parameters >>]**, and then select **[All RM Parameters]** or **[Select Desired RM Parameters]**. When you select **[Select Desired RM Parameters]**, 6 parameters at maximum can be selected.
- Print out all parameters for a reference loop by selecting your desired reference loop and then selecting **[Record]**.

25.7 Zeroing the RM Module

A zero calibration is carried out automatically every time when the patient monitor is switched on or the RM module connected, and then a zero calibration will automatically be triggered at a specific interval. Then, a zero calibration is triggered every 5 minutes. You can also start a manual zero calibration when there is a drift in the zero: in the **[RM Setup]** menu, select **[Zero RM]**.

25.8 Calibrating the Flow Sensor

A calibration must be performed every time when the RM module is connected to the patient monitor or the flow sensor is connected.

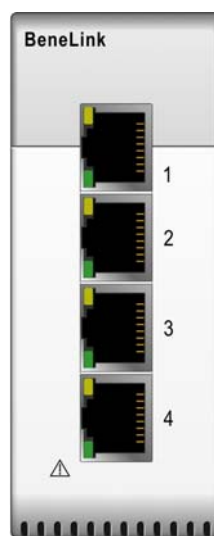
1. When calibration is needed, the RM parameter window displays **[Calibration Required]**. Select the parameter window to enter the **[Calibrate RM]** menu. You can also enter the **[Calibrate RM]** menu by selecting **[Calibrate RM >>]** from the **[RM Setup]** menu.
2. Select **[Sensor Type]** and then choose **[Infant One-time]**, **[Disposable]** or **[Reusable]** according to the sensor used.
3. Enter the positive and negative factor provided on the flow sensor and select **[Calibrate]**.
4. After the calibration is completed successfully, the last calibration time and the message **[Calibration Completed!]** are displayed. Otherwise, the message **[Calibration Failed!]** is displayed.

FOR YOUR NOTES

26 Device Integration

26.1 Introduction

BeneLink module is intended for connecting external devices, such as ventilators and anesthesia machines, to the BeneView patient monitor. It allows the information (patient data, alarms, etc.) from the external device to be displayed, saved, recorded, printed, or calculated through a BeneView patient monitor. If the patient monitor is connected with the CMS or gateway, information from the external device can also be transmitted to the CMS or gateway.



26.2 Safety Information

WARNING

- **Devices of the same category can not be connected to the BeneLink module simultaneously.**
 - **A patient monitor supports one BeneLink module only.**
 - **The signal labels used on the BeneView patient monitor may be different from those given on the external device. For details please see the description of parameters and alarms in corresponding sections of this chapter.**
 - **The alarms from the external device may be advanced or delayed before transmission to the BeneView patient monitor.**
 - **There can be differences between the alarm priorities displayed on your BeneView patient monitors and the priorities displayed on the external devices interfaced through BeneLink. Please see the list of Output Signals corresponding with each external device for the alarm priorities used by your patient monitor.**
-

26.3 Supported Devices

Category	Model
Anesthesia Machine	Mindray Wato 20/30/55/65
	Maquet Flow-i
	Draeger Fabius G5/Fabius Trio/Fabius Plus/Primus
	GE Datex-Ohmeda Aestiva 7900/Aestiva 7100/Avance/Aisys
Ventilator	Newport E360
	Puritan Bennett 840
	Maquet Servo-i/Servo-s
	Draeger Evita 2
	Draeger Evita 4/ Evita2 dura/Evita XL
	Hamilton G5/C2 /Galileo

NOTE

- **BeneLink module may support more devices than those listed in the above table. Please connect us or our service personnel for the most recent information on the supported devices.**

26.4 Differences in Displayed Values

In certain cases, there may be differences between the numerics seen on the BeneView patient monitor and those seen on the external device. The table below lists some situations and possible reasons.


Situation	Possible Reasons
Some parameter values are displayed as invalid values on the BeneView patient monitor.	The patient monitor and the external device may have different parameter configuration or displaying range of values. If the patient monitor displays a parameter that is not configured in the external device or a parameter value from the external device exceeds the displaying range of the patient monitor, the corresponding parameter value is displayed on the patient monitor as an invalid value.
The patient monitor and the external device may display the parameter values with different numbers of places of decimals.	The patient monitor displays the parameter values from the external device based on its own display rules. Same parameter value is displayed differently when the patient monitor and external device adopt different numbers of places of decimals of the value for display.
Non-continuously measured values and continuously measured values have the same displaying mode in the BeneView patient monitor.	Non-continuously measured values are displayed on the patient monitor as latest measured values until a new measurement is performed on the external device.
Differences between the parameter values displayed on the BeneView patient monitor and those displayed on the external device.	Some parameter values are converted to different units during transmission to the patient monitor so that they can be used for calculations. Sometimes, values from the external device may be advanced or delayed before transmission to the BeneView patient monitor.

NOTE

- When the pressure units are converted among cmH₂O, hPa and mbar, the parameter value remain unchanged, for example, 1cmH₂O=1hPa=1mbar, which may differ from some external devices.

26.5 Connecting an External Device

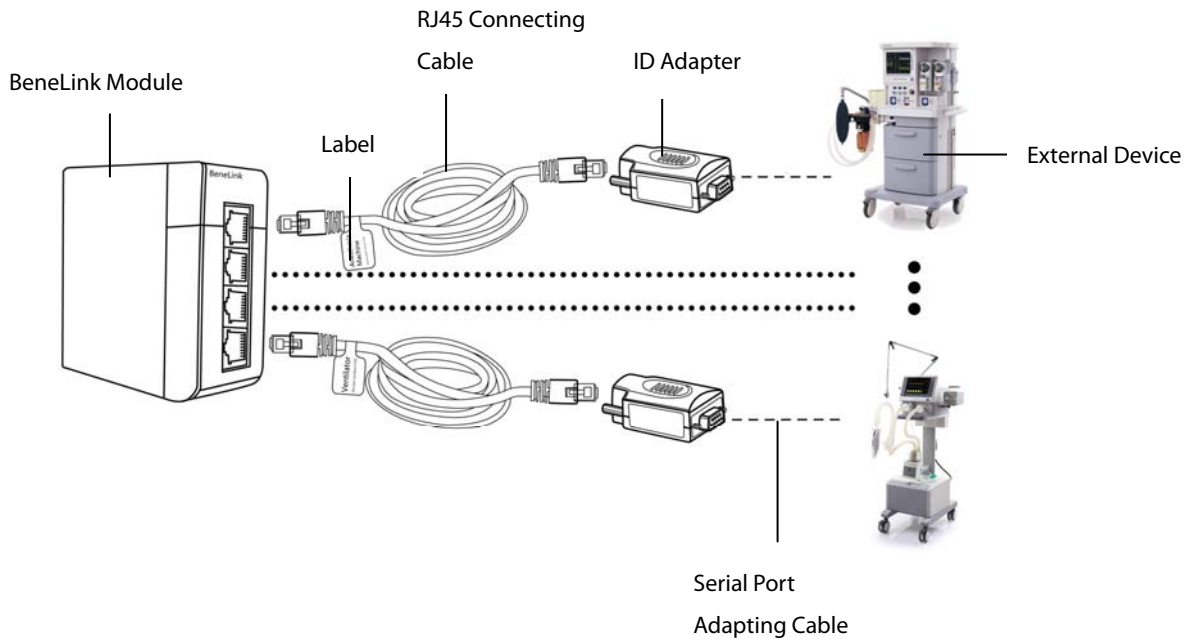The external device connects with the BeneLink module through an ID adapter, which supports only its matching device. Please refer to the following procedure to connect an external device:



1. Insert the BeneLink module into the module slot on the BeneView patient monitor.
2. Connect the ID adapter that matches the external device to the BeneLink module with an RJ45 connecting cable.
3. Plug the ID adapter into the RS232 port on the external device. Some external devices may have ports incompatible with the ID adapter. In this case, a serial port adapting cable is required. Please see the following table for the required adapting cable.
4. Stick a label indicating device name to the RJ45 connecting cable at the end nearby the BeneLink module. When the BeneLink module is connected to several external devices, you can tell the devices apart easily with these labels.
5. Switch the external device on.

After the external device is connected to the patient monitor, the indicating lights on both the ID adapter and the BeneLink module illuminate to show that the patient monitor communicates with the external device successfully.

The ID adapter has already been correctly configured before leaving the factory. If you want to re-configure the ID adapter, please select **[Main Menu]** → **[Maintenance>>]** → **[Factory Maintenance>>]** → enter the required password → **[Upgrade ID module>>]**, and follow this procedure:

1. Set **[Benelink Module Port]** to select which port the RJ45 connecting cable is connected to. You must connect the RJ45 connecting cable to the selected port when re-configuring the ID adapter. Otherwise, ID adapter re-configuration will fail.
2. Set **[ID]** to configure a new ID to the ID adapter.

External Device	ID for ID adapter	Type of serial port adapting cable
Mindray Wato 20/30/55/65	4D52B2AE	No need to use the adapting cable: the ID adapter can be plugged into the serial port of the external device directly.
Newport E360	4E50B1B0	Type B
Puritan Bennett 840	SNDF: 5042AFBE(recommended) SNDA: 5031AFCF(support less parameters than protocol SNDF)	No need to use the adapting cable. The ID adapter can be plugged into the serial port of the external device directly.
Maquet Flow-i	4D46B2BA	Type B
Maquet Servo-i/Servo-s	4D53B2AD	Type B
Draeger Evita 2 / Evita 2 dura / Evita 4/ Evita XL	4434BBCC	Type B
Hamilton G5 (protocol Block)	3542CABE	Type B
Hamilton G5 (protocol Polling)	3550CAB0	Type B
Hamilton C2	3270CD90	Type B
Hamilton Galileo	4750B8B0	Type B
GE Datex-Ohmeda Avance/Aisys	4F41B0BF	Type D
GE Datex-Ohmeda Aestiva 7100/7900	4F37B0C9	Type D
Drager Fabius GS/ Fabius Plus/ Fabius Trio	4446BBBA	Fabius GS: No need to use the adapting cable. The ID adapter can be plugged into the serial port of the external device directly. Fabius Plus: Type C Fabius Trio: Type C
Drager Primus	4450BBB0	Type C
Extended model	/	Type A


Serial port adapting cable	PN	Remark
Type A	009-001767-00	Male to female
Type B	009-001768-00	Male to male
Type C	009-001769-00	Male to male
Type D	009-002943-00	9-pin to 15-pin

 **WARNING**

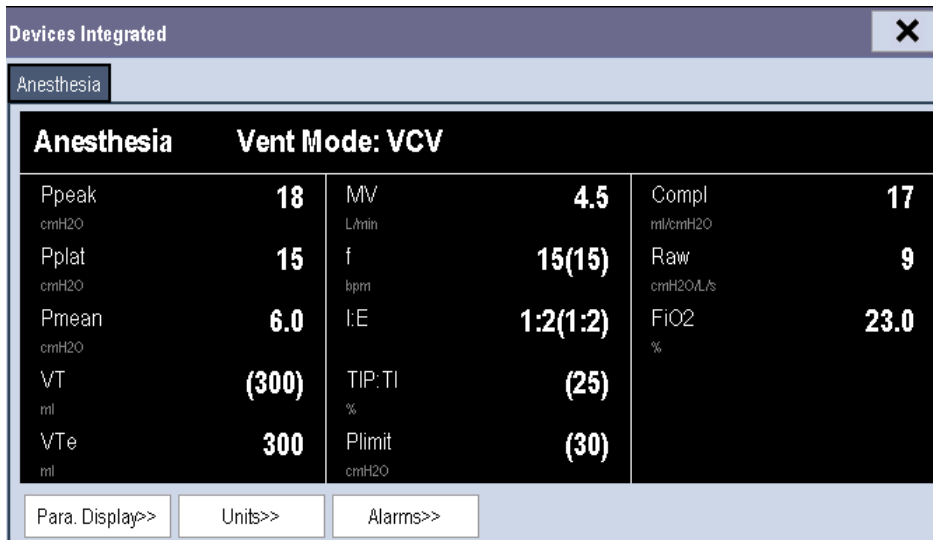- **First installation and debugging should be executed by our service personnel or authorized technician.**
 - **Please check the compatibility of the external device and the ID adapter before connection. Otherwise, unpredictable system failure may be resulted.**
 - **Ports on the BeneLink module are not normal network connectors. They are intended for connecting with the serial port of designated devices only. Do not connect them to public network interfaces.**
-

26.6 Devices Integrated Window

You can view the information of the external device in the [**Devices Integrated**] window, which provides the information of both individual devices and multi devices. In the individual device menu, you can select [**Para. Display>>**], [**Units>>**] or [**Alarms>>**] to set the parameters to be displayed or the parameter units, or view the alarm list.



Anesthesia		Vent Mode: VCV			
Ppeak cmH2O	18	MV L/min	4.5	Compl ml/cmH2O	17
Pplat cmH2O	15	f bpm	15(15)	Raw cmH2O/Ls	9
Pmean cmH2O	6.0	I:E	1:2(1:2)	FiO2 %	23.0
VT ml	(300)	TIP:TI %	(25)		
VTe ml	300	Plimit cmH2O	(30)		

The parameters in the [**Devices Integrated**] window are displayed in the order of priorities. In the case that the window can not display all the selected parameters, only parameters with higher priorities are displayed. Please refer to the following sections for parameter priorities.

For the parameter that is measured by the external device, the measurement displays directly after the parameter label. For the parameter that is controlled by the external device, its setting is enclosed in a parenthesis after the parameter label. For the parameter that can both be measured and controlled by the external device, both its measurement and setting are displayed after the parameter label, and the setting is also enclosed in a parenthesis. For example, PEEP 18 (20), in which PEEP is parameter label, 18 is the measurement, and (20) is the setting.

In the [**Devices Integrated**] window, you can select [**Multi Devices**] tab to view the parameter information of all the external devices interfaced currently. The displayed parameters are those selected in [**Para. Display**] menu of the individual device window. In the case that the patient monitor can not display all the selected parameters, only parameters of higher priorities are displayed.

26.7 System Functions of Patient Monitor

26.7.1 Alarms

The patient monitor does not display the realtime alarms from the external device. However, you can view current alarm list of the corresponding device by selecting [**Alarms>>**] in the individual device window. The alarm priority is defined by "*" before each alarm message. An alarm list can display up to 100 alarm messages.

26.7.2 Data Storage

The patient monitor can save and review the graphic trends, tabular trends, and alarm events of parameters from the external device. In [**Graphic Trends**] menu and [**Events**] menu, parameter from the external device is displayed in white. In [**Review**] menu, [**Trend Group**] menu, and [**Print Setup**] menu, a mark "+" is shown before each label of parameters from the external device. Please refer to the parameter list to see which parameters can be saved.

NOTE

- **Parameters from the external device are saved and displayed according to the time of the patient monitor.**
-

26.7.3 Recording and Printing

Information from the external device can be recorded and printed both in realtime and in graphic and tabular trends with BeneView patient monitor. Besides, the monitor can also record the frozen parameters of the external device.

26.8 Integrating the Anesthesia Machine

26.8.1 Wato 20/30/55/65

26.8.1.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
O ₂ %	Oxygen concentration	%	Yes
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	No
P _{peak}	Peak pressure	cmH ₂ O hPa mbar	Yes
P _{plat}	Plateau pressure	cmH ₂ O hPa mbar	Yes
P _{mean}	Mean pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
V _{Te}	Expiratory tidal volume	ml	Yes
MV	Minute volume	L/min	Yes
f _{tot}	Total breath rate	bpm	Yes
f	Breath rate	bpm	No
f _{SIMV}	Frequency of SIMV	bpm	No
Freq _{MIN}	Minimum breath frequency	bpm	No
I:E	Inspiratory time:Expiratory time ratio	/	No
TIP:TI	Percentage of inspiratory plateau time in inspiratory time	%	No
T _{slope}	Time for the pressure to rise to target pressure	s	No
T _{insp}	Time of inspiration	s	No
Trig Window	Trigger Window	%	No
P _{limit}	Pressure limit level	cmH ₂ O hPa mbar	No
P _{insp}	Pressure control level of inspiration	cmH ₂ O hPa mbar	No
P _{supp}	Pressure support level	cmH ₂ O hPa mbar	No
P-Trigger	Inspiratory trigger level (pressure trigger)	cmH ₂ O hPa mbar	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
Exp%	Inspiration termination level	%	No
Compl	Compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RAW	Airway resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
FiCO ₂	Fraction of inspired carbon dioxide	% mmHg kPa	Yes
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
EtO ₂	End-tidal O ₂	% mmHg kPa	Yes
FiN ₂ O	Fraction of inspired nitrous oxide	%	Yes
EtN ₂ O	End-tidal N ₂ O	%	Yes
FiDes	Inspired anesthetic agent	%	Yes
FiSev		%	Yes
FiEnf		%	Yes
Filso		%	Yes
FiHal		%	Yes
MAC		Minimum alveolar concentration	/
EtDes	End-tidal anesthetic agent	%	Yes
EtSev		%	Yes
EtEnf		%	Yes
Etlso		%	Yes
EtHal		%	Yes
N ₂ O Flow		N ₂ O flow	L/min
Air Flow	Air flow	L/min	No
O ₂ Flow	O ₂ flow	L/min	No
BIS	Bispectral index	/	Yes
SQI	Signal quality index	/	Yes
SR	Suppression ratio	/	Yes
EMG	Electromyograph	dB	Yes
SEF	Spectral edge frequency	Hz	Yes
TP	Total power	dB	Yes
BC	Burst count	/min	Yes

26.8.1.2 Output Signals—Alarms

BeneView		Wato
Priority	Label	Label
High	Apnea	Apnea Alarm
High	Volume Apnea > 2 min	Volume Apnea>2min
High	Paw Too High	Paw Too High
High	Paw Too Low	Paw Too Low
High	EtO ₂ Too High	EtO ₂ Too High
High	EtO ₂ Too Low	EtO ₂ Too Low
High	FiO ₂ Too High	FiO ₂ Too High
High	FiO ₂ Too Low	FiO ₂ Too Low
High	Drive Gas Pressure Low	Drive Gas Pressure Low
High	O ₂ Supply Failure	O ₂ Supply Failure
Mediate	FiO ₂ Too High	FiO ₂ Too High
Mediate	VTe Too High	TVe Too High
Mediate	VTe Too Low	TVe Too Low
Mediate	MV Too High	MV Too High
Mediate	MV Too Low	MV Too Low
Mediate	EtCO ₂ Too High	EtCO ₂ Too High
Mediate	EtCO ₂ Too Low	EtCO ₂ Too Low
Mediate	FiCO ₂ Too High	FiCO ₂ Too High
Mediate	FiCO ₂ Too Low	FiCO ₂ Too Low
Mediate	EtN ₂ O Too High	EtN ₂ O Too High
Mediate	EtN ₂ O Too Low	EtN ₂ O Too Low
Mediate	FiN ₂ O Too High	FiN ₂ O Too High
Mediate	FiN ₂ O Too Low	FiN ₂ O Too Low
Mediate	EtHal Too High	EtHal Too High
Mediate	EtHal Too Low	EtHal Too Low
Mediate	FiHal Too High	FiHal Too High
Mediate	FiHal Too Low	FiHal Too Low
Mediate	EtEnf Too High	EtEnf Too High
Mediate	EtEnf Too Low	EtEnf Too Low
Mediate	FiEnf Too High	FiEnf Too High
Mediate	FiEnf Too Low	FiEnf Too Low
Mediate	EtIso Too High	EtIso Too High
Mediate	EtIso Too Low	EtIso Too Low
Mediate	FiIso Too High	FiIso Too High
Mediate	FiIso Too Low	FiIso Too Low
Mediate	EtSev Too High	EtSev Too High
Mediate	EtSev Too Low	EtSev Too Low
Mediate	FiSev Too High	FiSev Too High
Mediate	FiSev Too Low	FiSev Too Low
Mediate	EtDes Too High	EtDes Too High
Mediate	EtDes Too Low	EtDes Too Low

BeneView		Wato
Priority	Label	Label
Mediate	FiDes Too High	FiDes Too High
Mediate	FiDes Too Low	FiDes Too Low
Mediate	BIS Too High	BIS Too High
Mediate	BIS Too Low	BIS Too Low
Mediate	Patient Circuit Leak	Patient Circuit Leak
Low	RR Too High	Rate Too High
Low	RR Too Low	Rate Too Low
Low	Pressure Limiting	Pressure Limiting
Low	O ₂ Sensor Unconnected	O ₂ Sensor Unconnected
Low	Battery in Use	Battery in Use
High	High Technical Alarm	Mechanical Ventilation Failure
		RT Clock Need Reset
		RT Clock Not Exist
		Keyboard Init Error
		Power System Comm Error
		Power System Comm Stop
		Power Supply Voltage Error
		Power Board High Temp
		Low Battery Voltage!
		System DOWN for battery depletion!
		Breathing Circuit Not Mounted
		Check Flow Sensors
		Ventilator Comm Error
		Ventilator Selftest Error
		Ventilator Hardware Error 01/02/03/04/05/06/07/08/09/10/11/12
		Auxi Ctrl Module Hardware Error 01/02/03/04/05
		Auxi Ctrl Module Comm Error
		Auxi Ctrl Module Comm Stop
		Flowmeter Hardware Error 01/02/03/04/05/06/07
		Flowmeter Cal. Data Error 01/02
		O ₂ -N ₂ O Ratio Error
		Flowmeter Comm Error
Flowmeter Comm Stop		
Device Fault, Ventilate Manually		
Paw < -10cmH ₂ O		
Mediate	Mediate Technical Alarm	Key Error
		IP Address Conflict
		Battery Undetected
		ACGO On
		O ₂ Flush Failure
		PEEP Valve Failure

BeneView		Wato
Priority	Label	Label
		Insp Valve Failure
		PEEP Safety Valve Failure
		Replace O ₂ sensor
		Pressure Monitoring Channel Failure
		Insp Reverse Flow
		Exp Reverse Flow
		TVe Below Control Range
		Ventilator Comm Stop
		Pressure Monitoring Channel Failure
		Volume Monitoring Disabled
		CO ₂ Canister Not Mounted
Low	Low Technical Alarm	Heating Module Failure
		3-way Valve Failure
		Flow Sensor Failure
		Calibrate Flow Sensor
		Calibrate O ₂ Sensor
		Calibrate PEEP Valve
		TV Comp Disabled
		TV Not Achieved
		Flowmeter Zero Failed
		N ₂ O Flow Too High
		O ₂ Flow Too High
		Air Flow Too High
		Pinsp Not Achieved
		TVe > TVi
		TV Delivery Too High
		Sensor Zero Failed
		Ventilator Init Error
Mediate	CO ₂ Module abnormal	CO ₂ Comm Stop
		CO ₂ Comm Error
		CO ₂ Sensor High Temp
		CO ₂ Sensor Low Temp
		CO ₂ High Airway Press.
		CO ₂ Low Airway Press.
		CO ₂ High Barometric
		CO ₂ Low Barometric
		CO ₂ Hardware Error
		CO ₂ Sampleline Occluded
		CO ₂ System Error
		CO ₂ No Watertrap
		EtCO ₂ Overrange
		FiCO ₂ Overrange

BeneView		Wato
Priority	Label	Label
		CO ₂ Zero Failed
		CO ₂ Cal. Failed
		CO ₂ Factory Cal. Invalid
		CO ₂ Check Airway
		CO ₂ No Sampleline
		CO ₂ Main Board Error
		CO ₂ Check Sensor or Main Board
		CO ₂ Replace Scrubber&Pump
		CO ₂ Replace Sensor
		CO ₂ 15V Overrange
		CO ₂ Init Error
		CO ₂ Selftest Error
		CO ₂ Temp Overrange
		CO ₂ Overrange
		CO ₂ Check Cal.
		CO ₂ Zero Error
		CO ₂ Sensor Error
		CO ₂ No Sensor
Mediate	AG Module abnormal	AG Hardware Error
		O ₂ Sensor Error
		AG Selftest Error
		AG Hardware Malfunction
		AG Init Error
		AG No Watertrap
		AG Change Watertrap
		AG Comm Stop
		AG Airway Occluded
		AG Comm Error
		AG Data Limit Error
		AG Zero Failed
		AG Cal. Failed
		AG Accuracy Error
		O ₂ Accuracy Unspecified
		N ₂ O Accuracy Unspecified
		CO ₂ Accuracy Unspecified
		Enf Accuracy Unspecified
		Iso Accuracy Unspecified
		Sev Accuracy Unspecified
		Hal Accuracy Unspecified
		Des Accuracy Unspecified
		Mixed anesthetic gas and MAC < 3
		Mixed anesthetic gas and MAC >= 3

BeneView		Wato
Priority	Label	Label
		EtCO ₂ Overage
		FiCO ₂ Overage
		EtN ₂ O Overage
		FiN ₂ O Overage
		EtHal Overage
		FiHal Overage
		EtEnf Overage
		FiEnf Overage
		EtIso Overage
		Filso Overage
		EtSev Overage
		FiSev Overage
		EtDes Overage
		FiDes Overage
		Mediate
BISx Disconnected		
BIS Comm Error		
BIS Overage		
SQI Overage		
SR Overage		
BIS High Imped.		
BIS Sensor Off		
BIS DSC Error		
BIS DSC Malf		
BIS No Cable		
BIS No Sensor		
BIS Wrong Sensor Type		
SQI<50%		
SQI<15%		
BIS Sensor Expired		
BIS Sensor Failure		
BIS Sensor Too Many Uses		
Disconnect/Reconnect BIS		
BIS Selftest Error		

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.8.2 Maquet Flow-i

26.8.2.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	No
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
VTi	Inspired tidal volume	ml	Yes
MV	Minute volume	L/min	Yes
MVe	Expiratory minute volume	L/min	Yes
MVi	Inspiratory minute volume	L/min	Yes
ftot	Total respiratory rate	bpm	Yes
f	Breath rate	bpm	No
I:E	Inspiratory time:Expiratory time ratio	/	No
TIP:TI	Percentage of inspiratory plateau time in inspiratory time	%	No
Rise Time%	rise time%	%	No
Tslope	Time for the pressure to rise to target pressure	s	No
Tinsp	Time of inspiration	s or %	No
Tapnea	Apnea time	s	No
PC above PEEP	PC above PEEP	cmH ₂ O hPa mbar	No
PS above PEEP	PS above PEEP	cmH ₂ O hPa mbar	No
P-Trigger	Inspiratory trigger level(pressure trigger)	cmH ₂ O hPa mbar	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
Insp Flow	Inspiratory flow	L/min	No
Exp Flow	Expiratory flow	L/min	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
Compl	Compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
FiCO ₂	Fraction of inspired carbon dioxide	% mmHg kPa	Yes
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
EtO ₂	End-tidal O ₂	% mmHg kPa	Yes
FiN ₂ O	Fraction of inspired nitrous oxide	%	Yes
EtN ₂ O	End-tidal N ₂ O	%	Yes
FiAA	Inspired anesthetic agent	%	Yes
EtAA	End-tidal anesthetic agent	%	Yes
FiAA 2nd	2nd Insp. Agent	%	Yes
EtAA 2nd	2nd Exp. Agent	%	Yes
MAC	Minimum alveolar concentration	/	Yes
PO ₂	Oxygen supply pressure	kPa	No
PN ₂ O	N ₂ O supply pressure	kPa	No
Pair	Air supply pressure	kPa	No
FG	Fresh gas flow	ml/min	No
Ti/Ttot	Duty cycle or ratio of inspiration time to total breathing cycle time (only during spontaneous breathing)	/	No

26.8.2.2 Output Signals—Alarms

BeneView		SynoVent
Priority	Label	Label
High	Apnea	Apnea
High	Paw Too High	Paw High
High	High Paw Sustained	High continuous pressure
Mediate	MV Too High	MV too high
Mediate	MV Too Low	MV too Low
Mediate	PEEP Too High	PEEP High
Mediate	PEEP Too Low	PEEP Low
Mediate	EtCO ₂ Too High	EtCO ₂ High

BeneView		SynoVent
Priority	Label	Label
Mediate	EtCO ₂ Too Low	EtCO ₂ Low
Mediate	FiCO ₂ Too High	FiCO ₂ High
Mediate	FiN ₂ O Too High	FiN ₂ O High
Mediate	EtIso Too High	EtISO High
Mediate	FiIso Too High	FiISO High
Mediate	FiIso Too Low	FiISO Low
Mediate	EtSev Too High	EtSEV High
Mediate	EtSev Too Low	EtSEV Low
Mediate	FiSev Too High	FiSEV High
Mediate	EtDes Too High	EtDES High
Mediate	EtDes Too Low	EtDES Low
Mediate	EtO ₂ Too High	EtO ₂ High
Mediate	EtO ₂ Too Low	EtO ₂ Low
Mediate	FiO ₂ Too High	FiO ₂ High
Mediate	FiO ₂ Too Low	FiO ₂ Low
Low	RR Too High	frequency high
Low	RR Too Low	frequency low
High	Circuit Occluded	Gas sampling tube Occlusion
High	High Technical Alarm	Mixture of Anesthesia agents
		Gas Supply
		Cross contamination of anesthetic Agents
		Vaporizer liquid level
		battery alarm
		patient Cassette remove
Mediate	Mediate Technical Alarm	patient Cassette exchange
		Gas Analyzer water trap
		Gas Analyzer water trap missing
Low	Battery in Use	internal communicaiton failure
		Battery operation

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.8.3 Draeger Fabius GS/Fabius Trio/Fabius Plus

26.8.3.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	No
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
Paw	Airway pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
VTe	Expiratory tidal volume	ml	Yes
MV	Minute volume	L/min	Yes
f	Breath rate	bpm	No
fspn	Spontaneous respiratory rate	bpm	Yes
I:E	Inspiratory time:Expiratory time ratio	/	No
TIP:TI	Percentage of inspiratory plateau time in inspiratory time	%	No
Tinsp	Time of inspiration	s	No
Pinsp	Pressure control level of inspiration	cmH ₂ O hPa mbar	No
Psupp	Pressure support level	cmH ₂ O hPa mbar	No
Pmax	Maximal breathing pressure	cmH ₂ O Mbar hPa	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
Insp Flow	Inspiration flow	L/min	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
Exp Flow	Expiratory flow	L/min	No
RRCO ₂	Respiratory rate of CO ₂	bpm	Yes
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
FiCO ₂	Fraction of inspired carbon dioxide	% mmHg kPa	Yes
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
FiN ₂ O	Fraction of inspired nitrous oxide	%	Yes
EtN ₂ O	End-tidal N ₂ O	%	Yes
FiDes	Inspired anesthetic agent	%	Yes
FiSev		%	Yes
FiEnf		%	Yes
Filso		%	Yes
FiHal		%	Yes
EtEnf		End-tidal anesthetic agent	%
EtDes	%		Yes
EtIso	%		Yes
EtSev	%		Yes
EtHal	%		Yes
FiAA	Inspired anesthetic agent		%
EtAA	End-tidal anesthetic agent	%	Yes
FiAA 2nd	2nd Insp. Agent	%	Yes
EtAA 2nd	2nd Exp. Agent	%	Yes
Insp. MAC	Inspired minimum alveolar concentration	/	No
Exp. MAC	Expired minimum alveolar concentration	/	No
MAC	Minimum alveolar concentration	/	Yes
ATMP	Barometric pressure	mmHg	No
HALlev	Anesthetic agent consumption	ml	No
ENFlev			
ISOlev			
DESlev			
SEVlev			
VO ₂	Oxygen consumption	ml/min	Yes
VO ₂ /m ²	Oxygen consumption per body surface area	ml/min/m ²	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
VO ₂ /kg	Oxygen consumption per body weight	ml/min/kg	No
VCO ₂	CO ₂ production	ml/min	No
EE	Energy expenditure	kcal/day	No
RQ	Respiratory quotient	/	No
PO ₂	Oxygen supply pressure	kPa	No
PN ₂ O	N ₂ O supply pressure	kPa	No
Pair	Air supply pressure	kPa	No
O ₂ cyl.	Oxygen cylinder pressure	kPa	No
O ₂ cyl.2nd	Secondary oxygen cylinder pressure	kPa	No
N ₂ O cyl.	N ₂ O cylinder pressure	kPa	No
air cyl.	Air cylinder pressure	kPa	No
FG	Fresh gas flow	ml/min	No
N ₂ O Flow	N ₂ O flow	L/min	No
Air Flow	Air flow	L/min	No
O ₂ Flow	O ₂ flow	L/min	No
Des flow	Anesthetic agent flow	ml/h	No
Enf flow			
Iso flow			
Hal flow			
Sev flow			
IBW	Ideal body weight	kg	No
BSA	Body surface area	m ²	No
BIS	Bispectral index	/	Yes
SQI	Signal quality index	/	Yes
SR	Suppression ratio	/	Yes
EMG	Electromyograph	dB	Yes
SEF	Spectral edge frequency	Hz	Yes
TP	Total power	dB	Yes
BC	Burst count	/min	Yes
SpO ₂	Arterial oxygen saturation from pulse oximetry	%	Yes
PR	Pulse rate	bpm	Yes

26.8.3.2 Output Signals—Alarms

BeneView		Fabius GS/Fabius Trio/Fabius Plus
Priority	Label	Label
High	Apnea	APNEA VENT
High	Volume Apnea > 2 min	APNEA VOL
High	Pressure Apnea	APNEA PRES
High	Paw Too High	PAW HIGH
High	Paw Too Low	PAW NEGATIVE
High	FiO ₂ Too Low	% O ₂ LOW
High	CONT PRES	CONT PRES
High	O ₂ Supply Failure	LO O ₂ SUPPLY
High	Check APL Valve	APL VALVE ?
High	No Fresh Gas	NO FRESHGAS
High	High Technical Alarm	VENT ERR
Mediate	FiO ₂ Too High	% O ₂ HIGH
Mediate	MV Too High	MIN VOL HIGH
Mediate	MV Too Low	MIN VOL LOW
Mediate	PEEP Too High	PEEP HIGH
Mediate	PRESS EXP High	PRESS EXP HI
Mediate	Check Expiration-Valve	EXP-VALVE ?
Mediate	Check Fresh Gas Supply	FRESH GAS ?
Mediate	Mediate Technical Alarm	BATTERY LOW
		PRESS ERR
		VOL ERR
Low	PRESSURE LIM	PRESSURE LIM
Low	Low Technical Alarm	SPEAKER FAIL
		POWER FAIL
		CAL % O ₂ ?
		% O ₂ ERR
		TIME LIMITED
		RS232COM ERR
		PORT 1 ERROR
		PORT 2 ERROR
THRESHOLD LO		

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.8.4 Draeger Primus

26.8.4.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
O ₂ %	Oxygen concentration	%	Yes
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	Yes
P _{peak}	Peak pressure	cmH ₂ O hPa mbar	Yes
P _{plat}	Plateau pressure	cmH ₂ O hPa mbar	Yes
P _{mean}	Mean pressure	cmH ₂ O hPa mbar	Yes
P _{aw}	Airway pressure	cmH ₂ O hPa mbar	Yes
V _T	Tidal volume	ml	No
V _{Ti}	Inspired tidal volume	ml	Yes
MV	Minute volume	L/min	Yes
MV _e	Expiratory minute volume	L/min	Yes
MVLEAK	Leakage minute volume	L/min	No
f _{tot}	Total respiratory rate	bpm	Yes
f	Breath rate	bpm	No
f _{mand}	Mandatory breathing frequency	bpm	No
f _{spn}	Spontaneous respiratory rate	bpm	Yes
FreqMIN	Minimum breath frequency	bpm	No
I:E	Inspiratory time:Expiratory time ratio	/	No
TIP:TI	Percentage of inspiratory plateau time in inspiratory time	%	No
T _{slope}	Time for the pressure to rise to target pressure	s	No
T _{insp}	Time of inspiration	s	No
P _{insp}	Pressure control level of inspiration	cmH ₂ O hPa mbar	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
Psupp	Pressure support level	cmH ₂ O hPa mbar	No
Pmax	Maximal breathing pressure	cmH ₂ O Mbar hPa	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
Compl	Compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RRCO ₂	Respiratory rate of CO ₂	bpm	Yes
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
FiCO ₂	Fraction of inspired carbon dioxide	% mmHg kPa	Yes
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
EtO ₂	End-tidal O ₂	% mmHg kPa	Yes
Δ O ₂	Difference between inspiratory and expiratory O ₂	% mmHg kPa	No
Tapnea	Apnea time	s	No
FiN ₂ O	Inspired anesthetic agent	%	Yes
FiIso		%	Yes
FiDes		%	Yes
FiEnf		%	Yes
FiSev		%	Yes
FiHal		%	Yes
EtN ₂ O		End-tidal anesthetic agent	%
EtEnf	%		Yes
EtDes	%		Yes
EtIso	%		Yes
EtSev	%		Yes
			%

BeneView		Unit	Is it saved in the trends?
Label	Description		
EtHal		%	Yes
FiAA	Inspired anesthetic agent	%	Yes
EtAA	End-tidal anesthetic agent	%	Yes
FiAA 2nd	2nd Insp. Agent	%	Yes
EtAA 2nd	2nd Exp. Agent	%	Yes
Insp. MAC	Inspired minimum alveolar concentration	/	No
Exp. MAC	Expired minimum alveolar concentration	/	No
MAC	Minimum alveolar concentration	/	Yes
HALlev	Anesthetic agent consumption	ml	No
ENFlev			
ISOLev			
DESlev			
SEVlev			
VO ₂	Oxygen consumption	ml/min	Yes
FG	Fresh gas flow	ml/min	No
N ₂ O Flow	N ₂ O flow	L/min	No
Air Flow	Air flow	L/min	No
O ₂ Flow	O ₂ flow	L/min	No
SpO ₂	Arterial oxygen saturation from pulse oximetry	%	Yes
PR	Pulse rate	bpm	Yes

26.8.4.2 Output Signals—Alarms

BeneView		Draeger Primus
Priority	Label	Label
High	Apnea	APNEA/APNEA VENT
High	Volume Apnea > 2 min	APNEA VOL
High	Pressure Apnea	APNEA PRES
High	Paw Too High	PAW HIGH
High	Paw Too Low	PAW NEGATIVE
High	FiO ₂ Too Low	% O ₂ LOW
High	CONT PRES	CONT PRES
High	CO ₂ Apnea	APNEA CO ₂
High	No Pulse	NO SPO ₂ PULS
High	PR Too Low	SPO ₂ PULS LO
High	SPO ₂ Too Low	SPO ₂ LOW

BeneView		Draeger Primus
Priority	Label	Label
High	O ₂ Supply Failure	O ₂ SUPPLY ?
High	No Fresh Gas	NO FRESHGAS
High	Circuit Occluded	CIRCLE OCCL
High	VENT DISC	VENT DISC
High	High Technical Alarm	VENT ERR
		INT.TMP.HIGH
		O ₂ CYL.DISCON
		CHK N ₂ O CYL
		NO N ₂ O DELIV
		NO O ₂ DELIV.
		NO AIR DELIV
		FG X-OVER ?
		VENT.UNLOCKD
		AW-TEMP HIGH
		NO N ₂ O
Mediate	FiO ₂ Too High	FI O ₂ HIGH
Mediate	VTe Too Low	TIDAL VOL. ?
Mediate	MV Too High	MIN VOL HIGH
Mediate	MV Too Low	MIN VOL LOW
Mediate	PEEP Too High	PEEP HIGH
Mediate	EtCO ₂ Too High	ET CO ₂ HIGH
Mediate	EtCO ₂ Too Low	ET CO ₂ Low
Mediate	FiCO ₂ Too High	INSP CO ₂ HIGH
Mediate	FiN ₂ O Too High	FI N ₂ O HIGH
Mediate	EtHal Too High	EXP. HAL HIGH
Mediate	FiHal Too High	% HAL HIGH
Mediate	FiHal Too Low	% HAL LOW
Mediate	EtEnf Too High	EXP. ENF HIGH
Mediate	FiEnf Too High	% ENF HIGH
Mediate	FiEnf Too Low	% ENF LOW
Mediate	EtIso Too High	EXP. ISO HIGH
Mediate	FiIso Too High	% ISO HIGH
Mediate	FiIso Too Low	% ISO LOW
Mediate	EtSev Too High	EXP. SEV HIGH
Mediate	FiSev Too High	% SEV HIGH
Mediate	FiSev Too Low	% SEV LOW
Mediate	EtDes Too High	EXP. DES HIGH

BeneView		Draeger Primus
Priority	Label	Label
Mediate	FiDes Too High	% DES HIGH
Mediate	FiDes Too Low	% DES LOW
Mediate	MAC Too Low	MAC LOW?
Mediate	PR Too High	SPO ₂ PULS HI
Mediate	SPO ₂ Too High	SPO ₂ HIGH
Mediate	Patient Circuit Leak	LEAKAGE
Mediate	Check Fresh Gas Supply	FRESH GAS ?
Mediate	Mediate Technical Alarm	POWER FAIL
		BATTERY LOW
		N ₂ O SUPPLY ?
		PRESSURE LIM
		MIXER INOP
		P MAX?
		SAFETY O ₂ ON
		FG.FLOW LIM.
		LOSS OF DATA
		HOSES MIXED?
		WRONG HOSES?
		% O ₂ ERR
		SET.CANCELED
		FG TOO HIGH
		FG ACTIVE
		FG AIR SENS?
FG O ₂ SENS?		
FG N ₂ O SENS?		
ABS. PRESENT?		
Mediate	AG Module abnormal	WATERTR. OLD?
		MIXED AGENT
		CO ₂ /AGT ERR
		N ₂ O ERR
		AGT ERR
		2nd AGENT
		FICO ₂ OFF
		CO ₂ LINE BLK
		CO ₂ ALRM OFF
Low	NO AIR	NO AIR
Low	NO O ₂ SUPPLY	NO O ₂ SUPPLY

BeneView		Draeger Primus
Priority	Label	Label
Low	Low Technical Alarm	FAN ERR
		PWR SPLY ERR
		PRESS ERR
		VOL ERR
		LO O ₂ SUPPLY
		CHK O ₂ CYL
		ID-FUNC-INOP
		HOSE OLD?
		HOSE MISSING
		COM VENT ERR
		APOLLO COM1?
		APOLLO COM2?
		O ₂ CYL OPEN
		N ₂ O CYL OPEN
		AIR CYL OPEN
		N ₂ O CYL.SENS?
		AIR CYL.SENS?
		O ₂ CYL.SENS?
		AIR CYL.?
		PRESS RELIEF
ABSORB. OLD?		
INSP VOL ERR		
Low	SpO ₂ Module abnormal	SPO ₂ SEN DISC
		SPO ₂ ALRM OF
		SPO ₂ ERR

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.8.5 GE Datex-Ohmeda Aestiva 7900/Aestiva 7100

26.8.5.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
VTe	Expiratory tidal volume	ml	Yes
MVe	Expiratory minute volume	L/min	Yes
O ₂ %	Oxygen concentration	%	Yes
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
Pmin	Minimum airway pressure	cmH ₂ O mbar hPa	No
VT	Tidal volume	ml	No
f	Breath rate	bpm	No
I:E	Percentage of inspiratory plateau time in inspiratory time	%	No
TIP:TI	Percentage of inspiratory plateau time in inspiratory time	%	No
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	No
Plimit	Pressure limit level	cmH ₂ O mbar hPa	No
Pinsp	Pressure control level of inspiration	cmH ₂ O mbar hPa	No

26.8.5.2 Output Signals—Alarms

BeneView		Datex-Ohmeda Aestiva 7900/Aestiva 7100
Priority	Label	Label
High	FiO ₂ Too Low	Low O ₂
High	Paw Too High	High Paw
High	Paw Too Low	Low Paw
High	High Paw Sustained	Sustained Paw (shutdown)
High	Volume Apnea > 2 min	Volume Apnea > 2 min
High	O ₂ Supply Failure	No O ₂ Pressure
High	No Fresh Gas	No Fresh Gas Flow
High	High Technical Alarm	Pinspired Not Achieved
		Inspiration Stopped
		+15V SIB Out-of-Range
		+15V Manifold Out-of-Range
		Display Voltage Out-of-Range
		Vaux_ref Out-of-Range
		Vext_ref Out-of-Range
		A/D Converter Failure
		CPU Failure
		Memory (EEPROM) Failure
		Memory (flash) Failure
		Memory (RAM) Failure
		Memory (video) Failure
		Bootup Memory Failure
		Software Watchdog Failure
		Hardware Watchdog Failure
		Internal Clock Too Fast
Internal Clock Too Slow		
CPU Internal Error		
Control Settings Input Has Failed		
Mediate	FiO ₂ Too High	High O ₂
Mediate	Sub-Atmospheric Paw	Sub-Atmospheric Paw
Mediate	MV Too Low	Low VE
Mediate	MV Too High	High VE
Mediate	VTe Too Low	Low Vte
Mediate	VTe Too High	High Vte
Mediate	Volume Apnea	Volume Apnea

BeneView		Datex-Ohmeda Aestiva 7900/Aestiva 7100
Priority	Label	Label
Mediate	Mediate Technical Alarm	No Pressure Mode/PEEP
		Inspiratory Overshoot
		Manifold Pressure Sensor Failure
		High Pressure Limit Reached (min sys)
		Inspiratory Reverse Flow
		Expiratory Reverse Flow
		Check Flow Sensors
		Flow Valve Failure
		Gas Inlet Valve Failure
		Bootup Gas Inlet Valve Failure
		Memory (redundant storage) Fail
		No Battery
		Low Battery Charge
		Low VE Limit Set
Low	Pressure Limiting	Sustained Paw
Low	Battery in Use	On Battery
Low	Low Technical Alarm	Check O ₂ Sensor
		O ₂ Calibration Error
		PEEP Not Achieved
		Vt Not Achieved
		No Inspiratory Flow Sensor
		No Expiratory Flow Sensor
		Insp Vt/Vte Mismatch
		Vdel Mismatch
		Bellows Empty
		'+V analog Failure
		'-V analog Failure
		Flow Sensor Cal Data Corrupt
		Low Battery
		Low Battery (shutdown)
		Battery Voltage Out Of Range
		Battery Current Out Of Range
		Circuit Auxiliary
		Auxiliary Breathing Circuit
Service Calibrations Due		

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.8.6 GE Datex-Ohmeda Avance/Aisys

26.8.6.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
Vte	Expiratory tidal volume	ml	Yes
MVe	Expiratory minute volume	L/min	Yes
ftot	Total respiratory rate	bpm	Yes
O ₂ %	Oxygen concentration	%	Yes
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
Pmin	Minimum airway pressure	cmH ₂ O mbar hPa	No
MVspn	Spontaneous breathed minute volume	L/min	Yes
fspn	Spontaneous respiratory rate	bpm	Yes
PEEPi	Intrinsic positive end-expiratory pressure	cmH ₂ O hPa mbar	No
Compl	Compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RAW	Airway resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
VTi	Inspired tidal volume	ml	Yes
MVi	Inspiratory minute volume	L/min	Yes
Paux Peak	Peak auxiliary pressure	cmH ₂ O hPa mbar	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
Paux Mean	Mean auxiliary pressure	cmH ₂ O hPa mbar	No
Paux Min	Minimum auxiliary pressure	cmH ₂ O hPa mbar	No
PEEPe	Extrinsic positive end-expiratory pressure	cmH ₂ O hPa mbar	No
PEEPtot	Total PEEP	cmH ₂ O hPa mbar	No
PEEPi time	Intrinsic PEEP age (elapsed time since last maneuver)	min	No
P0.1	100 ms occlusion pressure	cmH ₂ O hPa mbar	No
P0.1 time	P0.1 age (elapsed time since last maneuver)	min	No
ATMP	Barometric pressure	mmHg	No
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
EtO ₂	End-tidal O ₂	% mmHg kPa	Yes
Δ O ₂	Difference between inspiratory and expiratory O ₂	% mmHg kPa	No
FiCO ₂	Fraction of inspired carbon dioxide	% mmHg kPa	Yes
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
RRCO ₂	Respiratory rate of CO ₂	bpm	Yes
FiAA	Inspired anesthetic agent	%	Yes
EtAA	End-tidal anesthetic agent	%	Yes
FiAA 2nd	2nd Insp. Agent	%	Yes
EtAA 2nd	2nd Exp. Agent	%	Yes

BeneView		Unit	Is it saved in the trends?
Label	Description		
FiN ₂ O	Fraction of inspired nitrous oxide	%	Yes
EtN ₂ O	End-tidal N ₂ O	%	Yes
MAC	Minimum alveolar concentration	/	Yes
VO ₂	Oxygen consumption	ml/min	Yes
VO ₂ /m ²	Oxygen consumption per body surface area	ml/min/m ²	No
VO ₂ /kg	Oxygen consumption per body weight	ml/min/kg	No
VCO ₂	CO ₂ production	ml/min	No
EE	Energy expenditure	kcal/day	No
RQ	Respiratory quotient	/	No
PO ₂	oxygen supply pressure	kPa	No
PN ₂ O	N ₂ O supply pressure	kPa	No
Pair	air supply pressure	kPa	No
O ₂ cyl.	Oxygen cylinder pressure	kPa	No
O ₂ cyl.2nd	Secondary oxygen cylinder pressure	kPa	No
N ₂ O cyl.	N ₂ O cylinder pressure	kPa	No
air cyl.	Air cylinder pressure	kPa	No
Des flow	Anesthetic agent flow	ml/h	No
Enf flow			
Iso flow			
Hal flow			
Sev flow			
O ₂ Flow	O ₂ flow	L/min	No
N ₂ O Flow	N ₂ O flow	L/min	No
Air Flow	Air flow	L/min	No
T _{insp}	Time of inspiration	s	No
T _{exp}	Expiratory time	s	No
I:E	Inspiratory time:Expiratory time ratio	/	No
FRC	Fractional residual capacity	ml	No
VT	Tidal volume	ml	No
f	Breath rate	bpm	No
TIP:TI	Percentage of inspiratory plateau time in inspiratory time	%	No
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	No
Plimit	Pressure limit level	cmH ₂ O hPa mbar	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
Pinsp	Pressure control level of inspiration	cmH ₂ O hPa mbar	No
Psupp	Pressure support level	cmH ₂ O hPa mbar	No
Pmax	Maximal breathing pressure	cmH ₂ O Mbar hPa	No
Tapnea	Apnea time	s	No
IBW	Ideal body weight	Kg	No
BSA	Body surface area	m ²	No
Rise Time%	rise time%	%	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
P-Trigger	Inspiratory trigger level (pressure trigger)	cmH ₂ O hPa mbar	No
Tinsp	Time of inspiration	s or %	No
Tpause	Apnea Time	s or %	No

26.8.6.2 Output Signals—Alarms

BeneView		GE Datex-Ohmeda Avance/Aisys
Priority	Label	Label
High	Paw Too High	High Paw
High	Paw Too Low	Low Paw
High	High Paw Sustained	High Paw Sustained
High	Circuit Occluded	Circuit Occluded
High	Volume Apnea > 2 min	Volume Apnea > 2 min
High	O ₂ Supply Failure	No O ₂ Pressure
High	No Fresh Gas	No Fresh Gas Flow
High	EtO ₂ Too Low	Low etO ₂
High	EtO ₂ Too High	High etO ₂
High	FiO ₂ Too Low	Low FiO ₂
High	FiO ₂ Too High	High FiO ₂
High	CO ₂ Apnea	CO ₂ Apnea
High	High Technical Alarm	Pmax Reached

BeneView		GE Datex-Ohmeda Avance/Aisys
Priority	Label	Label
		Pinspired Not Achieved
		Other Priority Alarms (for high priority alarms not assigned a unique bit)
		No VO ₂ , High FiN ₂ O
		Low Drive Gas Pressure
		Low Battery Charge
		Low Battery (No AC)
		Control Settings Failure
		Standby ON (set when anesthesia system is not in therapy mode or when respiratory care ventilator is in standby)
		Therapy Computer Failure
		Monitoring Computer Failure
		Display Computer Failure
		System Error
		Mixer Failure
		Mixer Leak
		Mixer Control Failure
		Vent Failure
		Mechanical Ventilation Disabled
		Patient Detected (while in standby)
		High O ₂ Supply Pressure
		High Air Supply Pressure
Mediate	Sub-Atmospheric Paw	Sub-Atmospheric Paw
Mediate	MV Too Low	Low VE
Mediate	MV Too High	High VE
Mediate	VTe Too Low	Low Vte
Mediate	VTe Too High	High Vte
Mediate	Volume Apnea	Volume Apnea
Mediate	Patient Circuit Leak	Patient Circuit Leak
Mediate	RR Too High	Low RR
Mediate	RR Too Low	High RR
Mediate	EtCO ₂ Too Low	Low etCO ₂
Mediate	EtCO ₂ Too High	High etCO ₂
Mediate	FiCO ₂ Too High	High FiCO ₂
Mediate	EtAA Too Low	Low etAA
Mediate	EtAA Too High	High et AA

BeneView		GE Datex-Ohmeda Avance/Aisys
Priority	Label	Label
Mediate	FiAA Too Low	Low FiAA
Mediate	FiAA Too High	High FiAA
Mediate	AG Module abnormal	MGAS ANE_WARMING_UP (5-minute warming up)
		MGAS WARMING_UP (2-minute warming up)
		No VO ₂ , FiO ₂ > 85%
		Alternate O ₂ ON
		Air Only Mode
		MGAS Failure
		MGAS Outlet Occluded
		MGAS Filter Blocked
		MGAS Sample Line Blocked
		MGAS No Sample Line
		MGAS Replace Water Trap
		Module Not Compatible
		Vaporizer Cassette Failure
		Vaporizer Cassette Agent Level Low
		No Vaporizer Cassette
		Vaporizer Failure
		Vaporizer Leak
		AA Control Failure
AA Delivery Disabled		
Nebulizer Failure		
No Nebulizer		
Mediate	Mediate Technical Alarm	High Circuit O ₂
		Low Circuit O ₂
		No O ₂ Cell Sensor
		No Pressure Cntrl/PEEP
		Inspiration Stopped
		Inspiratory Reverse Flow
		Expiratory Reverse Flow
		Check Flow Sensors
		No Air Pressure
		No VO ₂ , Artifact
		No VO ₂ , High Bypass Flow
		No Battery

BeneView		GE Datex-Ohmeda Avance/Aisys
Priority	Label	Label
		Battery Failure
		Battery Charger Failure
		Non Circle Circuit Selected
		Expiratory Flow Sensed with Non Circle Circuit
		Verify Low VE Limit
		Fan Failure
		Heater Failure
		Power Supply Failure
		Display Failure
		Breathing System Failure
		Sensor Interface Board Failure
		ACGO Failure
		SCGO Failure
		Primary Audio Failure
		Backup Audio Failure
Low	Pressure Limiting	Sustained Paw
Low	PRESSURE LIM	Plimit Reached
Low	Battery in Use	Running On Battery (No AC)
Low	Low Technical Alarm	ASR on
		Replace O ₂ Cell
		O ₂ Cell Calibration Error
		PEEP Not Achieved
		Vt Not Achieved
		No Inspiratory Flow Sensor
		No Expiratory Flow Sensor
		Insp Vt/Vte Mismatch (VTE > Insp VT)
		Vdel Mismatch (System Leak)
		Bellows Empty
		No N ₂ O Pressure
		Memory (EEPROM) Failure
		Flow Sensor Cal Data Corrupt
Service Calibrations Due		

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.9 Integrating Ventilator

26.9.1 Newport E360

26.9.1.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	Yes
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
VTe	Expiratory tidal volume	ml	Yes
VTi	Inspiratory tidal volume	ml	Yes
MVspn	Spontaneous breathed minute volume	L/min	Yes
MVe	Expiratory minute volume	L/min	Yes
MVi	Inspiratory minute volume	L/min	Yes
ftot	Total respiratory rate	bpm	Yes
fspn	Spontaneous respiratory rate	bpm	Yes
f	Breath rate	bpm	No
I:E	Inspiratory time: Expiratory time ratio	/	No
Leak Comp	Leak compensation	%	No
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
Rstat	Static lung resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Rdyn	Dynamic lung resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes

BeneView		Unit	Is it saved in the trends?
Label	Description		
Cstat	Static compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
Cdyn	Dynamic compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RSBI	Rapid shallow breathing index	1/(min·L)	Yes
WOBimp	Imposed work of breathing	J/min	Yes
O ₂ Flow	O ₂ flow	L/min	No
Air Flow	Air flow	L/min	No
Insp.Flow	Inspiration flow	L/min	No
Exp. Flow	Expiratory flow	L/min	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
P-Trigger	Inspiratory trigger level (pressure trigger)	cmH ₂ O Mbar hPa	No
Psupp	Pressure support level	cmH ₂ O Mbar hPa	No
Plimit	Pressure limit level	cmH ₂ O mbar hPa	No
Tinsp	Time of inspiration	s	No
Pmax	Maximal breathing pressure	cmH ₂ O Mbar hPa	No
PEEP/CPAP	PEEP/CPAP	cmH ₂ O mbar hPa	No
PEEPtot	Total PEEP	cmH ₂ O hPa mbar	No

26.9.1.2 Output Signals—Alarms

BeneView		Newport E360
Priority	Label	Label
High	Paw Too High	High Paw
High	Paw Too Low	Low Paw
High	MV Too High	High Exhale MV
High	MV Too Low	Low Exhale MV
High	Apnea	Apnea Alarm
High	FiO ₂ Too High	FiO ₂ High
High	FiO ₂ Too Low	FiO ₂ Low
High	VT Not Achieved	Volume Target Not Met
High	Low Baseline	Low Baseline
High	High Baseline	High Baseline
High	Sustained Hbline	Sustained Hbline
High	Air Supply Pressure Low	Air Supply Loss
High	O ₂ Supply Pressure Low	O ₂ Supply Loss
High	Check Flow Sensors	Flow Sensor Error
High	Patient Disconnected	Patient Disconnect
High	Power Failure	Power Failure
High	Tinsp too Short	Insp Time too Short
Mediate	RR Too High	Resp. Rate Alarm
Mediate	O ₂ and air supply	Air & O ₂ Supply Loss
Mediate	O ₂ Sensor Unconnected	FiO ₂ Sensor Disconnected
Low	Battery in Use	Battery in Use
Low	Tinsp too Long	Insp Time too Long
High	High Technical Alarm	Device Alert
		No O ₂ Power-Up
		Control EEPROM Failure
		Low Battery
		Transducer Error
		Control RAM Failed
		Control ROM Failed
		Control CPU Failed
		Monitor RAM Failed
		Monitor ROM Failed
		Monitor CPU Failed
		Dual RAM Failed
		Monitor Tasks Failed
		Control Processor Failed
		Mon Internal System Failed
		Control Tasks Failed
Monitor Processor Failed		
Ctrol Internal System Failed		

BeneView		Newport E360
Priority	Label	Label
		Fan Failure
		Air Flow Sensor EEPROM Failure
		O ₂ Flow Sensor EEPROM Failure
		Air Servo Valve Leak
		O ₂ Servo Valve Leak
Mediate	Mediate Technical Alarm	Flow Sensor Cal Failed
		FiO ₂ Sensor Bad
		O ₂ Sensor Cal Failed
		External Battery
		Check Flow Sensor Board
		NO TEST
Low	Low Technical Alarm	I:E Ratio Inverse violation
		Plimit<Pbase
		Psupport+Pbase>60cmH ₂ O
		Pbase>Low Paw
		Tidal Volume Out of Range
		Flow Out of Range
		Ti Out of Range
		Rate Out of Range
		Psupport Out of Range
		Plimit Out of Range
		PEEP/CPAP Out of Range
		Flow Trigger Out of Range
		CPM Blinking
		EXH. VALVE CAL. Failed: Prox < 1
		EXH. VALVE CAL. Failed: Prox > 0.5
		EXH. VALVE CAL. Failed: Prox Low
EXH. VALVE CAL. Failed: Flow < 1		
LEAK TEST Leak Test Failed		

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.9.2 Puritan Bennett 840

26.9.2.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
O ₂ %	Oxygen concentration	%	Yes
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	Yes
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
Paw	Airway pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
VT _e	Expiratory tidal volume	ml	Yes
VT _i	Inspiratory tidal volume	ml	Yes
VT _e spn	Spontaneous expiratory tidal volume	ml	Yes
VT _{apnea}	Apnea tidal volume	ml	No
MV _{spn}	Spontaneous breathed minute volume	L/min	Yes
MV _e	Expiratory minute volume	L/min	Yes
ftot	Total respiratory rate	bpm	Yes
f _{apnea}	Breath rate for apnea ventilation	bpm	No
f	Breath rate	bpm	No
I:E	Inspiratory time: Expiratory time ratio	/	No
MVLEAK	Leakage minute volume	L/min	No
Leak Comp	Leak compensation	%	No
Rstat	Static lung resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Rdyn	Dynamic lung resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Cstat	Static compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes

BeneView		Unit	Is it saved in the trends?
Label	Description		
Cdyn	Dynamic compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RSBI	Rapid shallow breathing index	1/(min·L)	Yes
WOB	Work of breathing	J/L	Yes
Base Flow	Base Flow	L/min	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
P-Trigger	Inspiratory trigger level (pressure trigger)	cmH ₂ O Mbar hPa	No
Psupp	Pressure support level	cmH ₂ O Mbar hPa	No
Tplat	Plateau time	s	No
Rise Time%	Rise time	%	No
PEEP/CPAP	PEEP/CPAP	cmH ₂ O Mbar hPa	No
NIF	Negative inspiratory force	cmH ₂ O hPa mbar	No
P0.1	100 ms occlusion pressure	cmH ₂ O hPa mbar	No
PEEPi	Intrinsic positive end-expiratory pressure	cmH ₂ O hPa mbar	No
PEEPtot	Total PEEP	cmH ₂ O hPa mbar	No
Peak Flow	Peak flow	L/min	No
Tapnea	Apnea interval	s	No
IBW	Ideal body weight	kg	No
Ti max	Maximum inspiration time	s	No
Tube ID	Tube ID	mm	No

26.9.2.2 Output Signals—Alarms

BeneView		Puritan Bennett 840
Priority	Label	Label
High	Paw Too High	High Inspiratory Pressure
High	MV Too High	High Exhaled minute Volume
High	MV Too Low	low exhaled minute volume
High	Apnea	Apnea
High	FiO ₂ Too Low	Low O ₂ %
High	Ppeak Too Low	Low Ppeak
High	Air Supply Pressure Low	No Air Supply
High	O ₂ Supply Pressure Low	No O ₂ Supply
High	Airway Obstructed?	Severe Occlusion
High	Patient Disconnected	Circuit Disconnect
High	Power Failure	Loss of Power
Mediate	VTe Too High	High Exhaled Tidal Volume
Mediate	RR Too High	High frot
Mediate	VTe Too Low	Low Exhaled Mandatory Tidal Volume Alarm
Mediate	EtO ₂ Too High	High O ₂ Percent
High	High Technical Alarm	Compressor Inoperative
		Compliance Limited VT
		Procedure Error
		PAV Startup Too Long
		PAV R&C Not Assessed
		Volume Not Delivered
		Volume Not Delivered
Low	Low Technical Alarm	Inoperative Battery
		AC Power Loss
		Low Battery
Low	Tinsp too Long	Inspiration Too Long

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.9.3 Maquet Servo-i/Servo-s

26.9.3.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
O ₂ %	Oxygen concentration	%	Yes
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	Yes
P _{peak}	Peak pressure	cmH ₂ O hPa mbar	Yes
P _{plat}	Plateau pressure	cmH ₂ O hPa mbar	Yes
P _{mean}	Mean pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
V _{Te}	Expiratory tidal volume	ml	Yes
V _{Ti}	Inspired tidal volume	ml	Yes
MV	Minute volume	L/min	Yes
MV _{spn}	Spontaneous breathed minute volume	L/min	Yes
MV _e	Expiratory minute volume	L/min	Yes
MV _i	Inspiratory minute volume	L/min	Yes
f _{tot}	Total respiratory rate	bpm	Yes
f _{spn}	Spontaneous respiratory rate	bpm	Yes
f _{CMV}	CMV frequency	bpm	No
f _{SIMV}	Frequency of SIMV	bpm	No
f	Breath rate	bpm	No
I:E	Inspiratory time:Expiratory time ratio	/	No
Leak Comp	Leak compensation	%	No
C _{stat}	Static compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
C _{dyn}	Dynamic compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RSBI	Rapid shallow breathing index	1/(min·L)	Yes
WOB	Work of breathing	J/L	Yes
Exp. Flow	Expiratory flow	L/min	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
P-Trigger	Inspiratory trigger level (pressure trigger)	cmH ₂ O Mbar hPa	No
T _{insp}	Time of inspiration	s	No
T _{pause}	Apnea Time	s or %	No
Rise Time%	rise time%	%	No
P _{high}	Upper pressure level	cmH ₂ O mbar hPa	No
P _{low}	Lower pressure level	cmH ₂ O mbar hPa	No
T _{high}	Time for the upper pressure level	s	No
TPEEP	Time at PEEP level in Bi-Vent	s	No
Exp%	Inspiration termination level	%	No
PC above PEEP	PC above PEEP	cmH ₂ O mbar hPa	No
PS above PEEP	PS above PEEP	cmH ₂ O mbar hPa	No
PEEP/CPAP	PEEP/CPAP	cmH ₂ O mbar hPa	No
R _i	Inspiratory resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
R _e	Expiratory resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
PO ₂	oxygen supply pressure	kPa	No
P _{air}	air supply pressure	kPa	No
P0.1	100 ms occlusion pressure	cmH ₂ O hPa mbar	No
PEEP _{tot}	Total PEEP	cmH ₂ O hPa mbar	No
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
VCO ₂	CO ₂ production	ml/min	No
VT _{CO₂}	CO ₂ tidal elimination	ml	No

26.9.3.2 Output Signals—Alarms

BeneView		Puritan Bennett 840
Priority	Label	Label
High	Paw Too High	Airway pressure alarm Upper pressure limit exceeded
High	MV Too High	Exp.Minute volume too high
High	MV Too Low	Exp.Minute volume too low
High	Apnea	Apnea alarm
High	FiO ₂ Too High	O ₂ conc.too high
High	FiO ₂ Too Low	O ₂ conc.too low
High	PEEP Too Low	PEEP Low
High	No Gas Supply Pressure	Gas supply alarm
High	O ₂ cell disconnect	O ₂ cell disconnect
Mediate	EtCO ₂ Too High	EtCO ₂ conc.too high
Mediate	EtCO ₂ Too Low	EtCO ₂ conc.too low
Mediate	RR Too Low	Breath frequency Low
Mediate	RR Too High	Breath frequency High
Mediate	PEEP Too High	PEEP High
Low	Check tubing	Check tubing
High	High Technical Alarm	Breathing system uP Module error
		Inspiratory control uP Module error
		Monitoring System uP Module error
		Battery alarm
		Power Failure
		Mains Failure
		O ₂ potentiometer error
		CMV potentiometer error
		Range Switch error
		Mode Switch error
		Barometer error
		High continuous pressure
		Overrange
		Computer Interface Emulator hardware error
		NIV,Leakage out of range
NIV,Time in waiting position exceeds 2 min regulation pressure limited		
Mediate	Mediate Technical Alarm	Panel Interface uP Module error
		Exp.flow & CO ₂ linearization uP Module error
		Alarm buff
		CI Battery Voltage
		Pneumatic-Edi out of synch
		Edi activity low
		No Edi signal detected
Unsuccessful manual gas change alarm		

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.9.4 Draeger Evita 2

26.9.4.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
O ₂ %	Oxygen concentration	%	Yes
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	Yes
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
VTe	Expiratory tidal volume	ml	Yes
MV	Minute volume	L/min	Yes
MVspn	Spontaneous breathed minute volume	L/min	Yes
ftot	Total respiratory rate	bpm	Yes
fspn	Spontaneous breathing frequency	bpm	Yes
fSIMV	Frequency of SIMV	bpm	No
I:E	Inspiratory time:Expiratory time ratio	/	No
△int.PEEP	Intermittent PEEP	cmH ₂ O hPa mbar	No
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
Rdyn	Dynamic lung resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Cdyn	Dynamic compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes

BeneView		Unit	Is it saved in the trends?
Label	Description		
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
P-Trigger	Inspiratory trigger level (pressure trigger)	cmH ₂ O Mbar hPa	No
Phigh	Upper pressure level	cmH ₂ O mbar hPa	No
Plow	Lower pressure level	cmH ₂ O mbar hPa	No
Thigh	Time for the upper pressure level	s	No
Tlow	Time for the lower pressure level	s	No
Pmax	Maximum airway rressure	cmH ₂ O mbar hPa	No
Pmin	Minimum airway rressure	cmH ₂ O Mbar hPa	No
Vtrap	Trapped volume	ml	No
T	Inspiratory breathing gas temperature	°C °F	No
P0.1	100 ms occlusion pressure	cmH ₂ O hPa mbar	No
PEEPi	Intrinsic positive end-expiratory pressure	cmH ₂ O hPa mbar	No
EtCO ₂	End-tidal carbon dioxide	% kPa mmHg	Yes
Flow	Flow	L/min	No
Tapnea	Apnea Time	s	No
ASB ramp	ASB ramp	s	No
PASB	Assisted spontaneous breathing	cmH ₂ O hPa mbar	No
Vds	Dead space	ml	No
VCO ₂	CO ₂ production	ml/min	No

26.9.4.2 Output Signals—Alarms

BeneView		Puritan Bennett 840
Priority	Label	Label
High	Paw Too High	PAW HIGH
High	Paw Too Low	PAW LOW
High	MV Too High	MIN VOL HIGH
High	MV Too Low	MIN VOL LOW
High	Apnea	APNEA EVITA
High	FiO ₂ Too High	% O ₂ HIGH
High	FiO ₂ Too Low	% O ₂ LOW
High	AW-TEMP HIGH	AW-TEMP HIGH
High	PEEP Too High	PEEP HIGH
High	ASB>4s	ASB > 4 SEC
High	Air Supply Pressure Low	AIR SUPPLY ?
High	Check Flow Sensors	FLOW SENSOR?
High	EXP-VALVE?	EXP-VALVE ?
High	CLEAN CO ₂	CLEAN CO ₂
Mediate	EtCO ₂ Too High	ET CO ₂ HIGH
Mediate	EtCO ₂ Too Low	ET CO ₂ LOW
Mediate	VOL INCONST	VOL INCONST
Mediate	RR Too High	RESP RATE HI
High	High Technical Alarm	VOL ERR
		PRESS ERR
		AW-TEMP INOP
		AW-TEMP SENS
		CO ₂ NOT CAL
		% O ₂ ERR
		EVITA ERR
		COOLING INOP
Low	Low Technical Alarm	CYCLE FAILED
		CO ₂ ERR
		CO ₂ SENS?
		MIXER INOP
		SYNCHRO INOP

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.9.5 Draeger Evita 4/ Evita2 dura /Evita XL

26.9.5.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
O ₂ %	Oxygen concentration	%	Yes
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	Yes
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
VT	Tidal volume	ml	No
Vte	Expiratory tidal volume	ml	Yes
Vtapnea	Apnea tidal volume	ml	No
MV	Minute volume	L/min	Yes
MVspn	Spontaneous breathed minute volume	L/min	Yes
Ftot	Total respiratory rate	bpm	Yes
Fspn	Spontaneous breathing frequency	bpm	Yes
Fapnea	Breath rate for apnea ventilation	bpm	No
F	Breath rate	bpm	No
I:E	Inspiratory time: Expiratory time ratio	/	No
△int.PEEP	Intermittent PEEP	cmH ₂ O hPa mbar	No
FiO ₂	Fractional concentration of O ₂ in inspired gas	% mmHg kPa	Yes
Rdyn	Dynamic lung resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Cdyn	Dynamic compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RSBI	Rapid shallow breathing index	1/(min·L)	Yes
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
Tinsp	Time of inspiration	s	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
Pinsp	Pressure control level of inspiration	cmH ₂ O mbar hPa	No
Papnea	Apnea pressure	cmH ₂ O mbar hPa	No
Phigh	Upper pressure level	cmH ₂ O mbar hPa	No
Plow	Lower pressure level	cmH ₂ O mbar hPa	No
Thigh	Time for the upper pressure level	s	No
Tlow	Time for the lower pressure level	s	No
Pmax	Maximum airway pressure	cmH ₂ O mbar hPa	No
Pmin	Minimum airway pressure	cmH ₂ O mbar hPa	No
Vtrap	Trapped Volume	ml	No
T	Inspiratory breathing gas temperature	°C °F	No
NIF	Negative inspiratory force	cmH ₂ O hPa mbar	No
P0.1	100 ms occlusion pressure	cmH ₂ O hPa mbar	No
PEEPi	Intrinsic positive end-expiratory pressure	cmH ₂ O hPa mbar	No
EtCO ₂	End-tidal carbon dioxide	% kPa mmHg	Yes
Flow	Flow	L/min	No
Ext.Flow	External flow	L/min	No
Tapnea	Apnea time	s	No
ASB ramp	ASB ramp	s	No
PASB	Assisted spontaneous breathing	cmH ₂ O hPa mbar	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
FlowAssist	Flow assist	mbar.s/L cmH ₂ O.s/L hPa.s/L	No
Vol.Assist	Volume assist	mbar/L cmH ₂ O/L hPa/L	No
Tdisconnect	Delay time for "Airway pressure lower alarm limit"	s	No
Vds	Dead space	ml	No
VCO ₂	CO ₂ production	ml/min	No
ATC	Automatic tube Compensation	%	No
Tube ID	Tube ID	mm	No
PR	Pulse rate	bpm	Yes
SpO ₂	Arterial oxygen saturation from pulse oximetry	%	Yes

26.9.5.2 Output Signals—Alarms

BeneView		Puritan Bennett 840
Priority	Label	Label
High	Paw Too High	PAW HIGH
High	Paw Too Low	PAW LOW
High	MV Too High	MIN VOL HIGH
High	MV Too Low	MIN VOL LOW
High	Apnea	APNEA EVITA
High	FiO ₂ Too High	% O ₂ HIGH
High	FiO ₂ Too Low	% O ₂ LOW
High	AW-TEMP HIGH	AW-TEMP HI
High	PEEP Too High	PEEP HIGH
High	ASB>4s	ASB > 4 SEC
High	No Pulse	NO SPO ₂ PULS
High	PR Too Low	SPO ₂ PULS LO
High	SpO ₂ Too Low	SPO ₂ LOW
High	PR Too High	SPO ₂ PULS HI
High	SpO ₂ Too High	SPO ₂ HIGH
High	Air Supply Pressure Low	AIR SUPPLY ?
High	O ₂ Supply Pressure Low	LO O ₂ SUPPLY
High	Airway Obstructed?	TUBE OBSTRUC
High	Check Flow Sensors	FLOW SENSOR?
High	EXP-VALVE?	EXP-VALVE ?
High	CLEAN CO ₂	CLEAN CO ₂
Mediate	VTe Too High	TIDVOL HI

BeneView		Puritan Bennett 840
Priority	Label	Label
Mediate	EtCO ₂ Too High	ET CO ₂ HIGH
Mediate	EtCO ₂ Too Low	ET CO ₂ LOW
Mediate	VOL INCONST	VOL INCONST
Mediate	RR Too High	RESP RATE HI
Low	ASB > 1.5s	ASB > 1,5 SEC
Low	PPS-TI > 1.5s	PPS-TI > 1,5S
Low	ASB > Tinsp	ASB > TINSP
Low	Battery in Use	BATTERY ON
High	High Technical Alarm	VOL ERR
		PRESS ERR
		AW-TEMP INOP
		% O ₂ ERR
		EVITA ERR
		CYCLE FAILED
		N-VOL ERR
		NEO FLOW ?
		CO ₂ ZERO CAL
		SPO ₂ SEN DISC
		SPO ₂ ERR
		BATTERY ERR
		FAN ERR
		AIR PRESS HI
		HI O ₂ SUPPLY
LOSS OF DATA		
REM.PAD-ERR		
PEEP V ERR		
Mediate	Mediate Technical Alarm	BATT. < 2MIN
		CHECK EVITA
		EVITA STDBY
		AMB PRESS ?
		NEBULIZ OFF
		ERR MULTIPCB
Low	Low Technical Alarm	CO ₂ ERR
		CO ₂ SENSOR ?
		MIXER INOP
		SYNCHRO INOP
		INSPHOLD END
		EXSPHOLD END

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.9.6 Hamilton G5

26.9.6.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
f	Breath rate	bpm	No
VT	Tidal volume	ml	No
TPause	Apnea Time	s or %	No
P-Trigger	Inspiratory trigger level (pressure trigger)	cmH ₂ O hPa mbar	No
PEEP/CPAP	PEEP/CPAP	cmH ₂ O mbar hPa	No
Plow	Lower pressure level	cmH ₂ O mbar hPa	No
Psupp	Pressure support level	cmH ₂ O Mbar hPa	No
MV	Minute volume	L/min	Yes
Plimit	Pressure limit level	cmH ₂ O hPa mbar	No
Pinsp	Pressure control level of inspiration	cmH ₂ O hPa mbar	No
Phigh	Upper pressure level	cmH ₂ O mbar hPa	No
F-Trigger	Inspiratory trigger level (flow trigger)	L/min	No
I:E	Inspiratory time:Expiratory time ratio	/	No
Peak Flow	Peak flow	L/min	No
Exp%	Inspiration termination level	%	No
Ramp	Ramp	ms	No
IBW	Ideal body weight	kg	No
%MinVol	Percentage of minute volume to be delivered	%	No
Tlow	Time for the lower pressure level	s	No
Thigh	Time for the upper pressure level	s	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
Ti max	Maximum inspiration time	s	No
Tip	Inspiratory pause time	s	No
tube ID	Tube ID	mm	No
TRC	Tube resistance compensation	/	No
base flow	Base Flow	L/min	No
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
PEEP	Positive end-expiratory pressure	cmH ₂ O hPa mbar	Yes
Pmin	Minimum airway pressure	cmH ₂ O mbar hPa	No
PEEPi	Intrinsic positive end-expiratory pressure	cmH ₂ O hPa mbar	No
P0.1	100 ms occlusion pressure	cmH ₂ O hPa mbar	No
PTP	Pressure time product	cmH ₂ O.s mbar.s hPa.s	No
Insp.Flow	Inspiration flow	L/min	No
Exp. Flow	Expiratory flow	L/min	No
Vti	Inspired tidal volume	ml	Yes
Vte	Expiratory tidal volume	ml	Yes
VTe spn	Spontaneous expiratory tidal volume	ml	Yes
MVspn	Spontaneous breathed minute volume	L/min	Yes
ftot	Total respiratory rate	bpm	Yes
fspn	Spontaneous respiratory rate	bpm	Yes
Texp	Expiratory time	s	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
I:E	Inspiratory time:Expiratory time ratio	/	No
Ri	Inspiratory resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Re	Expiratory resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Cstat	Static compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
RCexp	Expiratory time constant	s	No
RCinsp	Inspiratory time constant	s	No
RSBI	Rapid shallow breathing index	1/(min·L)	Yes
O ₂ %	Oxygen concentration	%	Yes
WOB	Work of breathing	J/L	Yes
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
VCO ₂	CO ₂ production	ml/min	No
PR	Pulse rate	bpm	Yes
SpO ₂	Arterial oxygen saturation from pulse oximetry	%	Yes
fCMV	CMV frequency	bpm	No
fSIMV	Frequency of SIMV	bpm	No
%Tinsp	Time of inspiration	%	No
Tinsp	Time of inspiration	s	No

26.9.6.2 Output Signals—Alarms

BeneView		Hamilton G5
Priority	Label	Label
High	PawToo High	High pressure
High	PawToo Low	Low pressure
High	FiO ₂ Too High	High oxygen
High	FiO ₂ Too Low	Low oxygen
High	Apnea	Apnea
High	SpO ₂ Too Low	SpO ₂ too low

BeneView		Hamilton G5
Priority	Label	Label
High	SpO ₂ Too High	SpO ₂ too high
High	Patient Disconnected	Disconnection Patient or, Disconnection on patient side
High	Air Supply Pressure Low	Air supply failed
High	O ₂ Supply Pressure Low	Oxygen supply failed
High	O ₂ cell disconnect	O ₂ cell missing
High	O ₂ cell cal. Needed	O ₂ cell calibration needed
High	Power Failure	Loss of mains power
High	Check Flow Sensors	Check Flow Sensor type
High	No Gas Supply Pressure	All gas supplies failed
High	Disconnection ventilator side	Disconnection Ventilator or, Disconnection on ventilator side
High	Loss of PEEP	Loss of PEEP
High	MV Too Low	Low minute volume
High	MV Too High	High minute volume
High	High Technical Alarm	Fail to Cycle
		Wrong Flow Sensor type
		O ₂ cell defective
		Disconnction
		Low internal pressure
		High pressure during sigh
		Pressure not released
		Exhalation obstructed
		TF5514:Check loudspeaker
		Internal battery empty
		Ventilator unit connection lost
		Check internal battery
Mediate	RR Too High	High frequency
Mediate	RR Too Low	Low frequency
Mediate	EtCO ₂ Too High	High PetCO ₂
Mediate	EtCO ₂ Too Low	Low PetCO ₂
Mediate	O ₂ and air supply	Oxygen + air supplies failed
Mediate	O ₂ and heliox supply	Oxygen + heliox supplies failed
Mediate	Mediate Technical Alarm	Gas Supply
		High leak
		Low tidal volume
		High tidal volume

BeneView		Hamilton G5
Priority	Label	Label
		Turn the Flow Sensor
		APV init failed
		Internal battery low
		Panel connection lost
		Heliox supply failed
		SpO ₂ : sensor error(left slot)
		SpO ₂ : sensor error(right slot)
		SpO ₂ : no sensor ((left slot)
		SpO ₂ : no sensor ((left slot)
		SpO ₂ : patient disconnected ((left slot)
		SpO ₂ : patient disconnected ((right slot)
		SpO ₂ : light interference ((left slot)
		SpO ₂ : light interference ((right slot)
		SpO ₂ : poor signal ((left slot)
		SpO ₂ : poor signal ((right slot)
		Large change in FiO ₂
		Recruitment maneuver in process
		Brightness test alarm
		AERONEB disconnected
		Cuff disconnection
		Air + heliox supplies failed
		Oxygenation adjustment OFF (no SpO ₂)
		Ventilation adjustment OFF (no PetCO ₂)
		No hemodynamic status available
		High HLI
		MV oszillation
		FiO ₂ oszillation
		PEEP oszillation
		Cuff high pressure
		FiO ₂ set to 100% due to low saturation
Low	Low Technical Alarm	Operator
		General Alarm
		Volume too low for nebulizer
		ASV: Check high pressure limit
		APV: Check high pressure limit
		Pressure low limit reached
		Check %MinVol

BeneView		Hamilton G5
Priority	Label	Label
		Check Body Wt
		ASV: Cannot meet target
		Check PEEP/high pressure limit
		Check PEEP/Pcontrol
		Check PEEP/Psupport
		Check P-ramp
		Check trigger
		Check %TI
		Check pause
		Check I:E
		Check Vt
		Check rate
		Check peak flow
		Check TI
		Check Flow Pattern
		Flow sensor calibration needed
		Expiratory valve calibration needed
		Apnea ventilation ended
		Maximum leak compensation
		Low ExpMinVol alarm off
		CO ₂ sensor calibration needed
		Check CO ₂ airway adapter
		CO ₂ sensor disconnected
		CO ₂ sensor over temperature
		CO ₂ sensor faulty
		External battery empty
		Sensor simulation active
		IRV
		Cuff leak
		IntelliCuff not found
		Check VThigh limit
		AERONEB module disconnected
		Oxygenation adjustment OFF (no SpO ₂)
		Ventilation adjustment OFF (no PetCO ₂)
		Check CO ₂ sampling line
		Check INTELLiVENT PEEP limit setting
		Set low limit of ExpMinVol alarm

BeneView		Hamilton G5
Priority	Label	Label
		Recruitment in Progress
		Oxygenation Controller on Limit
		Ventilation Controller on Limit
		SBT conditions fulfilled
		SBT in progress

NOTE

- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

26.9.7 Hamilton C2 /Galileo

26.9.7.1 Output Signals—Parameters

BeneView		Unit	Is it saved in the trends?
Label	Description		
fCMV	CMV frequency	bpm	No
fSIMV	Frequency of SIMV	bpm	No
VT	Tidal volume	ml	No
%Tinsp	Time of inspiration	%	No
Tpause	Apnea Time	s or %	No
P-Trigger	Inspiratory trigger level(pressure trigger)	cmH ₂ O hPa mbar	No
PEEP/CPAP	PEEP/CPAP	cmH ₂ O mbar hPa	No
Psupp	Pressure support level	cmH ₂ O Mbar hPa	No
O ₂ %	Oxygen concentration	%	Yes
MV	Minute volume	L/min	Yes
I: E	Inspiratory time:Expiratory time ratio	/	No
Peak Flow	Peak flow	L/min	No
Exp%	Inspiration termination level	%	No
Ramp	Ramp	ms	No
IBW	Ideal body weight	kg	No
%MinVol	Percentage of minute volume to be delivered	%	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
EtCO ₂	End-tidal carbon dioxide	% mmHg kPa	Yes
SpO ₂	Arterial oxygen saturation from pulse oximetry	%	Yes
PR	Pulse rate	bpm	Yes
Vti	Inspired tidal volume	ml	Yes
Vte	Expiratory tidal volume	ml	Yes
ftot	Total respiratory rate	bpm	Yes
fspn	Spontaneous respiratory rate	bpm	Yes
Ppeak	Peak pressure	cmH ₂ O hPa mbar	Yes
Pmean	Mean pressure	cmH ₂ O hPa mbar	Yes
Pplat	Plateau pressure	cmH ₂ O hPa mbar	Yes
Texp	Expiratory time	s	No
Ri	Inspiratory resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Re	Expiratory resistance	cmH ₂ O/L/s hPa/L/s mbar/L/s	Yes
Cstat	Static compliance	ml/cmH ₂ O ml/hPa ml/mbar	Yes
Insp.Flow	Inspiration flow	L/min	No
VTe spn	Spontaneous expiratory tidal volume	ml	Yes
PEEPi	Intrinsic positive end-expiratory pressure	cmH ₂ O hPa mbar	No
Pmin	Minimum airway pressure	cmH ₂ O mbar hPa	No
Tinsp	Time of inspiration	s	No

BeneView		Unit	Is it saved in the trends?
Label	Description		
P0.1	100 ms occlusion pressure	cmH ₂ O hPa mbar	No
Exp. Flow	Expiratory flow	L/min	No
RCexp	Expiratory time constant	s	No
RCinsp	Inspiratory time constant	s	No
WOB	Work of breathing	J/L	Yes
PTP	Pressure time product	cmH ₂ O.s mbar.s hPa.s	No
Pinsp	Pressure control level of inspiration	cmH ₂ O mbar hPa	No

26.9.7.2 Output Signals—Alarms

BeneView		Hamilton C2 /Galileo
Priority	Label	Label
High	Paw Too High	High Pressue
High	Patient Disconnected	Disconnection Patient
High	Apnea	Apnea
High	Disconnection ventilator side	Disconnection Ventilator
High	Loss of PEEP	Loss of PEEP
High	MV Too Low	Low Min Vol
High	MV Too High	High Min Vol
High	High Technical Alarm	Fail to Cycle
Mediate	RR Too High	High Rate
Mediate	Mediate Technical Alarm	Gas Supply
Low	Low Technical Alarm	Operator
		General Alarm


NOTE

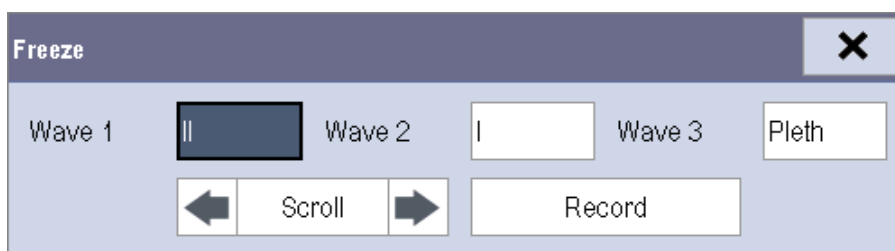
- The above alarm messages are derived from the open protocol of corresponding external device. For more information about these alarms, please see the Instructions for Use matching the device.

27 Freezing Waveforms

During patient monitoring, the freeze feature allows you to freeze the currently displayed waveforms on the screen so that you can have a close examination of the patient's status. Besides, you can select any frozen waveform for recording.

27.1 Freezing Waveforms

1. To freeze waveforms, select the  hardkey on the monitor's front.
2. The system closes the displayed menu (if any), and opens the **[Freeze]** menu.





3. All displayed waveforms are frozen, i.e. the waveforms stop being refreshed or scrolling.

The freeze feature exerts no effect on the split-screen view of minitrends, oxyCRG and other patients.

27.2 Viewing Frozen Waveforms




To view the frozen waveforms, you can either:

- Select the **[Scroll]** button and then rotate the Knob clockwise or counter-clockwise, or
- Directly select the  or  beside the **[Scroll]** button using a mouse or through the touchscreen.

The frozen waveforms will scroll left or right accordingly. And meanwhile, at the lower right corner of the bottommost waveform, there is an upward arrow. The freeze time is displayed below the arrow and the initial frozen time is **[0 s]**. With the waveforms scrolling, the freeze time changes at intervals of 1 second. This change will be applied for all waveforms on the screen.

27.3 Unfreezing Waveforms

To unfreeze the frozen waveforms, you can either:

- Select the  button at the upper right corner of the [**Freeze**] menu,
- Select the  hardkey on the monitor's front, or
- Perform any other action that causes the screen to be readjusted or opens a menu, such as plugging in or out a module, pressing the  hardkey, etc.

27.4 Recording Frozen Waveforms

1. In the [**Freeze**] menu, select, in turn, [**Wave 1**], [**Wave 2**] and [**Wave 3**] and then select your desired waveforms.
2. Select the [**Record**] button. The selected waveforms and all numerics at the frozen time are printed out by the recorder.

28 Review

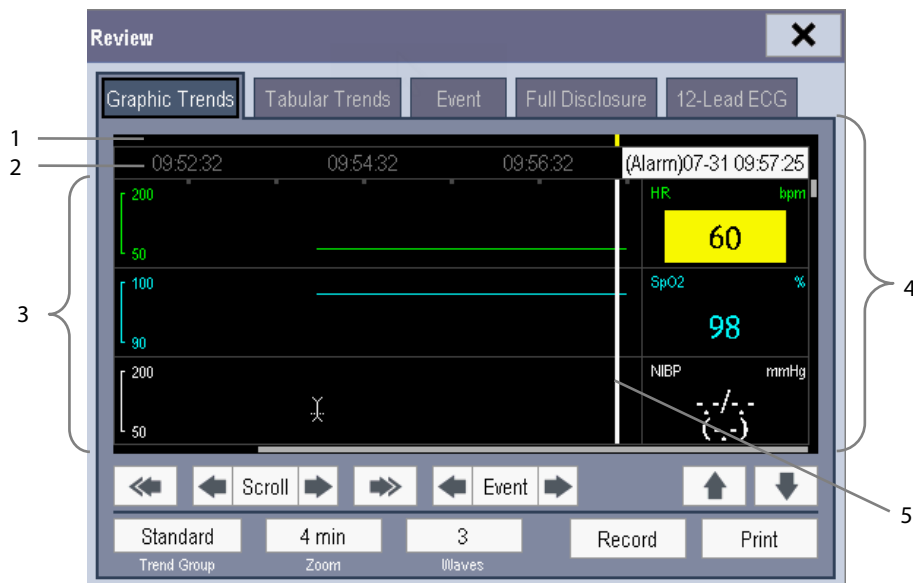
28.1 Accessing Respective Review Windows

1. Select the **[Review]** QuickKey, or **[Main Menu]**→**[Review >>]**.
2. Select **[Graphic Trends]**, **[Tabular Trends]**, **[Events]**, **[Full Disclosure]** or **[12-lead ECG]** to access their respective review windows.

For details about reviewing interpretation of resting 12-lead ECG results, refer to the chapter **8 Monitoring ECG**.

28.2 Reviewing Graphic Trends

In the **[Review]** menu, select **[Graphic Trends]** to access the following window.







1. Event mark area
2. Time axis
3. Graphic trends area
4. Parameter area
5. Cursor



Events are marked with colors in the event mark area. Red represents high level alarm event. Yellow represents medium/low level alarm event. Green represents manual event.

In this review window:

- Select **[Trend Group]** and you can select a trend group for viewing in the popup menu. If **[Custom 1]** or **[Custom 2]** is selected, you can further select **[Define Trend Group]**. Then you can select the parameters for viewing in the popup menu.
- You can set the time length of the review window by selecting **[Zoom]**.
- You can set the number of waves displayed in one page by selecting **[Waves]**.

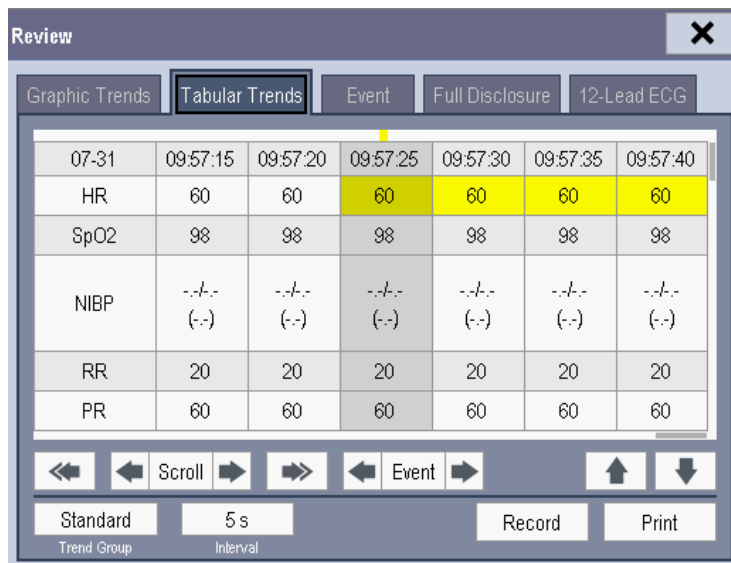
- To browse the graphic trends, you can either:
 - ◆ Select  or  beside **[Scroll]** to move the cursor one step to the left or right to navigate through the graphic trends, or
 - ◆ Select  or  to move the cursor one page to the left or right to navigate through the graphic trends.

A time indicating your current position is displayed above the parameter area. Numeric measurement values corresponding to the cursor location change as the cursor is moved. The measurement value that triggered high level alarm has red background. The one that triggered medium/low level alarm has yellow background.

- By selecting  or  beside **[Event]**, you can position the cursor to different event time.
- By selecting the **[Record]** button, you can print out the currently displayed graphic trends by the recorder.
- By selecting the **[Print]** button, you can set and print out the graphic trends report by the printer. For how to set the graphic trends report, please refer to the **Print** chapter.

28.3 Reviewing Tabular Trends

In the **[Review]** menu, select **[Tabular Trends]** to access the following window.







	07:31	09:57:15	09:57:20	09:57:25	09:57:30	09:57:35	09:57:40
HR		60	60	60	60	60	60
SpO2		98	98	98	98	98	98
NIBP		--/-- (--)	--/-- (--)	--/-- (--)	--/-- (--)	--/-- (--)	--/-- (--)
RR		20	20	20	20	20	20
PR		60	60	60	60	60	60

Standard 5 s Record Print
Trend Group Interval



Events are marked with colors in window's top area. Red represents high level alarm event. Yellow represents medium/low level alarm event. Green represents manual event.

In this review window:

- Select **[Trend Group]** and you can select a trend group for viewing in the popup menu. If **[Custom 1]** or **[Custom 2]** is selected, you can further select **[Define Trend Group]**. Then you can select the parameters for viewing in the popup menu.

- You can change the resolution of the trend data by selecting [**Interval**] and then selecting the appropriate setting:
 - ◆ [**5 s**] or [**30 s**]: select to view up to 4 hours of tabular trends at 5- or 30-second resolution.
 - ◆ [**1 min**], [**5 min**], [**10 min**], [**15 min**], [**30 min**], [**1 h**], [**2 h**] or [**3 h**]: select to view up to 120 hours of tabular trends at your selected resolution.
 - ◆ [**NIBP**]: select to view the tabular trends when NIBP measurements were acquired.
- To browse the tabular trends, you can either:
 - ◆ Select  or  beside [**Scroll**] to drag the scrollbar left or right to navigate through the trend database, or
 - ◆ Select  or  to scroll left or right to navigate through the trend database.

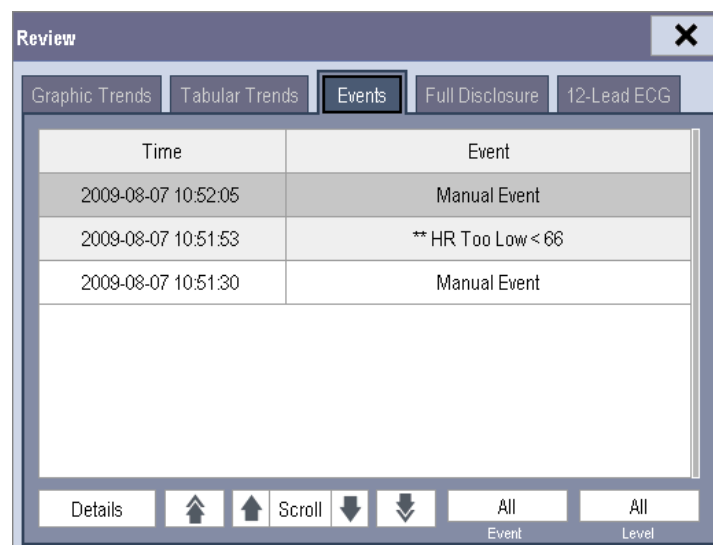
The measurement value that triggered high level alarm has red background. The one that triggered medium/low level alarm has yellow background.

- By selecting  or  beside [**Event**], you can position the cursor to different event time.
- By selecting the [**Record**] button, you can access the [**Record Setup**] menu and set the start and end time of the tabular trends you want to record. This feature is not available when reviewing a history patient. By further selecting [**Record**], you can print out the currently displayed tabular trends by the recorder.
- By selecting the [**Print**] button, you can set and print out the tabular trends report by the printer. For how to set the tabular trends report, please refer to the *Print* chapter.

28.4 Reviewing Events

In the [**Review**] menu, select [**Events**] to access the following window.


The events that can be reviewed include parameter alarm events, arrhythmia alarm events and manual events. When an event occurs, all the measurement numerics at the event trigger time and related waveforms 4 seconds respectively before and after the event trigger time are stored.



In this window:

- You can view the desired events by selecting **[Event]**.
- You can view the desired events according to the level by selecting **[Level]**.

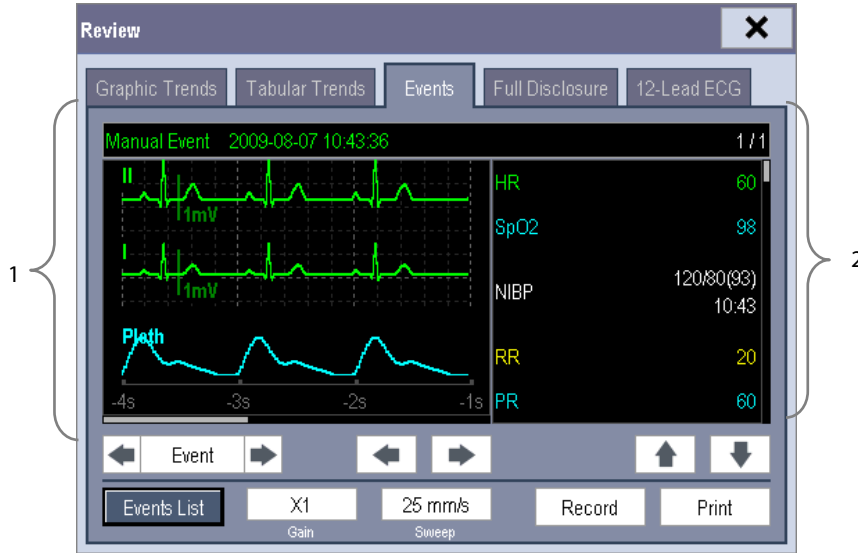After selecting the desired event, you can select **[Details]** to access the following window. In this window, the waveform area displays the waveforms related to the event, and the parameter area displays the parameter values happened at the event trigger time.




1. Waveform area 2. Parameter area

In this window:

- You can select or to navigate through the waveforms.
- You can select or beside the **[Event]** button to switch between events.
- You can set the desired **[Gain]** for ECG waveform.
- You can set the desired **[Sweep]**.
- By selecting the **[Record]** button, you can print out the currently displayed alarm events by the recorder.
- By selecting the **[Print]** button, you can print out the currently displayed alarm events by the printer.
- By selecting the **[Events List]** button, you can view the events list.

28.5 Reviewing Waveforms

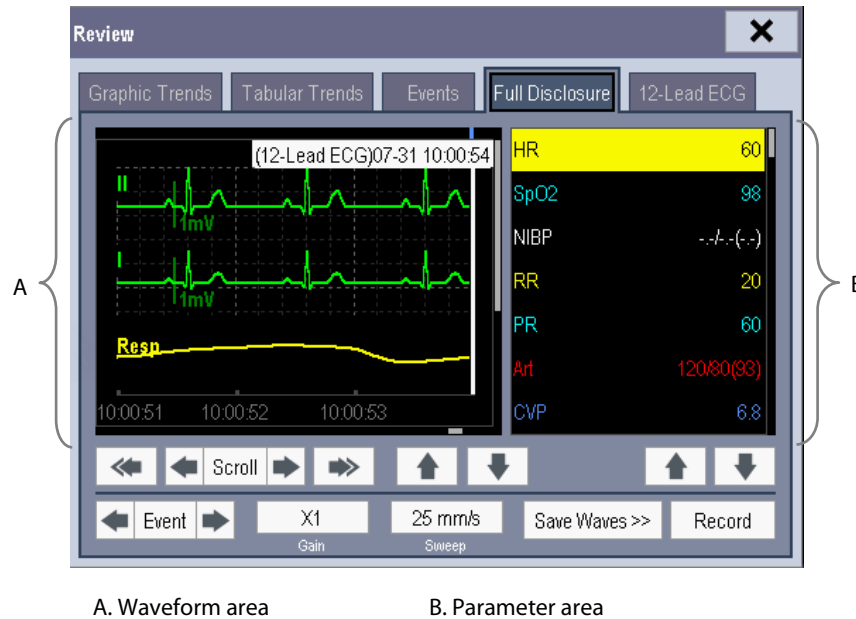In the **[Review]** menu, select **[Full Disclosure]** to access the following window.



In this review window:

- To review full-disclosure waveforms, you need to save waveforms first. Select **[Save Waves >>]** and then select the parameters whose waveforms you want to view. To save full-disclosure waveform, your monitor must be equipped with a CF storage card.
- To view the waveforms, you can either:
 - ◆ Select or beside the **[Scroll]** button to move the cursor one step left or right to navigate through the waveforms, or
 - ◆ Select or to move the cursor one page left or right to navigate through the waveforms.

A time indicating your current position is displayed at the top of the waveform area. Numeric measurement values corresponding to the cursor location are displayed in the parameter area, and change as the cursor is moved.

- You can change the ECG wave gain by selecting **[Gain]** and then selecting the appropriate setting.
- You can change the waveform sweep speed by selecting **[Sweep]** and then selecting the appropriate setting.
- By selecting the **[Record]** button, you can print out the first three waveforms and measurement numerics by the recorder.
- By selecting or beside the **[Event]** button, you can position the cursor between events.

FOR YOUR NOTES

29 Calculations

29.1 Introduction

The calculation feature is available with your patient monitor. The calculated values, which are not directly measured, are computed based on the values you provide.

You can perform the following calculations:

- Dose calculations
- Oxygenation calculations
- Ventilation calculations
- Hemodynamic calculations
- Renal calculations

To perform a calculation, select **[Main Menu]** → **[Calc >>]**, or the **[Calculations]** QuickKey and then select the calculation you want to perform.

NOTE

- **The calculation feature is independent of other monitoring functions and can be therefore used for patients being monitored by other monitors. Any operation in a calculation window does not affect the patient monitoring by the local patient monitor.**
-



WARNING

- **After the calculation is finished, verify the entered values are correct and the calculated values are appropriate. We assume no responsibility for any consequences caused by wrong entries and improper operations.**
-

29.2 Dose Calculations

29.2.1 Performing Calculations

To perform a dose calculation:

1. Select [**Main Menu**]→[**Calculations >>**]→[**Dose >>**], or select [**Calculations**] QuickKey→[**Dose >>**].
2. Select, in turn, [**Patient Cat.**] and [**Drug Name**] and then select the appropriate settings. The dose calculation program has a library of commonly used drugs, of which Drug A through Drug E are for those not specified in this library.
 - ◆ Drug A, B, C, D, E
 - ◆ Aminophylline
 - ◆ Dobutamine
 - ◆ Dopamine
 - ◆ Epinephrine
 - ◆ Heparin
 - ◆ Isuprel
 - ◆ Lidocaine
 - ◆ Nipride
 - ◆ Nitroglycerin
 - ◆ Pitocin
3. The system gives a set of default values when the above steps are finished. However, these values cannot be used as the calculated values. The user must enter values following the doctor's instructions, and then the calculated values can only be used
4. Enter the patient's weight.
5. Enter other values.
6. Verify if the calculated values are correct.

29.2.2 Selecting the Proper Drug Unit

Each drug has its fixed unit or unit series. Among a unit series, one unit may change to another automatically depending on the entered value.

The units for each drug are as follows:

- Drug A, B, C, Aminophylline, Dobutamine, Dopamine, Epinephrine, Isuprel, Lidocaine, Nipride and Nitroglycerin use the unit series: g, mg and mcg.
- Drug D, Heparin and Pitocin use the unit series: Unit, KU (kilo units) and MU (million units).
- Drug E uses the unit: mEq (milli-equivalents).

You must select the proper drug name (A, B, C, D or E) according to the units when you define a drug not listed in this library.

NOTE

- For neonate patients, [**Drip Rate**] and [**Drop Size**] are disabled.
-

29.2.3 Titration Table





To open the titration table, select [**Titration Table >>**] in the [**Dose Calculation**] window after the dose calculation is finished.

In the titration table, when you change:

- [**Reference**]
- [**Interval**]
- [**Dose Type**]

The titrated values change accordingly.

You can also:

- Select  or , or  or  beside the vertical scrollbar to view more values.
- Select [**Record**] to print out the currently displayed titrated values by the recorder.

29.3 Oxygenation Calculations

29.3.1 Performing Calculations

To perform an oxygenation calculation:

1. Select [**Main Menu**]→[**Calculations >>**]→[**Oxygenation >>**], or select [**Calculations**] QuickKey→[**Oxygenation >>**].
2. Enter values for calculation.
3. Select the [**Calculate**] button. The system performs a calculation per the current settings and displays the calculated values.
 - ◆ If a calculated value is outside the range, its background will highlight in yellow. You can select [**Range**] to view its normal range in the unit field.
 - ◆ Invalid values are displayed as [---].

In the [**Oxygenation Calculation**] window, you can:

- Change the pressure unit, Hb unit and oxygen content unit by selecting [**Press. Unit**], [**Hb Unit**] and [**OxyCont Unit**] and then selecting the appropriate settings. The changes take effect automatically.
- Trigger a recording by selecting the [**Record**] button. The currently displayed oxygenation calculations are printed out by the recorder.
- Review the previously performed calculations by selecting [**Review**].

29.3.2 Entered Parameters

Abbreviation	Unit	Full spelling
C.O.	L/min	cardiac output
FiO ₂	%	percentage fraction of inspired oxygen
PaO ₂	mmHg	partial pressure of oxygen in the arteries
PaCO ₂	mmHg	partial pressure of carbon dioxide in the arteries
SaO ₂	%	arterial oxygen saturation
PvO ₂	mmHg	partial pressure of oxygen in venous blood
SvO ₂	%	venous oxygen saturation
Hb	g/L	hemoglobin
CaO ₂	ml/L	arterial oxygen content
CvO ₂	ml/L	venous oxygen content
VO ₂	ml/min	oxygen consumption
RQ	None	respiratory quotient
ATMP	mmHg	atmospheric pressure
Height	cm	height
Weight	kg	weight

29.3.3 Calculated Parameters

Abbreviation	Unit	Full spelling
BSA	m ²	body surface area
VO ₂ calc	ml/min	oxygen consumption
C(a-v)O ₂	ml/L	arteriovenous oxygen content difference
O ₂ ER	%	oxygen extraction ratio
DO ₂	ml/min	oxygen transport
PAO ₂	mmHg	partial pressure of oxygen in the alveoli
AaDO ₂	mmHg	alveolar-arterial oxygen difference
CcO ₂	ml/L	capillary oxygen content
Qs/Qt	%	venous admixture
C.O. calc	L/min	calculated cardiac output

29.4 Ventilation Calculations

29.4.1 Performing Calculations

To perform a ventilation calculation:

1. Select [**Main Menu**]→[**Calculations >>**]→[**Ventilation >>**], or select [**Calculations**] QuickKey→[**Ventilation >>**].
2. Enter values for calculation. If the patient monitor is connected to an anesthesia machine or a ventilator, the system automatically loads the supported parameter values to the [**Ventilation Calculation**] window.
3. Select the [**Calculate**] button. The system performs a calculation per the current settings and displays the calculated values.
 - ◆ If a calculated value is outside the range, its background will highlight in yellow. You can select [**Range**] to view its normal range in the unit field.
 - ◆ Invalid values are displayed as [---].

In the [**Ventilation Calculation**] window, you can:

- Change the pressure unit by selecting [**Press. Unit**] and then selecting the appropriate setting. Corresponding pressure values shall convert and update automatically.
- Trigger a recording by selecting the [**Record**] button. The currently displayed ventilation calculations are printed out by the recorder.
- Review the previously performed calculations by selecting [**Review**].

29.4.2 Entered Parameters

Abbreviation	Unit	Full spelling
FiO ₂	%	percentage fraction of inspired oxygen
RR	rpm	respiration rate
PeCO ₂	mmHg	partial pressure of mixed expiratory CO ₂
PaCO ₂	mmHg	partial pressure of carbon dioxide in the arteries
PaO ₂	mmHg	partial pressure of oxygen in the arteries
TV	ml	tidal volume
RQ	None	respiratory quotient
ATMP	mmHg	atmospheric pressure

29.4.3 Calculated Parameters

Abbreviation	Unit	Full spelling
PAO ₂	mmHg	partial pressure of oxygen in the alveoli
AaDO ₂	mmHg	alveolar-arterial oxygen difference
Pa/FiO ₂	mmHg	oxygenation ratio
a/AO ₂	%	arterial to alveolar oxygen ratio
MV	L/min	minute volume
Vd	ml	volume of physiological dead space
Vd/Vt	%	physiologic dead space in percent of tidal volume
VA	L/min	alveolar volume

29.5 Hemodynamic Calculations

29.5.1 Performing Calculations

To perform a hemodynamic calculation:

1. Select [**Main Menu**]→[**Calculations >>**]→[**Hemodynamic >>**], or select [**Calculations**] QuickKey→[**Hemodynamic >>**].
2. Enter values for calculation.
 - ◆ For a patient who is being monitored, [**HR**], [**Art mean**], [**PA mean**] and [**CVP**] are automatically taken from the currently measured values. If you just have performed C.O. measurements, [**C.O.**] is the average of multiple thermodilution measurements. [**Height**] and [**Weight**] are the patient's height and weight you have entered. If the monitor does not provide these values, their fields appear blank.
 - ◆ For a patient who is not being monitored, confirm the values you have entered.
3. Select the [**Calculate**] button. The system performs a calculation per the current settings and displays the calculated values.
 - ◆ If a calculated value is outside the range, its background will highlight in yellow. You can select [**Range**] to view its normal range in the unit field.
 - ◆ Invalid values are displayed as [---].

In the [**Hemodynamic Calculation**] window, you can:

- Trigger a recording by selecting the [**Record**] button. The currently displayed renal calculations are printed out by the recorder.
- Review the previously performed calculations by selecting [**Review**].

29.5.2 Entered Parameters

Abbreviation	Unit	Full spelling
C.O.	L/min	cardiac output
HR	bpm	heart rate
PAWP	mmHg	pulmonary artery wedge pressure
Art Mean	mmHg	artery mean pressure
PA Mean	mmHg	pulmonary artery mean pressure
CVP	mmHg	central venous pressure
EDV	ml	end-diastolic volume
Height	cm	height
Weight	kg	weight

29.5.3 Calculated Parameters

Abbreviation	Unit	Full spelling
BSA	m ²	body surface area
C.I.	L/min/m ²	cardiac index
SV	ml	stroke volume
SI	ml/m ²	stroke index
SVR	DS/cm ⁵	systemic vascular resistance
SVRI	DS·m ² /cm ⁵	systemic vascular resistance index
PVR	DS/cm ⁵	pulmonary vascular resistance
PVRI	DS·m ² /cm ⁵	pulmonary vascular resistance index
LCW	kg·m	left cardiac work
LCWI	kg·m/m ²	left cardiac work index
LVSW	g·m	left ventricular stroke work
LVSWI	g·m/m ²	left ventricular stroke work index
RCW	kg·m	right cardiac work
RCWI	kg·m/m ²	right cardiac work index
RVSW	g·m	right ventricular stroke work
RVSWI	g·m/m ²	right ventricular stroke work index
EF	%	ejection fraction

29.6 Renal Calculations

29.6.1 Performing Calculations

To perform a renal calculation:

1. Selecting [**Main Menu**]→[**Calculations >>**]→[**Renal >>**], or select [**Calculations**] QuickKey→[**Renal >>**].
2. Enter values for calculation.
3. Select the [**Calculate**] button. The system performs a calculation per the current settings and displays the calculated values.
 - ◆ If a calculated value is outside the range, its background will highlight in yellow. You can select [**Range**] to view its normal range in the unit field.
 - ◆ Invalid values are displayed as [---].

In the [**Renal Calculation**] window, you can:

- Trigger a recording by selecting the [**Record**] button. The currently displayed renal calculations are printed out by the recorder.
- Review the previously performed calculations by selecting [**Review**].

29.6.2 Entered Parameters

Abbreviation	Unit	Full spelling
URK	mmol/L	urine potassium
URNa	mmol/L	urinary sodium
Urine	ml/24h	urine
Posm	mOsm/ kgH ₂ O	plasma osmolality
Uosm	mOsm/ kgH ₂ O	urine osmolality
SerNa	mmol/L	serum sodium
Cr	μmol/L	creatinine
UCr	μmol/L	urine creatinine
BUN	mmol/L	blood urea nitrogen
Height	cm	height
Weight	kg	weight

29.6.3 Calculated Parameters




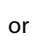
Abbreviation	Unit	Full spelling
URNaEx	mmol/24h	urine sodium excretion
URKEx	mmol/24h	urine potassium excretion
Na/K	%	sodium potassium ratio
CNa	ml/24h	clearance of sodium
Clcr	ml/min	creatinine clearance rate
FENa	%	fractional excretion of sodium
Cosm	ml/min	osmolar clearance
CH ₂ O	ml/h	free water clearance
U/P osm	None	urine to plasma osmolality ratio
BUN/Cr	None*	blood urea nitrogen creatinine ratio
U/Cr	None	urine-serum creatinine ratio

*: BUN/Cr is a ratio under the unit of mol.

29.7 Understanding the Review Window

With the review feature, you can review oxygenation, ventilation, hemodynamic and renal calculations. The review window for each calculation is similar. Take the hemodynamic calculations review window for example, you can access it by selecting **[Review]** in the **[Hemodynamic Calculation]** window.

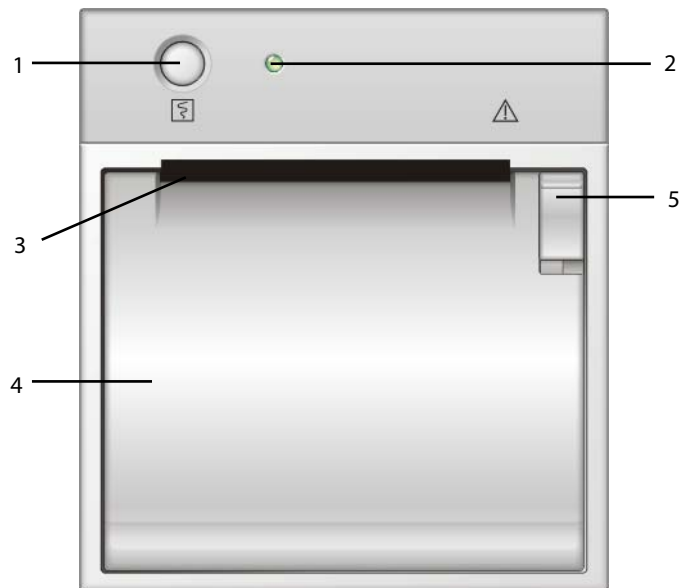
In this review window:

- You can select , ,  or  to view more values.
- The values that exceed the range are displayed in yellow background. The **[Unit]** field displays parameter units. If some parameter values are outside of their normal ranges, you can view their normal range in the **[Unit]** field by selecting **[Range]**.
- You can review an individual calculation by selecting its corresponding column and then selecting **[Original Calc]**. You can record the currently displayed calculations or perform another calculation in this window.

30 Recording

30.1 Using a Recorder

The thermal recorder records patient information, measurement numerics, up to three waveforms, etc.



1. Start/Stop key: press to start a recording or stop the current recording.
2. Indicator
 - ◆ On: when the recorder works correctly.
 - ◆ Off: when the monitor is switched off.
 - ◆ Flashes: if an error occurred to the recorder, e.g., the recorder runs out of paper.
3. Paper outlet
4. Recorder door
5. Latch

30.2 Overview of Recording Types

By the way recordings are triggered, the recordings can be classified into the following categories:

- Manually-triggered realtime recordings.
- Timed recordings.
- Alarm recordings triggered by an alarm limit violation or an arrhythmia event.
- Manually-triggered, task-related recordings.

The task-related recordings include:


- Frozen wave recording
- Graphic trends recording
- Tabular trends recording
- Events recording; parameter alarm recording, arrh. alarm recording, manual events recording
- Wave review recording
- Interpretation of resting 12-lead ECG result recording
- Titration table recording
- Hemodynamic calculations recording
- Oxygenation calculations recording
- Ventilation calculations recording
- Renal calculations recording
- oxyCRG recording
- C.O. curve recording
- PAWP recording
- Respiratory loops recording
- Monitor information recording

NOTE

- For details about alarm recording, refer to the chapter 7 *Alarms*.
 - For details about task-related recordings, refer to respective sections of this manual.
-

30.3 Starting and Stopping Recordings


To manually start a recording, you can either:

- Select the  hardkey on the front of either the patient monitor or the recorder module, or
- Select the **[Record]** button from the current menu or window.

Automatic recordings will be triggered in the following conditions:

- Timed recordings will start automatically at preset intervals.
- If both **[Alarm]** and **[Alm Rec]** for a measurement are set on, an alarm recording will be triggered automatically as alarms occur.

To manually stop a recording, you can either:

- Select the  hardkey again, or
- Select **[Clear All Tasks]** in the **[Record Setup]** menu.

Recordings stop automatically when:

- The runtime is over.
- The recorder runs out of paper.
- When the recorder has an alarm condition.

When a recording is stopped, the following markers will be added:

- Automatically stopped recording: print two columns of '*' at the end of the report.
- Manually or abnormally stopped recording: print one column of '*' at the end of the report.

30.4 Setting up the Recorder

30.4.1 Accessing the Record Setup Menu

By selecting [Main Menu]→[Record Setup >>], you can access the [Record Setup] menu.

30.4.2 Selecting Waveforms for Recording

The recorder can record up to 3 waveforms at a time. You can select, in turn, [Waveform 1], [Waveform 2] and [Waveform 3] in the [Record Setup] menu, and then select the waveforms you want. You can also turn off a waveform recording by selecting [Off]. These settings are intended for realtime and scheduled recordings.

30.4.3 Setting the Realtime Recording Length

After starting a realtime recording, the recording time depends on your monitor's settings. In the [Record Setup] menu, select [Length] and toggle between [8 s] and [Continuous].

- [8 s]: record 4-second waveforms respectively before and after current moment.
- [Continuous]: record the waveforms from the current moment until stopped manually.

30.4.4 Setting the Interval between Timed Recordings

Timed recordings start automatically at preset intervals. Each recording lasts 8 seconds. To set the interval between timed recordings: in the [Record Setup] menu, select [Interval] and then select the appropriate setting.

30.4.5 Changing the Recording Speed

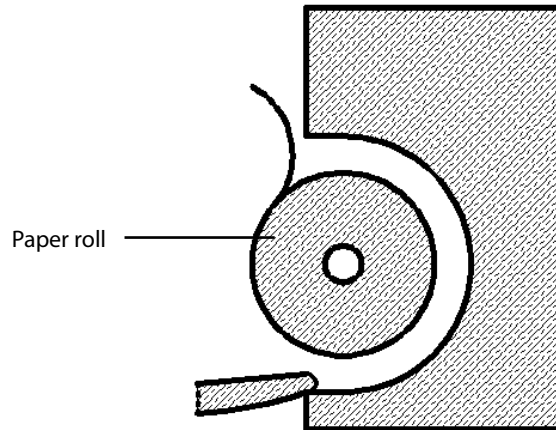
In the [Record Setup] menu, select [Paper Speed] and toggle between [25 mm/s] and [50 mm/s]. This setting is for all recordings containing waveforms.

30.4.6 Clearing Recording Tasks

In the [Record Setup] menu, select [Clear All Tasks]. All queued recording tasks are cleared and the current recording is stopped.

30.5 Loading Paper

1. Use the latch at the upper right of the recorder door to pull the door open.
2. Insert a new roll into the compartment as shown below.
3. Close the recorder door.
4. Check if paper is loaded correctly and the paper end is feeding from the top.



CAUTION

- **Use only specified thermal paper. Otherwise, it may cause damage to the recorder's printhead, the recorder may be unable to print, or poor print quality may result.**
 - **Never pull the recorder paper with force when a recording is in process. Otherwise, it may cause damage to the recorder.**
 - **Do not leave the recorder door open unless you reload paper or remove troubles.**
-

30.6 Removing Paper Jam

If the recorder works incorrectly or produces unusual sounds, check if there is a paper jam first. If a paper jam is detected, follow this procedure to remove it:

1. Open the recorder door.
2. Take out the paper and tear off the draped part.
3. Reload the paper and close the recorder door.

30.7 Cleaning the Recorder Printhead

If the recorder has been used for a long time, deposits of paper debris may collect on the printhead compromising the print quality and shortening the lifetime of the roller. Follow this procedure to clean the printhead:

1. Take measures against the static electricity such as Disposable Wrist Strap for the work.
2. Open the recorder door and take out the paper.
3. Gently wipe around the printhead using cotton swabs dampened with alcohol.
4. After the alcohol has completely been dried, reload the paper and close the recorder door.

 **CAUTION**

-
- **Do not use anything that may destroy the thermal element.**
 - **Do not add unnecessary force to the thermal head.**
-

FOR YOUR NOTES

31 Printing

31.1 Printer

The monitor can output patient reports via a connected printer. So far, the monitor supports the following printer:

- HP LaserJet 1505n
- HP LaserJet P2035n
- HP LaserJet P4015n

The specifications of the reports the monitor prints are:

- Paper: A4, Letter
- Resolution: 300 dpi

For more details about the printer, see the document accompanying the printer. With the upgrading of products, the monitor will support more printers and no prior notice will be given. If you have any doubt about the printer you have purchased, contact our company.

31.2 Connecting a printer

To print the reports or the trend data of a patient, you can choose either:

- the local printer

Connect the printer and the patient monitor directly with a network cable, and then start printing what you want,
or

- the Central Monitoring System

If your monitor is connected to a central monitoring system, it is recommended to use the central monitoring system for printing.

31.3 Setting Up the Printer

To set the printer's properties, select **[Main Menu]→[Print Setup >>]→[Printer Setup >>]**. In the **[Printer Setup]** menu, you can:

- Select a connected printer

Select **[Printer]** and then select a connected printer as the monitor's printer.

- Search for a printer

If your selected printer is not in the list or a new printer is added into the network, you can select the **[Search Printer]** to re-search for all printers in the network.

- Set up the paper

Select **[Paper Size]** and toggle between **[A4]** and **[Letter]**.

31.4 Starting Report Printouts

Reports	Contents	Procedures
ECG reports	ECG waveforms and relevant parameter values	Select [Main Menu]→[Print Setup >>]→[ECG Reports >>]→[Print]
Tabular trends	Depend on the selected parameter group, resolution and time period	Select [Main Menu]→[Print Setup >>]→[Tabular Trends Reports >>]→[Print] , or select [Main Menu]→[Review >>]→[Tabular Trends]→[Print]→[Print]
Graphic trends	Depend on the selected parameter group, resolution and time period	Select [Main Menu]→[Print Setup >>]→[Graphic Trends Reports >>]→[Print] , or select [Main Menu]→[Review >>]→[Graphic Trends]→[Print]→[Print]
Arrh. alarm review	ECG waveforms and relevant parameter values	Select [Print] in [Arrh. Events]
Parameter alarm review	Depend on the selected alarms	Select [Main Menu]→[Review >>]→[Alarms]→[Print]
Interpretation of resting 12-lead ECG	12-lead ECG waveforms and analysis results	Select [12-lead Analysis]→[Print] when a interpretation of resting 12-lead ECG is completed, or select [Main Menu]→[Review >>]→[12-lead Analysis]→[Print]
Realtime waves	Depend on the selected waveforms	Select [Main Menu]→[Print Setup >>]→[Realtime Reports >>]→[Print]

31.5 Stopping Reports Printouts

To stop report printouts, select **[Main Menu]→[Print Setup >>]→[Stop All Reports]**.

31.6 Setting Up Reports

31.6.1 Setting Up ECG Reports

You can print out ECG reports only under full-screen, half-screen or 12-lead monitoring screen. To set up ECG reports, select **[Main Menu]**→**[Print Setup >>]**→**[ECG Reports >>]**.

- **[Amplitude]**: set the amplitude of the ECG waveforms.
- **[Sweep]**: set the wave print speed.
- **[Auto Interval]**: If **[Auto Interval]** is set to **[On]**, the system will automatically adjust the space between waveforms to avoid overlapping.
- **[Gridlines]**: choose whether to show gridlines.
- **[12-Lead Format]**: If you select **[12×1]**, 12 waveforms will be printed on a paper from top to bottom. If you select **[6×2]**, 12 waveforms will be printed from left to right with 6 waveforms on each half part and a rhythm waveform will be printed at the bottommost. If you select **[3×4]**, 12 waveforms will be printed from left to right with 3 waveforms on each of the 4 columns and a rhythm waveform will be printed at the bottommost.

31.6.2 Setting Up Tabular Trends Reports

To set up tabular trends reports, select **[Main Menu]**→**[Print Setup >>]**→**[Tabular Trends Reports >>]**.

- **Start time**: You can set a time period whose trend data will be printed out by setting **[From]** and **[Back]**. For example, if you set **[From]** as 2007-4-2 10:00:00 and **[Back]** as **[2 h]**, the outputted data will be from 2007-4-2 08:00:00 to 2007-4-2 10:00:00. In addition, the **[Back]** can be set to either:
 - ◆ **[Auto]**: If **[Report Layout]** is set to **[Time Oriented]**, the report will be printed by time. If **[Report Layout]** is set to **[Parameter Oriented]**, the report will be printed by parameters.
 - ◆ **[All]**: If you select **[All]**, all trend data will be printed out. In this case, it is no need to set **[From]**.
- **[Interval]**: choose the resolution of the tabular trends printed on the report.
- **[Report Layout]**: If you select **[Time Oriented]**, the report will be printed by time. If you select **[Parameter Oriented]**, the report will be printed by parameters.
- **[Select Parameter >>]**: from the popup menu, you can:
 - ◆ **[Currently Displayed Trended Parameters]**: print the parameter trend data selected from the **[Tabular Trends]**.
 - ◆ **[Standard Parameter Group]**: select the standard parameter group for printing.
 - ◆ **[Custom]**: You can define a parameter group for printing from the parameters displayed in the low part of the menu.

31.6.3 Setting Up Graphic Trends Reports

To set up graphic trends reports, select **[Main Menu]**→**[Print Setup >>]**→**[Graphic Trends Reports >>]**. As setting up graphic trends reports is similar with tabular trends reports, you can refer to the **Setting Up Tabular Trend Reports** section for details.

31.6.4 Setting Up Realtime Reports

To set up realtime reports, select [**Main Menu**]→[**Print Setup >>**]→[**Realtime Reports >>**].

- [**Sweep**]: set the wave print speed.
- [**Select Wave >>**]: from the popup menu, you can:
 - ◆ [**Current**]: select the currently displayed waves for printing.
 - ◆ [**Select Wave**]: select the desired waves for printing.

31.7 End Case Reports

ECG reports, tabular trends reports, graphic trends reports, NIBP review reports and realtime reports can be set as end case reports. When you discharge a patient, the system will automatically print out all contents that are set as end case reports.

For example, to set ECG report as end case report:

1. select [**Main Menu**]→[**Print Setup >>**]→[**ECG Report >>**].
2. select [**End Case Report**]→[**Set as End Case Report**] and then select [**Ok**] from the popup dialog box.
3. set as described in the **31.6.1 Setting Up ECG Reports**.

31.8 Printer Statuses

31.8.1 Printer Out of Paper

When the printer runs out of paper, the print request will not be responded. If there are too many print jobs that are not responded, a printer error may occur. In these cases, you need to install paper and then re-send the print request. Restart the printer if necessary.

Therefore, you'd better ensure that there is enough paper in the printer before sending a print request.

31.8.2 Printer Status Messages

Printer Status Message	Possible causes and suggested action
Printer unavailable	The selected printer is not available. Check if the printer is switched on or correctly connected or installed with paper.

32 Other Functions

32.1 Marking Events

During patient monitoring, some events may exert effects on the patient and as a result change the waveforms or numerics displayed on the monitor. To help analysing the waveforms or numerics at that time, you can mark these events.

Select [**Main Menu**]→[**Mark Event >>**]. In the popup menu, you can select the waves to be stored when a manual event is triggered. You can select [**Trigger Manual Event**] from the [**Mark Event**] menu or the [**Manual Event**] QuickKey to trigger a manual event and store it at the same time.

When you are reviewing graphic trends, tabular trends or full-disclosure waveforms, the manual event symbol is displayed at the time the event is triggered.

32.2 Privacy Mode

Privacy mode is only available when a patient who is admitted at a patient monitor is also monitored by the central station.

To activate the privacy mode:

1. Select [**Main Menu**]→[**Screen Setup >>**].
2. Select [**Privacy Mode**] to activate the privacy mode.

The patient monitor behaves as follows as soon as the privacy mode is activated:

- The screen turns blank and [**Under monitoring. Press any key to exit the privacy mode.**] is displayed.
- Monitoring and data storing continue but patient data is only visible at the central station.
- Alarms can still be triggered. But all audible alarms are suppressed and the alarm light is deactivated at the patient monitor.
- All system sounds are suppressed, including heart beat tone, pulse tone, all prompt tones, etc.

WARNING

- **During privacy mode, all audible alarms are suppressed and the alarm light is deactivated at the patient monitor. Alarms sound only at the central station.**
-
-

To cancel the privacy mode, proceed as follows:

- Press any key.

The patient monitor exits the privacy mode automatically in one of the following situations:

- The patient monitor disconnects from central station.
- The alarm of **Battery Too Low** and **The monitor will quit soon. Please use AC power.** message appear.

32.3 Night Mode

To avoid disturbing the patient, night mode may be used.

To activate the night mode:

1. Select the **[Night Mode]** QuickKey or **[Main Menu]→[Screen Setup >>]→[Night Mode >>]**.
2. In the pop-up menu, set the desired brightness, alarm volume, QRS volume, key volume, NIBP end tone, or whether to stop NIBP measurement or not. When **[Stop NIBP]** is selected, all the NIBP measurements terminate after entering the night mode.
3. Select the **[Enter Night Mode]** button.

To cancel the night mode:

1. Select the **[Night Mode]** QuickKey or **[Main Menu]→[Screen Setup >>]→[Night Mode >>]**.
2. Select **[Ok]** in the popup.

WARNING

- **Before entering night mode, confirm the settings of brightness, alarm volume, QRS volume, and key volume. Pay attention to the potential risk when the setting value is a bit low.**
-

32.4 Analog Output

The patient monitor provides analog output signals to accessory equipment via the Micro-D connector on the rear of the monitor. To obtain analog output signals, connect the accessory equipment such as an oscillograph, etc. to the monitor and then follow this procedure:

1. Select **[Main Menu]** then **[Analog Output Setup>>]**.
2. Select **[Analog Out.]** and then select **[On]**.

NOTE

- **The analog output feature is seldom applied in clinical applications. You can contact your service personnel for more details.**
-

32.5 Transferring Data

You can transfer the patient data saved in the monitor to a PC via a crossover network cable or CF storage card, or within a LAN for data management, review or print.

32.5.1 Data Export System

You must install the data export system on the intended PC before performing the data transfer operation. Refer to the document accompanying the installation CD-ROM for installation instructions.

The data transfer feature supports patient management, data review, data format conversion, print, etc. in addition to data transfer. Refer to the help file of the system software for more details.

32.5.2 Transferring Data by Different Means

NOTE

- **Never enter the data transfer mode when the patient monitor is in normal operation or performs monitoring. You must re-start the patient monitor to exit the data transfer mode.**
-

Transfer data via a crossover network cable

Before transferring data using a crossover network cable, do as follows:

1. Connect one end of the crossover network cable to the patient monitor and the other end to the PC.
2. Set the IP address of the PC. This IP address must be in the same network segment with that of the patient monitor.
3. Make sure that the data export system is active on the PC.

Then, follow this procedure to transfer data:

1. Select **[Main Menu]**→**[Patient Data >>]**→**[Transfer Data]**.
2. Select **[Yes]** from the popup message box.
3. Input the IP address already set on the PC.
4. Select **[Start]** to start transferring data.

Transfer data within a LAN

Before transferring data within a LAN, do as follows:

1. Connect the patient monitor and the intended PC into the same LAN and acquire the PC's IP address.
2. Make sure that the data export system is active on the PC.

Follow the same procedure as via a crossover network cable to transfer data.

Transfer data via a CF storage card

1. Power off the patient monitor and remove the CF storage card from it. Refer to the **Basic Operations** section for details.
2. Run the data export system on the PC.
3. Insert the CF storage card into the card reader that connects the PC.
4. Perform the data transfer operation following the help file of the system software.

32.6 Nurse Call

The patient monitor also provides nurse call signals to a nurse call system connected to the monitor via the Nurse Call connector. To obtain nurse call signals, connect a nurse call system to the monitor and then follow this procedure:

1. Select [**Main Menu**]→[**Maintenance >>**]→[**User Maintenance >>**]→enter the required password.
2. Select [**Others >>**] to access the [**Others**] menu.
3. Select [**Nurse Call Setup >>**] to change the nurse call settings as follows:
 - Select [**Signal Type**] and toggle between [**Pulse**] and [**Continuous**].
 - ◆ [**Pulse**]: the nurse call signals are pulse signals and each pulse lasts 1 second. When multiple alarms occur simultaneously, only one pulse signal is outputted. If an alarm occurs but the previous one is not cleared yet, a new pulse signal will also be outputted.
 - ◆ [**Continuous**]: the nurse call signal lasts until the alarm ends, i.e. the duration of a nurse call signal equals to that of the alarm condition.
 - Select [**Contact Type**] and toggle between [**Normally Open**] and [**Normally Closed**].
 - ◆ [**Normally Open**]: select if your hospital's nurse call relay contact is normally open.
 - ◆ [**Normally Closed**]: select if your hospital's nurse call relay contact is normally closed.
 - Select [**Alm Lev**] and set the alarm level for nurse call-triggering alarms.
 - Select [**Alarm Cat.**] and then select the category to which the nurse call-triggering alarms belong.

Alarm conditions are indicated to nurses only when:

- The nurse call system is enabled,
- An alarm that meets your preset requirements occurs, and
- The monitor is not in the alarm paused or silence status.

NOTE

-
- If no setting is selected from [**Alm Lev**] or [**Alarm Cat.**], no nurse call signal will be triggered whatever alarms occur.
-



WARNING

-
- Do not rely exclusively on the nurse call system for alarm notification. Remember that the most reliable alarm notification combines audible and visual alarm indications with the patient's clinical condition.
-

32.7 Remote Display

This monitor enables remote display. It allows remote displays to be connected to the bedside monitor through network. The information coming from the monitor will be displayed on the remote display through the remote display driver so that clinical professionals can conveniently observe the patient's conditions from distance.

For details about remote display features, refer to the instructions for use accompanying the remote display driver.

NOTE

- **The contents displayed on the remote display are for convenient observance only and cannot be used for diagnostic interpretation.**
- **The user cannot operate the monitor through the remote display driver, namely, any operations performed through the remote display driver will not affect the monitor you observe.**

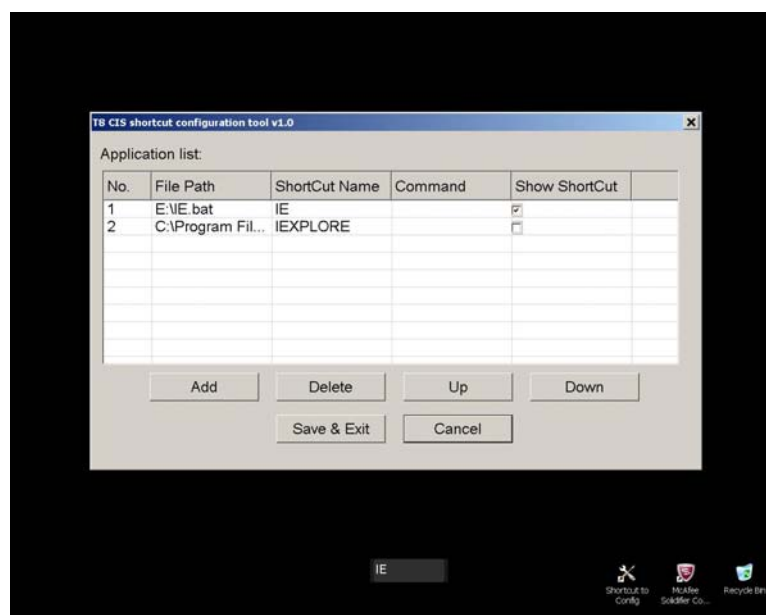
32.8 Built-in PC (only applicable to BeneView T8 patient monitor)

This monitor can be configured with Windows operating system. You can install and use the required PC application program on the monitor through Windows operating system.

32.8.1 Configuring Application Program ShortCuts

Select [CIS View] and PC ShortCuts area is displayed. Up to five PC application program ShortCuts can be displayed in this area. You can select from these ShortCuts to use the necessary software. To configure the ShortCuts,

1. Select [Main Menu]→[Maintenance]→[User Maintenance]→enter the required password→[CIS Maintenance].
2. To start the configuration tool, click "Config" on the desk or select [Start]→[My Computer]in the lower left corner of the desk. Run "Config.exe" under the path "C:\Program Files\Mindray".



Note

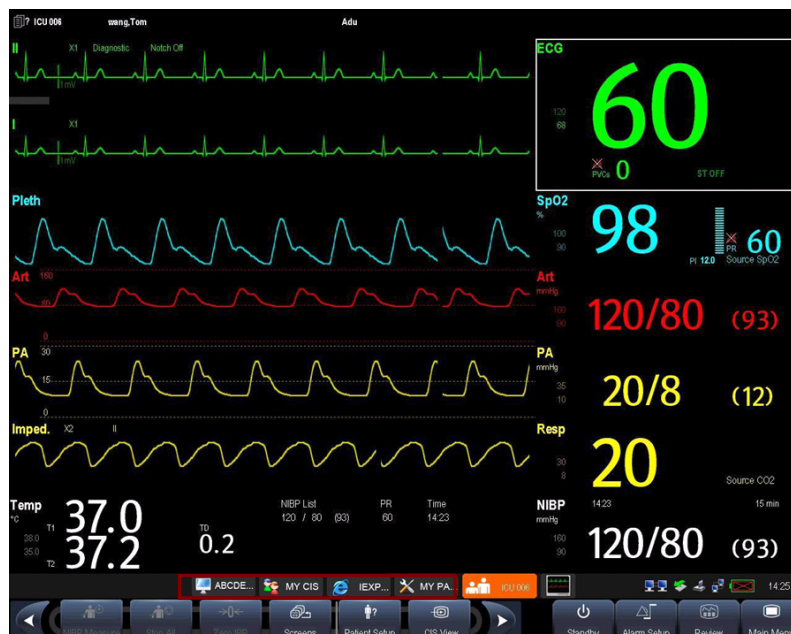
- The task bar is hidden automatically and is displayed when the mouse is placed at the bottom of the screen.

3. Select [**Add**] and select the application program to be added from the accessed dialog box. Then select [**Open**] to complete adding the application program.

You can select whether to display ShortCuts. [**Show ShortCut**] is ticked by default. If not ticked, the application program ShortCuts will not be displayed in the CIS ShortCuts area. Not selecting the checkbox usually occurs when application program is started up indirectly. In this case, add both startup program and started program into [**T8 CIS shortcut configuration tool**] and do not tick the started program. For example, if you want to start "iexplore.exe" application program to access "www.mindray.com" through "IE.bat" batch file, write parameters into the batch file. Then add "IE.bat" and "iexplore.exe" application programs into [**T8 CIS shortcut configuration tool**] and set "iexplore.exe" to unticked status. Finally, save the setting and exit.

4. Select [**Up**] or [**Down**] to change the display order of ShortCuts.
5. Select the cell under [**Shortcut Name**] to change the name of application program.
6. For the application program that can be started up together with parameter, select the cell under [**Command**] to configure a parameter of the application program. For example, if you add application program "iexplore.exe" into [**T8 CIS shortcut configuration tool**], set [**Command**] to "www.mindray.com". Then in the PC ShortCuts area, select the ShortCut of "iexplore.exe" and the system enters the website "www.mindray.com".
7. Select [**Save&Exit**] to finish ShortCut configuration.

Push [**Main Menu**] key on the monitor front panel to return to the main screen.



PC ShortCuts Area

32.8.2 Using PC Software


1. Select **[Main Menu]** and select **[CIS View]**. Or, select **[CIS View]** on the main screen directly. The ShortCuts of the PC software with which your monitor is configured will be displayed.
2. Select the ShortCut corresponding to the PC software you want to use to access the corresponding software screen. Only one PC software screen can be accessed at a time.

PC ShortCut area is automatically hidden while the PC software is running. It is automatically displayed when PC software display is minimized or turned off. You can adjust the size or display position of the window of application program via mouse.



-
- **All the waveforms and parameters on the monitor are hidden when PC software display is maximized. Pay attention to the risk arising from this operation.**
 - **Exit PC software or minimize PC software display when PC software is not in use.**
-

To hide PC software screen,

- Click  button in the upper right corner of the software screen.
- Click other area on the monitor screen.
- Push **[Main Menu]** key or **[Freeze]** key on the monitor front panel.

32.8.3 Switching over Users

The built-in PC supports two users. You can switch over between the two users through logout. The system default user is "CIS" and the login password is "MINDRAY". You can switch over the user to "Administrator" if necessary and the login password is "cisadmin".

32.8.4 Remote Login

If the monitor is networked via CIS connector, you can remotely log in to the monitor's built-in Windows system through the PC within the LAN to view the program running on the monitor and to carry out remote maintenance. Before exiting remote login, you need to restart the monitor's Windows system. To restart the monitor's Windows system, select **[Start]** → **[Run]** on the remote PC and then enter "Shutdown -f -r -t 0".

32.8.5 Using McAfee Solidifier

McAfee Solidifier is the default installation software of the built-in PC Windows system. McAfee Solidifier solidifies the executable files of the system, dynamic link library and batch files by way of dynamic white list. Executable files not included in the white list are held back so as to protect the system. You can update the application program or monitor Windows system via McAfee Solidifier.

Follow these steps to update an application program.

1. Enter update status

Before adding, updating or deleting an application program of the built-in PC, let McAfee Solidifier enter update status first. In this case, select "McAfee Solidifier" on the desk to enter command line dialog box and then enter command "sadmin bu".

Note

- **Before updating an application program, pay attention to anti-virus measures such as network anti-virus strategy and USB device virus scanning.**
-

2. Enter monitor status

After adding, updating or deleting an application program of the built-in PC, let McAfee Solidifier enter monitor status. In this case, select "McAfee Solidifier" on the desk to enter command line dialog box and then enter command "sadmin eu".

Other commonly used commands of McAfee Solidifier include:

- sadmin help: used to view the commonly used commands;
- sadmin status: used to view the status of McAfee.

32.9 Wireless Network

The patient monitors, each equipped with a wireless network card, constitute a wireless network via AP (access point). The designated service engineer or personnel shall be responsible for installing and configuring the wireless network for you and perform relative performance tests as well.

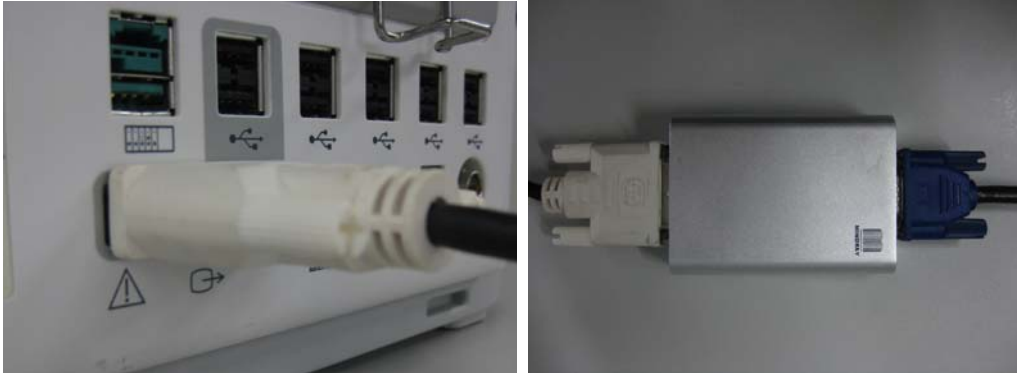
The radio device used in the monitor is in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC (Radio Equipment and Telecommunications Terminal Equipment Directive).

NOTE

- **The design, installation, reconstruction and maintenance of the wireless network's distribution shall be performed by authorized service personnel of our company.**
 - **The existence of obstacles (such as wall) will exert impact on data transferring or even cause network interruption.**
 - **The Central Monitoring System is capable of connecting up to 16 bedside monitors via the wireless network.**
-

32.10 Using DVI-VGA Adapter Box

The patient monitor can be connected with a VGA device via a DVI-VGA adapter box.



1. Connect the patient monitor's DVI output with DVI-VGA adapter box's DVI input.
2. Connect the DVI-VGA adapter box's VGA output with VGA device.


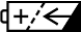
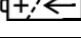
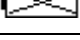
FOR YOUR NOTES

33 Batteries

33.1 Overview

The monitor is designed to operate on one or two Mindray LI23S002A rechargeable Lithium-ion battery whenever AC power supply is interrupted. The battery is charged whenever the patient monitor is connected to an AC power source regardless of whether or not the patient monitor is currently on. Since no external battery charger is supplied, the battery can only be charged inside the monitor so far. Whenever the AC power is interrupted during patient monitoring, the patient monitor will automatically run power from the internal batteries.

On-screen battery symbols indicate the battery status as follows:

	Indicates that batteries work correctly. The solid portion represents the current charge level of the batteries in proportion to its maximum charge level.
	Indicates that the batteries have low charge level and need to be charged.
	Indicates that the batteries are almost depleted and need to be charged immediately.
	Indicates that no batteries are installed or only one battery is installed to the BeneView T8 monitor.

The capacity of the internal battery is limited. If the battery capacity is too low, a technical alarm will be triggered and the [**Battery Too Low**] message displayed. At this moment, apply AC power to the patient monitor. Otherwise, the patient monitor will power off automatically before the battery is completely depleted.

NOTE

- **Take out the battery before the monitor is transported or will not be used for a long time.**
- **Use AC power supply when CIS is in use.**

WARNING

- **Keep the battery out of children's reach.**
 - **Use only specified batteries.**
-
-

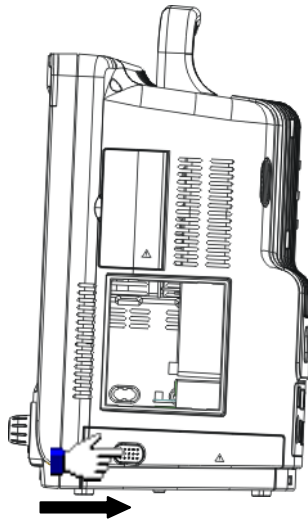
33.2 Installing or Replacing a Battery

BeneView T5

When the patient monitor uses two battery packs, one battery pack can be easily exchanged while the patient monitor operates from the other. If the patient monitor uses one battery pack, you should insert a new battery pack before the old one depletes.

To install or replace a battery, follow this procedure:

1. Push down the button on the battery door and then slide backward as indicated to open the battery door.



2. Push aside the latch latch fixing the battery and then remove the battery.
3. Place the new battery into the slot with its face up and its contact point inward.
4. If necessary, replace the other battery following the steps above.
5. Restore the latch to the original position and close the battery door.

NOTE

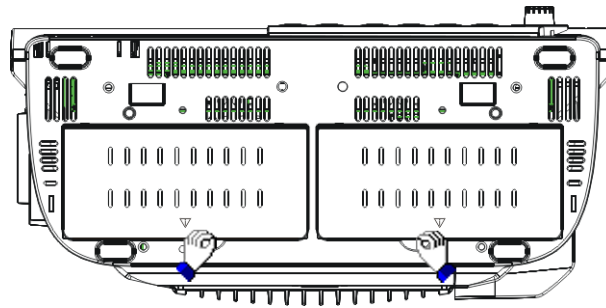
-
- **Using two batteries are recommended when SMR is connected.**
-

BeneView T8

The patient monitor uses two battery packs. If the two batteries have very different charge capacity, the message [**Diff. Battery Voltages**] is displayed. In this case, apply AC power to the patient monitor until the two batteries have approximately equal charge capacity or are both fully charged. You cannot use them before they have approximately equal charge capacity or are fully charged. In situations where no patient monitoring is performed or interrupting the patient monitoring is permitted, you can replace the batteries.

The patient monitor uses two batteries. You can install the batteries by following this procedure:

1. Turn off the patient monitor and disconnect the power cord and other cables.
2. Place the patient monitor with its face up.
3. Open the battery compartment door.



4. Place the batteries into the slots per the "+" and "-" indications.
5. Close the battery door and place the patient monitor upright.

33.3 Conditioning a Battery

A battery needs at least two conditioning cycles when it is put into use for the first time. A battery conditioning cycle is one complete, uninterrupted charge of the battery, followed by an uninterrupted discharge of the battery. Batteries should be conditioned regularly to maintain their useful life. Condition the batteries once when they are used or stored for two months, or when their run time becomes noticeably shorter.

To condition a battery, follow this procedure:

1. Disconnect the patient monitor from the patient and stop all monitoring and measuring procedures.
2. Insert the battery in need of conditioning into the battery slots of the patient monitor.
3. Apply AC power to the patient monitor and allow the battery to charge uninterruptedly for above 6 hours.
4. Remove AC power and allow the patient monitor to run from the battery until it shuts off.
5. Apply AC power again to the patient monitor and allow the battery to charge uninterruptedly for above 6 hours.
6. This battery is now conditioned and the patient monitor can be returned to service.

33.4 Checking a Battery

The performance of a rechargeable battery may deteriorate over time. To check the performance of a battery, follow this procedure:

1. Disconnect the patient monitor from the patient and stop all monitoring and measuring procedures.
2. Apply AC power to the patient monitor and allow the battery to charge uninterruptedly for above 6 hours.
3. Remove AC power and allow the patient monitor to run from the battery until it shuts off.
4. The operating time of the battery reflects its performance directly.

If the operating time of a battery is noticeably shorter than that stated in the specifications, replace the battery or contact your service personnel.

NOTE

-
- **Life expectancy of a battery depends on how frequent and how long it is used. For a properly maintained and stored lithium-ion battery, its life expectancy is about 3 years. For more aggressive use models, life expectancy can be less. We recommend replacing lithium-ion batteries every 3 years.**
 - **The operating time depends on the configuration and operation. For example, monitoring NIBP repeatedly will also shorten the operating time of the batteries.**
-

33.5 Recycling a Battery

When a battery has visual signs of damage, or no longer holds a charge, it should be replaced. Remove the old battery from the patient monitor and recycle it properly. To dispose of a battery, follow local laws for proper disposal.

WARNING

- **Do not disassemble batteries, or put them into fire, or cause them to short circuit. They may ignite, explode, or leak, causing personal injury.**
-

34 Care and Cleaning

Use only the substances approved by us and methods listed in this chapter to clean or disinfect your equipment. Warranty does not cover damage caused by unapproved substances or methods.

We make no claims regarding the efficacy of the listed chemicals or methods as a means for controlling infection. For the method to control infection, consult your hospital's Infection Control Officer or Epidemiologist.

In this chapter we only describe cleaning and disinfection of the main unit. For the cleaning and disinfection of other reusable accessories, refer to instructions for use of corresponding accessories.

34.1 General Points

Keep your equipment and accessories free of dust and dirt. To avoid damage to the equipment, follow these rules:

- Always dilute according to the manufacturer's instructions or use the lowest possible concentration.
- Do not immerse part of the equipment in liquid.
- Do not pour liquid onto the equipment or accessories.
- Do not allow liquid to enter the case.
- Never use abrasive materials (such as steel wool or silver polish), or erosive cleaners (such as acetone or acetone-based cleaners).



WARNING

- **Be sure to shut down the system and disconnect all power cables from the outlets before cleaning the equipment.**
-
-



CAUTION

- **If you spill liquid on the equipment or accessories, contact us or your service personnel.**
-
-

NOTE

- **To clean or disinfect reusable accessories, refer to the instructions delivered with the accessories.**
-
-

34.2 Cleaning

Your equipment should be cleaned on a regular basis. If there is heavy pollution or lots of dust and sand in your place, the equipment should be cleaned more frequently. Before cleaning the equipment, consult your hospital's regulations for cleaning the equipment.

Recommended cleaning agents are:

- sodium hypochlorite bleach (diluted)
- Hydrogen peroxide (3%)
- Ethanol (70%)
- Isopropanol (70%)

To clean your equipment, follow these rules:

1. Shut down the patient monitor and disconnect it from the power line.
2. Clean the display screen using a soft, clean cloth dampened with a glass cleaner.
3. Clean the exterior surface of the equipment using a soft cloth dampened with the cleaner.
4. Wipe off all the cleaning solution with a dry cloth after cleaning if necessary.
5. Dry your equipment in a ventilated, cool place.

34.3 Disinfecting

Disinfection may cause damage to the equipment and is therefore not recommended for this patient monitor unless otherwise indicated in your hospital's servicing schedule. Cleaning equipment before disinfecting is recommended.

The recommended disinfectants include: ethanol 70%, isopropanol 70%, Perform® classic concentrate OXY (KHSO₄ solution).



CAUTION

- **Never use EtO or formaldehyde for disinfection.**
-

35 Maintenance



WARNING

- **Failure on the part of the responsible individual hospital or institution employing the use of this equipment to implement a satisfactory maintenance schedule may cause undue equipment failure and possible health hazards.**
 - **The safety checks or maintenance involving any disassembly of the equipment should be performed by professional servicing personnel. Otherwise, undue equipment failure and possible health hazards could result.**
 - **If you discover a problem with any of the equipment, contact your service personnel or us.**
-

35.1 Safety Checks

Before every use, after your patient monitor has been used for 6 to 12 months, or whenever your patient monitor is repaired or upgraded, a thorough inspection should be performed by qualified service personnel to ensure the reliability.

Follow these guidelines when inspecting the equipment:

- Make sure that the environment and power supply meet the requirements.
- Inspect the equipment and its accessories for mechanical damage.
- Inspect all power cords for damage, and make sure that their insulation is in good condition.
- Make sure that only specified accessories are applied.
- Inspect if the alarm system functions correctly.
- Make sure that the recorder functions correctly and the recorder paper meets the requirements.
- Make sure that the batteries meet the performance requirements.
- Make sure that the patient monitor is in good working condition.
- Make sure that the grounding resistance and leakage current meet the requirement.

In case of any damage or abnormality, do not use the patient monitor. Contact the hospital's biomedical engineers or your service personnel immediately.

35.2 Maintenance and Testing Schedule

The following maintenance and tests, except for visual inspection, power on test, touchscreen calibration, and battery check, shall be carried out by the service personnel only. Contact your service personnel if any maintenance is required. Make sure to clean and disinfect the equipment before any test and maintenance.

Check/Maintenance Item		Recommended Frequency
Preventative Maintenance Tests		
Visual inspection		1. When first installed or reinstalled.
NIBP test	Pressure check	1. If the user suspects that the measurement is incorrect. 2. Following any repairs or replacement of relevant module. 3. At least once a year.
	Leakage test	
	Calibration	
Sidestream and Microstream CO ₂ tests	Leakage test	
	Performance test	
	Calibration	
AG tests	Leakage test	
	Performance test	
	Calibration	
Performance Tests		
ECG test and calibration	Performance test	1. If the user suspects that the measurement is incorrect. 2. Following any repairs or replacement of relevant module. 3. At least once every two years. Note: At least once a year is recommended for NIBP, CO ₂ and AG.
	Calibration	
Resp performance test		
SpO ₂ test		
NIBP test and calibration	Pressure check	
	Leakage test	
	Calibration	
Temp test		
IBP test and calibration	Performance test	
	Pressure calibration	
C.O. test		
Mainstream CO ₂ test and calibration		
Sidestream and Microstream CO ₂ tests and calibration	Leakage test	
	Performance test	
	Calibration	
AG test	Leakage test	

Check/Maintenance Item		Recommended Frequency
	Performance test	
	Calibration	
ICG test		
BIS test		
RM test		
CCO/SvO ₂ test	Interconnecting function	
	Output calibration	
PiCCO test		
ScvO ₂ test		
Nurse call relay performance test		
Analog output performance test		
Electrical Safety Tests		
Electrical safety tests		At least once every two years.
Other Tests		
Power on test		1. When first installed or reinstalled. 2. Following any maintenance or the replacement of any main unit parts.
Touchscreen calibration		1. When the touchscreen appears abnormal. 2. After the touchscreen is replaced.
Recorder check		Following any repair or replacement of the recorder.
Network print test		1. When first installed. 2. Whenever the printer is serviced or replaced.
Device integration check		1. When first installed. 2. Following any repair or replacement of the external device.
Battery check	Functionality test	1. When first installed. 2. Whenever a battery is replaced.
	Performance test	Once a year or if the battery run time reduced significantly.

35.3 Checking Monitor and Module Information

To view the information about system start time, selftest, etc., select [**Main Menu**]→[**Maintenance >>**]→[**Monitor Information >>**]. You can print out the information for the convenience of troubleshooting. The information will not be saved during shut down.

You can also view the information about the monitor configuration and system software version by selecting [**Main Menu**]→[**Maintenance >>**]→[**Software Version >>**].

35.4 Calibrating ECG

The ECG signal may be inaccurate due to hardware or software problems. As a result, the ECG wave amplitude becomes greater or smaller. In that case, you need to calibrate the ECG module.

1. Select the ECG parameter window or waveform area→[**Filter**]→[**Diagnostic**].
2. Select [**Main Menu**]→[**Maintenance >>**]→[**Calibrate ECG**]. A square wave appears on the screen and the message [**ECG Calibrating**] is displayed.
3. Compare the amplitude of the square wave with the wave scale. The difference should be within 5%.
4. After the calibration is completed, select [**Stop Calibrating ECG**]

You can print out the square wave and wave scale and then measure the difference between them if necessary. If the difference exceeds 5%, contact your service personnel.

35.5 NIBP Leakage Test

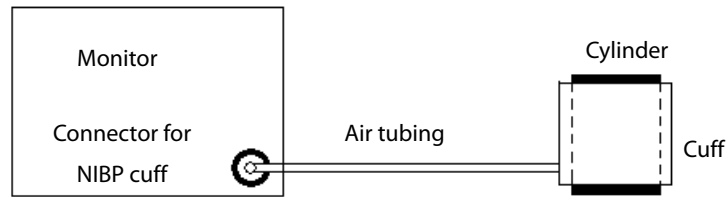
The NIBP leakage test checks the integrity of the system and of the valve. It is required at least once every two years or when you doubt the measured NIBP. If the test failed, corresponding prompt messages will be given. If no message is displayed, it means no leakage is detected.

Tools required:

- An adult cuff
- An air tubing
- A correct sized cylinder

Follow this procedure to perform the leakage test:

1. Set the patient category to [**Adu**].
2. Connect the cuff to the NIBP connector on the monitor.
3. Wrap the cuff around the cylinder as shown below.



4. Select [**Main Menu**]→[**Maintenance >>**]→[**NIBP Leakage Test**]. The NIBP display shows [**Leakage Testing...**].
5. After about 20 seconds, the monitor will automatically deflate. This means the test is completed.
6. If the message [**NIBP Pneumatic Leak**] is displayed, it indicates that the NIBP airway may have leakages. Check the tubing and connections for leakages. If you ensure that the tubing and connections are all correct, perform a leakage test again.

If the problem persists, contact your service personnel.

NOTE

-
- **The leakage test is intended for use to simply determine whether there are leakages in the NIBP airway. It is not the same as that specified in the EN 1060-3 standard.**
-

35.6 NIBP Accuracy Test

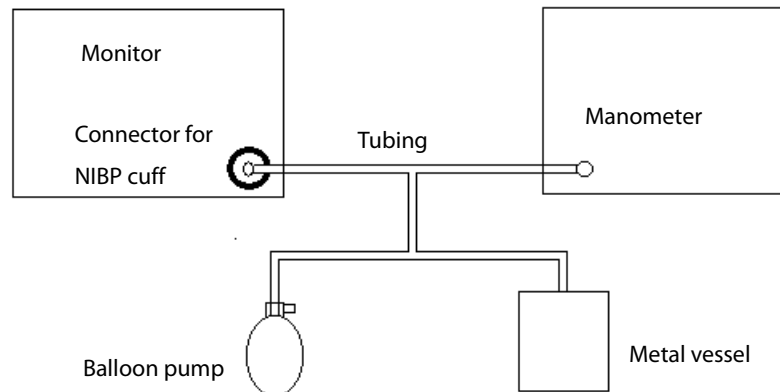
The NIBP accuracy test is required at least once every two years or when you doubt the measured NIBP.

Tools required:

- T-piece connector
- Appropriating tubing
- Balloon pump
- Metal Vessel (volume 500±25 ml)
- Reference manometer (calibrated with accuracy higher than 1 mmHg)

Follow this procedure to perform the accuracy test:

1. Connect the equipment as shown.



2. Before inflation, the reading of the manometer should be 0. If not, disconnect the airway and reconnect it until the readings is 0.
3. Select **[Main Menu]**→**[Maintenance >>]**→**[NIBP Accuracy Test]**.
4. Compare the manometer values with the displayed values. The difference between the manometer and displayed values should be no greater than 3 mmHg.
5. Raise the pressure in the metal vessel to 50 mmHg with the balloon pump. Repeat step 3 and 4.
6. Raise the pressure in the metal vessel to 200 mmHg with the balloon pump. Repeat step 3 and 4.

If the difference between the manometer and displayed values is greater than 3 mmHg, contact your service personnel.

35.7 Calibrating NIBP

NIBP is not user-calibrated. Cuff-pressure transducers must be verified and calibrated once every two years by a qualified service professional. Contact your service personnel when a calibration is necessary.

35.8 Calibrating CO₂

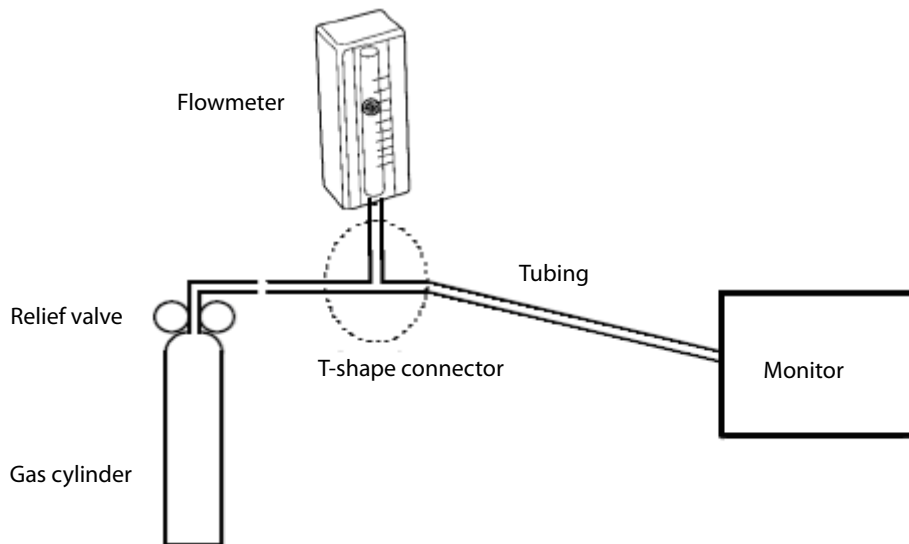
For sidestream and microstream CO₂ modules, a calibration is needed every year or when the measured values have a great deviation. For maintream CO₂ module, no calibration is needed. Calibration for sidestream CO₂ module can be performed only when the sidestream module enters the full accuracy mode.

Tools required:

- A steel gas cylinder with 6±0.05% CO₂ and balance gas N₂
- T-shape connector
- Tubing

Follow this procedure to perform a calibration:

1. Make sure that the sidestream or microstream CO₂ module has been warmed up or started up.
2. Check the airway for leakage and perform a leakage test as well to make sure the airway has no leakage.
3. Select [**Main Menu**] → [**Maintenance >>**] → [**User Maintenance >>**] → enter the required password → [**Maintain CO₂ >>**] → [**Calibrate CO₂ >>**].
4. In the [**Calibrate CO₂**] menu, select [**Zero**].
5. After the zero calibration is finished successfully, connect the equipment as follows:



6. Turn on and adjust the relief valve to make the flowmeter reads within 10-50mL/min and keeps stable as well.
7. In the [**Calibrate CO₂**] menu, enter the vented CO₂ concentration in the [**CO₂**] field.
8. In the [**Calibrate CO₂**] menu, the measured CO₂ concentration is displayed. After the measured CO₂ concentration becomes stable, select [**Calibrate CO₂**] to calibrate the CO₂ module.
9. If the calibration is finished successfully, the message [**Calibration Completed!**] is displayed in the [**Calibrate CO₂**] menu. If the calibration failed, the message [**Calibration Failed!**] is displayed. In this case, perform another calibration.

35.9 Calibrating AG

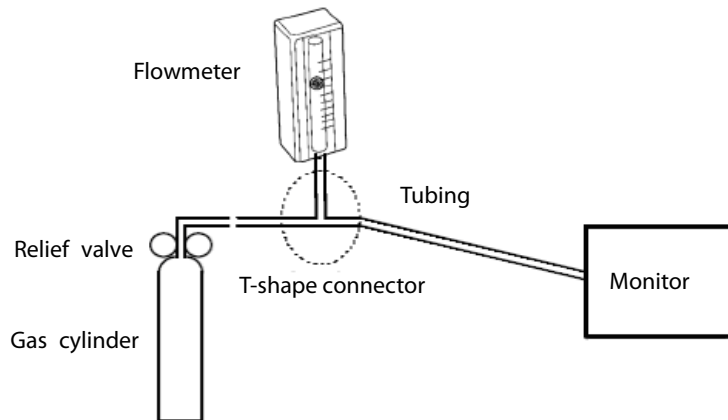
Calibrate the AG module every year or when the measured value has a great deviation.

Tools required:

- Gas bottle, with a certain standard gas or mixture gas. Gas concentration should meet the following requirements: AA>1.5%, CO₂>1.5%, N₂O>40%, O₂>40%, of which AA represents an anesthetic agent. $a/c \leq 0.01$ (a is the gas absolute concentration accuracy; c is the gas concentration)
- T-shape connector
- Tubing

Follow this procedure to perform a calibration:

1. Select [**Main Menu**]→[**Maintenance >>**]→[**User Maintenance >>**]→enter the required password→[**Calibrate AG >>**].
2. Check the airway and make sure that there are no occlusions or leaks.
 - ◆ Vent the tubing to the air and check if the [**Current FlowRate**] and [**Set FlowRate**] are approximately the same. If the deviation is great, it indicates that there is an occlusion in the tubing. Check the tubing for an occlusion.
 - ◆ Perform a leakage test to make sure that the airway has no leakage.
3. Connect the test system as follows:





4. Open the relief valve and vent a certain standard gas or gas mixture. Adjust the relief valve to make the flowmeter reads within 10-50mL/min and keeps stable as well.
5. In the [**Calibrate AG**] menu, the concentration and flowrate of each measured gas are displayed
 - ◆ If the difference between the measured gas concentration and the actual one is very small, a calibration is not needed.
 - ◆ If the difference is great, you should perform a calibration. Select [**Calibrate >>**] to enter the calibrate menu.
6. Enter the vented gas concentration. If you use only one gas for calibration, set other gases' concentration to 0.
7. Select [**Calibrate**] to start calibration.

8. If the calibration is finished successfully, the message [**Calibration Completed!**] is displayed. If the calibration failed, the message [**Calibration Failed!**] is displayed. Perform another calibration.

 **CAUTION**

-
- **If the O₂ module has been transported for long distance, calibrate it when installing the monitor.**
-

35.10 Calibrating the Touchscreen

1. Select [**Main Menu**]→[**Maintenance >>**]→[**User Maintenance >>**]→enter the required password→[**Cal. Touchscreen**].
2.  will, in turn, appear at different positions of the screen.
3. Select each  as it appears on the screen.
4. After the calibration is completed, the message [**Screen Calibration Completed!**] is displayed. Select [**Ok**] to confirm the completion of the calibration.

35.11 Electrical Safety Tests

Refer to **E Electrical Safety Inspection**.

35.12 Setting up IP Address

1. Select [**Main Menu**]→[**Maintenance >>**]→[**User Maintenance >>**]→enter the required password and then select [**Network Setup >>**] from the popup menu.
2. If your monitor is equipped with a wireless AP, you can set [**Network Type**] to [**WLAN**] in the network setup menu. Otherwise, the default setting is [**LAN**].
3. Set [**IP Address**].

If the patient monitor is connected to a CMS, its IP address should be set up. The user should not change the patient monitor's IP address randomly. If you want to know details about IP address setup, contact the technical personnel in charge of the CMS.

35.13 Entering/Exiting Demo Mode

To enter the Demo mode:

1. Select [**Main Menu**]→[**Maintenance >>**].
2. Select [**Demo >>**]. Enter the required password and then select [**Ok**].

To exit the Demo mode:

1. Select [**Main Menu**]→[**Maintenance >>**].
2. Select [**Exit Demo**] and then select [**Ok**].
3. The patient monitor exits the Demo mode.



WARNING

- **The Demo mode is for demonstration purpose only. To avoid that the simulated data are mistaken for the monitored patient's data, you must not change into Demo mode during monitoring. Otherwise, improper patient monitoring and delayed treatment could result.**
-
-

36 Accessories

WARNING

- Use accessories specified in this chapter. Using other accessories may cause damage to the patient monitor or not meet the claimed specifications.
- Single-use accessories are not designed to be reused. Reuse may cause a risk of contamination and affect the measurement accuracy.
- Check the accessories and their packages for any sign of damage. Do not use them if any damage is detected.

36.1 ECG Accessories

ECG Electrodes

Model	Quantity	Patient Category	Part No.
210	10 pieces	Adult	0010-10-12304
2245	50 pieces	Pediatric	9000-10-07469
2258-3	3 pieces	Neonate	900E-10-04880

12-Pin Trunk Cables

Leadwire supported	Compatible with	Type	Patient Category	Part No.
3-leadwire	AHA, IEC	Defibrillator-proof	Pediatric, neonate	0010-30-42720
3-leadwire	AHA, IEC	ESU-proof		0010-30-42724
3/5-leadwire	AHA, IEC	Defibrillator-proof	Adult, pediatric	0010-30-42719
3/5-leadwire	AHA, IEC	ESU-proof		0010-30-42723
10-leadwire	AHA	Defibrillator-proof		0010-30-42721
10-leadwire	IEC	Defibrillator-proof		0010-30-42722

Cable Sets

3-Electrode Cable Sets						
Type	Compatible with	Model	Patient Category	Part No.	Length	Remark
Clip	IEC	EL6304A	Adult, pediatric	0010-30-42732	1m	Long
		EL6306A	Neonate	0010-30-42897	1m	Long
	AHA	EL6303A	Adult, pediatric	0010-30-42731	1m	Long
		EL6305A	Neonate	0010-30-42896	1m	Long
Snap	IEC	EL6302B	Adult, pediatric	0010-30-42733	1m	Long
	AHA	EL6301B	Adult, pediatric	0010-30-42734	1m	Long

5-Electrode Cable Sets						
Type	Compatible with	Model	Patient Category	Part No.	Length	Remark
Clip	IEC	EL6502A	Adult, pediatric	0010-30-42728	0.6m	/
		EL6504A		0010-30-42730	1m to 1.4m	Long
	AHA	EL6501A		0010-30-42727	0.6m	/
		EL6503A		0010-30-42729	1m to 1.4m	Long
Snap	IEC	EL6502B		0010-30-42736	1.4m for F and N; 1m for others	Long
	AHA	EL6501B		0010-30-42735	1.4m for RL and LL; 1m for others	Long

10-Electrode Cable Sets						
Type	Compatible with	Model	Patient Category	Part No.	Length	Remark
Clip	IEC	EL6802A	Adult, pediatric	0010-30-42903	0.8m	Limb
		EL6804A		0010-30-42905	0.6m	Chest
	AHA	EL6801A		0010-30-42902	0.8m	Limb
		EL6803A		0010-30-42904	0.6m	Chest
Snap	IEC	EL6802B	Adult, pediatric	0010-30-42907	0.8m	Limb
		EL6804B		0010-30-42909	0.6m	Chest
	AHA	EL6801B		0010-30-42906	0.8m	Limb
		EL6803B		0010-30-42908	0.6m	Chest

36.2 SpO₂ Accessories

Extension Cable

Module type	Remarks	Part No.
Mindray SpO ₂ Module	/	0010-20-42710
Masimo SpO ₂ Module	8 pins, purple connector	040-000332-00
	7 pins, white connector	0010-30-42738
Nellcor SpO ₂ Module		0010-20-42712

SpO₂ Sensors

The SpO₂ sensor material that patients or other staff will come into contact with have undertaken the bio-compatibility test and is verified to be in compliance with ISO 10993-1.

Mindray SpO ₂ Module			
Type	Model	Patient Category	Part No.
Disposable	MAX-A	Adult (>30Kg)	0010-10-12202
	MAX-P	Pediatric (10 to 50Kg)	0010-10-12203
	MAX-I	Infant (3 to 20Kg)	0010-10-12204
	MAX-N	Neonate (<3Kg), Adult (>40Kg)	0010-10-12205
Single patient use	520A	Adult	520A-30-64101
	520P	Pediatric	520P-30-64201
	520I	Infant	520I-30-64301
	520N	Neonate	520N-30-64401
Reusable	DS-100A	Adult	9000-10-05161
	OXI-P/I	Pediatric, infant	9000-10-07308
	OXI-A/N	Adult, neonate	9000-10-07336
	518B	Adult, pediatric, neonate (Multi-sites)	518B-30-72107
	512E	Adult (Finger type)	512E-30-90390
	512F		512F-30-28263
	512G	Pediatric (Finger type)	512G-30-90607
	512H		512H-30-79061

Masimo SpO ₂ Module			
Type	Model	Patient Category	Part No.
Disposable	LNCS-NeoPt-L	Pediatric, neonate	0010-10-42626
	LNCS-Neo-L	Neonate	0010-10-42627
	LNCS-Inf-L	Infant	0010-10-42628
	LNCS-Pdt	Pediatric	0010-10-42629
	LNCS-Adt	Adult	0010-10-42630
Reusable	LNCS DC-I	Adult	0010-10-42600
	LNCS-DCIP	Pediatric	0010-10-42634
	LNCS YI	Adult, pediatric, neonate	0010-10-43016

Nellcor SpO ₂ Module			
Type	Model	Patient Category	Part No.
Disposable	MAX-A	Adult (>30Kg)	0010-10-12202
	MAX-P	Pediatric (10 to 50Kg)	0010-10-12203
	MAX-I	Infant (3 to 20Kg)	0010-10-12204
	MAX-N	Neonate (<3Kg), Adult (>40Kg)	0010-10-12205
Reusable	DS-100A	Adult	9000-10-05161
	OXI-P/I	Pediatric, infant	9000-10-07308
	OXI-A/N	Adult, neonate	9000-10-07336

- Wavelength emitted by the sensors intended for Mindray SpO₂ module: 512D: red light: 660 nm, infrared light: 940 nm; other SpO₂ sensors: red light: 660 nm, infrared light: 905 nm.
- Wavelength emitted by the sensors intended for Masimo SpO₂ module: red light: 660 nm, infrared light: 940 nm.
- Wavelength emitted by the sensors intended for Nellcor SpO₂ module: red light: 660 nm, infrared light: 890 nm.
- The maximum photic output consumption of the sensor is less than 18 mW.

The information about the wavelength range and maximum photic output consumption can be especially useful to clinicians, for example, clinicians performing photodynamic therapy.

36.3 NIBP Accessories

Tubing

Type	Patient Category	Part No.
Reusable	Adult, pediatric	6200-30-09688
	Neonate	6200-30-11560

Reusable Cuff

Model	Patient Category	Measurement Site	Limb Circumference (cm)	Bladder Width (cm)	Part No.
CM1201	Infant	Arm	10 to 19	9.2	0010-30-12157
CM1202	Pediatric		18 to 26	12.2	0010-30-12158
CM1203	Adult		24 to 35	15.1	0010-30-12159
CM1204	Large adult		33 to 47	18.3	0010-30-12160
CM1205	Thigh	Thigh	46 to 66	22.5	0010-30-12161
CM1300	Small infant	Arm	7 to 13	5.8	040-000968-00
CM1301	Infant		10 to 19	9.2	040-000973-00
CM1302	Pediatric		18 to 26	12.2	040-000978-00
CM1303	Adult		24 to 35	15.1	040-000983-00
CM1304	Large adult		33 to 47	18.3	040-000988-00
CM1305	Adult	Thigh	46 to 66	22.5	040-000993-00

Single-Patient Cuff

Model	Patient Category	Measurement Site	Limb Circumference (cm)	Bladder Width (cm)	Part No.
CM1500A	Neonate	Arm	3.1 to 5.7	2.2	001B-30-70692
CM1500B			4.3 to 8.0	2.9	001B-30-70693
CM1500C			5.8 to 10.9	3.8	001B-30-70694
CM1500D			7.1 to 13.1	4.8	001B-30-70695
CM1501	Infant		10 to 19	7.2	001B-30-70697
CM1502	Pediatric		18 to 26	9.8	001B-30-70698
CM1503	Adult		25 to 35	13.1	001B-30-70699
CM1504	Large adult		33 to 47	16.5	001B-30-70700
CM1505	Adult	Thigh	46 to 66	20.5	001B-30-70701

Disposable Cuff

Model	Patient Category	Measurement Site	Limb Circumference (cm)	Bladder Width (cm)	Part No.
M1872A	Neonate	Arm	7.1 to 13.1	5.1	900E-10-04873
M1870A			5.8 to 10.9	4.3	900E-10-04874
M1868A			4.3 to 8.0	3.2	900E-10-04875
M1866A			3.1 to 5.7	2.5	900E-10-04876

36.4 Temp Accessories

Extension Cable

Type	Model	Temp probe	Part No.
Reusable	MR420B	MR411, MR412	0011-30-37391

Temp Probes

Type	Model	Patient Category	Measurement Site	Part No.
Reusable	MR401B	Adult	Esophageal/Rectal	0011-30-37392
	MR403B		Skin	0011-30-37393
	MR402B	Pediatric, neonate	Esophageal/Rectal	0011-30-37394
	MR404B		Skin	0011-30-37395
Disposable	MR411	Adult, pediatric, neonate	Esophageal/Rectal	0011-30-37398
	MR412		Skin	0011-30-37397

36.5 IBP/ICP Accessories

Accessories Kit No.	Components	Part No.
6800-30-50876 (Hospira)	IM2201 12Pin IBP Cable	001C-30-70759
	Disposable Transducer	0010-10-42638
	Steady Rest for IBP Transducer and Clamp	M90-000133---
	Steady Rest for IBP Transducer and Clamp	M90-000134---
6800-30-50877 (BD)	IM2202 12Pin IBP Cable	001C-30-70757
	Disposable Pressure Transducer	6000-10-02107
	Transducer/Manifold Mount	0010-10-12156
ICP		
Model	Material	Part No.
Gaeltec TYPE.S13	12Pin ICP cable	0010-30-42742
Gaeltec ICT/B	Intracranial Pressure Transducer	0010-10-12151

It is proved through tests that the following accessories are compatible with the patient monitor. Only the accessories proceeded by "*" are available from our company. If you want to purchase other accessories, contact respective manufacturers and make sure if these accessories are approved for sale in local.

Manufacturer	Accessories
Smith Medical (Medex)	MX961Z14 Logical Cable, to be used in connection with the Adapter Cable (0010-20-42795) MX960 Reusable Transducer Kit MX9605A Logical 84in(213cm) Single Monitoring Kit MX960 Logical Transducer Mounting Plate MX261 Logical Clamp For Transducer Bracket MX262 Logical Clamp For 2 Transducer Mount Plates (More Logical Clamps are available from Medex. For detailed information, contact Medex.)
Braun	IBP Reusable Cable (REF: 5203511), to be used in connection with the Adapter Cable (0010-20-42795) Combitrans Monitoring Set (contact Braun for detailed information) Combitrans Attachment Plate Holder (REF:5215800) Combitrans Attachment Plate (contact Braun for detailed information)
Memscap	*Truck cable (0010-21-43082) SP844 Physiological Pressure Transducer 844-26 Monitoring Line Set 84X-49 Mounting Bracket
Utah	Reusable Blood Pressure Monitor Interface Cable (REF: 650-206) Deltran Disposable Pressure Transducer System (More Deltran sensors are available from Utah. For detailed information, contact Utah.) Pole Mount Unit (ERF: 650-150) Deltran Three Slot Organizer, Attaches to I.V. Pole Mount (REF: 650-100) Deltran Four Slot Organizer, Attaches to I.V. Pole Mount (REF: 650-105)

Edwards	<p>* IBP Truwave Reusable Cable (0010-21-12179)</p> <p>Pressure Monitoring Kit With Truwave Disposable Pressure Transducer.</p> <p>(More Truwave sensors are available from Edwards. For detailed information, contact Edwards.)</p> <p>DTSC IV Pole Clamp for Model DTH4 Backplate Holder</p> <p>DTH4 Disposable Holder for DPT</p>
---------	--

36.6 C.O. Accessories

Model	Material	Part No.
COC-001-SL	12Pin C.O. cable.	0010-30-42743
SP4042	TI Sensor	6000-10-02079
SP5045	TI Sensor Housing	6000-10-02080
MX387	12CC Control Syringe W/1CC Stop W/Rotator	6000-10-02081
131HF7	Dilution Hose	6000-10-02183
9850A	Cable kit with TI Sensor	0012-00-1519

36.7 CCO/SvO₂ Accessories

Material	PN
CCO/SvO ₂ cable	009-000259-00

36.8 CO₂ Accessories

Sidestream CO₂ module

Material	Patient Category	Remark	Part No.
DRYLINE Watertrap	Adult, pediatric	Reusable	9200-10-10530
DRYLINE Watertrap	Neonate		9200-10-10574
Sampling Line, Adult 2.5m	Adult, pediatric	Disposable	9200-10-10533
Sampling Line, Neonate, 2.5m	Neonate		9200-10-10555
Adult Nasal CO ₂ Sample Cannula	Adult		M02A-10-25937
Pediatric Nasal CO ₂ Sample Cannula	Pediatric		M02A-10-25938
Infant Nasal CO ₂ Sample Cannula	Infant		M02B-10-64509
DRYLINE Airway Adapter	Adult, pediatric	Straight	9000-10-07486

Microstream CO₂ Module

Disposable Airway Sampling Line			
Model	Patient Category	Remark	Part No.
XS04620	Adult, pediatric	/	0010-10-42560
XS04624		Humidified	0010-10-42561
007768		Long	0010-10-42563
007737		Long, humidified	0010-10-42564
006324	Infant, Neonate	Humidified	0010-10-42562
007738		Long, humidified	0010-10-42565

Disposable Nasal Sampling Line			
Model	Patient Category	Remark	Part No.
009818	Adult, intermediate	/	0010-10-42566
009822		Plus O ₂	0010-10-42568
009826		Long, plus O ₂	0010-10-42570
008174	Adult	/	0010-10-42577
008177		Humidified	0010-10-42572
008180		Humidified, plus O ₂	0010-10-42575
007266	Pediatric	/	0010-10-42567
008175		/	0010-10-42578
008178		Humidified	0010-10-42573
008181		Humidified, plus O ₂	0010-10-42576
007269		Plus O ₂	0010-10-42569
007743		Long, plus O ₂	0010-10-42571
008179	Infant, Neonate	Humidified	0010-10-42574

Mainstream CO₂ Module

Material	Model	Patient Category	Remark	Part No.
Airway adapter	6063	Adult	Disposable	0010-10-42662
	6421		Disposable, with mouthpiece	0010-10-42663
	6312	Neonate	Disposable	0010-10-42664
Mask	9960STD	Adult	/	0010-10-42670
	9960LGE		Adult large	0010-10-42671
	9960PED	Pediatric	/	0010-10-42669
Cable management straps	/	/	/	0010-10-42667
Sensor holding clips	/	/	/	0010-10-42668
Sensor	/	Adult, pediatric, neonate	Reusable	6800-30-50760

36.9 AG Accessories

Material	Patient Category	Remark	Part No.
Watertrap	Adult, pediatric	Reusable	9200-10-10530
	Neonate		9200-10-10574
Sampling line	Adult, pediatric	Disposable	9200-10-10533
	Neonate		9200-10-10555
Airway adapter	Adult, pediatric, neonate	Disposable, straight	9000-10-07486
	Adult, pediatric, neonate	Disposable, elbow	9000-10-07487

36.10 ICG Accessories

Material	Model	Part No.
BioZ tect ICG sensor	BZ-1550-50	0010-10-43258
BioZ Dx Patient Cable	5550	0010-10-42676
BioZ Dx Lead Wire Array	5561	0010-10-43259
BioZ Dx Patient Cable	5551	040-000543-00
BioZ Dx Lead Wire Array	5562	040-000544-00

36.11 BIS Accessories

Material	Patient Category	Part No.
BIS Cable	Adult, pediatric	6800-30-50761
BISx4 Cable	Adult, pediatric	115-005707-00

*If you need to purchase BIS Quatro, Pediatric, SRS, and CLIC sensors, please contact Covidien.

36.12 RM Accessories

Material	Patient Category	Remark	Part No.
Flow sensor	Adult, pediatric	Reusable	0010-30-42678
	Adult, pediatric	Disposable	0010-30-42679
	Infant	Disposable	0010-30-42680
RM connector	/	/	6800-20-50328

36.13 PiCCO Accessories

Material	Model	Part No.	Remark
12Pin IBP Y Cable	IM2203	040-000815-00	/
12Pin PiCCO Cable	CO7701	040-000816-00	/
2Pin Injectate Temperature Sensor Cable	040-000436-00	040-000817-00	/
Arterial Thermodilution Catheter	PV2015L20	/	Contact, germfree
	PV2013L07	/	Contact, germfree
PiCCO Monitoring Kits	PV8115	/	Contact, germfree

36.14 ScvO₂ Accessories

Material	Part No.	Remark
8Pin ScvO ₂ Module and Cable	115-008191-00	/
CeVOX Probe	/	Contact, germfree
	/	Contact, germfree

36.15 BeneLink Accessories

Material	Part No.
ID Adapter	115-008545-00
Serial port adapting cable, type A	009-001767-00
Serial port adapting cable, type B	009-001768-00
Serial port adapting cable, type C	009-001769-00
Serial port adapting cable, type D	009-002943-00
RJ45 connecting cable	009-001770-00

36.16 Others

Material	Part No.
Lithium battery	M05-010002-06
	022-000008-00
Power cord (India)	0000-10-10903
Domestic power cord (America)	DA8K-10-14452
Three-wire power cord (Britain)	DA8K-10-14453
Three-wire power cord (Europe)	DA8K-10-14454
Grounding cable	1000-21-00122
Defibrillator synchronization cable	6800-20-50781
Nurse call cable	8000-21-10361
Satellite module rack wall mount bracket	0010-30-42867
Keyboard wall mount bracket	0010-30-42868
Main unit wall mount bracket	0010-30-42955
Display wall mount bracket	0010-30-42956
Roll stand	0010-30-42943
Trolley-Mount Bracket	0010-30-42944
DVI-VGA adapter box	115-004861-00

A Product Specifications

NOTE

- For the specifications of BeneView T1, refer to *BeneView T1 Operating Manual*.

A.1 Monitor Safety Specifications

A.1.1 Classifications

The patient monitor is classified, according to IEC60601-1:

Components	Type of protection against electrical shock	Degree of protection against electrical shock	Degree of protection against harmful ingress of water	Degree of protection against hazards of explosion	Mode of operation
Main unit	I	Not marked	Ordinary	Not suitable	Continuous
Secondary display		Not marked			
MPM	NA	CF(*)			
IBP module					
SpO ₂ module					
C.O. module					
PiCCO module					
BIS module					
AG module		BF(*)			
CO ₂ module					
ICG module					
RM module					
ScvO ₂ module					
BeneLink module					
SMR	Not marked				
CCO/SvO ₂ module	Not marked				

- I: Class I equipment
- BF: Type BF applied part. (*Defibrillator-proof protection against electric shock.)
- CF: Type CF applied part. (*Defibrillator-proof protection against electric shock.)
- NA: Not applicable
- Ordinary: Ordinary equipment (enclosed equipment without protection against ingress of water)
- Not suitable: Equipment not suitable for use in the presence of a flammable anesthetic mixture with air with oxygen or nitrous oxide.

A.1.2 Environmental Specifications

Main unit, MPM, SpO₂ module, IBP module, C.O. module, Recorder, CCO/SvO₂ module, BIS module, BeneLink module		
Item	Operating conditions	Storage conditions
Temperature (°C)	0 to 40	-20 to 60
Relative humidity (noncondensing)	15% to 95%	10% to 95%
Barometric (mmHg)	427.5 to 805.5	120 to 805.5

Microstream CO₂ module		
Item	Operating conditions	Storage conditions
Temperature (°C)	0 to 40	-20 to 60
Relative humidity (noncondensing)	15% to 95%	10% to 95%
Barometric (mmHg)	430 to 790	430 to 790

Sidestream CO₂ module		
Item	Operating conditions	Storage conditions
Temperature (°C)	5 to 40	-20 to 60
Relative humidity (noncondensing)	15% to 95%	10% to 95%
Barometric (mmHg)	430 to 790	430 to 790

Mainstream CO₂ module		
Item	Operating conditions	Storage conditions
Temperature (°C)	0 to 40	-20 to 60
Relative humidity (noncondensing)	10% to 90%	10% to 90%
Barometric (mmHg)	427.5 to 805.5	400 to 805.5

AG module		
Item	Operating conditions	Storage conditions
Temperature (°C)	10 to 40	-20 to 60
Relative humidity (noncondensing)	15% to 95%	10% to 95%
Barometric (mmHg)	525 to 805.5	525 to 805.5

RM module		
Item	Operating conditions	Storage conditions
Temperature (°C)	5 to 40	-20 to 60
Relative humidity (noncondensing)	15% to 95%	10% to 95%
Barometric (mmHg)	427.5 to 805.5	120 to 805.5

ICG module		
Item	Operating conditions	Storage conditions
Temperature (°C)	10 to 40	0 to 50
Relative humidity (noncondensing)	15% to 95%	15% to 95%
Barometric (mmHg)	427.5 to 805.5	120 to 805.5

PiCCO module		
Item	Operating conditions	Storage conditions
Temperature (°C)	10 to 40	-20 to 60
Relative humidity (noncondensing)	15% to 75%	10% to 90%
Barometric (mmHg)	427.5 to 805.5	120 to 805.5

ScvO₂ module		
Item	Operating conditions	Storage conditions
Temperature (°C)	10 to 40	-20 to 60
Relative humidity (noncondensing)	15% to 75%	10% to 90%
Barometric (mmHg)	427.5 to 805.5	120 to 805.5

A.1.3 Power requirements

Line voltage	100 to 240 VAC
Current	BeneView T5: 2.5 to 1.4 A BeneView T8: 2.8 to 1.6 A
Frequency	50/60 Hz
Fuse	BeneView T5: Time-lag 250V T3.15A BeneView T8: Time-lag 250V T4A

A.2 Physical Specifications

Components	Weight	Size	Equipment type
Main unit (BeneView T5)	<6.6 kg	297×336×187 mm	Without modules, batteries, and recorder
Main unit (BeneView T8)	<9.9 kg	400×370×193 mm	Without modules, batteries, and recorder
SMR	<1.8 kg	142×402×151 mm	With no module inserted
MPM	<0.63 kg	136.5×80.5×102 mm	
SpO ₂ module	<0.26 kg	136.5×80.5×102 mm	
IBP module	<0.25 kg	136.5×40×102 mm	
C.O. module	<0.25 kg	136.5×40×102 mm	
Sidestream CO ₂ module	<0.48 kg	136.5×80.5×102 mm	
Microstream CO ₂ module	<0.37 kg	136.5×40×102 mm	
Mainstream CO ₂ module	<0.50 kg	136.5×40×102 mm	
M-type AG module	<1.75 kg	136.5×121×102 mm	With O ₂ and BIS modules
A-type AG module	<1.75 kg	136.5×121×102 mm	With O ₂ and BIS modules
ICG module	<0.35 kg	136.5×40×102 mm	
BIS module	<0.25 kg	136.5×40×102 mm	
RM module	<0.27 kg	136.5×40×102 mm	
CCO/SvO ₂ module	<0.25 kg	136.5×40×102 mm	
PiCCO Module	<0.28 kg	136.5×40×102 mm	
ScvO ₂ Module	<0.26 kg	136.5×40×102 mm	
BeneLink Module	<0.35kg	136.5×40×102 mm	

A.3 Hardware Specifications

A.3.1 Display

Host display	
Screen type	Color TFT LCD
Screen Size (diagonal)	12.1"(BeneView T5); 17"(BeneView T8)
Resolution	800×600 pixels(BeneView T5); 1280×1024 pixels(BeneView T8)
External display	
Screen type	Medical-grade TFT LCD
Screen Size	15", 17" or above
Resolution	800×600 pixels or above (BeneView T5); 1024×768 pixels or above (BeneView T8)
EMC	MPR II, CISPR 11B
Third certificate	UL, C-UL, TUV, CE, FCC

A.3.2 Recorder

Method	Thermal dot array
Horizontal resolution	16 dots/mm (25 mm/s paper speed)
Vertical resolution	8 dots/mm
Paper width	50 mm
Paper length	20 m
Paper speed	25 mm/s, 50 mm/s
Number of waveform channels	Maximum 3

A.3.3 Battery

Size	147.5×60.4×23.8 mm
Weight	350 g
Number of batteries	1 or 2 (BeneView T5); 2 (BeneView T8)
Battery Type	Chargeable Lithium-Ion
Voltage	11.1 VDC
Capacity	4500 mAh
Run time	BeneView T5: 330 minutes when powered by two new fully-charged batteries (25°C, ECG, SpO ₂ , Auto NIBP measurements at intervals of 15 minutes) BeneView T8: 120 minutes when powered by two new fully-charged batteries (25°C, ECG, SpO ₂ , Auto NIBP measurements at intervals of 15 minutes)
Charge time	nearly 5.5 h to 90% nearly 6 h to 100%
Shutdown delay	at least 5 min (after a low battery alarm first occurs)

A.3.4 LEDs

Alarm lamp	1 (two color coded: yellow and red)
Technical alarm lamp	1 (blue)
Power on LED	1 (green)
AC power LED	1 (green)
Battery LED	1 (green)

A.3.5 Audio Indicator

Speaker	Give alarm tones (45 to 85 dB), key tones, QRS tones; support PITCH TONE and multi-level tone modulation; alarm tones comply with IEC60601-1-8.
---------	---

A.3.6 Monitor Interface Specifications

Power	1 AC power input connector
Wired network	BeneView T5: 1 RJ45 connector, 100 Base-TX, IEEE 802.3 BeneView T8: 2 RJ45 connector, 100 Base-TX, IEEE 802.3
USB	BeneView T5: 4 connectors, USB 1.1; BeneView T8: up to 10 connectors, USB 1.1
SMR connector	1 connector, not standard USB
CF	50-pin CF revision 2.0 connector
Video interface	1 connector, standard DVI-D
Nurse call	1 connector, standard BNC
Equipotential Grounding Terminal	1
Micro-D connector	1 connector, It outputs ECG, IBP and defibrillator synchronization signals simultaneously
CIS connector (BeneView T5)	1 connector, for connecting the CIS box.

A.3.7 Outputs

Auxiliary Output	
Standard	Meets the requirements of IEC60601-1 for short-circuit protection and leakage current
ECG Analog Output	
Bandwidth (-3dB; reference frequency: 10Hz)	Diagnostic mode: 0.05 to 150 Hz Monitor mode: 0.5 to 40 Hz Surgical mode: 1 to 20 Hz ST mode: 0.05 to 40 Hz
QRS delay	≤25 ms (in diagnostic mode, and with Paced off)
Sensitivity	1V/mV ±5%
PACE rejection/enhancement	Pace enhancement Signal amplitude: $V_{oh} \geq 2.5V$ Pulse width: 10ms±5% Signal rising and falling time: ≤100µs
IBP Analog Output	
Bandwidth (-3dB; reference frequency:1Hz)	DC to 50 Hz
Max transmission delay	30 ms (with Notch off)
Sensitivity	1 V/100 mmHg ±5%
Nurse Call Signal	
Output mode	Relay
Electrical requirements	≤60W, ≤2A, ≤36VDC, ≤25VAC
Isolation voltage	1500 VAC
Contact type	Normally open or normally contact (optional)
Defib Sync Pulse	
Output impedance	≤100Ω
Max time delay	35 ms (R-wave peak to leading edge of pulse)
Amplitude	High level: 3.5 to 5 V, providing a maximum of 10 mA output current; Low level: < 0.5 V, receiving a maximum of 5 mA input current.
Pulse width	100 ms ±10%
Rising and falling time	≤1 ms
Digital video output (DVI-D connector)	
Video signals	Single Link TMDS
DDC signals	Signals 12C compliant
Alarm output (Network connector)	
Alarm delay time from BeneView patient monitor to remote equipment	The alarm delay time from the patient monitor to remote equipment is ≤2 seconds, measured at the BeneView signal output connector.

A.4 Data Storage

Trends	Trends: 120 hours, at 1 min resolution Mid-length trends: 8 hours, at 5 s resolution Minitrends: 1 hour, at 1 s resolution
Parameter alarms	100 alarms and manual events and related parameter waveforms. The waveform recording length can be 8s.
Arrh. events	100 arrhythmia events and relate waveforms and parameters. The waveform recording length can be 8s.
NIBP measurements	1000 sets
Interpretation of resting 12-lead ECG results	20 sets
Full-disclosure waveforms	48 hours at maximum. The specific storage time depends on the waveforms stored and the number of stored waveforms.

A.5 Wireless Network

Standards	IEEE 802.11g, Wi-Fi compatible						
Frequency range	2.412 to 2.462GHz						
Operating channel	China	America	Canada	Europe	Spain	France	Japan
	1 to 11				10, 11		2
	For other country, please refer to your local law.						
Safe distance	a circle centering AP with the radius of 10 m						

A.6 Measurement Specifications

The adjustable range of alarm limits is the same with the measurement range of signals unless otherwise specified.

A.6.1 ECG

ECG	
Standards	Meet standards of EC11, EC13, EN60601-2-27/IEC60601-2-27 and IEC60601-2-25
Lead set	3-lead: I, II, III 5-lead: I, II, III, aVR, aVL, aVF, V 12-lead: I, II, III, aVR, aVL, aVF, V1 to V6
ECG standard	AHA, IEC
Display sensitivity	1.25 mm/mV (X0.125), 2.5 mm/mV (X0.25), 5 mm/mV (X0.5), 10 mm/mV (X1), 20 mm/mV (X2), 40 mm/mV (X4), Auto
Sweep speed	6.25 mm/s, 12.5 mm/s, 25 mm/s, 50 mm/s
Bandwidth (-3dB)	Diagnostic mode: 0.05 to 150 Hz Monitor mode: 0.5 to 40 Hz Surgical mode: 1 to 20 Hz ST mode: 0.05 to 40 Hz
Common mode rejection ratio (with Notch off)	Diagnostic mode: >90 dB Monitor mode: >105 dB Surgical mode: >105 dB ST mode: >105 dB(with Notch on)
Notch	50/60 Hz Monitor and surgical mode: Notch turns on automatically. Diagnostic mode: Notch is turned on/off manually
Differential input impedance	≥5 MΩ
Input signal range	±8 mV (peak-to-peak value)
Accuracy of reappearing input signal	Use A and D methods based on EC11 to determine system total error and frequency response.
Electrode offset potential tolerance	±500 mV
Lead-off detection current	Measuring electrode: <0.1 μA Drive electrode: <1 μA
Input offset current	≤0.1 μA
Baseline recovery time	<5 s (after defibrillation)
Patient leakage current	<10 uA
Calibration signal	1mV (peak-to-peak value)
ESU protection	Cut mode: 300 W Coagulate mode: 100 W Recovery time: ≤10 s In compliance with the requirements in clause 4.2.9.14 of ANSI/AAMI EC 13:2002
ESU noise suppression	Based on the test method in clause 5.2.9.14 of EC 13, use ECG lead wires which are in compliance with AAMI. Compared with ECG baseline, the noise of peak to peak value ≤2 mV.
Pace Pulse	

Pace pulse markers	Pace pulses meeting the following conditions are labelled with a PACE marker: Amplitude: ± 2 to ± 700 mV Width: 0.1 to 2 ms Rise time: 10 to 100 μ s	
Pace pulse rejection	When tested in accordance with the ANSI/AAMI EC13-2002: Sections 4.1.4.1 and 4.1.4.3, the heart rate meter rejects all pulses meeting the following conditions.	
	Amplitude:	± 2 to ± 700 mV
	Width:	0.1 to 2 ms
Pacer pulse detector rejection of fast ECG signals	Rise time:	10 to 100 μ s
Pacer pulse detector rejection of fast ECG signals	20V/s RTI when measured in accordance with ANSI/AAMI EC13-2002 Section 4.1.4.3.	

Mindray algorithm

HR		
Measurement range	3-, 5-, and 12-lead ECG	Neonate: 15 to 350 bpm Pediatric: 15 to 350 bpm Adult: 15 to 300 bpm
Resolution	1 bpm	
Accuracy	3-, 5-, and 12-lead ECG: ± 1 bpm or $\pm 1\%$, whichever is greater.	
Sensitivity	200 μ V (lead II)	
HR averaging method	In compliance with the requirements in Clause 4.1.2.1 d) of ANSI/AAMI EC13-2002, the following method is used: If the last 3 consecutive RR intervals are greater than 1200 ms, the 4 most recent RR intervals are averaged to compute the HR. Otherwise, heart rate is computed by subtracting the maximum and minimum ones from the most recent 12 RR intervals and then averaging them. The HR value displayed on the monitor screen is updated every second.	
Response to irregular rhythm	In compliance with the requirements in Clause 4.1.2.1 e) of ANSI/AAMI EC13-2002, the heart rate after 20 seconds of stabilization is displayed as follows: Ventricular bigeminy (3a): -80 ± 1 bpm Slow alternating ventricular bigeminy (3b): -60 ± 1 bpm Rapid alternating ventricular bigeminy (3c): -120 ± 1 bpm Bidirectional systoles (3d): -90 ± 2 bpm	
Response time to heart rate change	Meets the requirements of ANSI/AAMI EC13-2002: Section 4.1.2.1 f). From 80 to 120 bpm: less than 11 s From 80 to 40 bpm: less than 11 s	
Time to alarm for tachycardia (not available in USA)	Meets the requirements of ANSI/AAMI EC13-2002: section 4.1.2.1 g).	
	Waveform	
	4ah - range:	11 s
	4a - range:	11 s
	4ad - range:	11 s
	Waveform 4bh - range:	11 s
	4b - range:	11 s
4bd - range:	11 s	

Tall T-wave rejection capability	When the test is performed based on part 4.1.2.1 c) of ANSI/AAMI EC 13-2002, the heart rate meter will reject all 100 ms QRS complexes with less than 1.2 mV of amplitude, and T waves with T-wave interval of 180 ms and those with Q-T interval of 350 ms.
Arrhythmia Analysis Classifications (Not available in USA)	Asystole, VFib/VTac, Vtac, Vent. Brady, Extreme Tachy, Extreme Brady, PVC, Couplet, Bigeminy, Trigeminy, R on T, VT>2, PVCs, Tachy, Brady, Missed Beats, Vent. Rhythm, PNP, PNC, Multif. PVC, Nonsus. Vtac, Pause, Irr. Rhythm
ST Segment Analysis (Not available in USA)	
Measurement range	-2.0 to 2.0 mV
Accuracy	-0.8 to 0.8 mV: ± 0.02 mV or $\pm 10\%$, whichever is greater. Beyond this range: Not specified.
Refreshing rate	10 s

Mortara algorithm

Only the differences from the Mindray algorithm are listed.

HR	
HR averaging method	In compliance with the requirements in Clause 4.1.2.1 d) of ANSI/AAMI EC13-2002, the following method is used: Heart rate is computed by averaging the most recent 16 RR intervals, unless the HR by averaging the most recent 4 heart beats is less than or equals to 48. The HR value displayed on the monitor screen is updated every second.
Time to alarm for tachycardia	Meets the requirements of ANSI/AAMI EC13-2002: section 4.1.2.1 g). Waveform 4ah – range: 11 s 4a – range: 11 s 4ad – range: 11 s 4bh – range: 11 s 4b – range: 11 s 4bd – range: 11 s
Arrhythmia Analysis Classifications	Asystole, Vfib, Vtac, Vent. Rhythm, Couplet, VT>2, Bigeminy, Trigeminy, R on T, Multif. PVC, Irr. Rhythm, Tachy, Brady, Missed Beats, PNP, PNC
ST Segment Analysis	
Refreshing rate	per 16 heartbeats

A.6.2 Resp

Technique	Trans-thoracic impedance	
Lead	Options are lead I and II. The default is lead II.	
Respiration excitation waveform	<300 μ A RMS, 62.8 kHz (\pm 10%)	
Respiration impedance range	0.3 to 5 Ω	
Baseline impedance range	200 to 2500 Ω (using an ECG cable with 1k Ω resistance)	
Differential input impedance	>2.5 M Ω	
Bandwidth	0.2 to 2 Hz (-3 dB)	
Sweep speed	6.25 mm/s, 12.5 mm/s or 25 mm/s	
Respiration Rate		
Measurement range	Adult:	0 to 120 rpm
	Pediatric, neonate:	0 to 150 rpm
Resolution	1 rpm	
Accuracy	7 to 150 rpm:	\pm 2 rpm or \pm 2%, whichever is greater
	0 to 6 rpm:	Not specified.
Apnea alarm time	10 s, 15 s, 20 s, 25 s, 30 s, 35 s, 40 s	
Alarm limit	Range (rpm)	Step (rpm)
RR High	Adult, pediatric: (low limit + 2) to 100 Neonate: (low limit + 2) to 150	1
RR Low	0 to (high limit - 2)	

A.6.3 SpO₂

Alarm limit	Range (%)	Step (%)
SpO ₂ High	(low limit + 2) to 100	1
SpO ₂ Low	Desat to (high limit - 2)	
Desat	Mindray, Masimo: 0 to (high limit - 2) Nellcor: 20 to (high limit - 2)	

Mindray SpO₂ Module

Standards	Meet standards of ISO9919
*Measurement accuracy verification: The SpO ₂ accuracy has been verified in human experiments by comparing with arterial blood sample reference measured with a CO-oximeter. Pulse oximeter measurement are statistically distributed and about two-thirds of the measurements are expected to come within the specified accuracy range compared to CO-oximeter measurements.	
Measurement range	0 to 100%
Resolution	1%
Accuracy	70 to 100%: \pm 2% (measured without motion in adult/pediatric mode) 70 to 100%: \pm 3% (measured without motion in neonate mode) 70 to 100%: \pm 3% (measured with motion) 0% to 69%: Not specified.

*Studies were performed to validate the accuracy of Pulse Oximeter with neonatal SpO₂ sensors by contrast with a CO-Oximeter. Some neonates aged from 1 day to 30 days with a gestation age of 22 weeks to full term were involved in this study. The statistical analysis of data of this study shows the accuracy (Arms) is within the stated accuracy specification. Please see the following table.

Sensor type	Totally neonates	Data	Arms
518B	97 (51 male & 46 female)	200 pairs	2.38%
520N	122 (65 male & 57 female)	200 pairs	2.88%
The Pulse Oximeter with neonatal SpO ₂ sensors was also validated on adult subjects.			
Refreshing rate	1 s		
SpO ₂ averaging time	7 s (When the sensitivity is set to High) 9 s (When the sensitivity is set to Medium) 11 s (When the sensitivity is set to Low)		

Masimo SpO₂ Module

SpO ₂	
Measurement range	1 to 100%
Resolution	1%
Accuracy	70 to 100%: ±2% (measured without motion in adult/pediatric mode) 70 to 100%: ±3% (measured without motion in neonate mode) 70 to 100%: ±3% (measured with motion) 0% to 69%: Not specified.
Refreshing rate	1 s
SpO ₂ averaging time	2-4 s, 4-6 s, 8 s, 10 s, 12 s, 14 s, 16 s
Low perfusion conditions	Pulse amplitude: >0.02% Light penetration: >5%
Low perfusion SpO ₂ accuracy	±2%

Nellcor SpO₂ Module

Measurement range	0 to 100%
Resolution	1%
Accuracy	70 to 100%: ±2% (adult/pediatric) 70 to 100%: ±3% (neonate) 0% to 69%: Not specified.
*: When the SpO ₂ sensor is applied for neonatal patients as indicated, the specified accuracy range is increased by ±1%, to compensate for the theoretical effect on oximeter measurements of fetal hemoglobin in neonatal blood.	

A.6.4 PR

Alarm limit	Range (bpm)	Step (bpm)
PR High	(low limit +2) to 300	1
PR Low	15 to (high limit-2)	

PR from Mindray SpO₂ Module

Measurement range	20 to 254 bpm
Resolution	1 bpm
Accuracy	±3 bpm (measured without motion) ±5 bpm (measured with motion)
Refreshing rate	1 s
SPO ₂ averaging time	7 s (when sensitivity is set to High) 9 s (when sensitivity is set to Medium) 11 s (when sensitivity is set to Low)

PR from Masimo SpO₂ Module

Measurement range	25 to 240 bpm
Resolution	1 bpm
Accuracy	±3 bpm (measured without motion) ±5 bpm (measured with motion)
Refreshing rate	1 s
SPO ₂ averaging time	2-4 s, 4-6 s, 8 s, 10 s, 12 s, 14 s, 16 s
Low perfusion conditions	Pulse amplitude: >0.02% Light penetration: >5%
Low perfusion PR accuracy	±3 bpm

PR from Nellcor SpO₂ Module

Measurement range	20 to 300 bpm
Resolution	1 bpm
Accuracy	20 to 250 bpm: ±3 bpm 251 to 300 bpm, not specified
Refreshing rate	1 s

PR from IBP Module

Measurement range	25 to 350 bpm
Resolution	1 bpm
Accuracy	±1 bpm or ±1%, whichever is greater
Refreshing rate	1 s

A.6.5 NIBP

Standards	Meet standards of EN60601-2-30/IEC60601-2-30, EN1060-1, EN1060-3,EN1060-4 and SP10			
Technique	Oscillometry			
Mode of operation	Manual, Auto and STAT			
Auto mode repetition intervals	1, 2, 2.5, 3, 5, 10, 15, 20, 30, 60, 90, 120, 180, 240 or 480 min			
STAT mode cycle time	5 min			
Max measurement time	Adult, pediatric: 180 s Neonate: 90 s			
Heart rate range	40 to 240 bpm			
Measurement ranges (mmHg)		Adult	Pediatric	Neonate
	Systolic:	40 to 270	40 to 200	40 to 135
	Diastolic:	10 to 210	10 to 150	10 to 100
	Mean:	20 to 230	20 to 165	20 to 110
Accuracy	Max mean error: ± 5 mmHg Max standard deviation: 8 mmHg			
Resolution	1mmHg			
Initial cuff inflation pressure range (mmHg)	Adult: 80 to 280 Pediatric: 80 to 210 Neonate: 60 to 140			
Default initial cuff inflation pressure (mmHg)	Adult: 160 Pediatric: 140 Neonate: 90			
Software overpressure protection	Adult: 297 ± 3 mmHg Pediatric: 240 ± 3 mmHg Neonate: 147 ± 3 mmHg			
PR				
Measurement range	40 to 240 bpm			
Resolution	1 bpm			
Accuracy	± 3 bpm or $\pm 3\%$, whichever is greater			
Alarm limit	Range (mmHg)			Step (mmHg)
Sys High	Adult: (low limit+5) to 270 Pediatric: (low limit+5) to 200 Neonate: (low limit+5) to 135			5
Sys Low	40 to (high limit-5)			
Mean High	Adult: (low limit+5) to 230 Pediatric: (low limit+5) to 165 Neonate: (low limit+5) to 110			
Mean Low	20 to (high limit-5)			
Dia High	Adult: (low limit+5) to 210 Pediatric: (low limit+5) to 150 Neonate: (low limit+5) to 100			
Dia Low	10 to (high limit-5)			

*Measurement accuracy verification: In adult and pediatric modes, the blood pressure measurements measured with this device are in compliance with the American National Standard for Electronic or Automated Sphygmomanometers (ANSI/AAMI SP10-1992) in terms of mean error and standard deviation by comparing with intra-arterial or auscultatory measurements (depending on the configuration) in a typical patient population. For auscultatory reference, the 5th Korotkoff sound was used to determine the diastolic pressure.

In neonatal mode, the blood pressure measurements measured with this device are in compliance with the American National Standard for Electronic or Automated Sphygmomanometers (ANSI/AAMI SP10-1992 and AAMI/ANSI SP10A-1996) in terms of mean error and standard deviation by comparing with intra-arterial measurements (depending on the configuration) in a typical patient population.

A.6.6 Temp

Standards	Meet standard of EN12470-4	
Technique	Thermal resistance	
Measurement range	0 to 50 °C (32 to 122 °F)	
Resolution	0.1 °C	
Accuracy	±0.1 °C or ±0.2 °F (without probe)	
Refreshing rate	1 s	
Minimum time for accurate measurement	Body surface: <100 s Body cavity: <80 s	
Alarm limit	Range	Step
T1/T2 High	(low limit +1) to 50 °C (low limit +1.8) to 122 °F	0.1 °C 0.1 °F
T1/T2 Low	0 to (high limit -1) °C 32 to (high limit -1.8) °F	
TD High	0 to 50 °C 0 to 90 °F	

A.6.7 IBP

Standards	Meet standard of EN60601-2-34/IEC60601-2-34.
Technique	Direct invasive measurement
IBP	
Measurement range	-50 to 300 mmHg
Resolution	1 mmHg
Accuracy	±2% or ±1 mmHg, whichever is greater (without sensor)
Refreshing rate	1 s
Pressure transducer	
Excitement voltage	5 VDC, ±2%
Sensitivity	5 µV/V/mmHg
Impedance range	300 to 3000Ω
Volume displacement (ABBOTT)	<0.04 mm ³ /100 mmHg

Alarm limit	Range (mmHg)	Step (mmHg)
Sys High	(low limit + 2) to 300	1
Mean High		
Dia High		
Sys Low	-50 to (high limit - 2)	
Mean Low		
Dia Low		

A.6.8 C.O.

Measurement method	Thermodilution method	
Measurement range	C.O.:	0.1 to 20 L/min
	TB:	23 to 43 °C
	TI:	0 to 27 °C
Resolution	C.O.:	0.1 L/min
	TB, TI:	0.1 °C
Accuracy	C.O.:	±5% or ±0.1 L/min, whichever is greater
	TB, TI:	±0.1 °C (without sensor)
Repeatability	C.O.:	±2% or ±0.1 L/min, whichever is greater
Alarm range	TB:	23 to 43 °C
Alarm limit	Range	Step
TB High	(low limit + 1) to 43 °C (low limit + 1.8) to 109.4 °F	0.1 °C
TB Low	23 to (high limit - 1) °C 73.4 to (high limit - 1.8) °F	0.1 °F

A.6.9 CCO

Operating mode	Interfaces with Edwards Vigilance II® or Vigileo™ monitor
Measured parameter	Consistent with CCO-related parameters outputted by Vigilance II® or Vigileo™ monitor
Parameter alarm	Vigilance II®: CCO/CCI,EDV/EDVI,SVR/SVRI,SV/SVI,RVEF
	Vigileo™: CCO/CCI, SV/SVI, SVV
Signal Outputs	
Standard	Meets the requirements of EN 60601-1 for short-circuit protection and leakage current
Output impedance	1000Ω
Isolation voltage	1500 VAC
ECG Analog Output	
Bandwidth (-3dB; reference frequency: 10Hz)	ST mode: 0.05~40Hz Diagnostic mode: 0.05~150Hz Monitor mode: 0.5~40Hz Surgical mode: 1~20Hz

Sensitivity	2V/mV \pm 5%
MAP Analog Signal Output	
Output voltage	0 to 5V (0 to 500mmHg)
Output voltage error	\pm 5%
CVP Analog Signal Output	
Output voltage	0 to 5V (0 to 100mmHg)
Output voltage error	\pm 5%

CCO-related Parameters Outputted by Vigilance II® Monitor		
Name	Range	Resolution
CCO	1 to 20 L/min	0.1
CCI	0 to 20 L/min/m ²	0.1
CO	1 to 20 L/min	0.1
CI	0 to 20 L/min/m ²	0.1
EDV	40 to 800 ml	1
EDVI	20 to 400 ml/m ²	1
SVR	0 to 3000 DS/cm ⁵	1
SVRI	0 to 6000 DS·m ² /cm ⁵	1
SV	0 to 300 ml/beat	1
SVI	0 to 200 ml/beat/m ²	1
BT	25 to 45 °C	0.1
RVEF	10 to 60%	1
ESV	10 to 700 ml	1
ESVI	5 to 400 ml/m ²	1
HRavg	30 to 250 bpm	1
CVP	0 to 100 mmHg	1
MAP	0 to 500 mmHg	1

CCO-related Parameters Outputted by Vigileo™ Monitor		
Name	Range	Resolution
CCO	1 to 20 L/min	0.1
CCI	0 to 20 L/min/m ²	0.1
SVR	0 to 3000 DS/cm ⁵	1
SVRI	0 to 6000 DS·m ² /cm ⁵	1
SV	0 to 300 ml/beat	1
SVI	0 to 200 ml/beat/m ²	1
SVV	0 to 99%	0.1
CVP	0 to 100 mmHg	1

Alarm Limit	Range	Step
CCO High	(Low limit+0.1) to 20 L/min	0.1 L/min
CCO Low	0 to (high limit-0.1)L/min	
CCI High	(Low limit+0.1) to 20 L/min/m ²	0.1 L/min/m ²
CCI Low	0 to (high limit-0.1)L/min/m ²	

Alarm Limit	Range	Step
EDV High	(Low limit+10)to 800 ml	10 ml
EDV Low	0 to (high limit-10)ml	
EDVI High	(Low limit+10) to 400 ml/m ²	10 ml/m ²
EDVI Low	0 to (high limit-10)ml/m ²	
SVR High	(Low limit+20) to 5000 DS/cm ⁵ or (low limit+2) to 500 kPa-s/l	20 DS/cm ⁵ or 2 kPa-s/l
SVR Low	0 to (high limit-20)DS/cm ⁵ or 0 to (high limit-2)kPa-s/l	
SVRI High	(Low limit+50) to 9950 DS-m ² /cm ⁵ or(low limit+5) to 995 kPa-s-m ² /l	50 DS-m ² /cm ⁵ or 5 kPa-s-m ² /l
SVRI Low	0 to (high limit-50)DS-m ² /cm ⁵ or 0 to (high limit-5)kPa-s-m ² /l	
SV High	(Low limit+5) to 300 ml/b	5 ml/b
SV Low	0 to (high limit-5)ml/b	
SVI High	(Low limit+5) to 200 ml/b/m ²	5 ml/b/m ²
SVI Low	0 to (high limit-5)ml/b/m ²	
RVEF High	(Low limit+5) to 100 %	5 %
RVEF Low	0 to (High limit-5)%	

A.6.10 SvO₂

Operating mode	Interfaces with Edwards Vigilance II® or Vigileo™ monitor
Measured parameter	Consistent with CCO-related parameters outputted by Vigilance II® or Vigileo™ monitor
Parameter alarm	SvO ₂ , ScvO ₂
Signal Output	
Standard	Meets the requirements of EN 60601-1 for short-circuit protection and leakage current
Output impedance	1000Ω
Isolation voltage	1500 VAC
SpO₂ Analog Signal Output	
Output voltage	0 to 10V (0 to 100%)
Output voltage error	±5%

SvO ₂ -related Parameters Outputted by Vigilance II® Monitor		
Name	Measurement Range	Resolution
SaO ₂	40 to 100%	1
VO ₂	0 to 999 ml/min	1
O ₂ El	0.0 to 99.9%	0.1
SNR	-10 to +20 dB	0.1
DO ₂	0 to 2000 ml/min	1

SvO ₂	0 to 99%	1
ScvO ₂	0 to 99%	1
SQI	1 to 4	1

A.6.11 PiCCO

Measured parameters	Measurement range	Coefficient of variation
CCO	0.25 l/min to 25.0 l/min	≤2%
C.O.	0.25 l/min to 25.0 l/min	≤2%
GEDV	40ml to 4800 ml	≤3%
SV	1ml to 250 ml	≤2%
EVLW	10ml to 5000 ml	≤6%
ITBV	50ml to 6000 ml	≤3%

* Coefficient of variation is measured using synthetic and/or database wave forms (laboratory testing). Coefficient of variation= SD/mean error.

A.6.12 ScvO₂

Measured parameters	Measurement range	Measurement accuracy
ScvO ₂	0 to 99%	50% to 80%: ±3% Other ranges: Not specified.

A.6.13 CO₂

Measurement mode	Sidestream, microstream, mainstream	
Technique	Infrared absorption	
Apnea time	10 s, 15 s, 20 s, 25 s, 30 s, 35 s, 40 s	
Alarm limit	Range	Step
EtCO ₂ High	(low limit + 2) to 99 mmHg	1 mmHg
EtCO ₂ Low	1 to (high limit - 2)mmHg	
FiCO ₂ High	1 to 99 mmHg	
awRR High	Adult, pediatric: (low limit + 2) to 100 rpm Neonate: (low limit + 2) to 150 rpm	1 rpm
awRR Low	0 to (high limit - 2) rpm	

Sidestream CO₂ Module

Standard	Meet standard of ISO 21647
CO ₂ Measurement range	0 to 99 mmHg
Accuracy*	0 to 40 mmHg: ±2 mmHg 41 to 76 mmHg: ±5% of the reading 77 to 99 mmHg: ±10% of the reading
Accuracy drift	Meet the requirement for measurement accuracy within 6 hours
Resolution	1 mmHg

Sample flowrate	Adult: 70 ml/min, 100 ml/min, 120 ml/min, 150 ml/min Pediatric, neonate: 70 ml/min, 100 ml/min	
Sample flowrate tolerance	15% or 15 ml/min, whichever is greater.	
Warm-up time	<1 min, enter the iso accuracy mode After 1 min, enters the full accuracy mode,	
Response time	Measured with a neonatal watertrap and a 2.5-meter neonatal sampling line: <3.5 s @ 100 ml/min <4 s @ 70 ml/min Measured with a adult watertrap and a 2.5-meter adult sampling line: <4.5 s @ 150 ml/min <5.5 s @ 120 ml/min <5.5 s @ 100 ml/min <7 s @ 70 ml/min	
Gas sampling delay time	Measured with a neonatal watertrap and a 2.5-meter neonatal sampling line: <3 s @ 100 ml/min <3.5 s @ 70 ml/min Measured with a adult watertrap and a 2.5-meter adult sampling line: <4 s @ 150 ml/min <5 s @ 120 ml/min <5 s @ 100 ml/min <6.5 s @ 70 ml/min	
awRR measurement range	0 to 120 rpm	
awRR measurement precision	±2 rpm	
Effect of interference gases on CO₂ measurements		
Gas	Concentration (%)	Quantitive effect*
N ₂ O	≤60	±1 mmHg
Hal	≤4	
Sev	≤5	
Iso	≤5	
Enf	≤5	
Des	≤15	±2 mmHg
*: means an extra error should be added in case of gas interference when CO ₂ measurements are performed between 0-40mmHg.		

Microstream CO₂ Module

Standard	Meet standard of ISO 21647
CO ₂ Measurement range	0 to 99 mmHg
Accuracy*	0 to 38 mmHg: ±2 mmHg 39 to 99 mmHg: ±5% of the reading+0.08% of (the reading-38)
Accuracy drift	Meet the requirement for measurement accuracy within 6 hours
<p>* Accuracy applies for respiration rate up to 80 rpm. For respiration rate above 80 rpm, the accuracy is 4 mmHg or ±12% of the reading, whichever is greater. for EtCO₂ exceeding 18 mmHg. For respiration rate above 60 rpm, the above accuracy can be achieved by using the CapnoLine H Set for Infant/Neonatal. In the presence of interfering gases, the above accuracy is maintained to within 4%.</p>	
Resolution	1 mmHg
Sample flow rate	50 ^{-7.5} ₊₁₅ ml/min
Initialization time	30 s (typical)
Response time	2.9 s (typical) (The response time is the sum of the rise time and the delay time when using a FilterLine of standard length) Rise time: <190 ms (10% to 90%) Delay time: 2.7 s (typical)
awRR measurement range	0 to 150 rpm
awRR measurement accuracy	0 to 70 rpm: ±1 rpm 71 to 120 rpm: ±2 rpm 121 to 150 rpm: ±3 rpm

Mainstream CO₂ Module

Standard	Meet standard of ISO 21647
CO ₂ Measurement range	0 to 150 mmHg
Accuracy	0 to 40 mmHg: ±2 mmHg 41 to 70 mmHg: ±5% of the reading 71 to 100 mmHg: ±8% of the reading 101 to 150 mmHg: ±10% of the reading
Accuracy drift	Meet the requirement for measurement accuracy within 6 hours
Resolution	1 mmHg
Response time	<60 ms
awRR measurement range	0 to 150 rpm
awRR measurement accuracy	1 rpm

A.6.14 tcGas

Operating mode	Interfaces with TCM CombiM or TCM TOSCA monitor	
Parameters	Measurement range	Measurement accuracy
tcpCO ₂	5 to 200 mmHg (0.7 to 26.7 kPa)	TOSCA Sensor 92, tc Sensor 54: 1 % CO ₂ : better than 1 mmHg (0.13 kPa) 10 % CO ₂ : better than 1 mmHg (0.13 kPa) 33 % CO ₂ : better than 3 mmHg (0.4 kPa) tc Sensor 84: 1 % CO ₂ : better than 1 mmHg (0.13 kPa) 10 % CO ₂ : better than 1 mmHg (0.13 kPa) 33 % CO ₂ : better than 5 mmHg (0.67 kPa)
tcpO ₂	0 to 800 mmHg (0.0 to 99.9 kPa)	tc Sensor 84: 0 % O ₂ : better than 1 mmHg (0.13 kPa) 21 % O ₂ : better than 3 mmHg (0.4 kPa) 50 % O ₂ : better than 5 mmHg (0.67 kPa) 90 % O ₂ : better than 25 mmHg (3.33 kPa)
SpO ₂	0 to 100 %	70 % to 100 %: ± 3 %
PR	25 bpm to 240 bpm	± 3 bpm
Power	0 to 1000 mW	± 20 % of reading

A.6.15 AG

Standards	Meet standard of ISO 21647	
Technique	Infrared absorption	
Warm-up time	Iso accuracy mode:	45 s
	Full accuracy mode:	10 min
Sample flow rate	Adult, pediatric:	120, 150, 200 ml/min
	Neonate:	70, 90, 120 ml/min
	Accuracy:	±10 ml/min or ±10%, whichever is greater
Measurement range	CO ₂ :	0 to 30%
	O ₂ :	0 to 100%
	N ₂ O:	0 to 100%
	Des:	0 to 30%
	Sev:	0 to 30%
	Enf:	0 to 30%
	Iso:	0 to 30%
	Hal:	0 to 30%
	awRR:	2 to 100 rpm
Resolution	CO ₂ :	1 mmHg
	awRR:	1 rpm

Iso accuracy	CO ₂ : ±0.3% _{ABS} N ₂ O: ±(8% _{REL} +2% _{ABS}) Other anesthetic gases: 8% _{REL}		
Full accuracy	Gases	Range (% _{REL})	Accuracy (% _{ABS})
	CO ₂	0 to 1	±0.1
		1 to 5	±0.2
		5 to 7	±0.3
		7 to 10	±0.5
		>10	Not specified
	N ₂ O	0 to 20	±2
		20 to 100	±3
	O ₂	0 to 25	±1
25 to 80		±2	
80 to 100		±3	
Des	0 to 1	±0.15	
	1 to 5	±0.2	
	5 to 10	±0.4	
	10 to 15	±0.6	
	15 to 18	±1	
>18	Not specified		
Sev	0 to 1	±0.15	
	1 to 5	±0.2	
	5 to 8	±0.4	
	>8	Not specified	
Enf, Iso, Hal	0 to 1	±0.15	
	1 to 5	±0.2	
	>5	Not specified	
awRR	2 to 60 rpm	±1 rpm	
	>60 rpm	Not specified	
Accuracy drift	Meet the requirement for measurement accuracy within 6 hours		
Apnea alarm time	10 s, 15 s, 20 s, 25 s, 30 s, 35 s, 40 s		
Refreshing rate	1 s		
Rise time (10 % ~ 90%)	gas sample flow rate 120ml/min, using the DRYLINE™ watertrap and neonatal DRYLINE™ sampling line (2.5m):		
	CO ₂	≤250 ms (fall time: 200ms)	
	N ₂ O	≤250 ms	
	O ₂	≤600 ms	
	Hal, Iso, Sev, Des	≤300 ms	
	Enf	≤350 ms	
	gas sample flow rate 200ml/min, using the DRYLINE™ water trap and adult DRYLINE™ sampling line (2.5m):		
	CO ₂	≤250 ms (fall time: 200 ms)	
	N ₂ O	≤250 ms	
	O ₂	≤500 ms	
	Hal, Iso, Sev, Des	≤300 ms	

	Enf	≤350 ms
Delay time	<4 s	
Anesthetic agent limit	Primary anesthetic agent In full accuracy mode: 0.15%, In ISO accuracy mode: 0.4%	
	Second anesthetic agent: In full accuracy mode: 0.3% or 5% REL (10% in ISO accuracy mode) of primary agent if primary agent is greater than 10% In ISO accuracy mode: 0.5%	

Effect of interference gases on AG measurements

Gas	Concentration(%)	Quantitive effect(%ABS)3)			
		CO ₂	N ₂ O	Agent 1)	O ₂
CO ₂	/	/	0.1	0	0.2
N ₂ O	/	0.1	/	0.1	0.2
Agent ^{1) 2)}	/	0.1	0.1	0.1	1
Xenon	<100%	0.1	0	0	0.5
Helium	<50%	0.1	0	0	0.5
Ethanol	<0.1%	0	0	0	0.5
Acetone	<1%	0.1	0.1	0	0.5
Methane	<1%	0.1	0.1	0	0.5
Saturated Isopropanol vapour	/	0.1	0	0	0.5
Metered dose inhaler propellants,	/	Unspecified	Unspecified	Unspecified	0.5

1) Agent represents one of Des, Iso, Enf, Sev, and Hal.

2) Multiple agent interference on CO₂, N₂O and O₂ is typically the same as single agent interference.

3) For CO₂, N₂O and Agents, maximum interference from each gas at concentrations within specified accuracy ranges for each gas. The total interference of all gases is never larger than 5%_{REL}.

Alarm limit	Range	Step
EtCO ₂ High	(low limit + 2) to 228 mmHg	1 mmHg
EtCO ₂ Low	0 to (high limit - 2)mmHg	
FiCO ₂ High	0 to 228 mmHg	
FiCO ₂ Low	0 to (high limit - 2)mmHg	
awRR High	Adult, pediatric: (low limit + 2) to 100 rpm Neonate: (low limit + 2) to 150 rpm	1 rpm
awRR Low	0 to (high limit - 2)rpm	
EtO ₂ High	(low limit + 0.3) to 100 %	0.1%
EtO ₂ Low	18 to (high limit - 0.3)%	
FiO ₂ High	(low limit + 0.3) to 100 %	
FiO ₂ Low	18 to (high limit - 0.3)%	
EtN ₂ O High	(low limit + 2) to 100 %	1%
EtN ₂ O Low	0 to (high limit - 2)%	
FiN ₂ O High	(low limit + 2) to 100 %	

FiN ₂ O Low	0 to (high limit - 2)%	
EtHal/Enf/Iso High	(low limit + 0.2) to 5.0 %	0.1%
EtHal/Enf/Iso Low	0 to (high limit - 0.2)%	
FiHal/Enf/Iso High	(low limit + 0.2) to 5.0 %	
FiHal/Enf/Iso Low	0 to (high limit - 0.2)%	
EtSev High	(low limit + 0.2) to 8.0 %	0.1%
EtSev Low	0 to (high limit - 0.2)%	
FiSev High	(low limit + 0.2) to 8.0 %	
FiSev Low	0 to (high limit - 0.2)%	
EtDes High	(low limit + 0.2) to 18.0 %	0.1%
EtDes Low	0 to (high limit - 0.2)%	
FiDes High	(low limit + 0.2) to 18.0 %	
FiDes Low	0 to (high limit - 0.2)%	

A.6.16 ICG

Technique	Thoracic electrical bioimpedance (TEB); Z MARC® algorithm(used for C.O. calculation)	
Measurement range	SV: 5 to 250 ml/beat HR: 44 to 185 bpm C.O.: 1.4 to 15 L/min	
Accuracy	SV: Not specified. HR: ±2 bpm C.O.: Not specified.	
Alarm limit	Range	Step
C.I. High	(low limit + 1.0) to 15.0 L/min/m ²	0.1 L/min/m ²
C.I. Low	0.0 to (high limit - 1.0)L/min/m ²	
TFC High	(low limit + 1) to 150 /kΩ	1 /kΩ
TFC Low	10 to (high limit - 1)/kΩ	

A.6.17 BIS

Standards	Meet standard of IEC 60601-2-26	
Technique	Bispectral index	
Measured parameters	EEG BIS, BIS L, BIS R: 0 to 100	
Calculated parameters	SQI, SQI L, SQI R:0 to 100 % EMG, EMG L, EMG R:0 to 100 dB SR, SR L, SR R:0 to 100% SEF, SEF L, SEF R:0.5 to 30.0 Hz TP, TP L, TP R:40 to 100 dB BC, BC L, BC R:0 to 30 sBIS L, sBIS R:0 to 10.0 sEMG L, sEMG R:0 to 10.0 ASYM:0 to 100%	
Impedance range	0 to 999 kΩ	
Sweep speed	6.25 mm/s, 12.5 mm/s, 25 mm/s or 50 mm/s	
Input impedance	>5 MΩ	
Noise (RTI)	<0.3 μV (0.25 to 50 Hz)	
Input signal range	±1 mV	
EEG bandwidth	0.25 to 100 Hz	
Patient leakage current	<10 μA	
Alarm limit	Range	Step
BIS High	(low limit + 2) to 100	1
BIS Low	0 to (high limit – 2)	

A.6.18 NMT

Operating mode	Interfaces with TOF-Watch® SX monitor	
Parameters	Measurement range	
TOF-Ratio	1%~160%	
TOF-Count	0~4	
Single	0%~160%	
PTC	0~15	
Tskin	20.0°C~41.5°C	

A.6.19 RM

Technique	Flow sensor	
Frequency response	≥30 Hz	
Dead space	≤11 ml	
Flow		
Measurement range	Adult/pediatric*:	± (2 to 120) L/min
	Infant:	± (0.5 to 30) L/min
Accuracy	Adult/pediatric*:	1.5 L/min or ±10% of the reading, whichever is greater
	Infant:	0.5 L/min or ±10% of the reading, whichever is greater
Resolution	0.1 L/min	
Paw		
Measurement range	-20 to 120 cmH ₂ O	
Accuracy	±3%	
Resolution	0.1 cmH ₂ O	
MVe/MVi		
Measurement range	Adult/Pediatric*:	2 to 60 L/min
	Infant:	0.5 to 15 L/min
Accuracy	±10%×reading	
TVe/TVi		
Measurement range	Adult/Pediatric*:	100 to 1500 ml
	Infant:	20 to 500 ml
Resolution	1 ml	
Accuracy	Adult/pediatric*:	±10% or 15 ml, whichever is greater
	Infant:	±10% or 6 ml, whichever is greater
RR (RM)		
Measurement range	4 to 120 rpm	
Accuracy	4 to 99 rpm	±1 rpm
	100 to 120 rpm	±2 rpm

*Pediatric in this form does not include neonate and infant.

Calculated Parameters		
	Measurement range	Measurement accuracy
I:E	4:1 to 1:8	Not specified.
FEV1.0%	0 to 100%	Not specified.
Pmean	0 to 120 cmH ₂ O	±10%×reading
TV	20 to 1500 ml	Adult/pediatric: ±10% or ±25 ml, whichever is greater. Infant: ±10% or ±6 ml, whichever is greater.
MV	2 to 60 L	±10%×reading
PEEP	0 to 120 cmH ₂ O	Not specified.
PEF	2 to 120 L/min	±2L/min or ±10% of reading, whichever is greater.
PIF	2 to 120 L/min	±2L/min or ±10% of reading, whichever is greater.
PIP	0 to 120 cmH ₂ O	1cmH ₂ O or ±10% of reading, whichever is greater.

Pplat	0 to 120 cmH ₂ O	Not specified.
Compl	0 to 200 ml/cmH ₂ O	
RSBI	0 to 4095 rpm/L	

Alarm limit	Range	Step
RR High	Adult, pediatric: (low limit + 2) to 100 rpm Neonate: (low limit + 2) to 150 rpm	1 rpm
RR Low	0 to (high limit -2) rpm	
PEEP High	(low limit +1) to 120 cmH ₂ O	1 cmH ₂ O
PEEP Low	0 to (high limit -1) cmH ₂ O	
PIP High	(low limit +1) to 120 cmH ₂ O	1 cmH ₂ O
PIP Low	0 to (high limit -1) cmH ₂ O	
MVe High	Adult and pediatric: (low limit +1.0) to 60.0 L/min Infant:(low limit +1.0) to 15.0	0.5 L/min
MVe Low	Adult and pediatric: 2.0 to (high limit -1.0) Infant:0.5 to (high limit -1.0)	

SvO ₂ -related Parameters Outputted by Vigileo™ Monitor		
Name	Measurement Range	Resolution
SvO ₂	0 to 99%	1
ScvO ₂	0 to 99%	1
SQI	1 to 4	1

Alarm Limit	Range (%)	Step (%)
SvO ₂ / ScvO ₂ High	(Low limit+1)to 99	1
SvO ₂ / ScvO ₂ Low	0 to(high limit-1)	

B EMC


The product is in radio-interference protection class A in accordance with EN55011. The product complies with the requirement of standard EN60601-1-2 "Electromagnetic Compatibility – Medical Electrical Equipment".

Note

- **Using accessories, transducers and cables other than those specified may result in increased electromagnetic emission or decreased electromagnetic immunity of the patient monitoring equipment.**
 - **The device or its components should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, the device or its components should be observed to verify normal operation in the configuration in which it will be used.**
 - **The device needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided below.**
 - **Other devices may affect this monitor even though they meet the requirements of CISPR.**
 - **When the inputted signal is below the minimum amplitude provided in technical specifications, erroneous measurements could result.**
-

Guidance and Declaration - Electromagnetic Emissions		
The device is suitable for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.		
Emission tests	Compliance	Electromagnetic environment - guidance
Radio frequency (RF) emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class A	The device is suitable for use in all establishments other than domestic and those indirectly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC61000-3-2	Class A	
Voltage Fluctuations/Flicker Emissions IEC 61000-3-3	Complies	

Guidance and Declaration - Electromagnetic Immunity			
The device is suitable for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.			
Immunity test	IEC60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air (±6kV air for CF card)	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV I/O for input/output lines (>3 m)	±2 kV for power supply lines ±1 kV I/O for input/output lines (>3 m)	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV common mode	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % U_T (>95 % dip in U_T) for 0.5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles <5 % U_T (>95 % dip in U_T) for 5 s	<5 % U_T (>95 % dip in U_T) for 0.5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles <5 % U_T (>95 % dip in U_T) for 5 s	Mains power quality should be that of a typical commercial or hospital environment. If the user of our product requires continued operation during power mains interruptions, it is recommended that our product be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 HZ) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note: U_T is the AC mains voltage prior to application of the test level.			

Guidance and Declaration - Electromagnetic Immunity		
The device is suitable for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.		
Immunity test	IEC 60601 Test level	Compliance level
Conducted RF IEC61000-4-6	3 Vrms 150k to 80M Hz	3 Vrms (BIS, ICG: 1Vrms)
Radiated RF IEC61000-4-3	3V/m 80M to 2.5G Hz	3V/m (Resp: 1V/m)
Electromagnetic environment - guidance		
<p>Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended Separation Distance:</p> $d = 1.2\sqrt{P} \quad (\text{BIS, ICG: } d = 3.5\sqrt{P})$ $d = 1.2\sqrt{P} \quad (\text{Resp: } d = 3.5\sqrt{P}) \quad 80 \text{ to } 800 \text{ MHz}$ $d = 2.3\sqrt{P} \quad 800\text{M to } 2.5\text{GHz}$ <p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range^b Interference may occur in the vicinity of equipment marked with the following</p> <p>symbol:  .</p>		
<p>Note 1: From 80 MHz to 800 MHz, the higher frequency range applies.</p> <p>Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>		
<p>a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the device.</p> <p>b Over the frequency ranges 150kHz to 80MHz, field strengths should be less than 3V/m. For BIS and ICG monitoring, the field strength should be less than 1V/m</p>		

Recommended Separation Distances between Portable and Mobile RF Communications Equipment and The device

The device is suitable for use in an electromagnetic environment in which radiated RF disturbance are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communication equipment.

Rated Maximum Output power of Transmitter Watts (W)	Separation Distance Meters (m) Corresponding to Frequency of Transmitter		
	150k to 80MHz $d = 3.5\sqrt{P}$	80M to 800MHz $d = 3.5\sqrt{P}$	800M to 2.5GHz $d = \left[\frac{7}{3}\right]\sqrt{P}$
0.01	0.35	0.35	0.23
0.1	1.11	1.11	0.74
1	3.5	3.5	2.34
10	11.07	11.07	7.38
100	35	35	23.34

For transmitters at a maximum output power not listed above, the separation distance can be estimated using the equation in the corresponding column, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: From 80 MHz to 800 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

C Default Configurations

This chapter lists some of the most important factory default settings for each department in configuration management. You cannot change the factory default configuration itself. However, you can make changes to the settings from the factory default configuration and then save the changed configuration as a user configuration. The last column of the following tables is for your notes and review.

Note: In this chapter, O.M means the monitor's operating mode. Column C refers to the settings that can be changed in configuration management. Column M refers to the settings that can be changed in monitoring mode.

C.1 Parameters Configuration

C.1.1 ECG

ECG Setup

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Lead Set		*	Auto (if auto lead detection is available); 5-lead (if auto lead detection is not available)					
Alm Source	*	*	HR					
Alarm	*	*	On					
Alm Lev	*	*	Med					
HR/PR High	Adu		120					
	Ped	*	160					
	Neo	*	200					
HR/PR Low	Adu		50					
	Ped	*	75					
	Neo	*	100					
Sweep	*	*	25 mm/s					
Beat Vol	*	*	2		1			
Paced		*	No					
Notch Filter	*	*	Weak					
Gain	*	*	X1					
Filter	*	*	Monitor	Surgery	Monitor	Diagnostic		
ECG Display	*	*	Normal					
Pacemaker Rate		*	60					

ST Analysis

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
ST Analysis	*	*	Off				On	
Alarm	*	*	Off					
Alm Lev	*	*	Med					
ST-X High	*	*	when ST Unit is mV:			0.20		
			when ST Unit is mm:			2.0		
ST-X Low	*	*	when ST Unit is mV:			-0.20		
			when ST Unit is mm:			-2.0		
ISO			-80 ms					
J	*	*	48 ms					
ST			J + 60 ms					

X represents I, II, III, aVR, aVL, aVF, V, V1, V2, V3, V4, V5 or V6.

Arrh. Analysis

Item Name	Algorithm	O.M		Genral	OR	ICU	NICU	CCU	User Defaults
		C	M						
Arrhythmia Threshold Settings									
PVCs High	Mindray	*	*	Adu, Ped:	10				
				Neo:	N/A				
Tachy		*	*	Adu:	120				
				Ped:	160				
				Neo:	N/A				
Brady		*	*	Adu:	50				
				Ped:	75				
				Neo:	N/A				
Asys. Delay		*	*	Adu, Ped:	5				
				Neo:	N/A				
Vtac Rate	*	*	Adu, Ped:	130					
			Neo:	N/A					
Vtac PVCs	*	*	Adu, Ped:	6					
			Neo:	N/A					
Multif. PVC's Window	*	*	Adu, Ped:	15					
			Neo:	N/A					
Extreme Tachy	*	*	Adu:	160					
			Ped:	180					
			Neo:	N/A					
Extreme Brady	*	*	Adu:	35					
			Ped:	50					
			Neo:	N/A					
Vbrd Rate	*	*	Adu, Ped:	40					
			Neo:	N/A					

Item Name	Algorithm	O.M		Genral	OR	ICU	NICU	CCU	User Defaults
		C	M						
Vbrd PVCs	Mortara	*	*	Adu, Ped: Neo:	5 N/A				
Pause Time		*	*	Ad, Ped: Neo:	2 N/A				
PVCs High		*	*	Adu, Ped: Neo:	10 N/A				
Asys. Delay		*	*	Adu, Ped: Neo:	4 N/A				
Vtac Rate		*	*	Adu, Ped: Neo:	130 N/A				
Vtac PVCs		*	*	Adu, Ped: Neo:	6 N/A				
Multif. PVC's Window		*	*	Adu, Ped: Neo:	15 N/A				
Tachy		*	*	Adu: Ped: Neo:	120 160 N/A				
Brady		*	*	Adu: Ped: Neo:	50 75 N/A				
Arrhythmia Alarm Settings									
PVCs/min Alarm	Mindray	*	*	Off				On	
R on T Alarm		*	*	Off				On	
Nonsus. Vtac Alarm		*	*	Off				On	
Vent. Rhythm Alarm		*	*	Off				On	
Bigeminy Alarm		*	*	Off				On	
Trigeminy Alarm		*	*	Off				On	
Asystole Alarm		*	*	On					
VFib/VTac Alarm		*	*	On					
Vtac Alarm		*	*	On					
Vent. Brady Alarm		*	*	On					
Extreme Tachy Alarm		*	*	On					
Extreme Brady Alarm		*	*	On					
X Alarm		*	*	Off					
Asystole Alm Lev		*	*	High					
VFib/VTac Alm Lev		*	*	High					
Vtac Alm Lev		*	*	High					
Vent. Brady Alm Lev	*	*	High						

Item Name	Algorithm	O.M		Genral	OR	ICU	NICU	CCU	User Defaults
		C	M						
Extreme Tachy Alm Lev		*	*	High					
Extreme Brady Alm Lev		*	*	High					
VT>2 Alm Lev		*	*	Low					
Pause Alm Lev		*	*	Low					
Couplet Alm Lev		*	*	Prompt					
PVC Alm Lev		*	*	Prompt					
Irr. Rhythm Alm Lev		*	*	Prompt					
PNP Alm Lev		*	*	Prompt					
PNC Alm Lev		*	*	Prompt					
Missed Beats Alm Lev		*	*	Prompt					
X Alm Lev		*	*	Med					
X Alm Rec		*	*	Off					
PVCs/min Alarm		Mortara	*	*	Off				On
R on T Alarm	*		*	Off				On	
Vent. Rhythm Alarm	*		*	Off				On	
Bigeminy Alarm	*		*	Off				On	
Trigeminy Alarm	*		*	Off				On	
Asystole Alarm				On					
VFib Alarm				On					
VTac Alarm				On					
X Alarm	*		*	On					
Asystole Alm Lev	*		*	High					
VFib Alm Lev	*		*	High					
VTac Alm Lev	*		*	High					
VT>2 Alm Lev	*		*	Low					
Couplet Alm Lev	*		*	Prompt					
PVC Alm Lev	*		*	Prompt					
Irr. Rhythm Alm Lev	*		*	Prompt					
PNP Alm Lev	*		*	Prompt					
PNC Alm Lev	*		*	Prompt					
Missed Beats Alm Lev	*		*	Prompt					
X Alm Lev	*		*	Med					
X Alm Rec	*		*	Off					

X represents a certain arrhythmia event. Refer to the Specifications chapter for details. The X in "X Alm Lev" refers to all arrhythmia events except for those specially marked ones.

C.1.2 RESP

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Sweep	*	*	6.25 mm/s					
Lead	*	*	Adu, Ped:		Auto			
			Neo:		II			
Gain	*	*	X2					
RR High	*	*	Adu, Ped:		30			
			Neo:		100			
RR Low	*	*	Adu, Ped:		8			
			Neo:		30			
Apnea Delay	*	*	Adu, Ped:		20			
			Neo:		15			
Detection Mode	*	*	Auto					
RR Source		*	Auto					

C.1.3 PR

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults	
	C	M							
Alarm	*	*	On						
Alm Lev	*	*	Med						
HR/PR High	Adu	*	*	120					
	Ped			160					
	Neo			200					
HR/PR Low	Adu	*	*	50					
	Ped			75					
	Neo			100					
PR Source	*	*	SpO ₂						
Beat Vol	*	*	2		1				

C.1.4 SpO₂

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
SpO ₂ High	*	*	Adu, Ped: 100 Neo: 95					
SpO ₂ Low	*	*	90					
Desat Limit	*	*	80					
Sweep	*	*	25 mm/s					
NIBP Simul		*	Off					
Sensitivity (Mindray)	*	*	Med					
Sensitivity (Masimo)	*	*	Normal					
Averaging (Masimo)	*	*	8 s					
Sat-Seconds (Nellcor)	*	*	0 s					

C.1.5 ΔSpO₂

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	Off					
Alm Lev	*	*	Mediate					
ΔSpO ₂ High	*	*	10 %					

C.1.6 Temp

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
T1/T2 High (°C)	*	*	38.0					
T1/T2 Low (°C)	*	*	35.0					
TD High (°C)	*	*	2.0					

C.1.7 NIBP

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults	
	C	M							
Alarm	*	*	On						
Alm Lev	*	*	Med						
Interval	*	*	15 min	5 min	15 min	30 min	15 min		
NIBP End Tone		*	Off						
Cuff Press. (mmHg)	Adu		80						
	Ped	*	*	60					
	Neo		40						
Initial Pressure (mmHg)	Adu		160						
	Ped	*	*	140					
	Neo		90						
Alarm Limits									
NIBP-S High (mmHg)	Adu		160						
	Ped	*	*	120					
	Neo		90						
NIBP-S Low (mmHg)	Adu		90						
	Ped	*	*	70					
	Neo		40						
NIBP-M High (mmHg)	Adu		110						
	Ped	*	*	90					
	Neo		70						
NIBP-M Low (mmHg)	Adu		60						
	Ped	*	*	50					
	Neo		25						
NIBP-D High (mmHg)	Adu		90						
	Ped	*	*	70					
	Neo		60						
NIBP-D Low (mmHg)	Adu		50						
	Ped	*	*	40					
	Neo		20						

C.1.8 IBP

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
P1 Measure	*	*	All					
P2 Measure	*	*	All					
P3 Measure	*	*	Mean					
P4 Measure	*	*	Mean					
Sensitivity	*	*	Med					

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Sweep	*	*	25 mm/s					
Filter	*		12.5 Hz					
Art, Ao, UAP, BAP, FAP, LV, P1-P2 Arterial Pressure Alarm Limits								
IBP-S High (mmHg)	Adu			160				
	Ped	*	*	120				
	Neo			90				
IBP-S Low (mmHg)	Adu			90				
	Ped	*	*	70				
	Neo			55				
IBP-M High (mmHg)	Adu			110				
	Ped	*	*	90				
	Neo			70				
IBP-M Low (mmHg)	Adu			70				
	Ped	*	*	50				
	Neo			35				
IBP-D High (mmHg)	Adu			90				
	Ped	*	*	70				
	Neo			60				
IBP-D Low (mmHg)	Adu			50				
	Ped	*	*	40				
	Neo			20				
PA Alarm Limits								
PA-S High (mmHg)	Adu			35				
	Ped	*	*	60				
	Neo			60				
PA-S Low (mmHg)	Adu			10				
	Ped	*	*	24				
	Neo			24				
PA-M High (mmHg)	Adu			20				
	Ped	*	*	26				
	Neo			26				
PA-M Low (mmHg)	Adu			0				
	Ped	*	*	12				
	Neo			12				
PA-D High (mmHg)	Adu			16				
	Ped	*	*	4				
	Neo			4				
PA-D Low (mmHg)	Adu			0				
	Ped	*	*	-4				
	Neo			-4				
CVP, LAP, RAP, ICP, UVP, P3-P4 Venous Pressure Alarm Limits								
IBP-M High (mmHg)	Adu	*	*	10				
	Ped			4				

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
	Neo		4					
IBP-M Low (mmHg)	Adu		0					
	Ped	*	0					
	Neo	*	0					
Art, Ao, BAP, FAP, LV, P1-P2 Arterial Pressure Scale								
Scale (mmHg)	*	*	0-160					
PA Scale								
Scale (mmHg)	*	*	0-30					
CVP, LAP, RAP, ICP, UVP Scale								
Scale (mmHg)	*	*	0-20					
UAP, P3-P4 Venous Pressure Scale								
Scale (mmHg)	*	*	0-80					

C.1.9 C.O.

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
TB High (°C)	*	*	39.0					
TB Low (°C)	*	*	36.0					
Comp. Const	*	*	0.542					
Auto TI	*	*	Auto					
Manual TI (°C)	*	*	2.0					
Measuring mode	*	*	Manual					

C.1.10 CCO/SvO₂ Setup (Vigilance II)

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Primary Parameter	*	*	C.O./CCO					
Secondary Parameters	*	*	SVR, EDV, SV					
CCO High	*	*	14					
CCO Low	*	*	2					
CCI High	*	*	7					
CCI Low	*	*	1					
EDV High	*	*	300					
EDV Low	*	*	80					
EDVI High	*	*	150					
EDVI Low	*	*	60					
SVR High	*	*	1500 DS/cm5					

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
SVR Low	*	*	500 DS/cm5					
SVRI High	*	*	3000 DS-m2/ cm5					
SVRI Low	*	*	1000 DS-m2/ cm5					
RVEF High	*	*	50					
RVEF Low	*	*	0					
SV High	*	*	120					
SV Low	*	*	20					
SVI High	*	*	60					
SVI Low	*	*	10					
SvO ₂ High	*	*	99					
SvO ₂ Low	*	*	10					
ScvO ₂ High	*	*	99					
ScvO ₂ Low	*	*	10					

C.1.11 CCO/SvO₂ Setup (Vigileo)

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Primary Parameter	*	*	CCO					
Secondary Parameters	*	*	SV, SVR, SVV					
CCO High	*	*	14					
CCO Low	*	*	2					
CCI High	*	*	7					
CCI Low	*	*	1					
SV High	*	*	120					
SV Low	*	*	20					
SVI High	*	*	60					
SVI Low	*	*	10					
SVV High	*	*	30					
SVV Low	*	*	0					
SVR High	*	*	1500 DS/cm5					
SVR Low	*	*	500 DS/cm5					
SVRI High	*	*	3000 DS-m2/ cm5					
SVRI Low	*	*	1000 DS-m2/ cm5					
SvO ₂ High	*	*	99					
SvO ₂ Low	*	*	10					
ScvO ₂ High	*	*	99					
ScvO ₂ Low	*	*	10					

C.1.12 PiCCO

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Inj. Volume		*	Adu: 15ml Ped: 10ml					
Manual Input pCVP		*	Unselected					
pCVP Value		*	Empty					
Auto C.O.		*	On					
PiCCO Parameters								
Main Parameter	*	*	CCI, CCO					
Secondary Parameter	*	*	GEDI, ELWI, SVRI					
Spider	*	*	CCI, ScvO ₂ , pArt-M, SVRI, SVV					
pArt/pCVP Setup								
Scale (mmHg)	*	*	pArt: 0~160mmHg pCVP: 0~20mmHg					
Sweep	*	*	25 mm/s					

C.1.13 CO₂

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Operating Mode	*	*	Measure					
Sweep	*	*	6.25 mm/s					
Scale (mmHg)	*	*	50					
Apnea Delay	*	*	Adu, Ped: 20 Neo: 15					
RR Source		*	Auto					
Sidestream CO₂ Setup								
Flow Rate	*	*	Adu: 120 ml/min					
			,Ped: 100 ml/min					
			Neo: 70 ml/min					
BTPS Compen	*	*	Off					
N ₂ O Compen	*	*	0					
O ₂ Compen	*	*	21	100	21			
Des Compen	*	*	0					
Microstream CO₂ Setup								
BTPS Compen		*	Off					
Max Hold		*	20 s					
Auto Standby (min)		*	0					
Mainstream CO₂ Setup								
Max Hold	*	*	10 s					
O ₂ Compen	*	*	Off					

Balance Gas	*	*	Room Air	
AG Compen	*	*	0	
Alarm Limits				
EtCO ₂ High (mmHg)	*	*	Adu, Ped: 50 Neo: 45	
EtCO ₂ Low (mmHg)	*	*	Adu, Ped: 25 Neo: 30	
FiCO ₂ High (mmHg)	*	*	Adu, Ped, Neo: 4	
RR High	*	*	Adu, Ped: 30 Neo: 100	
RR Low	*	*	Adu, Ped: 8 Neo: 30	

C.1.14 tcGas

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm Sound	*	*	Off					
Change Secondary Parameters	*	*	SpO ₂ , PR, Power					

C.1.15 AG

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Sweep	*	*	6.25 mm/s					
O ₂ Compen	*	*	Off	On	Off			
Operating Mode	*	*	Measure					
Flow Rate	*	*	Adu, Ped: 120 ml/min Neo: 70 ml/min					
Auto Standby	*	*	Off					
Apnea Time	*	*	20 s					
RR Source		*	Auto					
CO₂ Setup								
Wave Type	*	*	Draw					
Scale	*	*	when unit is mmHg: 50 when unit is % or KPa: 7.0					
EtCO ₂ High (mmHg)	*	*	Adu, Ped: 50 Neo: 45					
EtCO ₂ Low (mmHg)	*	*	Adu, Ped: 25 Neo: 30					
FiCO ₂ High (mmHg)	*	*	4					

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
RR High	*	*	Adu, Ped: 30 Neo: 100					
RR Low	*	*	Adu, Ped: 8 Neo: 30					
Gas Setup								
Agent	*	*	AA					
N ₂ O Scale	*	*	50					
O ₂ Scale	*	*	when unit is mmHg: 400 when unit is % or KPa: 50					
AA Scale	*	*	9.0					
Hal/Enf/Iso Scale	*	*	2.5					
Des Scale	*	*	9.0					
Sev Scale	*	*	4.0					
EtO ₂ High	*	*	88					
EtO ₂ Low	*	*	18					
FiO ₂ High	*	*	Adu, Ped: 100 Neo: 90					
FiO ₂ Low	*	*	18					
EtN ₂ O High	*	*	55					
EtN ₂ O Low	*	*	0					
FiN ₂ O High	*	*	53					
FiN ₂ O Low	*	*	0					
EtHal/Enf/Iso High	*	*	3.0					
EtHal/Enf/Iso Low	*	*	0.0					
FiHal/Enf/Iso High	*	*	2.0					
FiHal/Enf/Iso Low	*	*	0.0					
EtSev High	*	*	6.0					
EtSev Low	*	*	0.0					
FiSev High	*	*	5.0					
FiSev Low	*	*	0.0					
EtDes High	*	*	8.0					
EtDes Low	*	*	0.0					
FiDes High	*	*	6.0					
FiDes Low	*	*	0.0					

C.1.16 ICG

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Averaging	*	*	30					
Update Rate	*	*	10					
Sweep	*	*	12.5 mm/s					
Secondary Parameters	*	*	C.O., SVR, TFC					
C.I. High	*	*	5.0					
C.I. Low	*	*	1.5					
TFC High	*	*	60					
TFC Low	*	*	10					

C.1.17 BIS

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Smoothing Rate	*	*	15 s		Adu: 15s Ped: 15s Neo: N/A	15 s		
Display	*	*	EEG		Adu: EEG Ped: EEG Neo: N/A	EEG		
Filters	*	*	On		Adu: On Ped: On Neo: N/A	On		
Scale	*	*	100 μ V		Adu: 100 μ V Ped: 100 μ V Neo: N/A	100 μ V		
Sweep	*	*	25mm/s		Adu: 25mm/s Ped: 25mm/s Neo: N/A	25mm/s		
Trend Length	*	*	60 min		Adu: 60 min Ped: 60 min Neo: N/A	60 min		
Secondary Parameters	*	*	SR, SEF		Adu, Ped: SR, SEF Neo: N/A	SR,SEF		

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Display	*	*	BIS Trend			Adu: BIS Trend Ped: BIS Trend Neo: N/A	BIS Trend	
EEG Waveforms	*	*	All			Adu: All Ped: All Neo: N/A	All	
Parameter 1	*	*	BIS L			Adu: BIS L Ped: BIS L Neo: N/A	BIS L	
Parameter 2	*	*	EMG			Adu: EMG Ped: EMG Neo: N/A	EMG	
BIS High	*	*	70					
BIS Low	*	*	70					

C.1.18 NMT

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm Sound	*	*	Off					

C.1.19 RM

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alarm	*	*	On					
Alm Lev	*	*	Med					
Apnea Delay	*	*	Adu, Ped: 20 s Neo: 15 s					
Sensor Type		*	Disposable					
TV/MV	*	*	TV					
Flow/Vol	*	*	Flow					
Sweep	*	*	6.25 mm/s					
RR Source		*	Auto					
Paw Scale	*	*	Adu, Ped: 40 Neo: N/A					
Flow Scale	*	*	Adu, Ped: 60 Neo: N/A					
Vol Scale	*	*	Adu, Ped: 1200 Neo: N/A					

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Display Loop		*	PV Loop					
Reference Loop		*	On					
RR High	*	**	Adu, Ped: 30 Neo: 100					
RR Low	*	*	Adu, Ped: 8 Neo: 30					
PEEP High	*	*	10					
PEEP Low	*	*	0					
PIP High	*	*	40					
PIP Low	*	*	1					
MVe High	*	*	30.0					
MVe Low	*	*	2.0					

C.2 Routine Configuration

C.2.1 Alarm

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Alm Volume	*	*	2	1	2			
Reminder Vol	*	*	Low					
Recording Length	*	*	16 s					
Apnea Delay	*	*	Adu, Ped: 20 s Neo: 15 s					
Alarm Delay	*	*	6 s					
ST Alarm Delay	*	*	30 s					

C.2.2 Screens

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Choose Screen	*	*	Normal Screen					
Display the ST segments on ECG screen	*	*	Unselected					
Select Wave Sequence for Normal Screen	1	*	*	ECG1				
	2			ECG2				
	3			SpO ₂ +PR				
	4			Any IBP				
	5			Any IBP				
	6			CO ₂				

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
	7		Paw					
	8		Flow/Vol					
	9		ICG					
	10		BIS					
	11		Resp					
Select Parameters for Big Numerics Screen	Parameter 1		ECG					
	Parameter 2	*	SpO ₂ +PR					
	Parameter 3		Resp					
	Parameter 4		NIBP					

Item Name		Select QuickKeys (BeneView T5)
O.M	C	*
	M	
General		NIBP Measure→Stop All→Zero IBP→Review→Standby→Screens→Patient Setup→Manual Event→Realtime Print→Volume Setup
OR		
ICU		
NICU		NIBP Measure→Stop All→oxyCRG→Review→Standby→Screens→Patient Setup→Manual Event→Realtime Print→Volume Setup
CCU		NIBP Measure→Stop All→Zero IBP→Review→Standby→Screens→Patient Setup→Manual Event→Realtime Print→Volume Setup
User Defaults		

Item Name		Select QuickKeys (BeneView T8)
O.M	C	*
	M	
General		NIBP Measure→Stop All→Zero IBP→Screens→Patient Setup→Manual Event→Realtime Print→Print Setup→Minitrends→Volume Setup→Load Configuration→Privacy Mode
OR		NIBP Measure→Stop All→Zero IBP→Screens→Patient Setup→Manual Event→Realtime Print→Print Setup→Minitrends→Volume Setup→Load Configuration→PAWP
ICU		NIBP Measure→Stop All→Zero IBP→Screens→Patient Setup→Manual Event→Realtime Print→Print Setup→Minitrends→Volume Setup→Load Configuration→Privacy Mode
NICU		NIBP Measure→Stop All→oxyCRG→Screens→Patient Setup→Manual Event→Realtime Print→Minitrends→Zero IBP→Volume Setup→Load Configuration→Privacy Mode
CCU		NIBP Measure→Stop All→Zero IBP→Screens→Patient Setup→Manual Event→Realtime Print→Print Setup→Minitrends→Volume Setup→Load Configuration→Privacy Mode
User Defaults		

C.2.3 Waveform

Item Name		O.M		General	OR	ICU	NICU	CCU	User Defaults
		C	M						
Parameter/ Wave Colour	ECG			Green					
	NIBP			White					
	SpO ₂			Cyan					
	SpO _{2b}			Purple					
	ΔSpO ₂			Yellow					
	PR			Cyan					
	TEMP			White					
	Art/Ao/UAP/FAP /BAP/LV/P1~P4 (arterial pressure)			Red					
	PA			Yellow					
	CVP/ICP/P1~P4 (venous pressure)			Blue					
	LAP			Purple					
	RAP			Orange					
	UVP			Cyan					
	CO ₂ /tcpCO ₂			Yellow					
	RESP		*	Yellow					
	AA			Yellow					
	N ₂ O			Blue					
	O ₂ /tcpO ₂			Green					
	Hal			Red					
	Enf			Orange					
	Iso			Purple					
	Des			Cyan					
	Sev			Yellow					
	C.O.			White					
	Paw			Blue					
	Flow/Vol								
	EEG L/BIS L Trend			Yellow					
	EEG R/BIS R Trend			Blue					
	ICG			Purple					
	SvO ₂			Cyan					
ScvO ₂			Purple						
CCO			Yellow						
NMT			White						

X represents a waveform label, such as ECG, RESP, CO₂ and so forth. The ECG waveform cannot be set off.

C.2.4 Review

Item Name		O.M		General	OR	ICU	NICU	CCU	User Defaults
		C	M						
Tabular Trends	Interval	*	*	30 min	5min	30 min			
	Trend Group	*	*	Standard					
Graphic Trends	Trend Group	*	*	Standard					
Minitrend Length			*	2 h					
Full Disclosure	Save Waves	*	*	Save ECG1 by default.					

C.2.5 Event

Item Name		O.M		General	OR	ICU	NICU	CCU	User Defaults
		C	M						
Waveform 1			*	II					
Waveform 2			*	I			Pleth	I	
Waveform 3			*	Pleth			Resp	Pleth	

C.2.6 Record

Item Name		O.M		General	OR	ICU	NICU	CCU	User Defaults
		C	M						
Length			*	8 s					
Interval			*	Off					
Paper Speed			*	25 mm/s					
Alm Rec	X		*	Off					

X represents a parameter label.

C.2.7 Print

Item Name		O.M		General	OR	ICU	NICU	CCU	User Defaults
		C	M						
Paper Size			*	A4					
ECG Reports	Amplitude		*	10 mm/mV					
	Sweep		*	25 mm/s					
	Auto Interval		*	Off					
	12-Lead Format		*	12X1					
Tabular Trends Reports	Set as End Case Report		*	Unselected					
	Back		*	Auto					
	Spacing		*	Auto					
	Report Layout		*	Parameter Oriented					
	Currently Displayed Trended Parameters		*	Selected					
	Standard Parameter Group		*	Unselected					
Graphic Trends Reports	Custom		*	Unselected					
	Set as End Case Report		*	Unselected					
	Back		*	Auto					

	Zoom		*	Auto	
Realtime Report	Set as End Case Report		*	Unselected	
	Sweep		*	Auto	
	Select Wave		*	Current	

C.2.8 Others

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Brightness		*	5					
Key Volume		*	2					

C.3 User Maintenance Items

Item Name	O.M		General	OR	ICU	NICU	CCU	User Defaults
	C	M						
Changing Bed No.		*	Protected					
Atmospheric Pressure		*	760 mmHg					
Height Unit		*	cm					
Weight Unit		*	kg					
ST Unit		*	mV					
Press. Unit		*	mmHg					
CVP Unit		*	cmH ₂ O					
CO ₂ Unit		*	mmHg					
O ₂ Unit		*	%					
Temp Unit		*	°C					
Network Type		*	LAN					
Latching Alarms	*	*	No					
Alarm Pause Time	*	*	2 min					
Minimum Alarm Volume	*	*	2	1	2			
Alarm Sound		*	ISO					
Reminder Tone		*	Off					
Reminder Interval		*	1 min					
ECGLeadOff Lev.		*	Low					
SpO ₂ SensorOff Lev.		*	Low					
Alarm Tone Interval		*	High Level Alarm: 10 s Med/Low Level Alarm: 20 s					
Lethal Arrh. OFF		*	Disable					
Extended Arrh.		*	Enable					
Silence Other Bed		*	On					
Wave Line		*	Mediate					
Sweep Mode	*	*	Refresh					
Auxiliary Output		*	Analog Out.					
Primary Button		*	Left					





Item Name		O.M		General	OR	ICU	NICU	CCU	User Defaults
		C	M						
ECG Standard			*	AHA					
Notch Freq.			*	50 Hz					
Data Transfer Method			*	Off					
Transferred Data Length			*	4 h					
Apply Module Settings			*	On					
SpO ₂ Tone			*	Mode 1					
Signal Type			*	Continuous					
Contact Type			*	Normally Closed					
Nurse Call	Signal Type		**	Continuous					
	Contact Type		*	Normally Closed					
	Alm Lev	*	*	High, Med, Low					
	Alarm Cat.	*	*	Phys., Tech.					

FOR YOUR NOTES

D Alarm Messages

This chapter lists only the most important physiological and technical alarm messages. Some messages appearing on your monitor may not be included.

In this chapter:

The “I” field indicates how alarm indications are cleared: “A” means all alarm indications are cleared after the  hardkey is pressed, “B” indicates alarm lamp flashing and alarm tones are cleared and the alarm messages change to prompt messages after the  hardkey is pressed, and “C” indicates alarm lamp flashing and alarm tones are cleared and  appears before the alarm message after the  hardkey is pressed.

The “L” field indicates the alarm level: H means high, M means medium and L means low. “*” means the alarm level is user-adjustable.

XX represents a measurement or parameter label, such as ECG, NIBP, HR, ST-I, PVCs, RR, SpO₂, PR, etc.

In the “Cause and Solution” column, corresponding solutions are given instructing you to troubleshoot problems. If the problem persists, contact your service personnel.

D.1 Physiological Alarm Messages

Measurement	Alarm messages	L	Cause and solution
XX	XX Too High	M*	XX value has risen above the high alarm limit or fallen below the low alarm limit. Check the patient’s condition and check if the patient category and alarm limit settings are correct.
	XX Too Low	M*	
ECG	ECG Weak Signal	H	The ECG signal is so weak that the monitor can’t perform ECG analysis. Check the patient’s condition and the ECG connections.
	Asystole	H	Arrhythmia has occurred to the patient. Check the patient’s condition and the ECG connections.
	VFib/VTac	H	
	Vtac	H	
	Vent. Brady	H	
	Extreme Tachy	H	
	Extreme Brady	H	
	R on T	M*	
	VT>2	M*	
	Couplet	M*	
	PVCs/min	M*	
	Bigeminy	M*	
	Trigeminy	M*	
Tachy	M*		

Measurement	Alarm messages	L	Cause and solution
	Brady	M*	
	Missed Beats	M*	
	Irr. Rhythm	M*	
	Vent. Rhythm	M*	
	Multif. PVC	M*	
	Nonsus. Vtac	M*	
	Pause	M*	
	PNP	M*	
	PNC	M*	The pacer appears abnormal. Check the pacer.
Resp	Resp Apnea	H	The respiration signal was so weak that the monitor cannot perform respiration analysis. Check the patient's condition and the Resp connections.
	Resp Artifact	H	The patient's heartbeat has interfered with his respiration. Check the patient's condition and the Resp connections.
SpO ₂	SpO ₂ Desat	H	The SpO ₂ or SpO ₂ b value has fallen below the desaturation alarm limit. Check the patient's condition and check if the alarm limit settings are correct.
	SpO ₂ b Desat		
	No Pulse	H	The pulse signal was so weak that the monitor cannot perform pulse analysis. Check the patient's condition, SpO ₂ sensor and measurement site.
CO ₂	CO ₂ Apnea	H	The patient stops breathing, or the respiration signal was so weak that the monitor cannot perform respiration analysis. Check the patient's condition and the RM connections.
AG	AG Apnea	H	
RM	RM Apnea	H	
AG	FiO ₂ Too Low	H	Check the patient's condition, the ventilated O ₂ content and the AG connections.
tcGas	+tcpCO ₂ Alarm	M*	Parameter value has risen above the high alarm limit or fallen below the low alarm limit. Check the patient's condition and check if the patient category and alarm limit settings are correct.
	+tcpO ₂ Alarm	M*	
	+SpO ₂ Alarm	M*	
	+PR Alarm	M*	
NMT	TOF Alarm	M*	TOF value has risen above the high alarm limit or fallen below the low alarm limit. Check the patient's condition and check if the patient category and alarm limit settings are correct.

D.2 Technical Alarm Messages

Measurement	Alarm message	L	I	Cause and solution
XX	XX SelfTest Err	H	C	An error occurred to the XX module, or there is a problem with the communications between the module and the monitor. Re-plug the module and restart the monitor, or plug the module into another monitor.
	XX Init Err	H	A	
	XX Init Err N(N is between 1 to 8)	H	A	
	XX Comm Err	H	A	
	XX Comm Stop	H	C	

Measurement	Alarm message	L	I	Cause and solution
	XX Limit Err	L	C	XX parameter limit is accidentally changed. Contact your service personnel.
	XX Overrange	L	C	The measured XX value is not within the specified range for XX measurement. Contact your service personnel.
MPM	MPM 12V Err	H	C	An error occurred to the power supply part of the MPM module. Contact your service personnel.
	MPM 5V Err	H	C	
ECG	ECG Lead Off	L*	B	The electrode has become detached from the patient or the lead wire has become disconnected from the adapter cable. Check the connections of the electrodes and leadwires.
	ECG YY Lead Off	L*	B	
	Note: YY represents the leadwires, V (V1, V2, V3, V4, V5, V6,), LL, LA, RA, as per AHA standard, or C (C1, C2, C3, C4, C5, C6), F, L and R as per IEC standard.			
	ECG Noisy	L	A	The ECG signal is noisy. Check for any possible sources of signal noise around the cable and electrode, and check the patient for great motion.
	ECG Artifact	L	A	Artifacts are detected on the ECG analysis lead and as a result heart rate cannot be calculated and Asystole, Vfib and Vtac cannot be analyzed. Check the connections of the electrodes and leadwires and check for any possible source of interference around the cable and electrode. Check the patient's condition and check the patient for great motion.
	ECG High Freq. Noise	L	A	High frequency signals are detected on the ECG analysis lead. Check for any possible source of interference around the cable and electrode.
	ECG Low Freq. Noise	L	A	Low frequency signals are detected on the ECG analysis lead. Check for any possible source of interference around the cable and electrode.
	ECG Amplitude Too Small	L	C	The ECG amplitude didn't reach the detected threshold. Check for any possible source of interference around the cable and electrode.
	ECG Config. Err	L	C	ECG configuration is wrongly downloaded. Check the downloaded configuration and re-download the correct configuration.
Resp	Resp Disturbed	L	A	The respiration circuit is disturbed. Restart the monitor.
Temp	Temp Cal. Err	H	C	A calibration failed. Restart the monitor.
	T1 Sensor Off	L	A	The Temp sensor has become detached from the patient or the module. Check the sensor connections.
	T2 Sensor Off	L	A	
SpO ₂	SpO ₂ Sensor Off	L*	B	The SpO ₂ sensor has become detached from the patient or the module, or there is a fault with the SpO ₂ sensor, or an unspecified SpO ₂ sensor has been used. Check the sensor application site and the sensor type, and make sure if the sensor is damaged. Reconnect the sensor or use a new sensor.
	SpO ₂ b Sensor Off			
	SpO ₂ Sensor Fault	L	C	
	SpO ₂ b Sensor Fault			
	SpO ₂ No Sensor	L	B	
	SpO ₂ b No Sensor			
	SpO ₂ Unknown Sensor	L	C	

Measurement	Alarm message	L	I	Cause and solution
	SpO ₂ b Unknown Sensor			
	SpO ₂ Sensor Incompatible	L	C	
	SpO ₂ b Sensor Incompatible			
	SpO ₂ Too Much Light	L	C	There is too much light on the SpO ₂ sensor. Move the sensor to a place with lower level of ambient light or cover the sensor to minimize the ambient light.
	SpO ₂ b Too Much Light			
	SpO ₂ Low Signal	L	C	The SpO ₂ signal is too low or too weak. Check the patient's condition and change the sensor application site. If the error persists, replace the sensor.
	SpO ₂ b Low Signal			
	SpO ₂ Weak Signal			
	SpO ₂ No Pulse	L	C	The SpO ₂ sensor failed to obtain pulse signal. Check the patient's condition and change the sensor application site. If the error persists, replace the sensor.
	SpO ₂ b No Pulse			
	SpO ₂ Interference	L	C	The SpO ₂ signal has been interfered. Check for any possible sources of signal noise around the sensor and check the patient for great motion.
	SpO ₂ b Interference			
	SpO ₂ Comm abnormal	H	A	An error occurred to the SpO ₂ measurement module, or there is a problem with the communications between the module and the monitor. Re-plug the module and restart the monitor, or plug the module into another monitor.
	SpO ₂ b Comm abnormal			
	SpO ₂ Board Fault	L	C	There is a problem with the SpO ₂ measurement board. Do not use the module and contact your service personnel.
SpO ₂ b Board Fault				
SpO ₂ b has been closed	H	C	Different types of SpO ₂ measurement modules are applied. Use the same type of SpO ₂ measurement modules.	
NIBP	NIBP Loose Cuff	L	A	The NIBP cuff is not properly connected, or there is a leak in the airway.
	NIBP Air Leak	L	A	
	NIBP Pneumatic Leak	L	A	Check the NIBP cuff and pump for leakages.
	NIBP Cuff Type Wrong	L	A	The cuff type applied mismatches the patient category. Verify the patient category and replace the cuff.
	NIBP Air Pressure Err	L	A	An error occurred to the air pressure. Verify that the monitor application site meets the environmental requirements and check if there is any source that affects the air pressure.
	NIBP Weak Signal	L	A	The patient's pulse is weak or the cuff is loose. Check the patient's condition and change the cuff application site. If the error persists, replace the cuff.
	NIBP Signal Saturated	L	A	The NIBP signal is saturated due to excess motion or other sources.
	NIBP Overrange	L	A	The measured NIBP value is not within the specified range.
	NIBP Excessive Motion	L	A	Check the patient's condition and reduce the patient motion.
	NIBP Cuff Overpress.	L	A	The NIBP airway may be occluded. Check the airway and measure again.
NIBP Equip Err	H	A	An error occurred during NIBP measurement and therefore	

Measurement	Alarm message	L	I	Cause and solution
	NIBP Timeout	L	A	the monitor cannot perform analysis correctly. Check the patient's condition and NIBP connections, or replace the cuff.
	NIBP Measure Failed	L	A	
	NIBP Illegally Reset	L	A	An illegal reset occurred during NIBP measurement. Check if the airway is occluded.
IBP	YY Sensor Off			Check the sensor connection and reconnect the sensor.
	YY No Pulse	L	A	The catheter may be occluded. Please flush the catheter.
	YY represents an IBP label.			
C.O.	TB Sensor Off	L	A	Check the sensor connection and reconnect the sensor.
PiCCO	Invalid/Faulty PiCCO catheter	L	C	Erroneous or invalid catheter is used. Please use the proper catheter.
	TB Sensor Off	L	A	Check the sensor connections.
	PiCCO Comm Abnormal	H	A	Abnormal communication occurred between the PiCCO module and the system. Remove/connect the module again or restart the machine. If the problem remains, contact the Customer Services Dept. for help.
	PiCCO Comm Err	H	A	Erroneous communication occurred between the PiCCO module and the system. Remove/connect the module again or restart the machine. If the problem remains, contact the Customer Services Dept. for help.
	PiCCO Init Err	H	A	An error occurred to the module during the power-on self-test. Remove/connect the module again or restart the machine. If the problem remains, contact the Customer Services Dept. for help.
	Inject Temp. Sensor Err	L	C	An error occurred to the injectate temperature sensor or the sensor cable. Check/replace the sensor or the sensor cable.
	PiCCO Comm Stop	H	A	Remove/connect the module again or restart the machine. If the problem remains, contact the Customer Services Dept. for help.
ScvO ₂	Optical Module Err	L	C	Check the module connection. Change a module if necessary.
	ScvO ₂ Signal Too High	L	C	Check the sensor and reposition the catheter, then recalibrate the sensor.
	ScvO ₂ Signal Too Low	L	C	
	ScvO ₂ Too Much Light	L	C	Check and reposition the catheter, then recalibrate the sensor. Avoid the backlight which is excessively strong.
	Optical Module Disconnected	L	A	Connect the optical module.
	ScvO ₂ Comm Abnormal	H	A	Remove/connect the module again or restart the machine. If the problem remains, contact the Customer Services Dept. for help.
	ScvO ₂ Comm Err	H	A	

Measurement	Alarm message	L	I	Cause and solution
	ScvO ₂ Init Err	H	A	Remove/connect the module again. If the problem remains, contact the Customer Services Dept. for help.
	Unsupported CeVOX version	H	A	The module version is not compatible with the system. Please contact the Customer Services Dept. for help.
	ScvO ₂ Comm Stop	H	A	Remove/connect the module again or restart the machine. If the problem remains, contact the Customer Services Dept. for help.
CO ₂	CO ₂ Sensor High Temp	L	C	Check, stop using or replace the sensor.
	CO ₂ Sensor Low Temp	L	C	Check, stop using or replace the sensor.
	CO ₂ Temp Overrange	L	C	The operating temperature of the CO ₂ module goes beyond the specified range. After it restores within the specified range, the module will restart automatically.
	CO ₂ Airway High Press.	L	C	An error occurred in the airway pressure. Check the patient connection and patient circuit, and then restart the monitor.
	CO ₂ Airway Low Press.	L	C	
	CO ₂ High Barometric Press.	L	C	Check the CO ₂ connections, make sure that the monitor application site meets the requirements, and check for special sources that affect the ambient pressure. Restart the monitor.
	CO ₂ Low Barometric Press.	L	C	
	CO ₂ FilterLine Occluded	L	C	The airway or watertrap was occluded. Check the airway and remove the occlusion.
	CO ₂ No Watertrap	L	B	Check the watertrap connections.
	CO ₂ Check Adapter	L	A	There is a problem with the airway adapter. Check, clean or replace the adapter.
	CO ₂ FilterLine Err	L	C	Check if there is a leak in the CO ₂ sample line or the CO ₂ sample line has been occluded.
	CO ₂ Zero Failed	L	A	Check the CO ₂ connections. After the sensor's temperature becomes stabilized, perform a zero calibration again.
	CO ₂ System Err	L	A	Re-plug the module or restart the monitor.
	CO ₂ Check Cal.	L	C	Perform a calibration.
	CO ₂ Check Airway	L	C	An error occurred to the airway.
	CO ₂ No Filterline	L	A	Make sure that the filterline is connected.
	CO ₂ No Sensor	L	A	Make sure that the sensor is connected.
	CO ₂ Main Board Err	H	C	There is a problem with the CO ₂ module. Re-plug the module or restart the monitor.
	CO ₂ Checking Sensor	L	C	
	CO ₂ Replace Scrubber&Pump	L	C	
CO ₂ 15V Overrange	H	C		
CO ₂ Hardware Err	H	C		
tcGas	TCM Low Battery	M	C	Connect the TCM monitor with AC mains.
	TCM Battery Depleted	H	C	TCM monitor has less than 5 minutes running time on battery. Connect the TCM monitor with AC mains immediately.

Measurement	Alarm message	L	I	Cause and solution
	TCM Temperature Too High	H	C	The temperature in TCM CPU is too high. Please shut down TCM monitor immediately.
	TCM Alert	L	C	A TCM technical alarm is presented. Please check the TCM monitor to identify the cause of alarm.
AG	AG No Watertrap	L	B	Check the connections of the watertrap and re-connect it.
	AG Change Watertrap	L	A	Wait until the change is completed.
	AG Watertrap Type Wrong	L	A	Make sure that a correct watertrap has been used.
	O ₂ Accuracy Unspecified	L	A	The measured value has exceeded the specified accuracy range.
	N ₂ O Accuracy Unspecified	L	A	
	CO ₂ Accuracy Unspecified	L	A	
	Enf Accuracy Unspecified	L	A	
	Iso Accuracy Unspecified	L	A	
	Sev Accuracy Unspecified	L	A	
	Hal Accuracy Unspecified	L	A	
	Des Accuracy Unspecified	L	A	
	awRR Accuracy Unspecified	L	A	
	AG Hardware Err	H	A	
	AG Airway Occluded	L	A	Check the airway and remove the occlusion.
AG Zero Failed	L	A	Re-plug the module or restart the monitor, and then perform a zero calibration again.	
RM	RM No Sensor	L	A	Check and reconnect the sensor.
	RM Sensor Reversed	L	C	
	RM Zero Failed	L	C	Perform a zero calibration again.
	RM Power Err	L	A	There is a problem with the power supply. Re-plug the module or restart the monitor.
BIS	BIS High Imped.	L	A	Check and reconnect the BIS sensor.
	BIS Sensor Off	L	A	
	BIS DSC Err	L	C	An error occurred to the DSC during receiving signals. Check the DSC.
	BIS DSC Malf	L	C	The DSC automatically shuts down as a result of malfunction. Check the DSC.
	BIS No Cable	L	A	Check the BIS cables.
	BISx Disconnected	L	A	Check the BISx module.
	BIS No Sensor	L	A	Check the BIS sensor.
	BIS Wrong Sensor Type	L	A	Check or replace the sensor.
	BIS Sensor Too Many Uses	L	A	Replace the sensor.
	SQI<50%	L	A	The SQI value is too low. Check the patient's condition and the sensor connections.
	SQI<15%	L	A	
	BIS Sensor Expired	L	A	Replace the sensor.
	BIS Sensor Fault	L	C	Re-attach or Replace BIS Sensor
Disconnect/Reconnect BIS	L	C	Re-plug the BIS Module	
ICG	ICG Low Quality Signal	L	A	Check and reconnect the sensor.
	ICG Left Neck Sensor Off	L	A	

Measurement	Alarm message	L	I	Cause and solution
	ICG Right Neck Sensor Off	L	A	
	ICG L. Thorax Sensor Off	L	A	
	ICG R. Thorax Sensor Off	L	A	
	ICG Sensor Off	L	A	
NMT	TWSX Low Battery	M	C	Replace the battery.
	TWSX Battery Depleted	H	C	Replace the battery.
	TWSX No Acceleration Sensor	L	B	Connect the acceleration sensor.
	TWSX No Temp Sensor	L	B	Connect the temperature sensor.
	TWSX No Stimulation Cable	L	B	Connect the stimulation cable.
	TWSX Bad Electrode Connection	L	B	Reattach the electrode.
	TWSX Technical Alarm	L	C	An NMT technical alarm is presented. Please check the TOF-Watch® SX monitor to identify the cause of alarm.
Power	12V Too High	H	C	There is a problem with the system power supply. Restart the monitor.
	12V Too Low	H	C	
	5V Too High	H	C	
	5V Too Low	H	C	
	3.3V Too High	H	C	
	3.3V Too Low	H	C	
	Battery Too Low	H	C	Connect the monitor to an AC power source and allow the batteries to charge.
	Different Battery Voltages	M	C	The two batteries have different charge capacity, or the batteries unspecified have been used, or there is a problem with the batteries. Make sure that correct batteries are used and the batteries are not damaged, or replace the batteries.
	Battery Power Overload	H	C	The power consumption of the equipment is too high. Power the monitor using an AC power source.
RT Clock Not Exist	H	C	Contact your service personnel.	
Recorder	Recorder Init Err N	L	A	Restart the monitor.
	N is within 1 to 8.			
	Recorder SelfTest Err	L	A	Stop the recording and restart the monitor.
	Recorder Comm Err	L	A	
	Recorder S. Comm Err	L	A	
	Recorder Unavailable	L	A	
	Recorder Vlt High	L	C	An error occurred to the system power supply. Restart the monitor.
	Recorder Vlt Low	L	C	
	Recorder Head Hot	L	C	The recorder has been working for too long time. Stop the recording and resume the recording till the recorder's printhead cools down.
Rec Paper Wrong Pos.	L	A	Re-load the recorder paper.	
System	System Watchdog Err	H	C	An error occurred to the system. Restart the monitor.
	System Software Err	H	C	

Measurement	Alarm message	L	I	Cause and solution
	System CMOS Full	H	C	
	System CMOS Err	H	C	
	System FPGA Err	H	C	
	System Err N	H	C	
	N is within 2 to 12.			

FOR YOUR NOTES

E Electrical Safety Inspection

The following electrical safety tests are recommended as part of a comprehensive preventive maintenance program. They are a proven means of detecting abnormalities that, if undetected, could prove dangerous to either the patient or the operator. Additional tests may be required according to local regulations.

All tests can be performed using commercially available safety analyzer test equipment. These procedures assume the use of a 601PROXL International Safety Analyzer or equivalent safety analyzer. Other popular testers complying with IEC 60601-1 used in Europe, such as Fluke, Metron, or Gerb, may require modifications to the procedure. Please follow the instructions of the analyzer manufacturer.

The electrical safety inspection should be periodically performed every two years. The safety analyzer also proves to be an excellent troubleshooting tool to detect abnormalities of line voltage and grounding, as well as total current loads.

E.1 Power Cord Plug

Test Item		Acceptance Criteria
The power plug	The power plug pins	No broken or bent pin. No discolored pins.
	The plug body	No physical damage to the plug body.
	The strain relief	No physical damage to the strain relief. No plug warmth for device in use.
	The power plug	No loose connections.
The power cord		No physical damage to the cord. No deterioration to the cord.
		For devices with detachable power cords, inspect the connection at the device.
		For devices with non-detachable power cords, inspect the strain relief at the device.

E.2 Device Enclosure and Accessories

E.2.1 Visual Inspection

Test Item	Acceptance Criteria
The enclosure and accessories	No physical damage to the enclosure and accessories.
	No physical damage to meters, switches, connectors, etc.
	No residue of fluid spillage (e.g., water, coffee, chemicals, etc.).
	No loose or missing parts (e.g., knobs, dials, terminals, etc.).

E.2.2 Contextual Inspection

Test Item	Acceptance Criteria
The enclosure and accessories	No unusual noises (e.g., a rattle inside the case).
	No unusual smells (e.g., burning or smoky smells, particularly from ventilation holes).
	No taped notes that may suggest device deficiencies or operator concerns.

E.3 Device Labelling

Check the labels provided by the manufacturer or the healthcare facilities are present and legible.

- Main unit label
- Integrated warning labels

E.4 Protective Earth Resistance

1. Plug the probes of the analyzer into the device's protective earth terminal and protective earth terminal of the AC power cord.
2. Test the earth resistance with a current of 25 A.
3. Verify the resistance is less than limits.

LIMITS

For all countries, $R = 0.2 \Omega$ Maximum

E.5 Earth Leakage Test

Run an Earth Leakage test on the device being tested before performing any other leakage tests.

The following outlet conditions apply when performing the Earth Leakage test:

- normal polarity(Normal Condition),
- reverse polarity(Normal Condition),
- normal polarity with open neutral(Single Fault Condition),
- reverse polarity with open neutral(Single Fault Condition)

LIMITS

For UL60601-1,

- ◆ 300 μA in Normal Condition
- ◆ 1000 μA in Single Fault Condition

For IEC60601-1,

- ◆ 500 μA in Normal Condition
- ◆ 1000 μA in Single Fault Condition

E.6 Patient Leakage Current

Patient leakage currents are measured between a selected applied part and mains earth. All measurements have a true RMS only

The following outlet conditions apply when performing the Patient Leakage Current test.

- normal polarity(Normal Condition);
- reverse polarity(Normal Condition),
- normal polarity with open neutral(Single Fault Condition);
- reverse polarity with open neutral(Single Fault Condition).
- normal polarity with open earth(Single Fault Condition);
- reverse polarity with open earth(Single Fault Condition).

LIMITS

For CF  applied parts

- ◆ 10 μA in Normal Condition
- ◆ 50 μA in Single Fault Condition

For BF  applied parts

- ◆ 100 μA in Normal Condition
- ◆ 500 μA in Single Fault Condition



E.7 Mains on Applied Part Leakage

The Mains on Applied Part test applies a test voltage, which is 110% of the mains voltage, through a limiting resistance, to selected applied part terminals. Current measurements are then taken between the selected applied part and earth. Measurements are taken with the test voltage (110% of mains) to applied parts in the normal and reverse polarity conditions

The following outlet conditions apply when performing the Mains on Applied Part test.

- Normal Polarity;
- Reversed Polarity

LIMITS

- For CF  applied parts: 50 μ A
- For BF  applied parts: 5000 μ A


E.8 Patient Auxiliary Current

Patient Auxiliary currents are measured between any selected Applied Part connector and the remaining Applied Part connectors. All measurements may have a true RMS only response.

The following outlet conditions apply when performing the Patient Auxiliary Current test.

- normal polarity(Normal Condition);
- reverse polarity(Normal Condition),
- normal polarity with open neutral(Single Fault Condition);
- reverse polarity with open neutral(Single Fault Condition).
- normal polarity with open earth(Single Fault Condition);
- reverse polarity with open earth(Single Fault Condition).

LIMITS

For CF  applied parts,

- ◆ 10 μ A in Normal Condition
- ◆ 50 μ A in Single Fault Condition

For BF  applied parts,

- ◆ 100 μ A in Normal Condition
- ◆ 500 μ A in Single Fault Condition

NOTE

-
- **Make sure the safety analyzer is authorized comply with requirement of IEC61010-1.**
 - **Follow the instructions of the analyzer manufacturer.**
-

F Symbols and Abbreviations

F.1 Symbols

μA	microampere
μV	microvolt
μs	Microsecond
A	ampere
Ah	ampere hour
bpm	beat per minute
bps	bit per second
$^{\circ}\text{C}$	centigrade
cc	cubic centimeter
cm	centimeter
dB	decibel
DS	dyne second
$^{\circ}\text{F}$	fahrenheit
g	gram
GHz	gigahertz
GTT	gutta
h	hour
Hz	hertz
in	inch
kg	kilogram
kPa	kilopascal
L	litre
lb	pound
m	meter
mAh	milliampere hour
Mb	mega byte
mcg	microgram
mEq	milli-equivalents
mg	milligram
min	minute
ml	milliliter
mm	millimeter
mmHg	millimeters of mercury
cmH ₂ O	centimeters of water
ms	millisecond
mV	millivolt
mW	milliwatt

MΩ	megaohm
nm	nanometer
rpm	breath per minute
s	second
V	volt
VA	volt ampere
Ω	ohm
W	watt
-	minus, negative
%	percent
/	per; divide; or
+	plus
=	equal to
<	less than
>	greater than
≤	less than or equal to
≥	greater than or equal to
±	plus or minus
×	multiply

F.2 Abbreviations

AaDO ₂	alveolar-arterial oxygen gradient
AAMI	Association for Advancement of Medical Instrumentation
AC	alternating current
ACI	acceleration index
Adu	adult
AG	anaesthesia gas
AHA	American Heart Association
air cyl.	Air cylinder pressure
Air Flow	air flow
ANSI	American National Standard Institute
Ao	aortic pressure
Art	arterial
ATMP	Barometric pressure
aVF	left foot augmented lead
aVL	left arm augmented lead
aVR	right arm augmented lead
awRR	airway respiratory rate
BAP	brachial arterial pressure
Base Flow	base flow
BC	burst count
BIS	bispectral index

BP	blood pressure
BPSK	binary phase shift keying
BSA	body surface area
BT	blood temperature
BTPS	body temperature and pressure, saturated
C.I.	cardiac index
CCI	Continuous Cardiac Index
Cdyn	dynamic compliance
CCO	Continuous Cardiac Output
CaO ₂	arterial oxygen content
CCO	continuous cardiac output
CCU	cardiac (coronary) care unit
CE	Conformité Européenne
CFI	cardiac function index
CIS	Clinical Information System
CISPR	International Special Committee on Radio Interference
CMOS	complementary metal oxide semiconductor
CMS	central monitoring system
C.O.	cardiac output
CO ₂	carbon dioxide
COHb	carboxyhemoglobin
Compl	compliance
CP	cardiopulmonary
CPI	cardiac power index
CPO	Cardiac Power Output
Cstat	static compliance
CVP	central venous pressure
DC	direct current
Des	desflurane
Dia	diastolic
DPI	dot per inch
dPmx	left ventricular contractility
DVI	digital video interface
DO ₂	oxygen delivery
DO ₂ I	oxygen delivery index
ECG	electrocardiograph
EDV	end-diastolic volume
EE	Energy expenditure
EEC	European Economic Community
EEG	electroencephalogram
EMC	electromagnetic compatibility
EMG	electromyography
EMI	electromagnetic interference
Enf	enflurane
ESU	electrosurgical unit
Et	end-tidal

EtAA	End-tidal anesthetic agent
EtAA 2nd	2nd Exp. Agent
EtDes	
EtEnf	
EtHal	end-tidal anesthetic agent
EtIso	
EtSev	
EtCO ₂	end-tidal carbon dioxide
EtN ₂ O	end-tidal nitrous oxide
EtO	ethylene oxide
EtO ₂	end-tidal oxygen
ELWI	extravascular lung water index
EVLW	extravascular lung water
Exp%	expiration termination level
Exp. Flow	expiratory flow
Exp. MAC	Expired minimum alveolar concentration
f	breath rate
FAP	femoral arterial pressure
fapnea	breath rate for apnea ventilation
FCC	Federal Communication Commission
fCMV	CMV frequency
FDA	Food and Drug Administration
FEV1.0%	first second forced expiratory volume ratio
FG	Fresh gas flow
Fi	fraction of inspired
FiAA	Inspired anesthetic agent
FiAA 2nd	2nd Insp. Agent
FiDes	
FiEnf	
FiHal	inspired anesthetic agent
Filso	
FiSev	
FiCO ₂	fraction of inspired carbon dioxide
FiN ₂ O	fraction of inspired nitrous oxide
FiO ₂	fraction of inspired oxygen
Flow	flow
fmand	mandatory breathing frequency
FPGA	field programmable gate array
FRC	Fractional residual capacity
FreqMIN	minimum breath frequency
fsigh	sigh rate
fSIMV	frequency of SIMV
fspn	spontaneous breathing frequency
ftot	total breath rate
F-Trigger	inspiratory trigger level (flow trigger)
FV	flow-volume

GEDV	global end diastolic volume
GEDI	global end diastolic volume index
GEF	global ejection fraction
Hal	halothane
Hct	haematocrit
Hb	hemoglobin
Hb-CO	carbon mono-oxide hemoglobin
HbO ₂	oxyhemoglobin
HR	heart rate
I:E	inspiratory-expiratory ratio
IBP	invasive brood pressure
IBW	ideal body weight
ICG	impedance cardiography
ICP	intracranial pressure
ICT/B	intracranial catheter tip pressure transducer
ICU	intensive care unit
ID	identification
I:E	inspiratory time: Expiratory time ratio
IEC	International Electrotechnical Commission
IEEE	Institute of Electrical and Electronic Engineers
Ins	inspired minimum
Insp.Flow	inspiration flow
Insp. MAC	Inspired minimum alveolar concentration
△int.PEEP	intermittent PEEP
IP	internet protocol
Iso	isoflurane
IT	injectate temperature
ITBI	Intrathoracic Blood Volume Index
ITBV	Intrathoracic Blood Volume
LA	left arm
LAP	left atrial pressure
Lat	lateral
LCD	liquid crystal display
LCW	left cardiac work
LCWI	left cardiac work index
Leak Comp	leak compensation
LED	light emitting diode
LL	left leg
LVD	low voltage directive
LVDS	low voltage differential signal
LVET	left ventricular ejection time
LVSW	left ventricular stroke work
LVSWI	left ventricular stroke work index
MAC	minimum alveolar concentration
Art mean	mean arterial pressure
MDD	Medical Device Directive

MetHb	methemoglobin
%MinVol	Percentage of minute volume to be delivered
MRI	magnetic resonance imaging
MV	minute volume
MVe	expiratory minute volume
MVi	inspiratory minute volume
MVLEAK	leakage minute volume
MVspn	spontaneous breathed minute volume
N/A	not applied
N ₂	nitrogen
N ₂ O	nitrous oxide
N ₂ O cyl.	N ₂ O cylinder pressure
N ₂ O Flow	N ₂ O flow
Neo	neonate
NIBP	noninvasive blood pressure
NIF	negative inspiratory force
NMT	neuromuscular transmission
O ₂	oxygen
Δ O ₂	Difference between inspiratory and expiratory O ₂
O ₂ %	oxygen concentration
O ₂ CI	oxygen consumption index
O ₂ cyl.	Oxygen cylinder pressure
O ₂ cyl.2nd	Secondary oxygen cylinder pressure
O ₂ Flow	O ₂ flow
O ₂ R	oxygen extraction ratio
OR	operating room
oxyCRG	oxygen cardio-respirogram
PA	pulmonary artery
Pair	Air supply pressure
Papnea	apnea pressure
pArt-D	diastolic artery pressure
pArt-M	mean artery pressure
pArt-S	systolic artery pressure
Paux Mean	Mean auxiliary pressure
Paux Min	Minimum auxiliary pressure
Paux Peak	Peak auxiliary pressure
Paw	airway pressure
PAWP	pulmonary artery wedge pressure
PD	photodetector
Peak Flow	peak flow
Ped	pediatric
PEEP	positive end expiratory pressure
PEEP/CPAP	PEEP/CPAP
PEEPe	Extrinsic positive end-expiratory pressure
PEEPi	intrinsic positive end-expiratory pressure
PEEPi time	Intrinsic PEEP age (elapsed time since last maneuver)

PEEP _{tot}	total PEEP
PEF	peak expiratory flow
PEP	pre-ejection period
P _{high}	upper pressure level
PIF	peak inspiratory flow
P _{insp}	pressure control level of inspiration
PIP	peak inspiratory pressure
Pleth	plethysmogram
P _{limit}	pressure limit level
P _{low}	lower pressure level
P _{max}	maximum airway rressure
P _{mean}	mean pressure
PN ₂ O	N ₂ O supply pressure
PO ₂	Oxygen supply pressure
P _{peak}	peak pressure
P _{plat}	plateau pressure
PPV	Pulse Pressure Variation
PR	pulse rate
P _{supp}	pressure support level
PTC	post tetanic count
PTP	Pressure time product
P-Trigger	inspiratory trigger level (pressure trigger)
PVC	premature ventricular contraction
PVR	pulmonary vascular resistance
PVRI	pulmonary vascular resistance index
PVPI	pulmonary vascular permeability index
PPV	pulse pressure variation
p _{Art}	artery pressure
p _{CVP}	central venous pressure
P0.1	100 ms occlusion pressure
P0.1 time	P0.1 age (elapsed time since last maneuver)
R	right
RA	right arm
RAM	random access memory
Ramp	Ramp
RAP	right atrial pressure
RAW	airway resistance
RC _{exp}	Expiratory time constant
RC _{insp}	Inspiratory time constant
R _{dyn}	dynamic lung resistance
Re	expiratory resistance
Rec	record, recording
Resp	respiration
RHb	reduced hemoglobin
Ri	inspiratory resistance
Rise Time%	rise time

RL	right leg
RM	respiratory mechanics
RQ	Respiratory quotient
RR	respiration rate
RSBI	rapid shallow breathing index
Rstat	static lung resistance
SaO ₂	arterial oxygen saturation
SEF	spectral edge frequency
Sev	sevoflurane
SFM	self-maintenance
SI	stroke index
SMR	satellite module rack
SpO ₂	arterial oxygen saturation from pulse oximetry
SQI	signal quality index
SR	suppression ratio
STR	systolic time ratio
SV	stroke volume
SVI	Stroke Volume Index
SVR	systemic vascular resistance
SVRI	systemic vascular resistance index
SVV	stroke volume variation
SvO ₂	mixed venous oxygen saturation
ScvO ₂	central venous oxygen saturation
Sync	synchronization
Sys	systolic pressure
Tapnea	apnea interval
Taxil	axillary temperature
TB	Blood Temperature
TD	temperature difference
Temp	temperature
Texp	Expiratory time
TFC	thoracic fluid content
TFI	thoracic fluid index
TFT	thin-film technology
Thigh	time for the upper pressure level
Ti max	maximum inspiration time
Tinsp	time of inspiration
Tip	Inspiratory pause time
TIP:TI	percentage of inspiratory plateau time in inspiratory time
Tlow	time for the lower pressure level
TOF	train of four
Toral	oral temperature
TP	total power
Tplat	plateau time
TRC	Tube resistance compensation
Trect	rectal temperature

Trigger	trigger sensitivity
Trig Window	trigger window
Trise	rise time
Tslope	time for the pressure to rise to target pressure
Tube ID	tube ID
UAP	umbilical arterial pressure
UPS	uninterruptible power supply
USB	universal serial bus
UVP	umbilical venous pressure
VAC	volts alternating current
VCO ₂	CO ₂ production
VEPT	volume of electrically participating tissue
VI	velocity index
VO ₂	oxygen consumption
VO ₂ /kg	Oxygen consumption per body weight
VO ₂ /m ²	Oxygen consumption per body surface area
VO ₂ I	oxygen consumption index
VT _e /TV _e	expiratory tidal volume
VT _i /TV _i	inspiratory tidal volume
VT	tidal volume
VT _{apnea}	apnea tidal volume
VT _e spn	spontaneous expiratory tidal volume
VT _{sigh}	sigh tidal volume
WLAN	wireless local area network
WOB	work of breathing
WOB _{imp}	imposed work of breathing

FOR YOUR NOTE

